

915 KF Ti-Touch



Manual

8.915.8005EN / v10 / 2024-10-30



Metrohm AG
CH-9100 Herisau
Switzerland
+41 71 353 85 85
info@metrohm.com
www.metrohm.com

915 KF Ti-Touch

Program version 5.915.0046

Manual

This documentation is protected by copyright. All rights reserved.

This documentation is an original document.

This documentation has been prepared with great care. However, errors can never be entirely ruled out. Please send comments regarding possible errors to the address above.

Disclaimer

Deficiencies arising from circumstances that are not the responsibility of Metrohm, such as improper storage or improper use, etc., are expressly excluded from the warranty. Unauthorized modifications to the product (e.g. conversions or attachments) exclude any liability on the part of the manufacturer for resulting damage and its consequences. Instructions and notes in the Metrohm product documentation must be strictly followed. Otherwise, Metrohm's liability is excluded.

Melody for the BEEP command: excerpt from "En Altfrentsche", with kind permission of the Laseyer Quartett, Appenzell.

Table of contents

1	Introduction	1
1.1	Device description	1
1.1.1	Titration and measuring modes	1
1.1.2	Connectors	2
1.2	Displaying accessories	2
1.3	About the documentation	3
1.3.1	Symbols and conventions	3
2	Safety	5
2.1	Intended use	5
2.2	Responsibility of the operator	5
2.3	Requirements for operating personnel	6
2.4	Safety instructions	6
2.4.1	General notes on safety	6
2.4.2	Electrical safety	6
2.4.3	Tubing and capillary connections	7
2.4.4	Flammable solvents and chemicals	8
2.4.5	Recycling and disposal	8
3	Overview of the instrument	9
3.1	Front of the instrument	9
3.2	Rear of the instrument	10
4	Installation	11
4.1	Setting up the instrument	11
4.1.1	Packaging	11
4.1.2	Checks	11
4.1.3	Location	11
4.2	Connecting the power supply unit	11
4.3	Mounting the KF titration cell	13
4.4	Connecting MSB devices	14
4.4.1	Connecting a dosing device	15
4.4.2	Connecting an additional stirrer or titration stand	16
4.4.3	Connecting a Remote Box	16
4.5	Connecting USB devices	17
4.5.1	General	17
4.5.2	Connecting a USB hub	17
4.5.3	Connecting a printer	17
4.5.4	Connecting a balance	18

4.5.5	Connecting a PC keyboard	19
4.5.6	Connecting a barcode reader	19
4.5.7	Connecting a Sample Processor	20
4.5.8	Connecting the 885 Compact Oven SC	21
4.6	Setting up the titration vessel	22
4.7	Connecting sensors	22
4.7.1	General	22
4.7.2	Connecting a polarizable electrode	22
4.7.3	Connecting the temperature sensor	23
4.8	Connecting the Ti-Touch to a network	24
5	Titration	26
5.1	Water determination according to Karl Fischer (KFT)	26
6	Operation	27
6.1	Switching the instrument on and off	27
6.2	Fundamentals of operation	28
6.2.1	Touch-sensitive screen	28
6.2.2	Display elements and controls	29
6.2.3	Status display	30
6.2.4	Entering text and numbers	31
7	System settings	33
7.1	General system settings	33
7.1.1	Selecting the dialog language	33
7.1.2	Setting the date, time and local time	34
7.2	System-specific dialog options	35
7.3	User administration	38
7.3.1	Editing the user configuration	40
7.3.2	Creating an identification profile	43
7.3.3	Defining login options	44
7.3.4	Password options	46
7.3.5	Modification options	48
7.3.6	Reasons	48
7.3.7	Audit Trail	49
7.4	Measured value display	51
7.5	Acoustic signals	51
8	Titration	52
8.1	Adding a new titrant	53
8.2	Editing titrant data	54
8.3	Monitoring the working life	56

8.4	Dosing unit	57
8.4.1	Parameters for preparing (PREP) and emptying (EMPTY)	58
8.4.2	Tubing parameters	60
8.4.3	Shift direction of the valve disk	63
8.5	Exchange unit	63
8.5.1	Parameters for the preparation (PREP)	65
8.5.2	Tubing parameters	66
8.6	GLP test for exchange unit and dosing unit	67
8.7	Titer determination options and data	69
8.7.1	Titer validity	69
8.7.2	Properties of the previous titer determinations	70
9	Reagents	72
9.1	Editing reagent data	72
9.2	Reagent monitoring	73
10	Sensors	76
10.1	Adding a new sensor	77
10.2	Editing the sensor data	77
10.3	Monitoring the working life	78
11	Device manager	80
11.1	Adding a new device	81
11.2	Configuring the device	81
11.3	Ti-Touch	82
11.3.1	E-mail	83
11.3.2	PC/LIMS report	84
11.3.3	Shared memory	86
11.3.4	TCP/IP settings	87
11.4	Metrohm control instruments	88
11.4.1	Properties – Measuring input	89
11.4.2	Properties – MSB connector	91
11.4.3	Properties – Peripheral devices	92
11.5	Sample Processor	92
11.5.1	Properties – Sample Processor	93
11.5.2	Properties – Tower	94
11.5.3	Properties – Swing Head	95
11.6	Sample racks	100
11.6.1	Editing rack data	102
11.6.2	Rack adjustment	107
11.7	885 Compact Oven SC	108
11.7.1	Properties - 885 Compact Oven SC	109
11.7.2	Properties - Tower	110

11.7.3	Properties - Oven module	111
11.8	Printer	112
11.8.1	PDF settings	113
11.8.2	Network printer	114
11.8.3	More options	116
11.9	Balance	116
11.10	USB/RS-232 adapter	118
11.11	PC keyboard	120
11.12	Barcode reader	121
12	File manager	124
12.1	Managing files	124
12.1.1	Copying a file	126
12.1.2	Renaming a file	126
12.1.3	File properties	127
12.2	External storage medium	128
12.3	Creating backups / Restoring data	130
12.3.1	Restoring data	130
13	GLP manager	132
13.1	Automatic system test	133
13.2	Test tools	133
13.3	GLP tests for measurement and titration	134
13.3.1	Parameter description	135
13.4	System validation	138
13.4.1	Parameter description	139
13.5	System monitoring	142
13.5.1	Service interval	142
13.5.2	Backup interval	142
14	Common variables	144
14.1	Editing common variables	145
14.2	Properties of common variables	146
14.3	Monitoring validity	147
14.4	Assigning a result automatically to a common variable	148
15	Templates	150
15.1	Sample data	150
15.1.1	Sample identification list	151
15.1.2	Sample assignment table	152

15.2	Custom result templates	154
15.2.1	Editing result templates	155
15.3	Input lines	158
15.3.1	Editing the input signal	159
15.4	Output lines	160
15.4.1	Editing the output signal	162
15.5	Report header	163
16	Methods	165
16.1	Creating a new method	165
16.2	Saving a method	166
16.3	Loading a method	167
16.4	Editing a method	168
16.4.1	Inserting a command	169
16.5	Method options	170
16.5.1	Start options	171
16.5.2	Stop options	173
16.5.3	Sample data	174
16.5.4	Method properties	178
16.5.5	Note	179
16.5.6	Saving a determination automatically	180
17	Electronic signatures	182
17.1	Signing methods/determinations electronically	182
17.2	Deleting electronic signatures	183
18	Control	184
19	Favorites	188
19.1	Creating favorites	189
19.1.1	Editing favorites	189
20	Sample data	192
20.1	Entering sample data in the main dialog	192
20.2	Requesting sample data at the start of the determination	193
21	Sample table	195
21.1	General	195
21.2	Editing sample data	198
21.3	Properties	200

28 Manual control	248
28.1 Opening and closing the manual control	249
28.2 Measuring	250
28.2.1 Parameter description	251
28.3 Dosing	252
28.3.1 Continuous dosing	254
28.3.2 Dosing fixed volumes	256
28.3.3 Preparing	259
28.3.4 Emptying	260
28.3.5 Filling	260
28.3.6 Replacing reagent	260
28.4 Stirring	261
28.5 Remote	263
28.6 Sample Processor	264
28.6.1 Moving the lift	266
28.6.2 Moving to a rack position	268
28.6.3 External positions	269
28.7 885 Compact Oven SC	273
28.7.1 Controlling the heating/gas flow	275
29 Parameters	277
29.1 Volumetric Karl Fischer titrations (KFT)	277
29.1.1 Start conditions	277
29.1.2 Control parameters	279
29.1.3 Titration parameters	282
29.1.4 Stop conditions	283
29.1.5 Conditioning	284
29.1.6 Cell	286
29.1.7 Control instrument	286
29.1.8 Sensor	286
29.1.9 Dosing device	288
29.1.10 Stirrer	288
29.2 Measurements (MEAS)	289
29.2.1 Measuring parameters	289
29.2.2 Control instrument	292
29.2.3 Sensor	292
29.2.4 Stirrer	294
29.3 Evaluations (EVAL)	295
29.3.1 Fixed endpoint evaluation (EVAL FIX-EP)	295
29.3.2 Minimum and maximum evaluation (EVAL MIN/MAX)	296
29.3.3 Rate evaluation (EVAL RATE)	298
29.4 Calculations	299
29.4.1 Calculations (CALC)	299
29.4.2 Calculations (CALC LIVE)	306

X ■■■■■■

31.5	Manual control	353
31.6	Titration stand with pump	353
31.7	Miscellaneous	353
31.8	Volumetric Karl Fischer Titration	354
32	Appendix	356
32.1	Dosing unit	356
32.1.1	Maximum dosing and filling rate	356
32.1.2	Default parameters for preparing (PREP) and emptying (EMPTY)	356
32.2	Exchange unit	357
32.2.1	Maximum dosing and filling rate	357
32.2.2	Default parameters for preparing (PREP)	357
32.3	Stirring rate	357
32.4	Balance	358
32.5	Remote interface	359
32.5.1	Pin assignment of the remote interface	359
32.6	Result variables as parameter setting	362
32.7	Using AuditTrailViewer	363
32.7.1	Installing AuditTrailViewer	363
32.7.2	Opening the Audit Trail	365
32.7.3	Contents of the Audit Trail	366
32.7.4	Filtering the Audit Trail	367
32.7.5	Exporting the Audit Trail	368
32.8	Diagnosis	368
32.8.1	LCD test	369
32.8.2	Formatting an external storage medium	369
32.8.3	Removing an external storage medium	370
32.8.4	Adjusting the touch screen	370
32.8.5	Testing the touch screen	371
32.8.6	Software update (loading program versions and language files)	372
32.8.7	Service	375
33	Technical specifications	376
33.1	Touch screen	376
33.2	Measuring inputs	377
33.2.1	Polarizer	377
33.2.2	Temperature	377
33.3	Built-in magnetic stirrer	378
33.4	Pump	378
33.5	Interfaces	378

Table of figures

Figure 1	Front 915 KF Ti-Touch	9
Figure 2	915 KF Ti-Touch rear	10
Figure 3	Connecting the power supply unit	12
Figure 4	MSB connections	14
Figure 5	MSB connector	15
Figure 6	Connecting a polarizable electrode	23
Figure 7	Connecting the temperature sensor	24
Figure 8	Connecting the Ti-Touch to a network	25
Figure 9	Reagent dosing for KFT	26
Figure 10	Dosing unit – port assignment	62
Figure 11	Exchange unit – tubing connections	67
Figure 12	Directory structure on the external storage medium	129
Figure 13	Live display "Preparing the dosing unit"	259
Figure 14	Live display "Preparing the exchange unit"	260
Figure 15	Evaluation of minimum and maximum	297
Figure 16	Rotational speed depending on the stirring rate	358
Figure 17	Connectors of the Remote Box	359
Figure 18	Pin assignment of remote socket and plug	359

1 Introduction

1.1 Device description

The 915 KF Ti-Touch is a compact titration system for volumetric Karl Fischer titration. This titrator combines in a single device the touch-sensitive color monitor for convenient and efficient operation, the titration unit, the magnetic stirrer and the integrated membrane pump for adding and aspirating solvents. The upper side of the housing offers space for the titrant and the titration cell. The titrator is standard-equipped for operation with an external dosing drive of the *800 Dosino* type with a dosing unit. You can, however, also use a *805 Dosimat* with an exchange unit. Thanks to its compact construction, you can use the 915 KF Ti-Touch in a small space as a stand-alone titrator.

You manage titrants, sensors, methods, etc. conveniently in the internal memory of the 915 KF Ti-Touch. You can also save your files externally through the USB connector, e.g. on a USB flash drive. On this storage medium you can not only store your methods and determinations, but also create a backup together with all of the data and settings of your system.

The integrated Ethernet connector is available to you should you wish to connect your 915 KF Ti-Touch to a network. The network connection offers you the following advantages:

- Saving data to a computer within the network
- Printing reports on a network printer
- Sending displayed messages as e-mails

1.1.1 Titration and measuring modes

The 915 KF Ti-Touch supports the following titration and measuring modes.

- **KFT**
Volumetric water content determination according to Karl Fischer.
Measuring modes:
 - **Ipol** (voltametric measurement with selectable polarization current)
 - **Upol** (amperometric measurement with selectable polarization voltage)



- **MEAS**

Measuring modes:

- **Ipol** (voltametric measurement with selectable polarization current)
- **Upol** (amperometric measurement with selectable polarization voltage)
- **T** (temperature measurement)

1.1.2 Connectors

The 915 KF Ti-Touch is equipped with the following connectors:

- **Power connection**

For connecting to the power grid using the power supply unit provided.

- **Two MSB connectors (Metrohm Serial Bus)**

For connecting dosing devices, one additional stirrer or a Remote Box.

- **USB connector**

For connecting peripheral devices (printer, PC keyboard, etc.), a USB flash drive, a Sample Processor or a USB hub.

- **Sensor connectors**

One connector each for:

- Polarizable electrodes
- Temperature sensor (Pt1000 or NTC)

- Ethernet connector

For connecting the Ti-Touch to a network.

- **Two connection nipples**

For connecting tubing for aspirating solvent and extracting the contents of the titration cell.

1.2 Displaying accessories

Up-to-date information on the scope of delivery and on optional accessories can be found on the Metrohm website.

1 Searching for a product on the website

- Go to <https://www.metrohm.com>.
- Click on .
- Enter the article number of the product (e.g. **2.1001.0010**) into the search field and press **[Enter]**.

The search result is displayed.

2 Displaying product information

- To display the products matching the search term, click on **Product models**.
- Click on the desired product.

Detailed information regarding the product is displayed.

3 Displaying accessories and downloading the accessories list

- To display the accessories, scroll down to **Accessories and more**.
 - The **scope of delivery** is displayed.
 - Click on **[Optional parts]** for the optional accessories.
- To download the accessories list, click on **[Download accessories PDF]** under **Accessories and more**.



NOTE

Metrohm recommends keeping the accessories list for reference purposes.

1.3 About the documentation



CAUTION

Please read through this documentation carefully before putting the device into operation. The documentation contains information and warnings which the user must follow in order to ensure safe operation of the device.

1.3.1 Symbols and conventions

The following symbols and formatting may appear in this documentation:

(5- 12)	Cross-reference to figure legend The first number refers to the figure number, the second to the instrument part in the figure.
1	Instruction step Perform the steps one after the other.
Method	Dialog text, parameter in the software
File ► New	Menu or menu item
[Continue]	Button or key
	WARNING This symbol draws attention to a possible life-threatening hazard or risk of injury.



WARNING

This symbol draws attention to a possible hazard due to electrical current.



WARNING

This symbol draws attention to a possible hazard due to heat or hot instrument parts.



WARNING

This symbol draws attention to a possible biological hazard.



WARNING

Warning of optical radiation



CAUTION

This symbol draws attention to possible damage to instruments or instrument parts.



NOTICE

This symbol highlights additional information and tips.

2 Safety

2.1 Intended use

The 915 KF Ti-Touch is designed for usage as a titrator in analytical laboratories. Its application field is volumetric Karl Fischer titration.

This device is suitable for processing chemicals and flammable samples. Therefore, the use of the device requires the user to have basic knowledge and experience in handling toxic and caustic substances. Knowledge of the application of the fire prevention measures prescribed for laboratories is also mandatory.

2.2 Responsibility of the operator

The operator must ensure that basic regulations on occupational safety and accident prevention in chemical laboratories are observed. The operator has the following responsibilities:

- Instruct personnel in the safe handling of the product.
- Train personnel in the use of the product according to the user documentation (e.g. install, operate, clean, eliminate faults).
- Train staff on basic occupational safety and accident prevention regulations.
- Provide personal protective equipment (e.g. protective glasses, gloves).
- Provide suitable tools and equipment to carry out the work safely.

The product may be used only when it is in perfect condition. The following measures are required to ensure the safe operation of the product:

- Check the condition of the product before use.
- Remedy defects and malfunctions immediately.
- Maintain and clean the product regularly.

Requirements for operating personnel

Only qualified personnel may operate the product. Qualified personnel are persons who meet the following requirements:

- Basic regulations on occupational safety and accident prevention for chemical laboratories are known and complied with.
- Knowledge of handling hazardous chemicals is present. Personnel have the ability to recognize and avoid potential dangers.
- Knowledge regarding the application of fire prevention measures for laboratories is available.
- Safety-relevant information is communicated and understood. The personnel can operate the product safely.
- The user documentation has been read and understood. The personnel operate the product according to the instructions in the user documentation.

Safety instructions

General notes on safety



WARNING

Operate this instrument only according to the information contained in this documentation.

This instrument left the factory in a flawless state in terms of technical safety. To maintain this state and ensure non-hazardous operation of the instrument, the following instructions must be observed carefully.

Electrical safety

The electrical safety when working with the instrument is ensured as part of the international standard IEC 61010.



WARNING

Only personnel qualified by Metrohm are authorized to carry out service work on electronic components.

**WARNING**

Never open the housing of the instrument. The instrument could be damaged by this. There is also a risk of serious injury if live components are touched.

There are no parts inside the housing which can be serviced or replaced by the user.

Supply voltage**WARNING**

An incorrect supply voltage can damage the instrument.

Only operate this instrument with a supply voltage specified for it (see rear panel of the instrument).

Protection against electrostatic charges**WARNING**

Electronic components are sensitive to electrostatic charges and can be destroyed by discharges.

Do not fail to pull the power cord out of the power socket before you set up or disconnect electrical plug connections at the rear of the instrument.

2.4.3 Tubing and capillary connections**CAUTION**

Leaks in tubing and capillary connections are a safety risk. Tighten all connections well by hand. Avoid applying excessive force to tubing connections. Damaged tubing ends lead to leakage. Appropriate tools can be used to loosen connections.

Check the connections regularly for leakage. If the instrument is used mainly in unattended operation, then weekly inspections are mandatory.



2.4.4 Flammable solvents and chemicals




WARNING

All relevant safety measures are to be observed when working with flammable solvents and chemicals.

- Set up the instrument in a well-ventilated location (e.g. fume cupboard).
- Keep all sources of flame far from the workplace.
- Clean up spilled liquids and solids immediately.
- Follow the safety instructions of the chemical manufacturer.

2.4.5 Recycling and disposal

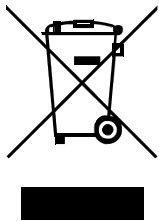

Properly dispose of chemicals and of the product to reduce negative effects on the environment and public health. Local authorities, waste disposal companies or dealers provide more detailed information on disposal. Observe the WEEE EU directive (WEEE = Waste Electrical and Electronic Equipment) for the proper disposal of waste electronic equipment within the European Union.

3 Overview of the instrument

3.1 Front of the instrument

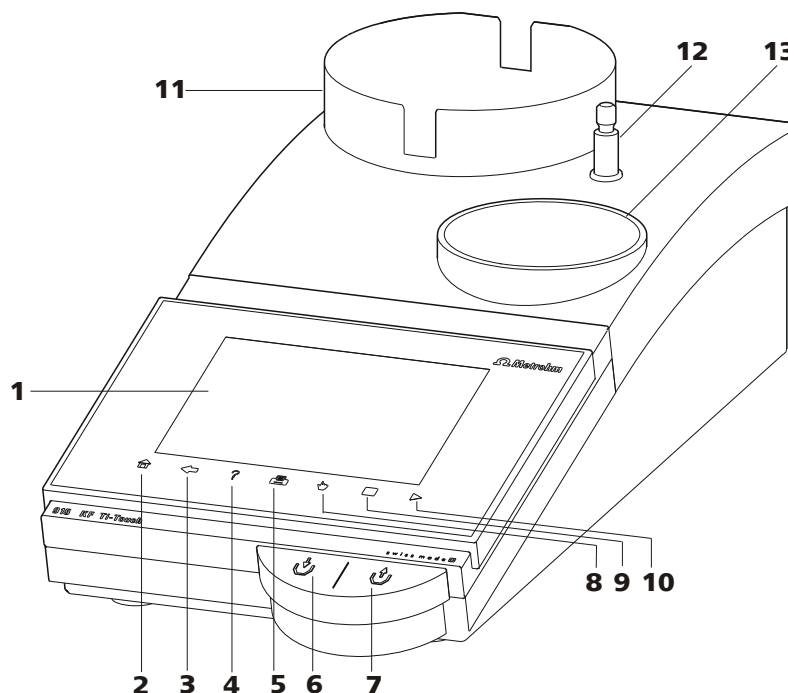


Figure 1 Front 915 KF Ti-Touch

1 Display

Touch screen.

3 Fixed key [Back]

Saves the entry and opens the next-higher dialog page.

5 Fixed key [Print]

Opens the print dialog.

7 Key

Pressing the key aspirates air out of the aspiration bottle. The vacuum in the aspiration bottle suctions the liquid out of the KF titration cell and into the aspiration bottle.

9 Fixed key [STOP]

Cancels the running determination.

2 Fixed key [Home]

Opens the main dialog.

4 Fixed key [Help]

Opens the online help for the dialog displayed.

6 Key

Pressing the key pumps air into the solvent bottle. The overpressure in the solvent bottle pushes solvent into the KF titration cell.

8 Fixed key [Manual]

Opens the manual control.

10 Fixed key [START]

Starts a determination.

- ## 11 Bottle holder

With holding clips, for reagent bottle.

- ### 13 Titration stand

With built-in magnetic stirrer and membrane pump for placement of the titration cell.

- ## 12 Support rod (lower part)

For mounting the support rod (upper part).

3.2 Rear of the instrument

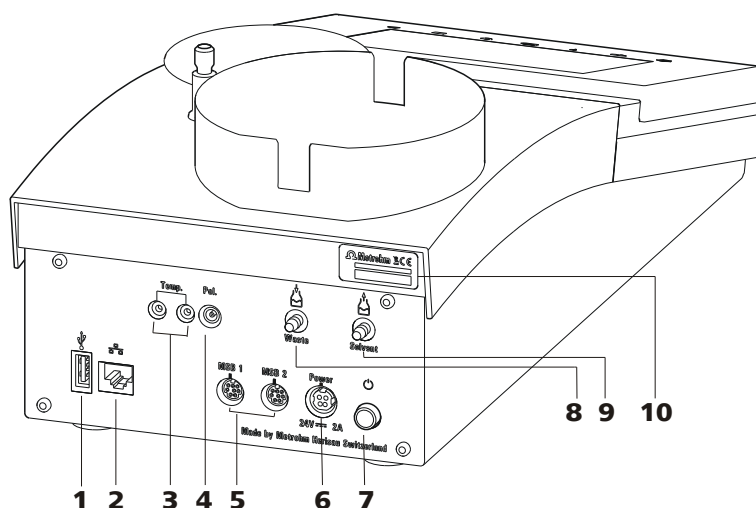


Figure 2 915 KF Ti-Touch rear

- 1 USB connector (type A)**

For connecting a printer, USB flash drive, USB hub, Sample Processor etc.

- ## 2 Ethernet connector (RJ-45)

For connecting to a network.

- ### 3 Temperature sensor connector (Temp.)

For connecting temperature sensors (Pt1000 or NTC). Two B sockets, 2 mm.

- #### 4 Electrode connector (Pol.)

For connecting polarizable electrodes, e.g. double Pt electrodes. Socket F.

- ## 5 MSB connector (MSB 1 and MSB 2)

Metrohm Serial Bus. For connecting external dosing devices, stirrers or a Remote Box.
Mini DIN, 8-pin.

- ## 6 Power socket (Power)

For connecting the external power supply unit.

- ## 7 Power switch

Switches the instrument on/off.

- ## 8 Connection nipple for PVC tubing

For aspirating the contents of the titration cell.

- ## 9 Connection nipple for PVC tubing

For aspirating solvent.

- ## 10 Type plate

Contains the serial number.

4 Installation

4.1 Setting up the instrument

4.1.1 Packaging

The instrument is supplied in protective packaging together with the separately packed accessories. Keep this packaging, as only this ensures safe transportation of the instrument.

4.1.2 Checks

Immediately after receipt, check whether the shipment has arrived complete and without damage by comparing it with the delivery note.

4.1.3 Location

The instrument has been developed for operation indoors and may not be used in explosive environments.

Place the instrument in a location of the laboratory which is suitable for operation and free of vibrations and which provides protection against corrosive atmosphere and contamination by chemicals.

The instrument should be protected against excessive temperature fluctuations and direct sunlight.

4.2 Connecting the power supply unit

The 915 KF Ti-Touch has an external power supply unit for a 24 V power supply (DC). This is connected to the power socket of the Ti-Touch.



WARNING

An incorrect supply voltage can damage the instrument.

Operate the instrument only with the supply voltage specified for it. Use the supplied power supply unit exclusively, see chapter *Technical specifications*.

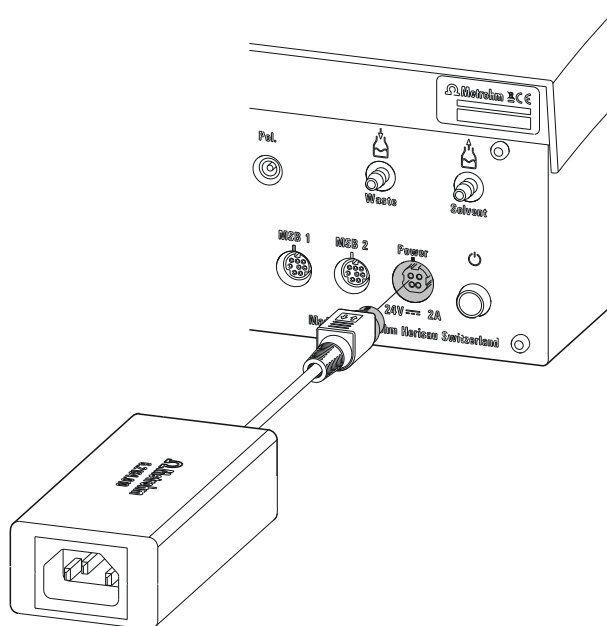


Figure 3 Connecting the power supply unit

Proceed as follows:

- 1 Connect the plug of the external power supply unit with the power socket of the Ti-Touch (see figure 3, page 12).



NOTE

The plug of the power supply unit is protected against accidental disconnection of the cable by means of a pull-out protection feature. If you wish to pull out the plug, you first need to pull back the outer plug sleeve marked with arrows.

- 2 Connect the power cord with the external power supply unit of the Ti-Touch and with the power grid.



CAUTION

Switch off the Ti-Touch correctly by pressing the power switch before you disconnect the electricity supply. If this is not done, then there is a danger of data loss.

4.3 Mounting the KF titration cell

The KF titration cell can be mounted with the aid of the upper part of KF titration cell.



NOTE

To prevent damage to the KF titration cell, the permissible height above the stirrer can be fixed with the aid of the clamping ring at the upper part of the KF titration cell.

Proceed as follows:

- 1** Push the clamping ring (6.2013.010) down as far as possible onto the mounted support rod (6.2016.050). Make sure while doing so that the indent is facing upwards.
- 2** Fix the 6.1414.030 vessel lid of the KF titration cell (with correctly inserted sealing ring from the 6.1244.040 sealing set) to the support rod.
- 3** Keep the locking lever pressed down and push the vessel lid down.
- 4** Fasten the 6.1415.220 (or 6.1415.250) KF titration cell with a 6.1903.020 (or 6.1903.030) stirring bar inside on the vessel lid. Fold back the holding bracket upwards while doing so. The markings on the vessel lid and on the plastic ring must be aligned above one another. Afterwards, press the holding bracket downwards in order to fix the KF titration cell. The levers of the holding bracket must enclose the pins of the plastic ring on the KF titration cell in order to ensure a secure hold.
- 5** Adjust the height of the KF titration cell by pressing the locking lever. It should almost touch the surface of the stirrer.
- 6** Readjust the clamping ring to the KF titration cell position and screw tight to fix the position of the KF titration cell. Make sure that the nose of the electrode holder is placed in the opening of the clamping ring.

If the height of the KF titration cell has been adjusted correctly, the entire cell can be raised and swiveled as required by pressing the locking lever.

4.4 Connecting MSB devices

In order to connect MSB devices, e.g. dosing device or Remote Box, the Ti-Touch has two connectors at what is referred to as the *Metrohm Serial Bus* (MSB). Various peripheral devices can be connected in sequence (daisy chain) at a single MSB connector (8-pin Mini DIN socket) and be controlled simultaneously by the Ti-Touch. In addition to the connection cable, stirrers and the Remote Box are each equipped with their own MSB socket for this purpose.



NOTE

The additional stirrer must not be connected to the MSB 1!

The following figure provides an overview of the instruments that can be connected to an MSB socket, along with a number of different cabling variations.

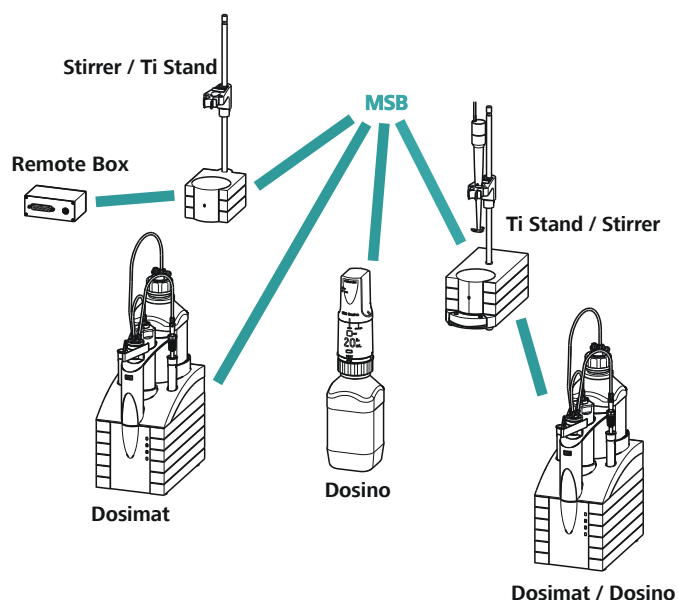


Figure 4 MSB connections



- No additional stirrer may be connected to the MSB 1!
- Only one device of the same type may be used at a single MSB connector at one time.
- When making the connection, take care to ensure that the flat part of the MSB plug marked with arrows is pointing in the direction of the marking on the MSB connector (*see figure 5, page 15*).

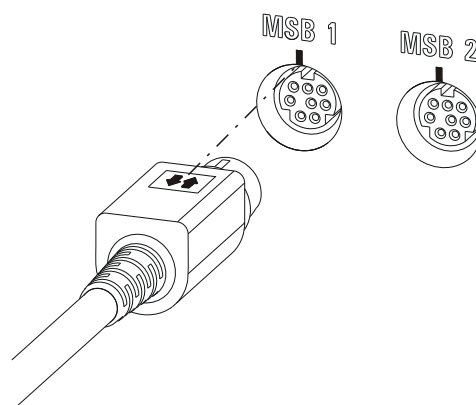


Figure 5 MSB connector



Switch off the Ti-Touch before you plug in MSB devices. When it is switched on, the Ti-Touch automatically recognizes which device is connected to which MSB connector. The connected MSB devices are entered automatically in the device manager.

MSB connections can be extended with the 6.2151.010 cable. The maximum connection length permitted is 6 m.

4.4.1 Connecting a dosing device

You can connect two dosing devices to the Ti-Touch.

The types of dosing devices that are supported are:

- 800 Dosino
- 805 Dosimat

Proceed as follows:

- 1** Switch off the Ti-Touch.

- 2 Connect the dosing device connection cable to an MSB connector (2-**5**) on the rear side of the Ti-Touch.
- 3 Switch on the Ti-Touch.

4.4.2 Connecting an additional stirrer or titration stand

In addition to the built-in magnetic stirrer, you can also use the magnetic stirrers 801 Stirrer, 803 Ti Stand or the 804 Ti Stand with the propeller stirrer 802 Stirrer.

Proceed as follows:

- 1 Switch off the Ti-Touch.
- 2 Connect the connection cable of the magnetic stirrer or of the titration stand to MSB 2 (2-**5**) on the rear of the Ti-Touch.
- 3 Switch on the Ti-Touch.

4.4.3 Connecting a Remote Box

Instruments that are controlled via remote lines and/or that send control signals via remote lines can be connected via the 6.2148.010 Remote Box. In addition to Metrohm, other instrument manufacturers also use similar connectors that make it possible to connect different instruments together. These interfaces are also frequently given the designations "TTL Logic", "I/O Control" or "Relay Control" and they generally have a signal level of 5 volts.

Control signals are understood to be electrical line statuses or brief (> 200 ms) electrical pulses which display the operating status of a device or which trigger or report an event. Sequences on a variety of instruments can thus be coordinated in a single complex automation system. However, no exchange of data is possible.

Proceed as follows:

- 1 Switch off the Ti-Touch.
- 2 Connect the Remote Box connection cable to an MSB connector (2-**5**) on the rear side of the Ti-Touch.
- 3 Switch on the Ti-Touch.

You can connect an 885 Compact Oven SC. The Remote Box also has an MSB socket at which a further MSB device, e.g. a dosing device, can be connected.

You will find precise information concerning the pin assignment of the interface on the Remote Box in Appendix (*see chapter 32.5, page 359*).

4.5 Connecting USB devices

4.5.1 General

The 915 KF Ti-Touch has a USB connector (type A socket) for peripheral devices with USB interface and for Sample Processors. If you wish to connect more than one device to the USB, you can use a commercially available USB hub.



NOTE

We recommend that the Ti-Touch be switched off while you set up or disconnect connections between the devices.

4.5.2 Connecting a USB hub

Use a USB hub with its own power supply.

Connect the USB hub as follows:

- 1 With the help of the 6.2151.030 cable (length 0.6 m) or the 6.2151.020 cable (length 1.8 m), connect the USB connector of the Ti-Touch (Type A) with the USB connector of the hub (Type B, see manual for the USB hub).

The USB hub is recognized automatically.

4.5.3 Connecting a printer

Printers connected to the 915 KF Ti-Touch must meet the following requirements:

- Printer languages: HP-PCL, HP-PCL-GUI, Canon BJL Commands or Epson ESC P/2
- Paper size: A4 or Letter.

Connect the printer as follows:

- 1 With the aid of the 6.2151.020 cable, connect the USB connector of the Ti-Touch (type A) with the USB connector of the printer (type B, see manual for the printer).

- 2 Configure the printer in the device manager of the Ti-Touch (see *chapter 11.8, page 112*).

4.5.4 Connecting a balance

If you want to connect a balance to the Ti-Touch you need a USB/RS-232 adapter (6.2148.050).

The following table offers an overview of the balances that you can use together with the Ti-Touch and of which cable you will need for connection to the RS-232 interface:

Balance	Cables
AND ER, FR, FX with RS-232 interface (OP-03)	6.2125.020 + 6.2125.010
Mettler AB, AG, PR (LC-RS9)	In the scope of delivery for the balance
Mettler AM, PM, PE with interface option 016 or Mettler AJ, PJ with interface option 018	6.2146.020 + 6.2125.010 Also from Mettler: ME 47473 adapter and either ME 42500 hand switch or ME 46278 foot switch
Mettler AT	6.2146.020 + 6.2125.010 Also from Mettler: ME 42500 hand switch or ME 46278 foot switch
Mettler AX, MX, UMX, PG, AB-S, PB-S, XP, XS	6.2134.120
Mettler AE with interface option 011 or 012	6.2125.020 + 6.2125.010 Also from Mettler: ME 42500 hand switch or ME 46278 foot switch
Ohaus Voyager, Explorer, Analytical Plus	Cable AS017-09 from Ohaus
Precisa balances with RS-232-C interface	6.2125.080 + 6.2125.010
Sartorius MP8, MC, LA, Genius, Cubis	6.2134.060
Shimadzu BX, BW	6.2125.080 + 6.2125.010

Connect the balance as follows:

- 1** Connect the USB plug of the USB/RS-232 adapter with the USB connector of the Ti-Touch (Type A).
The USB/RS-232 adapter will be recognized automatically and entered in the device manager of the Ti-Touch.
- 2** Connect the RS-232 interface of the USB/RS-232 adapter with the RS-232 interface of the balance (see table for cable).
- 3** Switch on the balance.
- 4** If necessary, switch on the RS-232 interface of the balance.
- 5** Configure the RS-232 interface of the USB/RS-232 adapter in the device manager of the Ti-Touch (*see chapter 11.10, page 118*).
- 6** Enter and configure the balance in the device manager of the Ti-Touch (*see chapter 11.9, page 116*).
- 7** Make sure that the parameters of the USB/RS-232 adapter configured in the device manager match those of the balance.

4.5.5 Connecting a PC keyboard

The PC keyboard is used as an aid for text and numerical input.

Connect the PC keyboard as follows:

- 1** Connect the USB plug of the keyboard with the USB connector of the Ti-Touch (Type A).
- 2** Enter and configure the keyboard in the device manager of the Ti-Touch (*see chapter 11.11, page 120*).

4.5.6 Connecting a barcode reader

The barcode reader is used as an aid for text and numerical input. You can connect a barcode reader with USB interface.

Connect the barcode reader as follows:

- 1** Connect the USB plug of the barcode reader with the USB connector of the Ti-Touch (Type A).

- 2 Enter and configure the barcode reader in the device manager of the (see chapter 11.12, page 121).

Settings on the barcode reader:

Program the barcode reader as follows (also see manual for the barcode reader):

- 1 Switch the barcode reader to programming mode.
- 2 Specify the desired layout for the keyboard (USA, Germany, France, Spain, German-speaking Switzerland).

This setting must match the setting in the device manager.
- 3 Make sure that the barcode reader is set in such a way that Ctrl characters (ASCII 00 to 31) can be sent.
- 4 Program the barcode reader in such a way that the ASCII character 02 (STX or Ctrl B) is sent as the first character. This first character is normally referred to as the "Preamble" or "Prefix Code".
- 5 Program the barcode reader in such a way that the ASCII character 04 (EOT or Ctrl D) is sent as the last character. This last character is normally referred to as the "Postamble", "Record Suffix" or "Postfix Code".
- 6 Exit the programming mode.

4.5.7 Connecting a Sample Processor

If you wish to integrate your Ti-Touch in an automation system, then you can connect the following Sample Processors to the USB connector:

- 814 USB Sample Processor
- 815 Robotic USB Sample Processor XL
- 810 Sample Processor (from firmware version 5.915.0040)

To connect a Sample Processor, you need the 6.2151.000 controller cable.

Connect the Sample Processor as follows:

- 1 Connect the Sample Processor to the power grid.
- 2 Connect the Sample Processor to the Ti-Touch with the controller cable.

The Sample Processor will be recognized automatically and entered in the device manager of the Ti-Touch.



NOTE

The plug on the controller cable is protected against accidental disconnection by means of a pull-out protection feature. If you wish to pull out the plug, you first need to pull back the outer plug sleeve marked with arrows.

- 3 Configure the Sample Processor in the device manager (*see chapter 11.5, page 92*).

4.5.8 Connecting the 885 Compact Oven SC



NOTE

In order for the 885 Compact Oven SC to be controlled with the 915 KF Ti-Touch, the following firmware versions must be installed on the instruments as a minimum requirement:

- 915 KF Ti-Touch 5.915.0045
- 885 Compact Oven SC 5.885.0013

The 885 Compact Oven SC is connected via the USB interface. The 6.2151.110 USB cable is required for this.

If the USB interface on the 915 KF Ti-Touch is also needed for other USB devices (USB flash drive, USB keyboard, USB/RS-232 adapter for balance), use a USB hub with its own power supply.

Connect the 885 Compact Oven SC as follows:

- 1 Switch off the 915 KF Ti-Touch.
- 2 Connect the 885 Compact Oven SC to the Ti-Touch using a USB cable.
- 3 Switch on the 885 Compact Oven SC.
- 4 Switch on the 915 KF Ti-Touch.



NOTE

When the 915 KF Ti-Touch is started, the 885 Compact Oven SC is added to the device manager and the settings from the device manager are transferred to the instrument.

The switch-on sequence must be strictly observed with the 885 Compact Oven SC.

- 5 Configure the 885 Compact Oven SC in the device manager (see *chapter 11.7, page 108*).

4.6 Setting up the titration vessel

The tutorial for the 915 KF Ti-Touch describes in detail how to set up the Karl Fischer titration cell and what you need to observe while doing so.

4.7 Connecting sensors

4.7.1 General

The measuring interface includes one measuring input (**Pol.**) for a polarizable electrode and one measuring input (**Temp.**) for a temperature sensor (Pt1000 or NTC).

4.7.2 Connecting a polarizable electrode

Connect the polarizable electrode as follows:

- 1** Plug the electrode plug into the **Pol.** socket of the Ti-Touch.

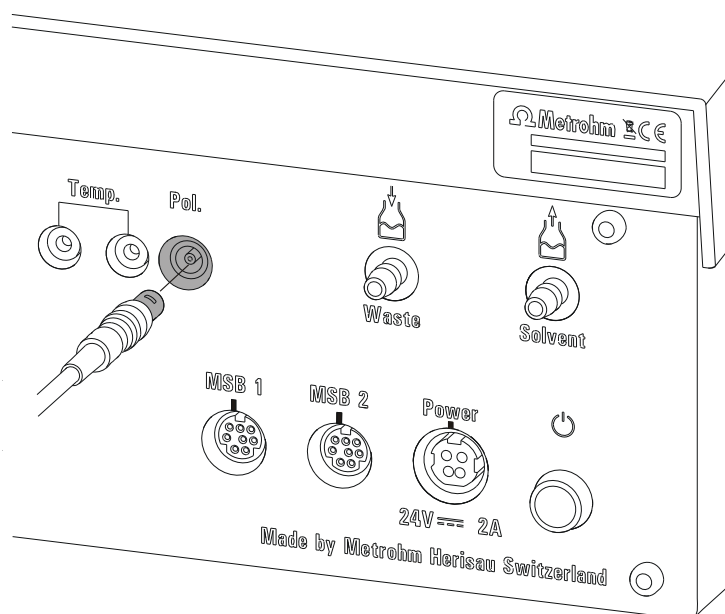


Figure 6 Connecting a polarizable electrode



NOTE

The electrode cable is protected against accidental disconnection of the cable by means of a pull-out protection. If you wish to pull out the plug again, you first need to pull back the outer plug sleeve.

4.7.3 Connecting the temperature sensor

A temperature sensor of the Pt1000 or NTC type can be connected to the **Temp.** connector.

Connect the temperature sensor as follows:

- 1 Plug the temperature sensor plugs into the **Temp.** sockets of the Ti-Touch.

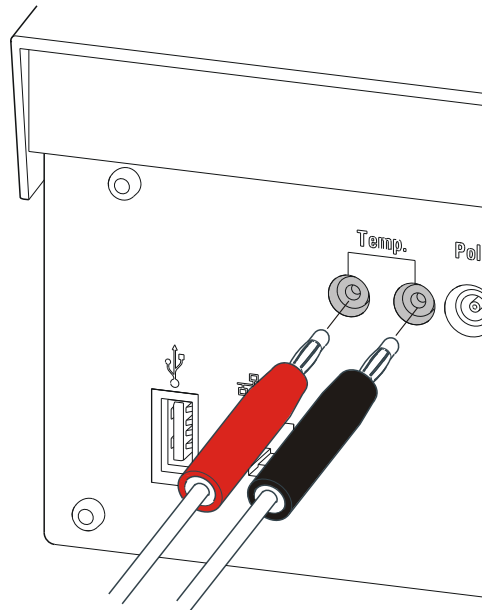


Figure 7 Connecting the temperature sensor




NOTE

The red plug must always be plugged into the red socket for the purpose of shielding against disruptions.

4.8 Connecting the Ti-Touch to a network

The 915 KF Ti-Touch has a network connection (Ethernet). This can be used to integrate your Ti-Touch in your network. You can, for example, store data on a PC within the network or print reports on a network printer. In *Chapter 11, page 80*, you will find information as to which settings are necessary for the network connection.

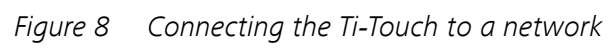

5 Titrations

5.1 Water determination according to Karl Fischer (KFT)

Karl Fischer Titration is a method for volumetric water determination. Conditioning is carried out automatically both before and after the actual titration. The reagent dosing is controlled in such a way that a predefined endpoint is reached as quickly and as accurately as possible. The volume steps and the rate of reagent dosing are regulated by the difference between the current measured value and the predefined endpoint. This means that titration is performed more slowly in the control range and that smaller volumes are added. The titration is stopped at the endpoint either drift-controlled or after a waiting time. The volume dosed until the endpoint is used to calculate the water content of the sample.

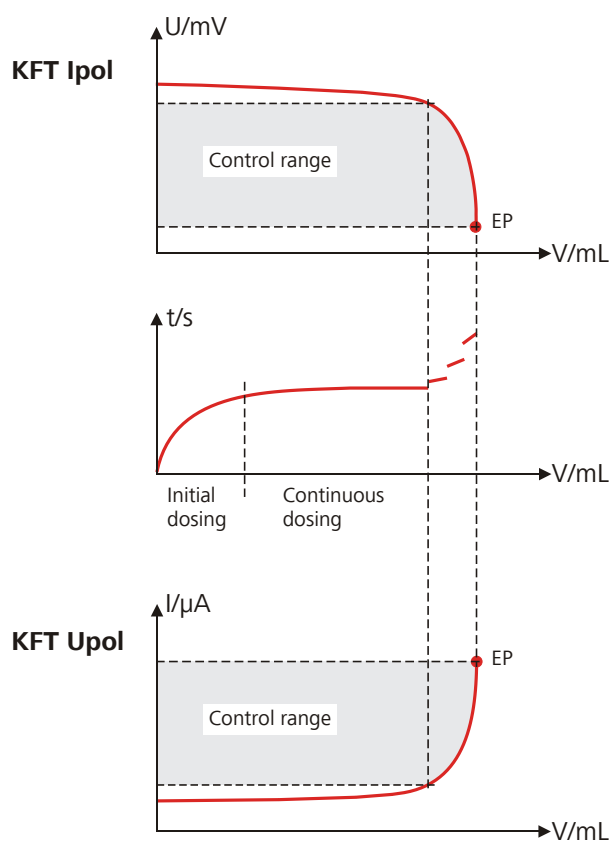


Figure 9 Reagent dosing for KFT

6 Operation

6.1 Switching the instrument on and off

Switching on the instrument



CAUTION

Peripheral devices (e.g., printers, 885 Compact Oven SC, etc) must be connected and switched on before you switch on the 915 KF Ti-Touch.



NOTE

English is set as the default dialog language when the instrument is switched on for the first time.


Proceed as follows:

- 1 ▪ Press the power switch on the left-hand side of the back panel of the 915 KF Ti-Touch.
The 915 KF Ti-Touch is initialized. A system test is performed. This process takes some time.



NOTE

If a buret unit is connected, then a request appears to carry out the **Prepare** function:

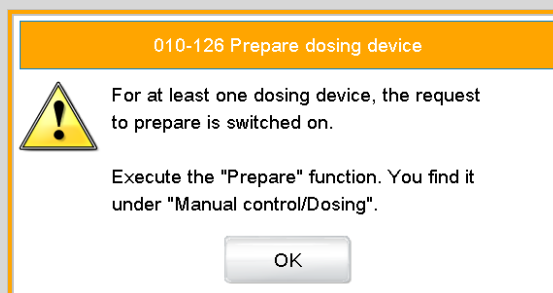

All tubings and the cylinder are rinsed with the **Prepare** function. The preparing of the buret unit is described in *chapter 28.3.3, page 259*.

- Confirm the message with **[OK]**.

The main dialog is displayed:

Switching off the instrument



CAUTION

The 915 KF Ti-Touch must be switched off by pressing the power switch on the rear of the instrument before the electricity supply is interrupted. If this is not done, then there is a danger of data loss.

Proceed as follows:

- 1 Press the power switch on the left-hand side of the back panel of the 915 KE Ti-Touch.

The current data is saved and the system is shut down. This process takes just a short time. At the same time, all other instruments connected to the 915 KF Ti-Touch via a USB cable are also being switched off (except the 885 Compact Oven SC).

6.2 Fundamentals of operation

6.2.1 Touch-sensitive screen

The entire 915 KF Ti-Touch user interface is touch-sensitive. Simply touch a few of the buttons on the interface to learn how a touch-sensitive screen reacts. You can always return to the main dialog by touching .

In order to enable an element on the 915 KF Ti-Touch user interface, just touch the screen with your fingertip, the eraser of a pencil or a stylus (special tool for operating instruments with touch-sensitive screens).




CAUTION

Never touch the touch screen with a pointed or sharp object such as a ballpoint pen.

In the default setting, the software is configured in such a way that an acoustic signal will be generated every time an active control is touched. This setting can be deactivated in the system settings (*see chapter 7.5, page 51*).

6.2.2 Display elements and controls

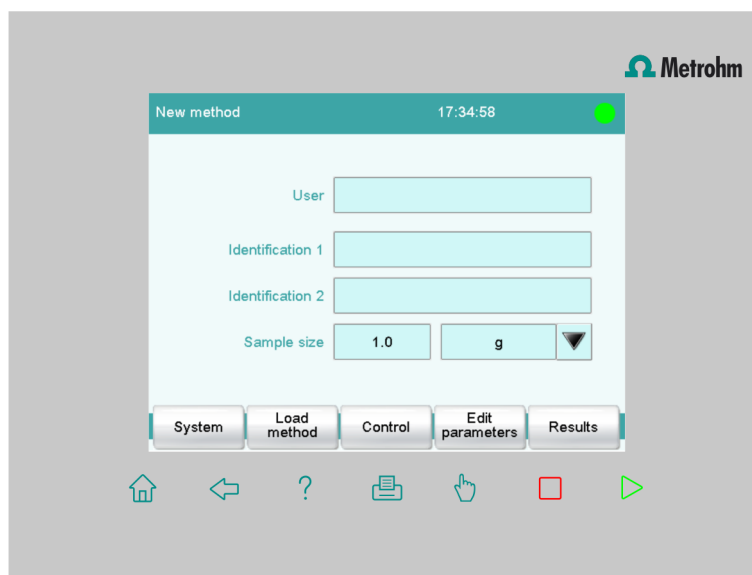

The following display elements and controls are available:

Table 1 Fixed keys which are always available

	[Home] always opens the main dialog.
	[Back] saves the entry and opens the next-higher dialog page.
	[Help] opens the online help for the dialog displayed.
	[Print] opens the printing dialog.
	[Manual] opens the manual control.
	[Stop] cancels the running determination.
	[Start] starts a determination.

The file name of the loaded method, the time and the system status are displayed in the main dialog in the **Title bar**.


In the other dialogs, the title bar shows the headings of the next upper level and of the displayed dialog. This is an aid for orientation during navigation through the user dialog.

Table 2 Screen elements

Buttons open a new dialog when they are tapped.



Inactive buttons with gray lettering indicate that the respective function is not available at the moment.

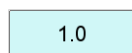

Input fields open an input dialog when tapped.



Tapping on the **selection symbol** opens a selection list.



A **check box** can also be activated or deactivated by tapping on it.

6.2.3 Status display

The current status of the system is displayed in the upper right-hand corner of the title bar.

Table 3 Status displays



The instrument is in normal status.



The working medium is being conditioned.



Conditioning has been paused.



The working medium is conditioned.



A method has been started.



A method has been paused.



An action has been started in manual control.

6.2.4 Entering text and numbers

In the editing dialog for text or numerical input, enter the individual characters by tapping in the input field. The following functions are available:

Text editor



Table 4 Editing functions

Editing function	Description
[OK]	The modification is applied and the editing dialog is exited.
[Cancel]	The editing dialog is exited without applying the modification.
[Delete entry]	The content of the input field is deleted completely.
[X]	The character in front of the cursor is deleted.
[⇐]	The cursor within the input field is shifted to the left by one character.
[⇒]	The cursor within the input field is shifted to the right by one character.
[a...z]	The lower-case letters are displayed. The label changes to [A...Z]. The upper-case letters are displayed again by tapping.
[0-9]	Numbers and mathematical characters are displayed.
[Special characters]	Special characters are displayed. You can use the button [More] to navigate through all available characters.

Number editor

Edit command / Stop conditions

Stop time

Input:
1 ... 999999

Default value:
off

7	8	9	off
4	5	6	R1 
1	2	3	
0	+/-	.	

Cancel Delete entry OK

Table 5 Editing functions

Editing function	Description
[OK]	The modification is applied and the editing dialog is exited.
[Cancel]	The editing dialog is exited without applying the modification.
[Delete entry]	The content of the input field is deleted completely.
[off]	If not only numbers but also special values (e.g., off) can be entered, then the corresponding buttons will be shown to the right of the numerical keypad.
[R1]	For many parameters, a result previously defined in the method can also be entered in place of a number (<i>see chapter 32.6, page 362</i>). You can select the result variable by touching [R1] .




NOTE

A commercially available USB keyboard can be connected to make text and numerical input easier. The key assignment is described in *chapter 11.11, page 120*.

7 System settings

Main dialog: **System ► System settings**

This chapter describes the various system settings and configurations.

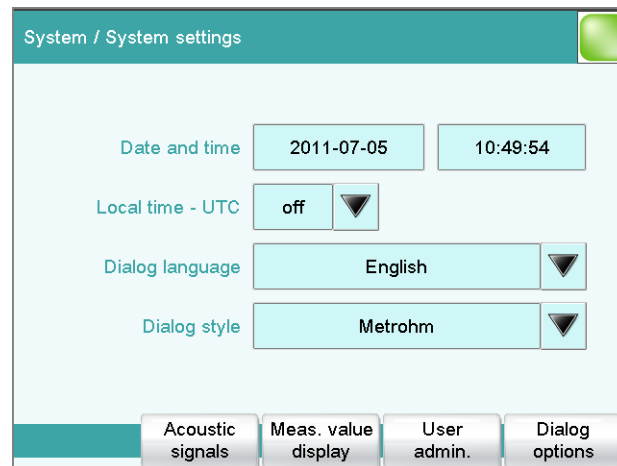

- Selecting the dialog language (see chapter 7.1.1, page 33).
- Setting the date and time (see chapter 7.1.2, page 34).
- Defining system-specific dialog options (see chapter 7.2, page 35).
- User administration (see chapter 7.3, page 38).
- Defining settings for the measured value display (see chapter 7.4, page 51).
- Configuring settings for acoustic signals (see chapter 7.5, page 51).

7.1 General system settings

7.1.1 Selecting the dialog language

The user interface is available in several languages. In addition to the two default dialog languages *English* and *German*, additional languages can be selected.

Proceed as follows to select the dialog language:

1 Opening the system settings

- In the main dialog, tap on **[System]**.
- Tap on **[System settings]**.

The **System / System settings** dialog is displayed.

2 Selecting the dialog language

- Tap on the **Dialog language** list box and select the desired language.

3 Saving the settings

Tap on the fixed keys [] or [].

The main dialog is displayed in the respective dialog language.

7.1.2 Setting the date, time and local time

The Ti-Touch displays the date and time in accordance with ISO standard 8601.

Proceed as follows to set the date and time:

1 Opening the system settings

- In the main dialog, tap on **[System]**.
- Tap on **[System settings]**.

The **System / System settings** dialog is displayed.

2 Entering the date

- Tap on the input field for the date.
The editor opens.
- Enter the current date in the format **YYYY-MM-DD** and confirm with **[OK]**.
The arrow keys **[⇐]** and **[⇒]** are used to move the cursor to the left or to the right by one character.

The input is saved and the editor is closed.

3 Entering the time

- Tap on the input field for the time.
The editor opens.
- Enter the current time in the format **hh:mm:ss** (24-hour format) and confirm with **[OK]**.
The arrow keys **[⇐]** and **[⇒]** are used to move the cursor to the left or to the right by one character.

The input is saved and the editor is closed.

- 4

Enter the local time
 - Tap on the list box **Local time - UTC** and select the difference from the UTC (Coordinated Universal Time).
The selection **off** means that the time is saved with no difference from the UTC.
- 5


Saving the settings

Tap on the fixed keys [↩] or [🏠].

The time settings are saved.

7.2 **System-specific dialog options**

Main dialog: **System ▶ System settings ▶ Dialog options**

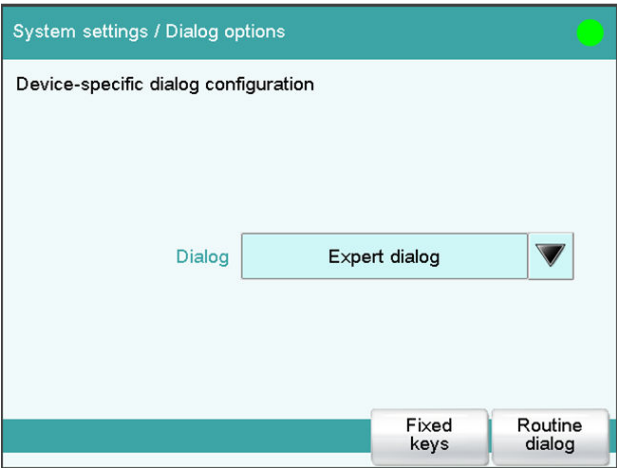

If you work without a login function, in this dialog you can define whether the system should generally be operated in expert mode or in routine mode. If you work with the login function activated, you must define this setting separately for each user (*see chapter 7.3.1, page 40*).

Dialog

Dialog mode in which the user may operate the system.

Selection	Expert dialog	Routine dialog
	Expert dialog All functions that are supported by the system are available.	Routine dialog The user dialog can be limited for routine operations. Only selected functions are available (<i>see Configuring the routine dialog, page 37</i>).

[Fixed keys]

Block unneeded fixed keys (see "Blocking unneeded fixed keys", page 36).

[Routine dialog]

Configure functions for the routine dialog (see "Configuring the routine dialog", page 37).

Selecting the dialog mode

Proceed as follows to change the dialog mode:

1 Select the dialog mode

Open the selection list **Dialog** and select either **Expert dialog** or **Routine dialog**.

2 Saving the settings

Tap on the fixed keys [] or [].

The setting will apply to all dialogs.



NOTE

If you have selected **Routine dialog** and if the routine dialog was configured in such a way that the dialog **System settings / Dialog options** is blocked, then you can switch back over to the expert dialog as follows:

- Operation without login function:
In the main dialog, enter **User = Metrohm**.
- Operation with login function:
A user who works with expert dialog must log in.

Blocking unneeded fixed keys

This following configurations apply for **both** dialog modes: routine dialog and expert dialog.

Blocking fixed keys

Proceed as follows to block unneeded fixed keys:

1 Display fixed keys which can be blocked

Tap on the **[Fixed keys]** button.

All fixed keys which can be blocked are displayed.

2 Deactivate fixed keys

Deactivate those fixed keys which are not permitted to be used.

3 Saving the settings

Tap on the fixed keys [↩] or [🏠].

Deactivated fixed keys cannot be used.

Configuring the routine dialog

A suitable **Standard configuration** has already been saved for routine operations:

- Methods can only be loaded, but not modified or created.
- Determinations cannot be recalculated.

You can readjust this standard configuration by disabling additional functions or re-enabling disabled functions.



NOTE


The configuration of the routine dialog applies for all routine users. You also have the option of defining user-specific routine settings. To do this, you must work with the login function activated and create an identification profile for each user (*see chapter 7.3.2, page 43*).

Proceed as follows to modify the configuration for the routine dialog:

1 Opening the dialog

Tap on the **[Routine dialog]** button.

The list of all buttons in the main dialog, in the manual control, etc. is displayed:

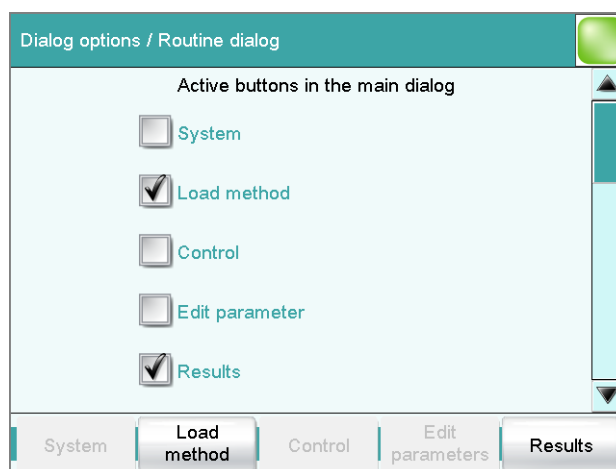

2 Deactivate buttons

Deactivate those buttons which are not permitted to be used.

All deactivated buttons will appear grayed-out, i.e. they are inactive.

3 Deactivate other functions

Many additional buttons and parameters can be disabled in the dialogs of **[System]**, **[Load method]**, **[Control]**, **[Edit parameters]** and **[Results]**. The corresponding option must be activated in order to enable these buttons.

4 Saving the settings

Tap on the fixed keys [] or [].

All deactivated functions will appear grayed-out, i.e. they are inactive.

7.3 User administration

Main dialog: **System** ► **System settings** ► **User admin.**

All of the functions of the user administration are described in this chapter. No user administration is mandatory for the simple operation of a titration system. If, however, you wish to be in compliance with *FDA Guidance 21 CFR Part 11*, then you must use the functions of the user administration. Additional information regarding the requirements of the FDA Guidance can be found in the document *Compliance Guide Touch Control*.



NOTE

If you work with the login function activated, then the user administration is accessible only for users with administrator rights. This means that you must ensure that at least two users have administrator rights so at least one of them will be available. Keep the access rights for a user with administrator rights in a safe place so that they are accessible in an emergency.

System settings / User administration		
User	Dialog	Status
Administrator	Expert dialog	active
ak	Routine dialog	active
ga	Routine dialog	active
jb	Expert dialog	active
wl	Expert dialog	active

Login options

Create ID profile

New

Delete

Edit

The following data is displayed in the user list for each user:

- Name
- Dialog mode in which the user may operate the system.
- Status

You can use the user list two different ways:

- Operation with login function:
If you work with the login function activated, i.e. if each user must log into the system before starting work, then only those users entered in the list can log in. The user who is currently logged in is shown in the main dialog.
- Operation without login function:
If you work without the login function, the users whose names are entered in the list can be selected in the main dialog or a user name can be entered. This makes it possible to document who has operated the titration system.

[Login options]

Define the settings for the login, Audit Trail, etc. (see chapter 7.3.3, page 44).

[Create ID profile]

Create an identification profile for the selected user on a storage medium (see chapter 7.3.2, page 43).

[New]

Adding a new user to the list (see chapter 7.3.1, page 40).

[Delete]

Delete the selected user from the list.



NOTE

Once users have worked with the login function and password protection, they can no longer be deleted, even if the password protection is disabled again. The status of these users must be set to **inactive** (requirement of FDA Guidance *21 CFR Part 11*).

The last user with administrator rights cannot be deleted.

[Edit]

Editing the data of the selected user (*see chapter 7.3.1, page 40*).

7.3.1 Editing the user configuration

User list: **User** ► **New** / **Edit**

User administration / Edit user

User

Full name

Dialog

Expert dialog

Status

active

Admin. rights

☐

Cancel

Favorites

Signature method

Signature determ.

User

The designation of the user is used for unambiguous identification, e.g. the company internal shorthand symbol or the personal number. The user name is printed out in all reports containing determination data and stored in the determination file. Each file always contains the name of the user who created it and the name of the last user to edit it.

Entry	24 characters maximum
Default value	empty

Full name

Complete name of the user.

Entry	24 characters maximum
Default value	empty

Dialog

Dialog mode in which the user may operate the system.

Selection	Expert dialog Routine dialog
Expert dialog	All functions that are supported by the system are available.
Routine dialog	The user dialog can be limited for routine operations. Only selected functions are available (<i>see Configuring the routine dialog, page 37</i>).

Status

Status of the user. Users can be deactivated. This function is useful, for instance, if the user is no longer authorized to operate the system or no longer works for the company.

Selection	active inactive
Default value	active
active	The user is authorized to operate the system.
inactive	The user is not authorized to operate the system and can no longer log in.

Admin. rights

on | off (Default value: **off**)

If this option is activated, then the user has administrator rights.



NOTE

The last user with administrator rights cannot be deleted anymore.

[Favorites]

Saving methods and sample tables as user-specific favorites (*see chapter 19, page 188*).

[Signature method]

Assign the permission to use and sign methods (*see "Edit user / Signature method" dialog, page 41*).

[Signature determ.]

Assign the permission to use and sign determinations (*see "Edit user / Signature determination" dialog, page 42*).

"Edit user / Signature method" dialog

User: **Edit ► Signature method**

Authorizations for using and signing methods can be defined in this dialog. These settings are only effective if you are working with the login function enabled and password protection.

Use only released methods

on | off (Default value: **off**)

If this option is activated, then the user may only start methods that have been released.

Review methods (signature level 1)

on | off (Default value: **off**)

If this option is activated, then the user may only sign methods which have the status **saved**. The method is assigned the status **reviewed**.

Release methods (signature level 2)

on | off (Default value: **off**)

If this option is activated, then the user may only sign methods which have the status **reviewed**. The method is assigned the status **released**.



NOTE

If the two options **Review methods** and **Release methods** are selected for a user, then this user can sign different methods at either level 1 or level 2, but not the same method at both level 1 and level 2.

Delete signatures

on | off (Default value: **off**)

If this option is activated, then the user may delete the signatures of a method which has been released. The method is assigned the status **saved**. The signatures of a method can only be deleted if it has the status **released**.

"Edit user / Signature determination" dialog

User: **Edit ► Signature determ.**

Authorizations for signing determinations can be defined in this dialog. These settings are only effective if you are working with the login function enabled and password protection.

Review determinations (signature level 1)

on | off (Default value: **off**)

If this option is activated, then the user may only sign determinations at the first level. The determination is assigned the status **reviewed**.

Release determinations (signature level 2)

on | off (Default value: **off**)

If this option is activated, then the user may only sign determinations at the second level. The determination is assigned the status **released**.

Delete signatures

on | off (Default value: **off**)

If this option is activated, then the user may delete the signatures of a determination which has been released. The signatures of a determination can only be deleted if it has the status **released**.

7.3.2 Creating an identification profile

User list: **User** ► **Create ID profile**

If you plan to carry out the login with an identification profile (*see chapter 7.3.3, page 44*), then you must first create an identification profile on a storage medium for each user. A check can then be made at the time of login as to whether or not the user does exist and whether or not he or she is working in the expert dialog or in the routine dialog. After a successful login the routine dialog settings stored on the card are loaded.



NOTE

In addition to the user name, the current routine dialog settings are also stored in this identification profile. This means you can define user-specific routine dialog settings for each user. However, you must configure them (see *"Configuring the routine dialog"*, page 37), **before** you create the identification profile.

Before you create the identification profile, check whether the desired routine dialog settings are active.

1 Plug in a storage medium

Plug in the USB storage medium.

2 Select a user

In the user list, select the user for whom the profile is to be created.

3 Create the identification profile

Tap on **[Create ID profile]**.

The user configuration and the current routine dialog settings are saved.

7.3.3 Defining login options

Main dialog: **System ► System settings ► User admin. ► Login options**

There are a number of different ways to log onto the system:

- Without login
- Login via user name
- Login via user name and password
- Login via identification profile
- Login via identification profile and password



NOTE


If you work with the login function activated, then the user administration is accessible only for users with administrator rights. This means that you must ensure that at least two users have administrator rights so at least one of them will be available. Keep the access rights for a user with administrator rights in a safe place so that they are accessible in an emergency.



NOTE

If you exit this dialog with [] or [], and if you have selected one of the login variants **Login via user name** or **Login via identification profile**, then the login dialog will open automatically and you must also log in to the system.

This means that you must make sure that you have first defined all the users and created the identification profiles before you activate the login function.

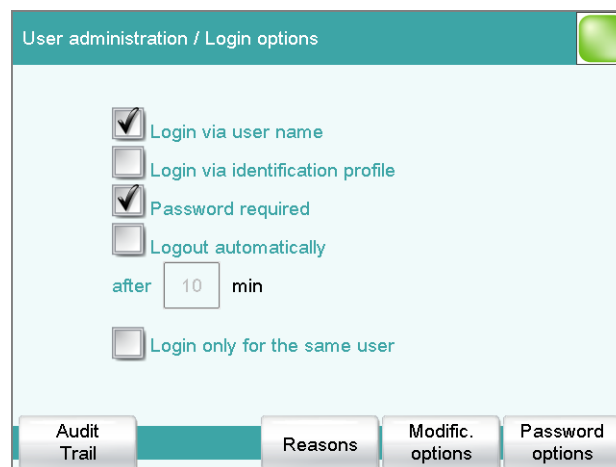

Password expires

on | off (Default value: **off**)

If this option is activated, then the user must define a new password after the time specified. A password that has already been used cannot be used again.

Input range	1 to 999 days
Default value	365 days

Forgotten password



CAUTION

If a user has forgotten his password, a new user name must be defined. The same user name can only be used again after a re-installation and re-creating the user list.

Proceed as follows:

1 Create a backup

Create a backup (see chapter 12.3, page 130).

2 Executing a system initialization

Switch off the Ti-Touch and execute a system initialization (see *chapter 30.1, page 349*).

3 Restore backed-up data

Use the function **Restore** to reload the data from the backup into your system (see chapter 12.3.1, page 130).

Deactivate the options **User list** and **System settings / User admin**.

4 Switch on the Ti-Touch

Switch the Ti-Touch back on again after a few seconds.

5 Restore the user list

Create the user list again and redefine the login options.

7.3.5 Modification options

Main dialog: **System ► System settings ► User admin. ► Login options ► Modific. options**

In the dialog **Login options / Modification options**, you can define for which actions a reason must be entered. These reasons are documented in the Audit Trail (*see chapter 7.3.7, page 49*) together with the modification. The reason for the last modification is shown in the properties of the method or determination.



NOTE

The reasons are only requested when working with activated login function and password.

Saving modified method

on | off (Default value: **off**)

If this option is activated, then a reason must be given when saving a method modification.

Recalculating determination

on | off (Default value: **off**)

If this option is activated, then a reason must be given when determinations are being recalculated.

7.3.6 Reasons

Main dialog: **System ► System settings ► User admin. ► Login options ► Reasons**

In the dialog **Login options / Reasons**, you can create a selection list containing reasons from which a selection can be made when signing and modifying methods and determinations. Some reasons have already been provided.



[Delete]

[Edit]

Change the designation of the selected reason.

Main dialog: **System** ► **System settings** ► **User admin.** ► **Login options** ► **Audit Trail**

You can also use the Audit Trail function specifically to record the data which are of interest to you.

You can view, filter and export the Audit Trail on a PC with the supplied software program *AuditTrailViewer*. Details concerning utilization of the *AuditTrailViewer* can be found in the Appendix (*chapter 32.7, page 363*).

If necessary, you can delete all of the entries in the Audit Trail with **[Delete Audit Trail]**. Do not fail to create a backup before doing so, however.

Security log

on | off (Default value: **off**)

If this option is activated, then the following events will be recorded: user login/logout, password changes, automatic deactivation of users and messages displayed during the login procedure.

User administration log

on | off (Default value: **off**)

If this option is activated, then all changes with respect to user administration will be recorded (modifying user list/user data, changing login options, etc.).

Method log

on | off (Default value: **off**)

If this option is activated, then all changes with respect to methods and determinations will be recorded (saving, deleting, renaming, copying and loading methods; editing method parameters; starting, stopping, pausing determinations).

Data log

on | off (Default value: **off**)

If this option is activated, then the following events will be recorded: changing sample data; settings in the **Control** dialog, recalculating determinations, changing statistics data, changing the titer/concentration of a titrant, changing the calibration data of a sensor, changing the value of a common variable.

System log

on | off (Default value: **off**)

If this option is activated, then the program start and the displayed messages of the following types will be recorded:

- : General warning messages
- : Error messages

- Data on exchange/dosing unit:
 - Parameters for the function **PREP**
 - Length and diameter of the tubings
 - Port assignment of the dosing unit
 - Cylinder volume
 - Serial number
 - etc.
- etc.

**NOTE**

If data is read out from the data chip, then a check is made whether the titrant list already contains a titrant of the same type with the identical serial number. If this is the case, then the older data set will **always** be overwritten by the new data set, no matter whether the data set in the titrant list or the data set on the data chip is the most recent one.

[New]

Adding a new titrant to the list (see chapter 8.1, page 53).

[Delete]

Delete the selected titrant from the list.

[Edit]

Editing the data of the selected titrant (see chapter 8.2, page 54).

8.1 Adding a new titrant

Before you can use a titrant, you must add it to the titrant list. To do this, use the button **[New]**.

- Exchange unit/dosing unit with data chip:
The exchange unit or dosing unit must be attached. All dosing devices on which non-configured exchange or dosing units have been detected are included in a selection list. Tapping on the button **[Edit]** opens the properties dialog, see following chapter.
- Exchange unit without data chip:
The properties dialog is opened after the dosing unit type has been selected, see following chapter.

8.2 Editing titrant data

Titrant list: **Titrants** ► **New / Edit**

Titration / Edit			
Titrant	Composite 5		▼
Concentration			▼
Comment			
Titer	4.955	mg/mL	▼
Date titer det. 2011-07-08 10:02:09			
Working unit		Dosing unit	Titer options

All of the data for the selected titrant is displayed in the dialog **Titrants / Edit**.

Titrant

The designation of the titrant is used for unambiguous identification.

Entry	24 characters maximum
Selection	Selection of frequently used titrants

Concentration

Concentration of the titrant.

Input range	–999999999 to 999999999
Default value	1.000

Unit of the concentration.

Entry	10 characters maximum
Selection	μmol/mL mmol/L mol/L g/L mg/L mg/mL μg/L ppm % mEq/L
Default value	mol/L

The concentration of the titrant is not relevant with Karl Fischer titration (KFT).

Comment

Entry	24 characters maximum
-------	-----------------------

Titer

Titer of the titrant.

Input range	–999999999 to 999999999
Default value	1.000
Unit of the titer.	
Entry	10 characters maximum
Default value	empty
Selection	μmol/mL mmol/L mol/L g/L mg/L mg/mL μg/L ppm % mEq/L

**NOTE**

If you modify the titer of the titrant in a loaded determination at a later date and would then like to recalculate the determination with the corrected value, then you must modify the value in the determination data under **View data / Titrant data** (see "Titrant data", page 215).

Date titer det.

Date and time of the last titer determination. For new titrants, the time the preparation was made is specified until after the first time a titer determination has been carried out.

[Working life]

Defining the working life of the titrant (see chapter 8.3, page 56).

[Dosing unit]

This button is only displayed for **Type = IDU**.

Defining the properties of the dosing unit used (see chapter 8.4, page 57).

[Exchange unit]

This button is only displayed for **Type = IEU** or **EU**.

Defining the properties of the exchange unit used (see chapter 8.5, page 63).

[Titer options]

Displaying the properties for titer determination (see chapter 8.7, page 69).

8.3 Monitoring the working life

Titrant: **Edit ► Working life**

In the dialog **Edit titrant / Working life**, you can define the time interval after which the titrant must be replaced. This is particularly important if your titrant has a limited working life. If you do not wish to monitor the working life, then you can enter only the date of manufacture for documentation purposes.

Preparation date

Date on which the reagent was manufactured or the bottle was opened.
For new titrants, the time the preparation was made will be specified.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the working life will be monitored.

Working life

If you define a time interval for the working life, then the **Expiry date** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Expiry date

If you define an expiry date, then the **Working life** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

8.4 Dosing unit

Titrant: **Edit ► Dosing unit**

Edit titrant / Dosing unit

Titrant: Titrant 5

Name

Order number 6.3032.150

Serial number 10170557

Cylinder volume 5 mL

Serial no. cyl.

Valve disk GLP test PREP param. Tubing param.

You can edit data for the dosing unit of the titrant in this dialog.

Name

Designation of the exchange or dosing unit.

Entry	24 characters maximum
Default value	empty

Order number

Order number of the exchange or dosing unit. It will be read out automatically on units with integrated data chips.

Entry	24 characters maximum
Default value	empty

Serial number

Serial number of the exchange or dosing unit. It will be read out automatically on units with integrated data chips.

Entry	8 digits maximum
-------	-------------------------

Cylinder volume

Cylinder volume of the dosing unit. It will be read out automatically on dosing units with integrated data chips.

Selection	2 5 10 20 50
Default value	20

Serial no. cyl.

Serial number of the dosing cylinder. It will be read out automatically on exchange or dosing units with integrated data chips. The number can be changed manually at any time, e.g. when a cylinder is replaced.

Entry	8 digits maximum
-------	-------------------------

[Valve disk]

Specifying the shift direction of the valve disk (see chapter 8.4.3, page 63).

[GLP test]

Defining the time interval for the GLP test (see chapter 8.6, page 67).

[PREP param.]

Entering the parameters for the preparation (see chapter 8.4.1, page 58).

[Tubing param.]

Entering the parameters for the connected tubing (see chapter 8.4.2, page 60).

8.4.1 Parameters for preparing (PREP) and emptying (EMPTY)

Titrant: **Edit** ▶ **Dosing unit** ▶ **PREP param.**

In the dialog **Dosing unit / PREP parameters**, you can adjust the parameters for the execution of the **Prepare** (command PREP) and **Empty** (command EMPTY) functions. The **Prepare** function is used to rinse the cylinder and tubing of the dosing unit and fill it air bubble-free. You should carry out this function before the first determination or once a day. The **EMPTY** function empties the cylinder and the tubings of the dosing unit.

Dosing port PREP/EMPTY

Dosing port through which the cylinder contents are ejected.

Selection	Dosing port 1 Dosing port 2 Fill port Special port
Default value	Dosing port 1

Dosing rate Dos. port 1

Rate used for the aspiration and ejection of the reagent via dosing port 1.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32.1, page 356*).

Dosing rate Dos. port 2

Rate used for the aspiration and ejection of the reagent via dosing port 2.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32.1, page 356*).

Dosing rate Fill port

Rate used for the aspiration and ejection of the reagent via fill port.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32.1, page 356*).

Dosing rate Spec.port

Rate used for the aspiration and ejection of the reagent via the special port.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32.1, page 356*).

Diameter

Input range	0.0 to 999.9 cm
Default value	0.0 cm

The setting **0.0** means that this tubing will neither be rinsed nor emptied.

Port

Port to be used as fill port for the **PREP** and **EMPTY** functions (see figure 10, page 62).

Selection	Port 1 Port 2 Port 3 Port 4
Default value	Port 2

Diameter

Input range	0.0 to 9.9 mm
Default value	2.0 mm

Port

Port to be used as special port for the **PREP** and **EMPTY** functions (see figure 10, page 62).

Selection	Port 1 Port 2 Port 3 Port 4
Default value	Port 4

Length

Length of the tubing.



Input range	0.0 to 999.9 cm
Default value	0.0 cm
The setting 0.0 means that this tubing will neither be rinsed nor emptied.	

Diameter

Diameter of the tubing.

Input range	0.0 to 9.9 mm
Default value	2.0 mm

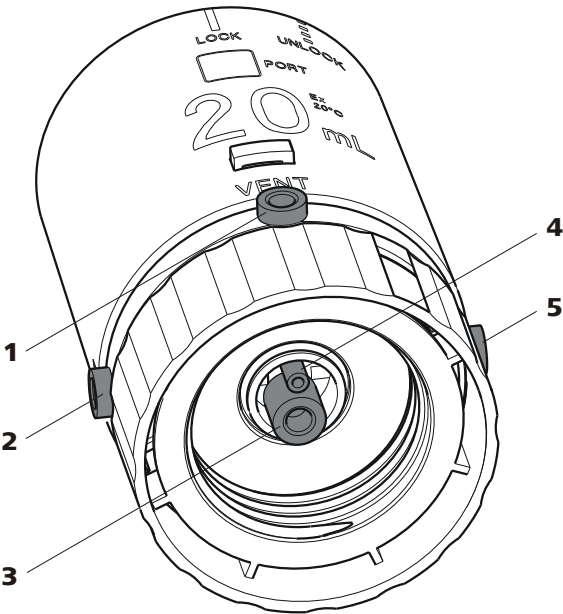


Figure 10 Dosing unit – port assignment

1 VENT This port is set up for the deaeration of the reagent bottle. An adsorber tube (filled with desiccant) is usually mounted here.	2 Port 1 The default definition of this port is dosing port 1.
3 Port 2 The default definition of this port is the filling port. A riser tube is usually mounted to it.	4 Port 4 Air is suctioned through this port during the Emptying function.
5 Port 3 The default definition of this port is dosing port 2.	

8.4.3 Shift direction of the valve disk

Titrant: **Edit** ▶ **Dosing unit** ▶ **Valve disk**

In this dialog, you can specify the shift direction of the valve disk.

Direction

Shift direction of the valve disc.

Selection	descending ascending automatic not over
Default value	automatic

descending

The ports are moved to in descending order.

ascending

The ports are moved to in ascending order.

automatic

The ports are moved to by the shortest path.

not over

Define a protected port.

Not over port

This parameter can only be edited with **Direction = not over**.

Define a protected port if the valve disc is not to be rotated over this port (useful with pipetting functions). The protected port can, however, be moved to directly.

Selection	1 2 3 4
Default value	4

8.5 Exchange unit

Titrant: **Edit** ► **Exchange unit**

Edit titrant / Exchange unit

Titrant: Composite 5

Name

Order number

6.3026.150

Serial number

10950001

Cylinder volume

5 mL

Serial no. cyl.

41424

GLP test

PREP param.

Tubing param.

You can edit data for the exchange unit of the titrant in this dialog.

Name

Designation of the exchange or dosing unit.

Entry	24 characters maximum
Default value	empty

Order number

Order number of the exchange or dosing unit. It will be read out automatically on units with integrated data chips.

Entry	24 characters maximum
Default value	empty

Serial number

Serial number of the exchange or dosing unit. It will be read out automatically on units with integrated data chips.

Entry	8 digits maximum
-------	------------------

Cylinder volume

Cylinder volume of the exchange unit. It will be read out automatically on exchange units with integrated data chips.

Selection	1 5 10 20 50
Default value	20

Serial no. cyl.

Serial number of the dosing cylinder. It will be read out automatically on exchange or dosing units with integrated data chips. The number can be changed manually at any time, e.g. when a cylinder is replaced.

Entry	8 digits maximum
-------	------------------

[GLP test]

Defining the time interval for the GLP test (see chapter 8.6, page 67).

[PREP param.]

Entering the parameters for the preparation (see chapter 8.5.1, page 65).

[Tubing param.]

Entering the parameters for the connected tubing (see chapter 8.5.2, page 66).

8.5.1 Parameters for the preparation (PREP)

Titrant: **Edit ► Exchange unit ► PREP param.**

In the dialog **Exchange unit / PREP parameters**, you can adjust the parameters for the execution of the **Prepare** function (command PREP). This function is used to rinse the cylinder and tubings of the exchange unit and fill it air bubble-free. You should carry out this function before the first determination or once a day.

Volume

Volume of titrant dosed during a rinsing cycle.

Input range	0.00000 to 99999.9 mL
Selection	Cylinder volume
Default value	Cylinder volume

Cylinder volume

The entire cylinder volume is being dosed.

Cycles

Number of rinsing cycles. We recommend carrying out at least two rinsing cycles in order to remove all air bubbles.

Selection	1 2 3 4 5 6 7 8 9
Default value	2

Dosing rate

The rate at which dosing takes place.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32, page 356*).

Filling rate

Rate at which the dosing cylinder is filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum filling rate depends on the cylinder volume (*see chapter 32, page 356*).

8.5.2 Tubing parameters

Titrant: **Edit** ▶ **Exchange unit** ▶ **Tubing param.**

You can enter the length and diameter of the connected tubings in the dialog **Exchange unit / Tubing parameters**. The values which have already been entered correspond to the dimensions of the supplied standard tubings.

Dosing tip

Tubing to the dosing tip (11-**2**).

Length

Length of the tubing.

Input range	0.0 to 999.9 cm
Default value	40.0 cm

Diameter

Diameter of the tubing.

Input range	0.0 to 9.9 mm
Default value	2.0 mm

Dosing cylinder

Tubing to the dosing cylinder (11-**3**).

Length

Length of the tubing.

Input range	0.0 to 999.9 cm
Default value	13.0 cm

Diameter

Diameter of the tubing.

Input range	0.0 to 9.9 mm
Default value	2.0 mm

Reagent bottle

Tubing to the reagent bottle (**11-1**).

Length

Length of the tubing.

Input range	0.0 to 999.9 cm
Default value	25.0 cm

Diameter

Diameter of the tubing.

Input range	0.0 to 9.9 mm
Default value	2.0 mm

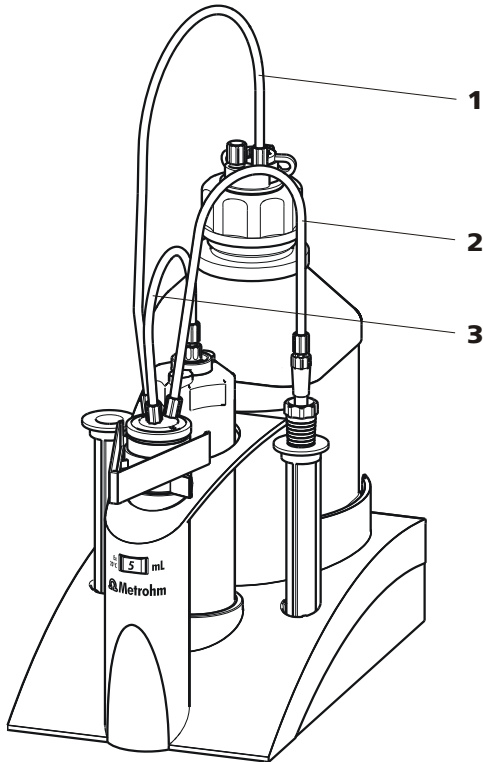


Figure 11 Exchange unit – tubing connections

- 1 Tubing to the reagent bottle
- 3 Tubing to the dosing cylinder

- 2 Tubing to the dosing tip

8.6 GLP test for exchange unit and dosing unit

Titrant: **Edit ► Dosing unit / Exchange unit ► GLP test**

In the dialog **Exchange unit / GLP test** or **Dosing unit / GLP test**, respectively, you can define the time interval after which a GLP test must be carried out again for the exchange unit or dosing unit.

GLP test date

Date on which the last GLP test was carried out.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the time interval after which a GLP test has to be carried out again will be monitored.

GLP test interval

If you define a time interval for the GLP test, then the date in **Next GLP test** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next GLP test

If you define a date for the next GLP test, then the **GLP test interval** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

8.7 Titer determination options and data

Titrant: **Edit ► Titer options**

Detailed information concerning the titer determination is displayed in the **Edit titrant / Titer options** dialog:

- **Titer method**
Method by which the titer was determined. If the titer was entered manually, then **manual** will be displayed.
- **User**
User who carried out the titer determination.
- **Statistical data**
The following information is also displayed for automatically assigned titers if the mean value of the results has been saved as the titer (*see "Save as titer", page 156*):
 - **n (titer det.)**
Number of titer determinations.
 - **s abs**
Absolute standard deviation
 - **s rel**
Relative standard deviation

[Validity]

Defining the time interval for the titer validity (*see chapter 8.7.1, page 69*).

[History]

Displaying information about the last ten titer determinations (*see chapter 8.7.2, page 70*).

8.7.1 Titer validity

Titrant: **Edit ► Titer options ► Validity**

In the dialog **Titer options / Validity**, you can define the time interval after which the titer must be determined again.

Date titer det.

Date and time of the last titer determination. For new titrants, the time the preparation was made is specified until after the first time a titer determination has been carried out.

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the titer validity will be monitored.

Titer validity

If you define a time interval for the validity of the titer, then the date in **Next titer determ.** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next titer determ.

If you define a date for the next titer determination, then the time interval for the **Titer validity** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

8.7.2 Properties of the previous titer determinations

Dialog "Titer options / History"

Titrant: **Edit** ▶ **Titer options** ▶ **History**

The date, time and titer of the last ten titer determinations are displayed in tabular form in the dialog **Titer options / History**. Titrers that were determined automatically will be displayed in green; manually entered titer values will be displayed in black with the designation **(m)**. You can delete these entries, e.g. if you have opened a new bottle.

**NOTE**

If multiple determinations have been carried out in order to determine the titer, then only one entry will be made in the history.

[Delete History]


Delete the entire history.

[Graph]

Open the diagram of the titer values, see following chapter.

Dialog "History / Graph"

Titrant: **Edit ► Titer options ► History ► Graph**

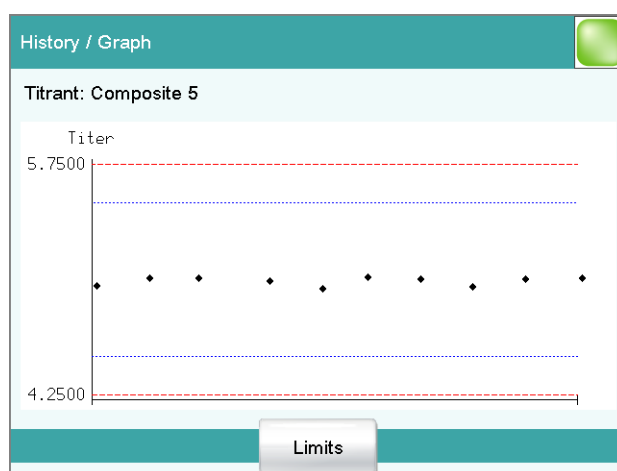

In this diagram, titer values are plotted against the date of the titer determination. You can define warning limits (blue dashed lines) and intervention limits (red dashed lines). These limits will not, however, be monitored.

[Limits]

Define warning and intervention limits.

9 Reagents

Main dialog: **System ► Reagents**

This chapter describes how you can create a list of the reagents used in the system.

System / Reagents	
Reagent	Type
Reagent A	Volumetric
Reagent B	Volumetric
Reagent C	Volumetric

New

Delete

Edit

The list of reagents gives the designation and type of each reagent.

[New]

Add a new reagent to the list, see following chapter.

[Delete]

Delete the selected reagent from the list.

[Edit]

Edit the data of the selected reagent, see following chapter.

9.1 Editing reagent data

Main dialog: **System** ► **Reagents** ► **New / Edit**

Reagent

The designation of the reagent is used for unambiguous identification.

Entry	24 characters maximum
-------	------------------------------

Comment

Entry	24 characters maximum
-------	------------------------------

[Reagent monitoring]

Set the parameters for the reagent monitoring, see following chapter.

9.2 Reagent monitoring

The conditions for the monitoring of the reagent are defined in the dialog **Edit reagent / Reagent monitoring**.

If one of the following values is reached, then the reagent must be replaced. The values are checked in the following cases:

- at the start of the determination.

- at the end of the determination.

Number of determ.

The number of determinations to be carried out with a certain amount of reagent depends on the type of sample and its amount.

Input range	1 to 999
Selection	off
Default value	off

Working life

Working life of the reagent.

Input range	1 to 999 days
Selection	off
Default value	off

Volume

Volume of titrant dosed.

Input range	1.0 to 999.9 mL
Selection	off
Default value	off

[Status]

Display the status overview of the current values of the reagent monitoring.

[Reagent replacement]

Edit the parameters for the reagent exchange.

"Reagent monitoring / Status" dialog

The current reagent monitoring values are displayed in this dialog.

[Reset]

Reset the values to zero.

Dialog "Reagent monitoring / Reagent replacement"

The parameters for the reagent replacement are defined in this dialog.

Reagent replacement

The reagent can either be exchanged manually or automatically.

Selection	manual auto
Default value	manual

manual

If a monitored parameter has reached the limit set, a message is being displayed. Then the reagent has to be exchanged manually.

auto

If a monitored parameter has reached the limit set, the method defined below is started automatically.

Memory

This parameter can only be edited with **Reagent replacement = auto**.

Memory location the method is loaded from. All memory locations are selectable, even if they are currently not accessible.

Selection	Internal memory External memory 1 External memory 2 Shared memory
Default value	Internal memory

Shared memory

Shared directory in the network.

Method

This parameter can only be edited with **Reagent replacement = auto**.

Method used for emptying the titration cell.

**NOTE**

Make sure that the memory is accessible.

Entry	32 characters maximum
Selection	Selection of stored methods

10 Sensors

Main dialog: **System ► Sensors**

This chapter describes how you can create a list of the sensors used in the system.

System / Sensors		
Sensor	Sensor type	Meas. input
Double Pt-wire electrode	Metal	
Metal electrode	Metal	
Temperature sensor	Temp.	

New

Delete

Edit

Two standard sensors are defined in the sensor list: **Metal electrode** and **Temperature sensor**. These sensors cannot be deleted or renamed. A maximum of 25 additional sensors can be added to these sensors.

For each sensor, the following data is displayed in the sensor list:

- Designation
- Type
 - **Metal:** Metal electrode
 - **Temp.:** Temperature sensor
 - **Other:** Other sensor

The following sensor data is stored in the list of sensors:

- Name
Each sensor in the system is identified by its unambiguous name.
- Working life
- etc.

[New]

Adding a new sensor to the list (see chapter 10.1, page 77).

[Delete]

Delete the selected sensor from the list.

[Edit]

Editing the data of the selected sensor (see chapter 10.2, page 77).

10.1 Adding a new sensor

Before you can use a sensor, you must add it to the sensor list. To do this, use the button **[New]**.

- The properties dialog is opened after the sensor type has been selected, see following chapter. The following sensor types can be selected:
 - Metal electrode (Pt electrode)
 - Other sensor
 - Temperature sensor

10.2 Editing the sensor data

Sensor list: **Sensor ► New / Edit**

All of the data for the selected sensor is displayed in the dialog **Sensors / Edit**.

Sensor

The designation of the sensor is used for unambiguous identification.

Entry	24 characters maximum
-------	------------------------------

Order number

Order number of the sensor.

Entry	24 characters maximum
Default value	empty

Serial number

Serial number of the sensor.

Display message

Document message

Cancel determination

The determination is stopped.

11 Device manager

Main dialog: **System ► Device manager**

This chapter describes how you can configure the Ti-Touch, a Sample Processor and the connected peripheral devices. A detailed description of the hardware installation of the Sample Processor can be found in the corresponding installation manual.

System / Device manager	
Device name	Device type
915 KF Ti-Touch	Ti-Touch operation
Ti-Touch	Ti-Touch titrator
Sample Processor	815 Sample Processor
Printer	Printer
Balance	Balance
USB/RS-232 adapter	USB/RS-232 adapter

New

Delete

Edit

The instrument list gives the name and type of every instrument. The **Ti-Touch titrator** is a component of the 915 KF Ti-Touch and therefore permanently entered as control instrument in the device manager. A connected Sample Processor with the peripheral devices (dosing devices, stirrers, etc.) connected to the MSB connector is entered in the list automatically. A USB/RS-232 adapter will also be automatically recognized and entered in the instrument list with default settings. A printer is entered in the instrument list in the default settings. You must enter a PC keyboard, a barcode reader or a balance in the instrument list yourself.

[New]

Add a new device to the list (see chapter 11.1, page 81).

[Delete]

Delete the selected device from the list.



NOTE

Devices which are recognized automatically cannot be deleted from the list while they are still connected.

[Edit]

Configure the selected device (see chapter 11.2, page 81).

11.1 Adding a new device

Of the following device types, you can enter one device each in the device list, even if it is not yet connected:

- Sample Processor
- Balance
- Barcode reader
- USB/RS-232 adapter
- PC keyboard

Proceed as follows:

1 Displaying the device selection

Tap on **[New]**.

2 Selecting the instrument

Tap on the button for the desired device.

The new device is entered in the list.

11.2 Configuring the device

Instrument list: **Instrument ► Edit**

The data stored for an instrument depends on the type of instrument. You can define an instrument name and a comment for each instrument. The instrument name of the Ti-Touch is printed out in the standard report header.

The description of the individual instruments can be found in the following specific chapters:

- Ti-Touch (*see chapter 11.3, page 82*)
- Metrohm control instruments (*see chapter 11.4, page 88*)
- Sample Processor (*see chapter 11.5, page 92*)
- Printer (*see chapter 11.8, page 112*)
- Balance (*see chapter 11.9, page 116*)
- PC keyboard (*see chapter 11.11, page 120*)
- USB/RS-232 adapter (*see chapter 11.10, page 118*)
- Barcode reader (*see chapter 11.12, page 121*)

11.3 Ti-Touch

Instrument list: **915 KF Ti-Touch** ► **Edit**

Device manager / Edit

Device type: Ti-Touch operation

Device name

915 KF Ti-Touch

Comment

Switch off display

off

min

Program version

5.915.0010

Serial number

395612

Control Remote Box

Not available

E-mail

PC/LIMS report

Shared memory

TCP/IP settings

More information

Device name

This designation is used for identification purposes when selecting control devices (command, manual control).

Entry	24 characters maximum
-------	-----------------------

Comment

Entry	24 characters maximum
-------	-----------------------

Switch off display

If this time interval has expired without the Ti-Touch having been operated, then the display will be switched off. Touching the display will switch it back on again at any time.

The function is no longer available with program version 5.915.0046 or higher.

Input range	1 to 999 min
Selection	off
Default value	off

Program version

Program version of the software.

Serial number

Shows the serial number of the device.

Control Remote Box

Shows to which control device and MSB connector the Remote Box is connected.

The "Control Remote Box" is the interface via which the system can be started and stopped externally. If multiple Remote Boxes are connected, then the one that is recognized first when the program starts will be used as the "Control Remote Box."

Selection	Name of the control device / Number of the MSB connector Not available
-----------	---

11.3.1 E-mail

915 KF Ti-Touch: **Edit** ► **E-mail**

The system allows you to send displayed messages as e-mails. The Ti-Touch must be connected to a network for this to function. The following types of messages can be sent:

- : general warning messages
- : error messages

Configuring e-mail dispatch

Proceed as follows so that messages can be sent as e-mails:

1 Activating the option

- In the instrument properties of the 915 KF Ti-Touch, tap on the **[E-mail]** button.
- In the **Edit device / E-mail** dialog, activate the **Send the following messages as e-mail:** option.

2 Configuring e-mail addresses

- Tap on the **[E-mail settings]** button.
The **E-mail / Settings** dialog is displayed.
- Enter the addresses of the mail server, the sender and the desired recipient.

Parameter description

Send the following messages as e-mail:

on | off (Default value: **off**)

If this parameter is activated, then messages with the following symbols will be sent as e-mails:

- : General warning messages

- : Error messages

Send only during running determination

on | off (Default value: **on**)

If this parameter is activated, then messages will be sent as e-mails only if a determination is running. Deactivate this parameter if messages are also to be sent in normal status.

Mail server

Address of the mail server for outgoing mail, e.g. mail.metrohm.ch. You can find the address of the mail server either in your e-mail program settings or obtain it from your IT department.

Entry	50 characters maximum
Default value	empty

Sender

E-mail address of the sender. This address must be formatted as an e-mail address, but need not necessarily correspond to an existing e-mail account, e.g. ti-touch@metrohm.com.

Entry	50 characters maximum
Default value	empty

Recipient

The messages will be sent to this e-mail address.

Entry	50 characters maximum
Default value	empty

11.3.2 PC/LIMS report

915 KF Ti-Touch: **Edit ► PC/LIMS report**

You can generate a machine-readable report with all of the important data concerning a determination, which is referred to as a PC/LIMS report. This report can be saved as a TXT file (according to ISO/IEC 8859-1) or as a UTF-8 file or sent to a terminal program or a LIMS via an RS-232 interface:

- manually with the  fixed key (see chapter 27, page 240).
- automatically at the end of a determination (see chapter 16.5.6, page 180).

The file name of the report file is constructed as follows: *PC_LIMS_Report-ID1-YYYYMMDD-hhmmss.txt*. A detailed description of the contents of the PC/LIMS report can be found in the *PC/LIMS Report Guide*.

Memory

Memory location where the PC/LIMS report is stored as a file. The report will be saved in the directory *pc_lims_report*. This directory will be created the first time a PC/LIMS report is generated. All three memory locations are listed as possible selections, even if they cannot be accessed at the moment.

Selection	off External memory 1 External memory 2 Shared memory
Default value	off

off

The report will not be saved as a file.

Shared memory

The report will be saved in a shared directory on the network. The shared directory is selected in the **Edit device / Shared memory** dialog (see chapter 11.3.3, page 86).

RS-232

The RS-232 interface via which the PC/LIMS report is sent. The interface parameters are adjusted in the **Edit device / Port parameters** dialog (see chapter 11.10, page 118).

Selection	off COM 1 COM 2
Default value	off

off

The report will not be sent via an RS-232 interface.

COM 2

This interface is inactive.

Coding

Format in which the PC/LIMS report is coded and stored.

Selection	ISO 8859-1 UTF-8
Default value	ISO 8859-1

ISO 8859-1

This format is recommended for all languages that use the extended ASCII code (e.g. German, English, Spanish, etc.).

UTF-8

This format is required for all languages that do not use the extended ASCII code (e.g. Russian, Chinese, Korean, etc.).

11.3.3 Shared memory

915 KF Ti-Touch: **Edit** ► **Shared memory**

If you have your Ti-Touch connected to your network, then you can define in this dialog a shared memory location on a computer within your network for the purpose of saving data (methods, determinations, etc.).



CAUTION

If the computer on which you share a memory location does not have the same subnet as the Ti-Touch, then a WINS server must be present. The computer must be entered in this server. In order for the communication to work, the SMB protocol version 1 (SMBv1) must be activated. Since this function is not standard anymore in Windows 10, it has to be activated manually. From firmware version 5.915.0044, the SMB protocol version 2 and 3 is also supported.

Computer

Host name of the computer on which a memory location is to be shared. If you are working with Windows, then you will find the host name of the computer as follows: In the Windows **command prompt** window, enter the command **ipconfig -all**. The host name is listed together with other parameters of the computer.

Do not under any circumstances enter an IP address in this input field.

Entry	max. 32 characters
Default value	empty

Share name

Share name of the shared memory location (file directory) on the above-specified computer. Please note that the Share name of a file directory often does not match the name of the file directory. The share name can be found on the Release tab in the Properties dialog of the shared file directory.

Entry	max. 32 characters
Default value	empty

Domain

Network domain in which the above-specified computer is located. If you are working with Windows, then you will find the domain name as follows: In the Windows **command window**, enter the command **ipconfig -all**. The domain is listed together with other parameters of the computer. Leave the field empty if the computer is not located within a domain.

User

Password

[Connect]

11.3.4 TCP/IP settings

If you have connected your Ti-Touch to your network, then you will have to define the network-relevant settings in this dialog. The Ti-Touch requires an IP address so that it is identifiable on the network. The Ti-Touch can acquire the IP address either dynamically from a DHCP server or you can enter the address directly.

on | off (Default value: **on**)

If this parameter is activated, then the Ti-Touch will obtain its IP address directly from a DHCP server. In this case, the remaining parameters can no longer be edited.

IP address for the Ti-Touch. IP addresses are 32-bit numbers and are written as sequences of four decimals, each separated by a period, e.g. "10.157.212.8".

915 KF Ti-Touch

Subnet mask

The net mask or subnet mask, together with the IP address, indicates to which network the device to be connected belongs. Subnet masks are 32-bit numbers and are written as sequences of four decimals, each separated by a period.

Entry	x.x.x.x
Default value	255.255.255.0

Default gateway

IP address for the standard gateway. A gateway sets up connections to several networks. It is located in the same subnet as the device to be configured.

Entry	x.x.x.x
Default value	192.168.10.10

MAC address

Unambiguous address of the network adaptor for identification within a network. This parameter cannot be edited.

NetBIOS name

Unambiguous designation of the device within the network. This parameter cannot be edited.


11.4 Metrohm control instruments

The following Metrohm devices can be used as control instruments:

- Ti-Touch Titrator (permanently entered as control instrument in the device manager)
- Sample Processors

The following specifications apply for all of the control instruments mentioned above.

Additional device configurations for the Sample Processor are described in a separate chapter (*see chapter 11.5, page 92*).

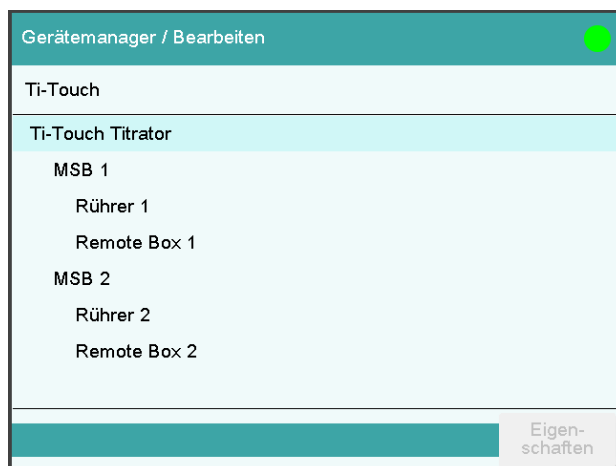

The dialog shows the control instrument with its measuring inputs, the MSB connectors and connected peripheral devices (dosing device, stirrer, Remote Box).

[Properties]


Open the properties dialog of the highlighted entry.

The following properties can be displayed and edited to a certain extent:

- Properties of the measuring inputs (*see chapter 11.4.1, page 89*)
- Properties of the MSB connectors (*see chapter 11.4.2, page 91*)
- Properties of the peripheral devices on the MSB connector (*see chapter 11.4.3, page 92*)

11.4.1 Properties – Measuring input

Instrument list: **Control instrument** ► **Edit** ► **Measuring input**
1 ► **Properties**

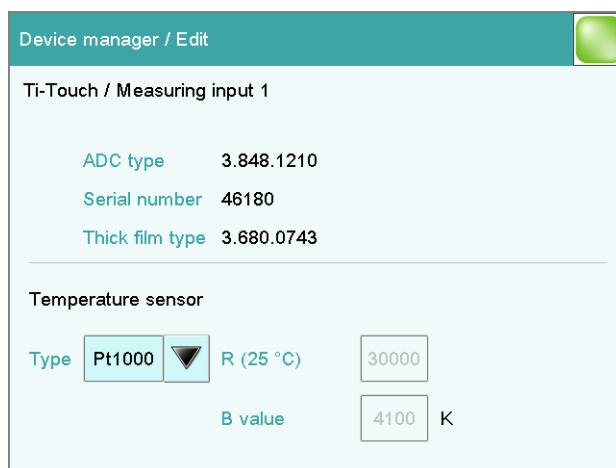

In this dialog, you define the type of temperature sensor you have connected to the selected measuring input.

11.4.2 Properties – MSB connector

Instrument list: **Control instrument ► Edit ► MSB connector
1/2 ► Properties**

Device manager / Edit

Ti-Touch / MSB 1

When would you like to be asked to prepare the dosing device on MSB 1?

☒ Switch on

☒ Attach an exchange or dosing unit

☐ Time interval h

In this dialog, you can define when the request to carry out the **Prepare** function for connected dosing devices is to be displayed. This setting applies for all dosing devices of the selected MSB connector. The manual control contains a description of how to prepare the exchange unit and/or the dosing unit (*see chapter 28.3.3, page 259*).

Switch on

on | off (Default value: **on**)

If this parameter is activated, then you will be requested to prepare the dosing device when the Ti-Touch is switched on.

Attach an exchange or dosing unit

on | off (Default value: **on**)

If this parameter is activated, then you will be requested to prepare the dosing device when the exchange/dosing unit is attached.

Time interval

on | off (Default value: **off**)

Activate this parameter if you wish to receive a regular request to prepare the dosing unit.

Input range	0.1 to 999.9 h
Default value	12.0 h


11.4.3 Properties – Peripheral devices

Instrument list: **Control instrument** ► **Edit** ► **Peripheral device** ► **Properties**

The following data is displayed for the connected dosing devices and stirrers:

- Type
The display includes the device number, e.g. "800" (Dosino).
- Serial number

11.5 Sample Processor

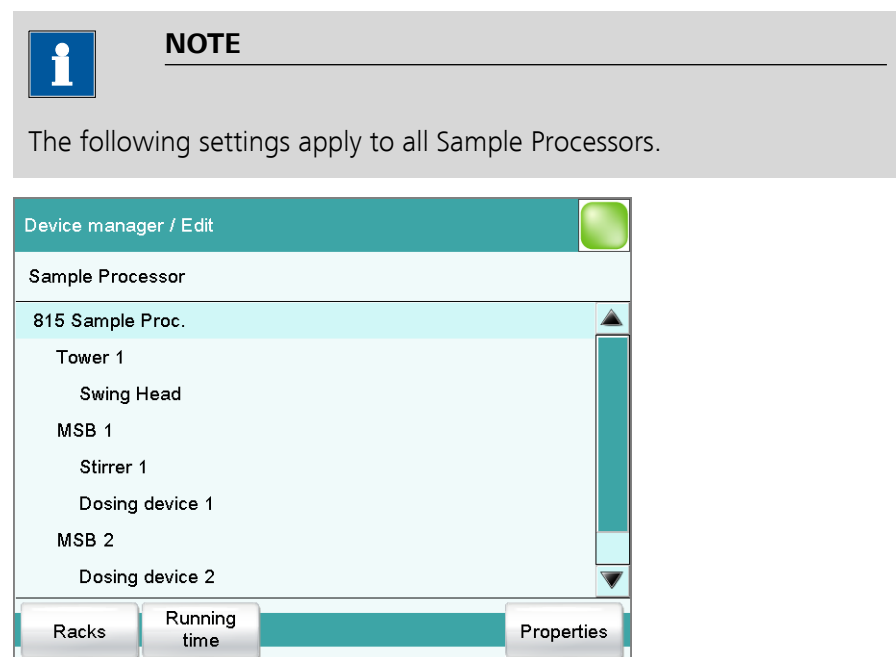

The dialog shows the Sample Processor with its tower (or towers), the MSB connectors and connected peripheral devices (dosing device, stirrer, Remote Box). If the control instrument is not connected, then the properties dialog of the control instrument will be displayed directly (*see chapter 11.5.1, page 93*).

The settings of the MSB connectors and the connected peripheral devices are described in *Chapter "Properties – MSB connector", page 91*.

[Racks]

Configure sample racks (see chapter 11.6, page 100).

[Running time]

Configure the running time meter.

[Properties]

11.5.1 Properties – Sample Processor

Device name

Entry	24 characters maximum
-------	-----------------------

Entry	24 characters maximum
-------	-----------------------

Program version of the instrument software.

Shows the serial number of the device.

Name of the currently attached sample rack.

Display the internal adjustment data (EEPROM data) of the Sample Processor.



■■■■■■■■ 93

11.5.2 Properties – Tower

Instrument list: **Sample Processor** ► **Edit** ► **Tower 1/2** ► **Properties**

Device manager / Edit

Sample Processor / Tower 1

Maximum stroke path	235	mm
Min. beaker radius	off	mm
Rotation increment	5.00	°
Axial distance	196.00	mm

NOTE

In the case of a Sample Processor with 2 towers, Tower 1 is the tower on the right and Tower 2 is the tower on the left, when viewed from the front.

Maximum stroke path

Lowest permissible lift position for automatic and manual operation.

This is a safety setting. A correct specification of this value can prevent the glass of an electrode from breaking, because a titration head can not be moved below the position specified.

Input range	1 to 235 mm
Default value	235 mm

Min. beaker radius

Smallest permissible radius of the sample vessels for automatic and manual operation.

This is a safety setting. It is intended to prevent a completely loaded, wide titration head from being moved into a narrow sample vessel.

Input range	1.0 to 100.0 mm (Increment: 0.1) If the lift is to be moved to the work position, a check is made to see whether the beaker radius defined for the attached rack (see " <i>Beaker radius samples</i> ", page 102) is less than this value.
Selection	off
Default value	off

No check takes place.

Amount by which the sample rack can be rotated forward or backward relative to its current position. This parameter is required in the command **MOVE** for the settings **Destination = rotate +** and **rotate -**.

[Ext. pos. 1-4]

Edit the settings that apply for all external positions (see "Properties – External positions 1-4", page 96).

[Robotic arm]

Edit the robotic arm settings (see *"Properties – Robotic arm"*, page 97).

[Adjustment data]

Display the overview of the EEPROM data of the Swing Head. This dialog cannot be edited.

[Edit]

Edit specific settings of the selected external position (see "Properties – External position", page 98).

Properties – External positions 1-4

Swing Head: **Properties** ► **Ext. pos. 1-4**

Properties / Edit external positions 1-4

Sample Processor / Tower 1 / Swing Head 786

External positions 1-4

Shift position	0	mm
Rinse position	0	mm
Swing increment	10.00	°

Properties which apply to all four external positions can be defined in the dialog **Properties / Edit external positions 1-4**.

Shift position

Lift position at which the robotic arm rotates to the external positions.

Input range	0 to 'Maximum stroke path' mm
Default value	0 mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (see " <i>Maximum stroke path</i> ", page 94).

Rinse position

Lift position used for rinsing.

Input range	0 to 'Maximum stroke path' mm
Default value	0 mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (see " <i>Maximum stroke path</i> ", page 94).

Swing increment

Amount by which the robotic arm can be swung relative to its current position. This parameter is used in the command **MOVE** for the settings **Destination = swing +** and **swing -**.

Input range	0.00 to 180.00 °
Default value	10.00 °

Properties – Robotic arm

Swing Head: **Properties** ► **Robotic arm**

Properties / Robotic arm

Robotic arm Tower 1

Swing offset

0.00

°

Maximum swing angle

84.00

°

Swing radius

110.00

mm

Rotation offset

0.00

°

Swing direction

-



Properties of the robotic arm can be defined in the dialog **Properties / Robotic arm**.

Swing offset

The swing offset is a physical angle offset of a specific robotic arm model. The required values can be found in the leaflet for the robotic arm.

Input range	-270.00 to 270.00 °
Default value	0.00 °

Maximum swing angle

Useable swing range of the robotic arm. Each robotic arm model displays a different value on the basis of its construction. The range can also be reduced if necessary. You will find the required values in the leaflet for the robotic arm.

Input range	0.00 to 330.00 °
Default value	60.00 °

Swing radius

Maximum swing radius of the robotic arm. The swing radius is dependent on the length of the robotic arm and is, together with the axial distance (see "Axial distance," page 95), the most important variable for precise movement to a rack position. You will find the required values in the leaflet for the robotic arm.

Input range	30.00 to 300.00 mm
Default value	110.00 mm

Rotation offset

The rotation offset is the offset from the center of the tower to the center of the robotic arm. This value does not usually need to be changed. If a Swing Head is to be mounted on the tower with a lateral offset, then this value can be determined by the service technician when the rack is adjusted.

Input range	-270.00 to 270.00 °
Default value	0.00 °

Swing direction

The swing direction of the robotic arm depends on its type. For a 2-tower model, the robotic arm must be defined as right-swinging on Tower 1 and as left-swinging on Tower 2.

Selection	+ -
Default value	-

- + Left-swinging.
- Right-swinging.

Properties – External position

Swing Head: **Properties** ► **Edit**

Properties which apply only to the selected external position can be defined in the dialog **Properties / Edit external position**.

Swing angle for the selected external position.

Work position for the selected external position.

11.6 Sample racks

Instrument list: **Sample Processor** ▶ **Edit** ▶ **Racks**

Edit device / Sample racks		
Rack	Positions	Code
6.2041.310	12	000001
6.2041.320	16	000010
6.2041.340	24	001000
6.2041.350	48	010000
6.2041.360	12	100000
6.2041.370	14	000011
6.2041.380	14	000101
6.2041.390	16	100001

Load

Copy

Delete

Edit

In the dialog **Edit device / Sample racks**, you will find the list of configured racks. The attached rack is displayed in green. New racks can be loaded or created and existing ones can be edited or deleted in this dialog window.

The following data is displayed in the list:

- **Rack**
Name of the sample rack. Metrohm standard racks are designated by their order numbers.
- **Positions**
Number of positions on the rack.
- **Code**
The rack code corresponds to the arrangement of magnets on the base of the rack and is read in by the Sample Processor in order to recognize the rack.

[Load]

Load a new sample rack (see "Loading the sample rack", page 101).

[Copy]

Create a new sample rack by copying an existing rack (see "Creating a new sample rack", page 102).

[Delete]

Delete the selected sample rack from the list.

[Edit]


Edit the data of the selected sample rack (see chapter 11.6.1, page 102).

Loading the sample rack

When you receive the file for a new sample rack from your Metrohm representative, you can easily import this file into your existing system. Proceed as follows:

1 Copying the file to an external storage medium

- Copy the file to a subdirectory of "Files".

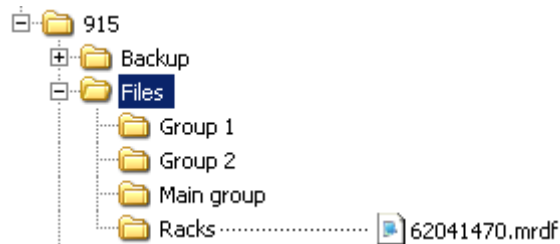

If this structure is not maintained, the new rack will not be found because the software directly accesses subdirectories of "Files".

- Plug in the external storage medium at the Ti-Touch.

2 Display the list of saved sample racks

- In the **Edit device / Sample racks** dialog, tap on the **[Load]** button.

The selection of file groups on the external storage medium is displayed. If only one group is available, then the list of the saved sample rack files will be displayed directly.

- Select the group with the desired sample rack.
- Tap on **[Show files]**.

The list with the saved sample rack files is opened.

3 Load the sample rack file

- Select the desired file.
- Tap on **[Load]**.

The new sample rack is now loaded and appears in the list of available racks.

Creating a new sample rack

Own sample racks can be simply and conveniently created by copying an existing sample rack. Proceed as follows:

1 Copy the existing rack

- In the dialog **Edit device / Sample racks**, select a sample rack which is to be used as a template.
- Tap on **[Copy]**.

The **Sample rack / Copy** dialog is displayed.

2 Enter the rack name and rack code

- In the field **New rack name**, enter a name for the new rack.
- In the field **New rack code**, **110000** is suggested. As a rule, this rack code is used for special racks. Apply this code or enter a new rack code and confirm with **[OK]**.

The new sample rack appears in the list of available racks.

11.6.1 Editing rack data

Sample rack list: **Rack** ► **Edit**

Sample rack / Edit rack data

Rack 6.2041.410, Code 001010, 142 Positions

Beaker radius samples mm

Beaker sensor 

1. Calibration pos. 

Rack offset °

Adjust rack Lift pos. tower 2 Lift pos. tower 1 Special beakers

You can edit the data of the selected rack in the dialog **Sample rack / Edit rack data**.

Beaker radius samples

Actual radius of the sample vessels at the general sample positions of the rack.

This beaker radius may not be less than the minimum beaker radius defined in the tower properties (see *"Min. beaker radius,"* page 94). If the

lift is to be moved to the work position, then these two values will be compared with one another.

Input range	1.0 to 100.0 mm
Selection	off
Default value	off

off

No check takes place.

Beaker sensor

Each time a sample position is to be moved to with the **MOVE** command, the beaker sensor checks whether a vessel is present. In the **MOVE** command, you define the action that takes place if the beaker sensor does not detect a vessel at the position being moved to.

Selection	Tower Robotic arm off
Default value	off

Robotic arm

A Swing Head with beaker sensor must be mounted. In addition, a suitable work position must be defined for the lift, so that the robotic arm touches the sample vessel. The work position is moved to for the purpose of beaker detection.

off

No check takes place.

Rack offset

The rack offset is a production-related tolerance value between the upper and lower sections of the rack. This value is determined when the rack is adjusted and is displayed here. It can be changed if necessary.

Input range	-10.00 to 10.00 °
Default value	0.00 °

[Adjust rack]

Adjust the rack (*see chapter 11.6.2, page 107*).

[Lift pos. tower 1]/[Lift pos. tower 2]

Define rack specific lift positions (*see "Lift positions Tower 1/2", page 103*).

[Special beakers]

Define rack positions as special beakers (*see "Editing special beakers", page 105*).

Lift positions Tower 1/2

Sample rack list: **Rack ► Edit ► Lift pos. tower 1 / Lift pos. tower 2**

Edit rack data / Lift positions tower 1

Rack 6.2041.410, Code 001010, 142 Positions

Work position	120	mm
Rinse position	0	mm
Shift position	0	mm
Special position	0	mm

You can define rack-specific lift positions in the dialog **Edit rack data / Lift positions tower 1/2**. These then apply for all rack positions except those that have been defined as special beakers.

These positions can be moved to directly under manual control and with the LIFT command. Only lift positions within the maximum stroke path can be entered. This is defined in the device properties of the tower.



NOTE

These lift positions can also be assigned directly in manual control after moving to the desired lift height (*see chapter 28.6.1, page 266*).

Work position

At this lift position the electrodes, stirrer and buret tips are optimally positioned for work.

Input range	0 to 'maximum stroke path' mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (see "<i>Maximum stroke path</i>", page 94).
-------------	---

Rinse position

Lift position used for rinsing.

Shift position

Input range	0 to 'maximum stroke path' mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (<i>see "Maximum stroke path", page 94</i>).
-------------	---

Input range	0 to 'maximum stroke path' mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (<i>see "Maximum stroke path", page 94</i>).
-------------	---

Sample rack list: **Rack** ► **Edit** ► **Special beakers**

■■■■■■■ 105

A maximum of 16 rack positions can be defined as special beakers for each sample rack. The dialog **Edit rack data / Special beakers** shows a list of all special beakers with their assigned rack positions.

[Edit]

Edit the data of the selected special beaker, see the following.

Special beaker / Edit

Special beaker 1

Rack position	142	
Work pos. tower 1	100	mm
Work pos. tower 2	0	mm
Beaker radius	off	mm
Beaker sensor	Tower	

You can edit the data of the selected special beaker in the dialog **Special beaker / Edit**.

Rack position

Number of rack position for selected special beaker. Any rack position can be defined as a special beaker. It is, however, preferable to set them at high rack positions in order to be able to begin sample series at rack position 1. Rack positions defined as special beakers can no longer be used as sample positions.

Input range	0 to 'maximum number of rack positions'
-------------	---

Work pos. tower 1/2

Work position for the selected special beaker. One specific work position can be defined for tower 1 and tower 2.

Input range	0 to 'maximum stroke path' mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (see "<i>Maximum stroke path</i>", page 94).
-------------	---

Beaker radius

Actual radius of the selected special beaker.

This beaker radius may not be less than the minimum beaker radius defined in the tower properties (see *"Min. beaker radius,"* page 94). If the

Input range	1.0 to 100.0 mm
Selection	off
Default value	off

No check takes place.

Beaker sensor

Selection	Tower Robotic arm off
Default value	off

A Swing Head with beaker sensor must be mounted. In addition, a suitable work position must be defined for the lift, so that the robotic arm touches the sample vessel. The work position is moved to for the purpose of beaker detection.

No check takes place.

11.6.2 Rack adjustment


If necessary, each sample rack can be finely adjusted, i.e. the rack offset is determined in the direction of rotation. The adjustment of a rack is usually not necessary. However, if the exact positioning of a robotic arm is necessary for a high degree of accuracy (e.g. for very small sample beakers), then a fine adjustment can be carried out.

Proceed as follows:

1 Open the dialog for rack adjustment

- The **Sample rack / Edit rack data** dialog is displayed.

The **Edit rack data / Adjust rack** dialog is displayed:

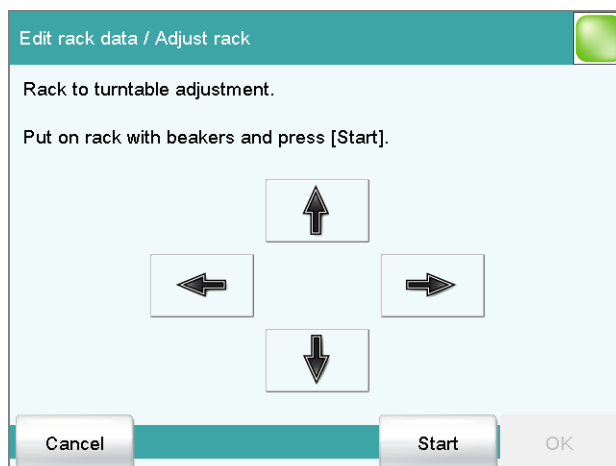

2 Carry out the adjustment

- Tap on the **[Start]** button (NOT on the **[▶]** fixed key).
The rack is being initialized. The rack then moves to position 1 and the lift is lowered to the work position.
- If necessary, the lift position can be corrected with the arrow keys **[↑]** and **[↓]**.
- Now use the arrow keys **[←]** and **[→]** to rotate the rack so that the robotic arm or the center of the titration head is positioned exactly above the center of rack position 1.
- Finish the adjustment with **[OK]**.

The lift is moved to the uppermost stop.


11.7 885 Compact Oven SC



NOTE

In order for the 885 Compact Oven SC to be controlled with the 915 KF Ti-Touch, at least the following firmware versions must be installed on the instruments:

- 915 KF Ti-Touch 5.915.0045
- 885 Compact Oven SC 5.885.0013

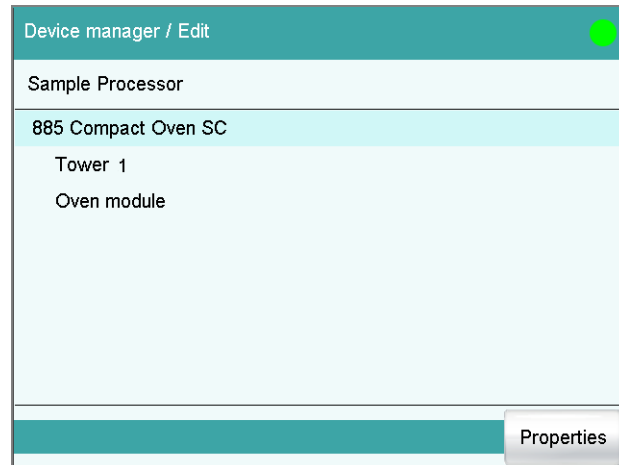

The dialog displays the 885 Compact Oven SC with its tower and the oven module. If the control instrument is not connected, then the properties dialog of the control instrument will be displayed directly (*see chapter 11.7.1, page 109*).

Properties

Open the Properties dialog of the highlighted entry.

11.7.1 Properties - 885 Compact Oven SC

Instrument list: **885 Compact Oven SC** ▶ **Edit** ▶ **Properties**

Device name

This designation is used for identification purposes when selecting control instruments (command, manual control).

Entry	24 characters maximum
-------	-----------------------

Comment

Entry	24 characters maximum
-------	-----------------------

Program version

Program version of the device software.

Serial number

Display of the serial number of the device.

11.7.3 Properties - Oven module

Instrument list: **885 Compact Oven SC** ► **Edit** ► **Oven module** ► **Properties**

Device manager / Edit

Sample Processor / Oven module

Program version

5.860.0011

Serial number

07134

Initial temperature

110

°C

Action when turning off the touch-device

☒

Turn off heater

☒

Turn off gas flow

Program version

Program version of the oven module.

Serial number

Display of the serial number of the oven module.

Initial temperature

Entering an initial temperature means that the oven will heat up to the set temperature when the 915 KF Ti-Touch is switched on.

Input range	50 to 250 °C (Increment: 1)
Default value	110 °C
Selection	off

off

The oven is not switched on.

Action when switching off the touch device.

Switch off heater

Setting for the behavior of the heater of the oven module when the Ti-Touch is switched off.

Selection	on off
Default value	on

Switch off the gas flow

Setting for the behavior of the gas flow of the oven module when the Ti-Touch is switched off.

off

The reports will not be printed out on paper.

Connector

Selection of the connection type for the printer.

Selection	USB Ethernet
Default value	USB

USB

Local printer at a USB interface.

Ethernet

Network printer.



NOTE

The **Ti-Touch** prints the reports with a fixed resolution of 300 dpi. If you are using a printer with a resolution of 360 dpi (or a multiple thereof, e.g. an Epson), then the text will be printed out somewhat smaller than with printers with a resolution of 300 dpi (or a multiple thereof, e.g. a Canon or HP).

11.8.1 PDF settings

Printer: **Edit** ► **PDF settings**

The settings for saving a report as a PDF file are defined in this dialog.

Memory

Memory location where the PDF file will be saved. The report will be saved in the directory *PDF_Report*. This directory is created when a PDF file is generated. The file name is appended with the date and time when the determination was started and the date and time when the PDF was generated (e.g.: **Identification 1**-20130320-141044-20130320-141220). The first date and time set indicates the start of the determination, the second date and time set indicates the generation of the PDF file. This prevents previously created PDF files from being overwritten each time a new file is created, e.g. when a result is recalculated and the report is generated again.

All three memory locations are listed as possible selections, even if they cannot be accessed at the moment.

Selection	off External memory 1 External memory 2 Shared memory
Default value	off

off

The report will not be saved as a PDF file.



Password

Password of the user configured on the computer.

Entry	max. 32 characters
Default value	empty

11.8.3 More options

Paper size

Selection of the paper size.



NOTE

The paper size selected will be ignored when the report is printed on a **CUSTOM NEOS** printer.

Selection	A4 (210 mm x 297 mm) Letter (216 mm x 279 mm)
Default value	A4 (210 mm x 297 mm)

Color

on | off (Default value: **off**)

If this option is activated, then the report will be printed out in color.


Spooler

on | off (Default value: **on**)

If this option is activated, then you can continue to work while the printer is printing. If this option is deactivated, then the Ti-Touch will be blocked while the printout is being produced.

11.9 Balance

List of devices: **Balance** ► **Edit**

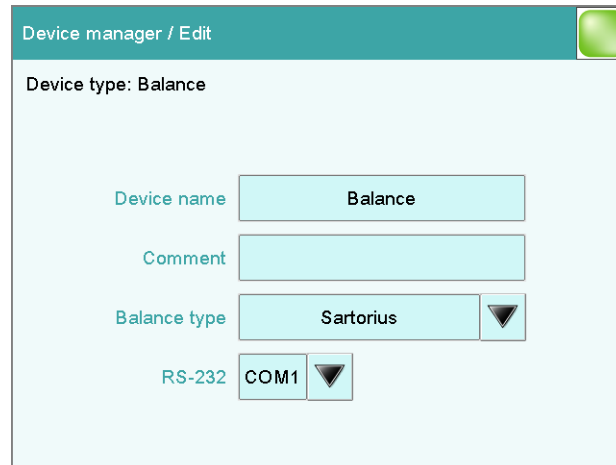

You define the type of balance and its connector in the dialog **Device manager / Edit**.

Device name

Here you can enter a device name of your choice.

Entry	24 characters maximum
-------	-----------------------

Comment

Entry	24 characters maximum
-------	-----------------------

Balance type

If you have connected a balance, then you have to specify the balance type here.

Selection	AND Mettler Mettler AT Mettler AX
	Ohaus Precisa Sartorius Shimadzu
Default value	Sartorius

RS-232

RS-232 interface to which the balance is connected. The interface parameters are adjusted in the **Edit device / Port parameters** dialog (see chapter 11.10, page 118).

Selection	COM 1 COM 2
Default value	COM 1

COM 2

The interface is disabled.

The following table indicates the balance type that needs to be selected for the balance model:

Balance	Balance type
AND	AND
Mettler AB, AE, AG, AM, AJ, PE, PM, PJ, PR, XP, XS	Mettler
Mettler AT	Mettler AT
Mettler AX, MX, UMX, PG, AB-S, PB-S	Mettler AX
Ohaus Voyager, Explorer, Analytical Plus	Ohaus
Precisa	Precisa
Sartorius	Sartorius
Shimadzu BX, BW	Shimadzu

11.10 USB/RS-232 adapter

Instrument list: **USB/RS-232 adapter** ▶ **Edit** ▶ **COM interface** ▶ **Edit**

As a rule, balances are equipped with a serial RS-232 interface. To connect a balance, you will require the 6.2148.050 cable. You can adjust the interface parameters in the dialog **Edit device / Port parameters**. These settings must match the settings on the connected device.



NOTE

If you make changes to parameter settings in this dialog, then you must switch the Ti-Touch off and back on in order for the changes to take effect.

Edit device / Port parameters

Baud rate

9600

Data bits

8

Parity

none

Stop bits

1

Handshake

Hardware (DTR/CTS)

Transfer rate in characters per second.

Number of data bits.

Type of parity testing.

Number of stop bits.

Type of data transfer protocol.



Function in the editing dialog	PC keyboard key
Enter respective character	Numbers, letters and special characters + Shift
[OK]	[Enter] / [↵]
[□]	[←] (backspace)
Line break (for multi-line text input)	[Enter] / [↵]
The cursor is moved up or down by one character.	[↑] [↓]
The cursor is moved to the left or to the right by one character.	[←] [→]
	[F1]
	[F2]
	[F3]
	[F4]
	[F5]
	[F6]
	[F7]



NOTE

The lettering of the USB keyboard may differ from above lettering, depending on the country-specific keyboard used.

11.12 Barcode reader

List of devices: **Barcode reader** ► **Edit**

Device manager / Edit

Device type: Barcode reader

Device name

Barcode reader

Comment

Barcode input target

Active input field

Keyboard layout

English US

A barcode reader can be connected to read in sample data or other texts. A connected barcode reader will be recognized automatically and entered in the list of devices with default settings.

You will hear an acoustic signal as confirmation that a character string has been transmitted by the barcode reader and accepted.



NOTE

The data will be read only if the system is in normal status, i.e. when no determination is running.

Device name

Here you can enter a device name of your choice.

Entry	24 characters maximum
-------	-----------------------

Comment

Entry	24 characters maximum
-------	-----------------------

Barcode input target

Selection of the input field for the character string read in by the barcode reader.

Selection	Active input field	Method	Identification 1	Identification 2	Sample size
-----------	--------------------	--------	------------------	------------------	-------------

Active input field

The character string is entered in the input field of the opened text- or number-input dialog.

Method

The character string is entered in the input field **Method**.

Identification 1

The character string is entered in the input field **Identification 1**.

The character string is entered in the input field **Identification 2**.

The character string is entered in the input field **Sample size**. Character strings containing characters other than numbers and decimal separators will be ignored.

Keyboard layout


Specify the country-specific keyboard layout for the emulation of the PC keyboard. This setting must match the setting on the barcode reader (see documentation for the barcode reader).

Selection	English US German DE French FR Spanish ES German CH
Default value	English US



NOTE

If you use an external storage medium with the FAT or FAT32 file system, then you can save a maximum of 999 files per group. If you find it necessary to store more than 999 files in a single group, then you must reformat the storage medium with the file system **ExFAT** (see chapter 32.8.2, page 369).

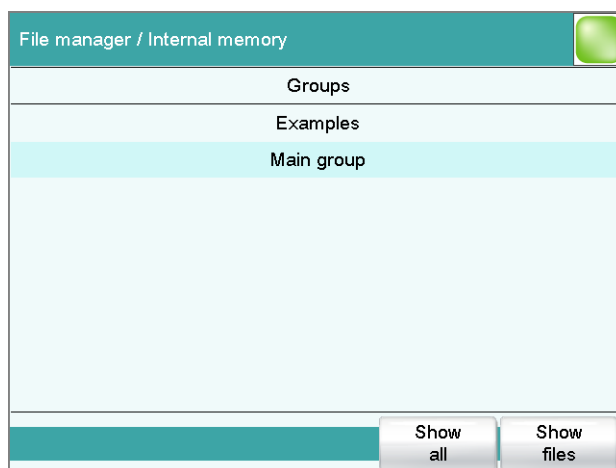


[Show all]

Display the list of all the files in the selected memory location.

[Show files]

Display the list of all the files in the selected group.

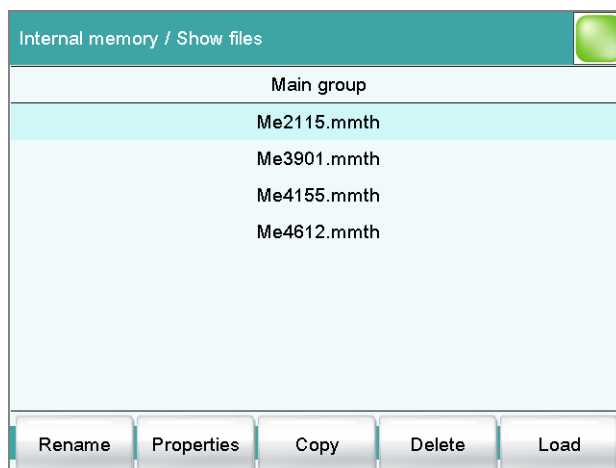

[Rename]

Rename the selected file (see chapter 12.1.2, page 126).

[Properties]

Display the file properties (see chapter 12.1.3, page 127).

[Copy]

Copy the selected file to a different memory location (see chapter 12.1.1, page 126).

[Delete]

Delete the selected file.

[Load]

Load the selected file.

12.1.1 Copying a file

Proceed as follows to copy a file:

1 Selecting the file

- Select the desired file.
- Tap on **[Copy]**.

The selection of memory locations is displayed.

2 Copy the file

- Tap on the desired memory location.
Only memory locations currently being accessible are active.

The file is copied. The file group is retained, i.e. the group will be created again in the new memory location if it does not already exist there.

12.1.2 Renaming a file

Proceed as follows to rename a file:

1 Selecting the file

- Select the desired file.
- Tap on **[Rename]**.

The **Show files / Rename** dialog is displayed.

2 Change the file name

- Tap on the input field **File name**.
The text editor opens.
- Enter a new file name (max. 32 characters) and confirm with **[OK]**.
- Tap on **[OK]**.

The file is saved under the new name.


12.1.3 File properties

File manager: **Internal memory / External memory 1 / External memory 2 ► Show all / Show files ► Properties**

Detailed information concerning the file is displayed in the **Show files / Properties** dialog. They can be edited to a certain extent.

- **File name**
Name of the file.
- **File group**
Entry or selection of the group in which the file will be saved. If you enter a new name, then the file will be moved into the new group.
- **Write protection**
If this option is activated, then the file cannot be saved, deleted, moved or renamed. This is only an internal write protection feature and is independent of the write protection of the operating system on your computer.
- **Created by**
User who saved the file initially.
- **Created on**
Date and time at which the file was saved for the first time.
- **Last saved by**
User who saved the file most recently.
- **Last saved on**
Date and time at which the file has been saved most recently. Additionally, the version of the file is displayed. The version number will be increased by one each time the file is saved.
- **Size**
File size.
- **Program version**
Program version of the software with which the file was last saved.

12.2 External storage medium

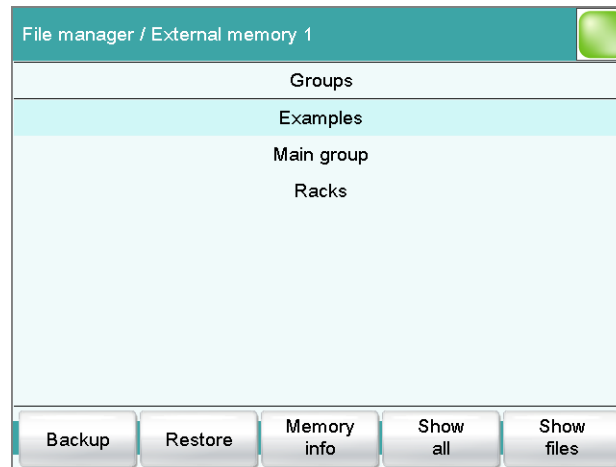

[Backup]

Creating a backup of all data and settings on this storage medium (see *chapter 12.3, page 130*).

[Restore]

Load the backup. This function is active only if a backup is available (see *chapter 12.3, page 130*).

[Memory info]

Display information on the storage medium, see following chapter.

Storage medium – Information

Detailed information about the storage medium is displayed in the dialog **External memory 1/2 / Memory info**:

- **Name**
Designation of the storage medium.
- **Write protection**
If this option is activated, then the file cannot be saved, deleted, moved or renamed. This is only an internal write protection feature and is independent of the write protection of the operating system on your computer.
- **Storage capacity**
Total capacity of the storage medium.
- **Used memory**
Storage capacity used on the storage medium.
- **Free memory**
Free storage capacity on the storage medium.

Directory structure

The directory structure appears as follows on the external storage medium:

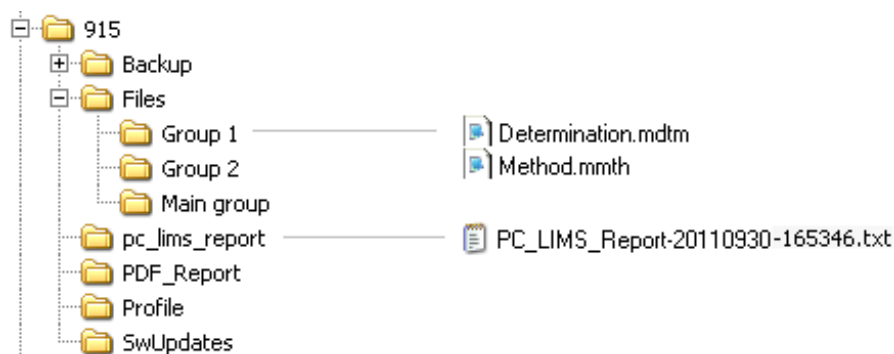


Figure 12 Directory structure on the external storage medium

Backup	All of the files of the backup are stored in this directory. The directory will be created the first time a backup is created.
Files	This directory contains all of the groups as well as the files contained in them (methods, determinations, etc.).
pc_lims_report	PC/LIMS reports are stored in this directory as TXT files. The directory will be created the first time a PC/LIMS report is printed out.
PDF_Report	PDF reports are stored in this directory. The directory will be created the first time a report is saved as a PDF file.
Profiles	The identification profile is stored in this directory. The directory will be created the first time an identification profile is generated.
SwUpdates	Files for software updates must be saved in this directory.

12.3 Creating backups / Restoring data

File manager: **External memory 1 / External memory 2**

You can use the **[Backup]** function to easily create a backup containing all the data and settings of your system. You should make a backup at regular intervals in order to avoid data loss.



NOTE

Only **one** backup can be created on a storage medium.

If a backup has already been stored on the medium, then it will be overwritten when this function is carried out again.

12.3.1 Restoring data

You can use the function **[Restore]** to restore either a complete backup or only certain data.



CAUTION

All of the methods in the internal memory will be deleted.

Backups are not backward compatible, i.e. backups of the latest version cannot be restored using earlier versions.

Proceed as follows to restore backed-up data:

1 Select the storage medium

- Connect the storage medium on which the backup is located.
- Select the storage medium under **System ► File manager**.

2 Select data

- Tap on **[Backup]**.
The selection of the data that can be restored separately is displayed (see the following).
- Deselect the data sets which are not to be restored.
- Tap on **[Load data]**.

The selected data blocks are restored.

3 Complete the restoration

- Switch the Ti-Touch off and on again.

Restorable data blocks

The following data blocks can be loaded individually:

- **Methods**
All of the methods stored in the internal memory.
- **Sample table**
Current sample table.
- **Result table**
Current result table.
- **Current determination data**
All the data for the current determination (including the method with which the determination was carried out).
- **Control**
Settings in the control dialog.
- **User list**
Settings for each user under **System settings / User administration**.
- **System settings / User admin.**
All of the system settings, including device-specific dialog configuration and dialog options for the command list and fixed keys, device-specific settings for the user administration (login options and password options).
- **Titration**
All titration with all their data.
- **Sensors**
All sensors with all their data.
- **Device data**
All of the devices configured in the device manager, with all of the data.
- **GLP data**
All data from the GLP manager.
Deactivate if the backup is to be loaded onto a different system.
- **Common variables**
All common variables.
- **Templates**
All templates for sample data, result calculations, etc.
- **Routine dialog settings**
Current routine dialog settings (see **Dialog options / Routine dialog**).
- **Sample racks**
The sample racks present in the properties of the USB sample processor.
- **Own subsequences**
All subsequences created and stored by the user.

[GLP tests HW/SW]

Document the GLP tests "Measurement" and "Titration" (*see chapter 13.3, page 134*).

[System validation]

Document the system validation (*see chapter 13.4, page 138*).

[Monitoring]

Defining the time intervals for system maintenance and backups (*see chapter 13.5, page 142*).

[View test data]

Display the data of the selected test.

13.1 Automatic system test

The system test is carried out automatically when the Ti-Touch is switched on. The result of each individual test is shown in green if no error has occurred. If a result is shown in red, then an error occurred when the respective test was carried out. Switch the Ti-Touch off and back on again. If the error still occurs please notify Metrohm Service.

Print report at system start

on | off (Default value: **off**)

If this parameter is activated, then the result of the system test will be printed out automatically when the Ti-Touch is switched on.

13.2 Test tools

Main dialog: **System ► GLP manager ► Test tools**

You can create a list of test tools to be used in the tests in the GLP manager.

The following test tools have already been defined:

- **767 Calibrated Reference:** Device for checking measuring inputs and electrode cable.
- **822 Titr.Curve Simulator:** Curve simulator for checking the hardware and software.
- **773 pH/mV Simulator:** Device for checking measuring inputs and electrode cable.
- **868 UR Generator:** Device for checking measuring inputs and electrode cable.

[New]

Add a new test tool to the list.

[Delete]

Delete the selected test tool from the list.

[Edit]

Change the designation of the selected test tool.

13.3 GLP tests for measurement and titration

Main dialog: **System ► GLP manager ► GLP tests HW/SW ► GLP test "Measurement" / GLP test "Titration"**

You can document the GLP tests for measurements and titrations in the dialog **GLP manager / GLP tests Hardware/Software** and its subdialogs. The following procedure describes the GLP test "Measurement," but also applies for the GLP test "Titration".

Proceed as follows:

1 Open the properties dialog

- In the **System / GLP manager** dialog, tap on the **[GLP tests HW/SW]** button and then on the **[GLP test "Measurement"]** button.

GLP tests HW/SW / GLP test "Measurement" 

Method	SOP_1231a	
User	chk	
Test date	2011-08-01	
Test result	Test OK	
Comment		

GLP test interval Hardware

2 Edit data

- Define the test method, the test results, etc. (see "Dialogs "GLP tests HW/SW / GLP test "Measurement"" and "GLP tests HW/SW / GLP test "Titration""", page 135).

3 Define the hardware used

- Tap on **[Hardware]**.

- The dialog **GLP tests HW/SW / GLP test "Measurement"** appears again.

4

- Tap on the fixed key [].

13.3.1 Parameter description

Dialogs "GLP tests HW/SW / GLP test "Measurement"" and "GLP tests HW/SW / GLP test "Titration""

Method

Method the GLP test has been carried out with.

Entry	32 characters maximum
Selection	Selection of methods stored in the internal memory

User

User who carried out the GLP test.

Entry	24 characters maximum
Selection	Selection of configured users

Test date

Date on which the GLP was carried out.

Format: YYYY:MM:DD

Test result

Result of the test.

Selection	Test OK Test not OK
Default value	Test OK

Comment

Entry	24 characters maximum
-------	-----------------------



[GLP test interval]

Defining the time interval for the GLP test (see "Dialogs "GLP test "Measurement" / Test interval" and "GLP test "Titration" / Test interval", page 136).

[Hardware]

Documenting the hardware with which the GLP test was carried out (see "Dialogs "GLP test "Measurement" / Hardware" and "GLP test "Titration" / Hardware", page 137).

Dialogs "GLP test "Measurement" / Test interval" and "GLP test "Titration" / Test interval"

Test date

Date on which the GLP was carried out.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the time interval after which a GLP test has to be carried out again will be monitored.

GLP test interval

If you define a time interval for the GLP test, then the date in **Next GLP test** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next GLP test

If you define a date for the next GLP test, then the **GLP test interval** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

Dialogs "GLP test "Measurement" / Hardware" and "GLP test "Titration" / Hardware"

Test tool

Selection of the test tool. Test tools are defined at **GLP manager ► Test tools**.

Selection	Selection of configured test tools
-----------	------------------------------------

Control device

Selection of the control device with which the GLP test has been carried out.

Selection	Selection of configured control devices
-----------	---

Measuring input

Selection of the measuring input used for the GLP test. The selection is not dependent on whether the control device has one or two measuring interfaces.

Selection	1 2
Default value	1

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

The dialog **GLP manager / System validation** appears again.

5 Insert a note

- Tap on **[Note (SOP)]**.
- Enter a brief description, e.g. a summary of the SOP (standard operating procedure) according to which the system validation was carried out (*see ""System validation / Note (SOP)" dialog", page 140*).
- Tap on the fixed key **[↩]**.

13.4.1 Parameter description

"GLP manager / System validation" dialog

Method

Method with which the system validation has been carried out.

Entry	32 characters maximum
Selection	Selection of methods stored in the internal memory

User

User who carried out the system validation.

Entry	24 characters maximum
Selection	Selection of configured users

Test date

Date on which the system validation was carried out.

Format: YYYY:MM:DD

Test result

Result of the test.

Selection	Test OK Test not OK
Default value	Test OK

Comment

Entry	24 characters maximum
-------	-----------------------

[Note (SOP)]

Enter a brief description, e.g. a summary of the SOP (standard operating procedure) according to which the system validation has been carried out.

[Validation interval]

Defining the time interval for the system validation (see "Dialog "System validation / Validation interval"", page 140).

[Test statistics]

Documenting the statistical data of the system validation (see *""System validation / Test statistics" dialog"*, page 141).

"System validation / Note (SOP)" dialog

In this dialog, you can enter a brief text, e.g. a summary of the SOP (standard operating procedure) according to which the system validation was carried out.

Dialog "System validation / Validation interval"

Last validation

Date on which the last system validation was carried out.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the time interval after which a system validation has to be carried out again will be monitored.

Validation interval

If you define a time interval for the system validation, then the date in **Next validation** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next validation

If you define a time interval for the next system validation, then the **Validation interval** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

In the determination data it will be documented, that the time interval has been expired.

The determination is stopped.

"System validation / Test statistics" dialog

In the dialog **System validation / Test statistics**, you can document the statistical data for the last system validation.

Number (n)

Number of determinations carried out.

Input range	1 to 9999999999
Default value	empty

Mean value

Mean value out of single results.

Input range	-999999999.00000 to 999999999.00000
Default value	empty

s abs

Absolute standard deviation of the results.

Input range	0.00000 to 999999999.00000
Default value	empty

s rel

Relative standard deviation of the results.

Input range	0.00000 to 100.00000 %
Default value	empty %

d rel

Systematic deviation of the results.

Input range	0.00000 to 100.00000 %
Default value	empty %

a sys

Systematic error.

Input range	0.00000 to 9999999999.00000
Default value	empty

13.5 System monitoring

13.5.1 Service interval

Main dialog: **System ► GLP manager ► Monitoring ► Service interval**

In the dialog **GLP manager / Service interval**, you can define the time interval for the maintenance of the system by the Metrohm Service department. The service interval is checked each time the system is started.

Last service

Date on which the last servicing was carried out.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the time interval after which system maintenance must be carried out again will be monitored.

Service interval

If you define a time interval for the system maintenance, then the date in **Next service** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next service

If you specify a date for the next system maintenance, then the **Service interval** will be tracked automatically.

Format: YYYY:MM:DD

13.5.2 Backup interval

Main dialog: **System ► GLP manager ► Monitoring ► Backup interval**

In the dialog **GLP manager / Backup interval**, you can define the time interval for backups. The backup interval is checked each time the system is started.

Last backup

Date on which the last backup was created.

Format: YYYY:MM:DD

14 Common variables

Main dialog: **System** ► **Common variables**

You can save 25 **method-independent variables**, or common variables. These variables can be used in future calculations (as variables **CV01...CV25**). Common variables are useful, e.g. for the following applications:

- Determination of a blank value which will be taken into account during the content determination of the sample.
- Determination of the content of a standard solution, which will be taken into account during the content determination of the sample.

System / Common variables		
CV	Name	Value
01	Blank value	0.0143 mL
02		
03		
04		
05	Density	0.986 g/mL
06		
07		
08		

Delete

Edit

The designation and the value (including the unit) are displayed for every common variable in the list.

[Delete]

Delete the selected common variable from the list.

[Edit]

Edit the data of the selected common variable, see following chapter.

14.1 Editing common variables

List of common variables: **Common Variable ► Edit**

The common variables can be modified as follows:

- Edit manually, see the following.
- Automatic assignment from the determination run. A calculation result must be configured accordingly for this purpose (*see chapter 14.4, page 148*).

Name

Designation of the common variable.

Entry	24 characters maximum
Default value	empty

Value

Value of the common variable.

Input range	–999999999 to 999999999
Default value	empty

Unit of the common variable.

Entry	10 characters maximum
Default value	empty
Selection	% mol/L mmol/L g/L mg/L mg/mL mg/100 g ppm g mg µg mL µL mg/ piece °C µg/min mL/min µL/min
Default value	%

14.3 Monitoring validity

Common variable: **Edit ► Validity**

In the dialog **Edit common variables / Validity**, you can define the time interval after which a new value must be assigned to the common variable.

Last assignment

Date on which the common variable was last assigned a value.

Format: YYYY:MM:DD

Monitoring

on | off (Default value: **off**)

If this parameter is activated, then the time interval after which the common variable must be assigned a new value will be monitored.

Validity

If you define a time interval for the validity of the common variable, then the date in **Next assignment** will be tracked automatically.

Input range	1 to 999 days
Default value	999 days

Next assignment

If you specify a date for the next assignment, then the time interval for the **Validity** will be tracked automatically.

Format: YYYY:MM:DD

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message

For all three options it is documented in the determination data (see dialog **More determination data / Messages**), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

14.4 Assigning a result automatically to a common variable



NOTE

This instruction is based on the assumption that the method contains a calculation command with a calculation.

Proceed as follows to assign a result to a common variable:

1 Opening the editing dialog of the result

- In the command list, select the command **CALC**.
- Tap on the **[Edit command]** button.
- Select the result whose value is to be assigned to a common variable and tap on **[Edit]**.

The editing dialog of the result is displayed.

2 Define result options

- Tap on the **[Result options]** button.

R1 Blank value

Variable for mean value

Save as titer

Variable

- Activate the parameter **Save as common variable**.
- Define **Variable** = CV01...CV25.

3

Tap on the fixed keys [←] or [🏠].

In the future, the calculated result will be assigned to the selected common variable (result name, value and unit).

- ### Use sample identification list

Activate this parameter so that the sample identifications defined in the list will be displayed as selections.

on | off (Default value: **off**)

If this parameter is activated, the following will no longer be possible:

- Assignment ident.**


Identification which is used as assignment identification to load the correct method.

Request assignment identification

If this parameter is activated, then the sample identification will be requested automatically at the start of a determination.

15.1.1 Sample identification list

Main dialog: **System ► Templates ► Sample data ► Sample identification list**

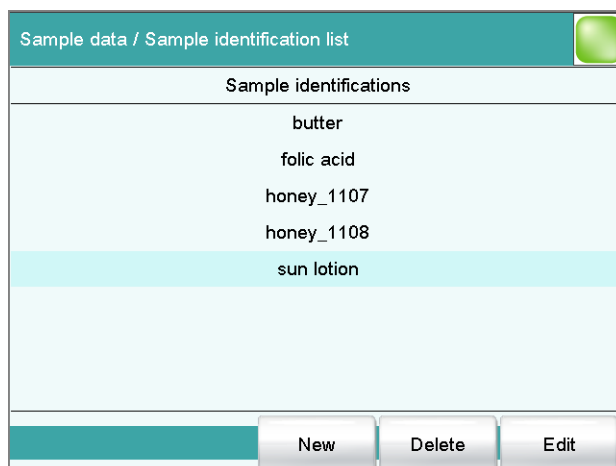

You can create a system-specific list with sample identifications in the dialog **Sample data / Sample identification list**. This list makes it easier for you to enter frequently used sample identifications, i.e. this list is available to you for entering them (main dialog, command **REQUEST**, etc.). It may be a good idea to define the unchanging part of the identification as a template and to add the variable part during sample data input.

[New]

Add a new sample identification to the list.

[Delete]

Delete the selected sample identification from the list.

[Edit]

Change the designation of the selected sample identification.

15.1.2 Sample assignment table

Main dialog: **System ► Templates ► Sample data ► Sample assignment table**

Sample data / Sample assignment table	
Sample ident.	Method
butter	Me2115
folic acid	Me4155
honey*	Me3901
sun lotion	Me4612

New
Delete
Edit

[New]

[Delete]

[Edit]

Editing the sample assignment

Identification of the sample.



You can place an `*` as a wildcard at the beginning or end of the character string. Doing so allows you to prefix or suffix a sequential number, for example, to the identification, which will be ignored when methods are being assigned.

Entry	24 characters maximum
Selection	Selection of defined identifications in the sample identification list

Memory location the method is loaded from. All memory locations are selectable, even if they are currently not accessible.

Selection	Internal memory External memory 1 External memory 2 Shared memory
Default value	Internal memory

Shared directory in the network.

Method that is loaded when a sample with the above-specified identification is processed. You can also enter a method that does not yet exist in the selected memory. When a determination is started there is a check whether the method is available.

15.2.1 Editing result templates

List of result templates: **Result template** ► **New / Edit**

Result name

The result name is the text which will be shown in the results display and in the report.

Entry	max. 24 characters
Default value	R

Calculation formula

Shows the calculation formula. A special editor is opened for the definition (see chapter 29.4.3, page 308).

Entry	max. 100 characters
Default value	empty

Decimal places

Number of decimal places used to display the result.

Input range	0 to 5
Default value	2

Result unit

The result unit is displayed and saved together with the result.

Entry	max. 10 characters
Selection	% mol/L mmol/L g/L mg/L mg/mL mg/100 g ppm g mg µg mL µL mg/ piece °C µg/min mL/min µL/min
Default value	%

[Note]

Entering a note on the calculation.

[Note for wizard]

Enter a note on the calculation. This note is displayed when loading the result template.

[Result options]

Defining additional settings for the calculation.

"Edit result template / Result options" dialog

In the dialog **Edit result template / Result options**, settings for how to process the calculated result are defined.

Variable for mean value

If the statistic calculations have been activated (see method options), the mean value of the single results will be saved as variable SMN1 to SMN9. As default value, always the first free variable is being displayed.

Selection	SMN1 ... SMN9	off
-----------	----------------------	------------

off

For the result not statistic calculations are carried out.

Save as titer

The result can be saved as titer for the selected titrant.



NOTE

The titer is assigned to the titrant defined in the last titration command before the calculation. Please note that the CALC command containing the titer assignment is inserted after the titration command with which the titer is determined.

Selection	Single value	Mean value	off
Default value	off		

Single value

The result of the current determination is saved as titer.

Mean value

If statistic calculations are carried out for the result, the current mean value of the determination series is saved.

off

The result is not saved as titer.

Save as common variable

on | off (Default value: **off**)

The calculated result can be saved as a method-independent variable, a so-called common variable. The result is then also available in other methods for calculations. In all cases, only the single value will be saved, even if the statistics function is enabled.

Variable

Selection of the common variable to which the result is being assigned.

Selection	CV01 ... CV25
-----------	---------------

on | off (Default value: **on**)

Save result in result table

on | off (Default value: **off**)

Precision

Setting, with which accuracy the result is used in additional calculations.

Selection	Round Truncate Full precision
Default value	Round
Round	The result is rounded to the defined number of decimal places (commercial rounding, in accordance with the US Pharmacopeia USP). If the digit at the first dropped decimal place is 1, 2, 3 or 4 , then it will be rounded off; if this digit is 5, 6, 7, 8 or 9 , then it will be rounded up. Negative digits will be rounded in accordance with their amount, i.e. away from zero.
Truncate	The result is cut to the number of decimal places defined.
Full precision	The result is used with full accuracy (floating point number either in "single precision" (32 bit) or in "double precision" (64 bit), according to the standard IEEE 754).

[New]

[Delete]

[Edit]

Edit the data of the selected input signal, see following chapter.

15.3.1 Editing the input signal

List of input signals: **Input signal** ► **New / Edit**

Signal name

Designation of the template.

Input signal

Entering the desired bit pattern:

Entering the bit pattern:

- 0 = line inactive
- 1 = line active
- * = retain line status

The input lines are always numbered from right to left, i.e. with the signal *****1 line 0 is expected to be active.



NOTE

We recommend masking lines that are of no interest or for which no defined condition can be predicted with an asterisk (*).

Entry	Bit pattern of exactly 8 characters
Default value	*****

Signal name	Output signal	Function
Meter Mode T	*****0010*	switches the 780/781 pH/Ion meter to temperature measurement and starts the measurement.
Meter Mode U	*****0011*	switches the 780/781 pH/Ion meter to voltage measurement and starts the measurement.
Start device*	*****p****p	starts Devices 1 and 2 (e.g. Titrino, Titrand, etc.) *) The signal is transmitted as a short pulse of 200 ms.
Start device1	*****p	starts Device 1 (e.g. Titrino, Titrand, etc.) *) The signal is transmitted as a short pulse of 200 ms.
Start device2	*****p*****	starts Device 2 (e.g. Titrino, Titrand, etc.) *) The signal is transmitted as a short pulse of 200 ms.
Start Dos*	****p*p*****	starts Dosimat at Devices 1 and 2 (Titrino via "activate"). The signal is transmitted as a short pulse of 200 ms.
Start Dos1	*****p*****	starts Dosimat at Device 1 (Titrino via "activate"). The signal is transmitted as a short pulse of 200 ms.
Start Dos2	****p*****	starts Dosimat at Device 2 (Titrino via "activate"). The signal is transmitted as a short pulse of 200 ms.

*) In the case of the 780/781 pH/Ion meter, a result report is triggered.

[New]

Add a new output signal to the list, see following chapter.

[Delete]

Delete the selected output signal from the list.

[Edit]

Edit the data of the selected output signal, see following chapter.

15.5 Report header

Main dialog: **System ► Templates ► Report header**

In the dialog **Templates / Report header**, you can create a system-specific report header that contains information concerning the laboratory, for example. The report header can contain a maximum of four lines of 46 characters each and will always be printed out in front of the standard report header. The standard report header contains the instrument type, serial number and program version of the software, instrument name and printing date.

Print logo

on | off (Default value: **off**)

If this parameter is activated, then the logo defined below will be printed out in the right margin.

Logo

Selection	Metrohm logo Custom logo
Default value	Metrohm logo

Custom logo

Select this setting if your logo is to be printed out in the report header. How to create a logo is described below.

Applying the custom logo

You can create a graphics file of your own and use it as a logo in the report header in place of the Metrohm logo.



NOTE

Observe the size of the graphics object. The maximum size of your graphics file is 64 kB. We recommend a maximum of 500 × 200 pixels.

Proceed as follows to create your own logo:

- Creating the logo**
 - Create a graphics file with the desired logo in an image processing program.
 - Save the graph as a JPG file under the file name "CustomImage.jpg".



2 Copying the file to an external storage medium

- Copy the file "CustomImage.jpg" to the directory "915". If this structure is not maintained, the file will not be found because the software directly accesses the "915" directory.
- Plug in the external storage medium at the Ti-Touch.

3 Loading the graphics file

- In the **Templates / Report header** dialog, tap on the **[Load]** button.
The file is copied to the internal memory.
- Switch the instrument off and back on again.

4 Applying the custom logo

- In the **Templates / Report header** dialog, activate the **Print logo** option.
- Define **Logo** = **Custom logo**.

Your logo will now be printed in the report header of all reports in the future.

16 Methods


16.1 Creating a new method

Proceed as follows to create a new method:

1 Opening the method table

- In the main dialog, tap on **[Load method]**.
- Tap on **[New method]**.

The method table with the stored templates opens:

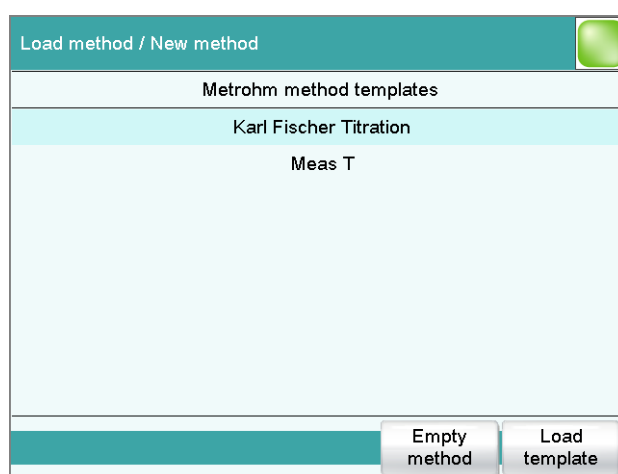

2 Loading the method

- Tap on **[Empty method]**.

or

- Select the desired template and tap on **[Load template]**.

The method is now loaded. **New method** is displayed in the main dialog in the title line.

If a new method has been created, then the individual parameters can be modified with **[Edit parameters]**.

- The method will be saved and the command sequence is displayed.

- Tap on the input field **File name**.
The text editor opens.

- The method will be saved and the command sequence is displayed.

Proceed as follows to load a method:

- In the main dialog, tap on **[Load method]**.
The selection of memory locations is displayed. This selection will be skipped if only the internal memory is accessible.
- Select the memory location where the desired method is saved.
The selection of file groups is displayed (*see chapter 12.1, page 124*). If only one group is available, then the method table will be displayed directly.
- Select the group with the desired method.
- Tap on **[Show files]**.

Internal memory / Show files

Main group

Me2115

Me3901

Me4155

Me4612

Delete

Load

- Select the desired method.

[Insert command]

Insert a new method command. It is inserted before the selected command.

**NOTE**

Not all commands are available for subsequences. Commands which cannot be inserted into subsequences are disabled.

[Delete command]

Delete the selected method command.

[Edit command]

Edit the selected method command.

16.4.1 Inserting a command

A method is comprised of individual commands. When you start a determination, the commands are executed one after the other.

The following constraints apply:


- Maximum of 99 commands
- Maximum of nine commands for titrations and measurements

Inserting a command

To insert a command in a method, proceed as follows:

1 Open the command overview

- In the main dialog, tap on **[Edit parameters]**.
- Tap on the **[Insert command]** button.

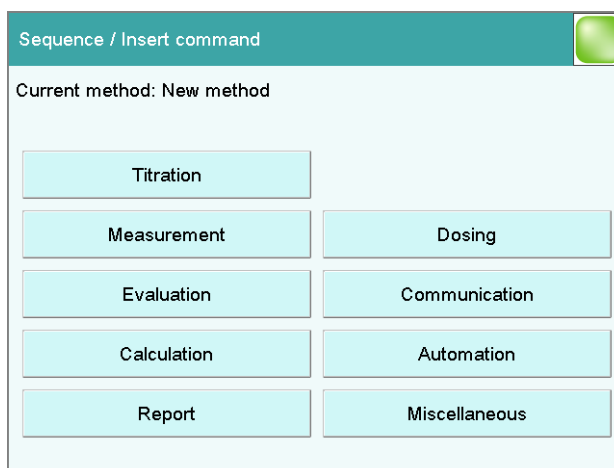

This dialog contains all of the commands which can be inserted into a method run, organized into thematic groups.

2 Select a command group

- Tap on the desired command group.

In the case of **[Report]**, the command is inserted directly into the command list. For all other command groups, the available commands will be displayed (titration / measuring commands, etc.).


3 Inserting a command

Tap on the desired command.

The dialog is closed and the command is displayed in the command list.

16.5 Method options

Main dialog: **Edit parameters** ► **Method options**

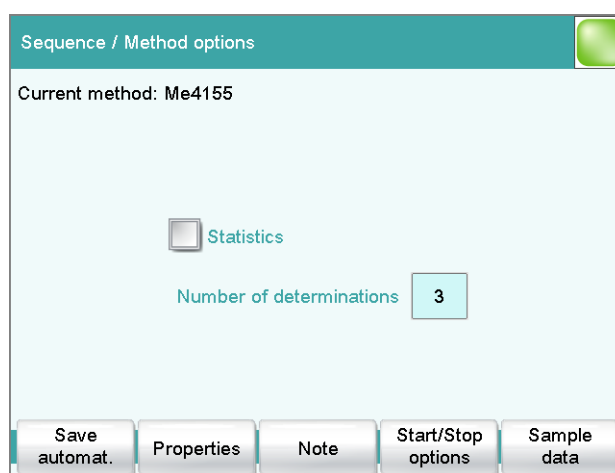

Method options are settings which apply to the method as a whole and not to an individual command.

The following settings are available:

- Activate/deactivate statistics calculations
- Save determination automatically
- Sign a method electronically
- Check the method when the determination starts
- Define a fixed sample size
- etc.

Statistics

on | off (Default value: **off**)

If this parameter is activated, then statistics calculations will be carried out for all of the defined results. The prerequisite here, however, is that a variable for the mean value is defined in the result options (see "Dialog "Edit calculation / Result options"", page 304).

Number of determinations

The number of determinations that are carried out for the statistics calculations.

If an additional determination has to be added to the determination series, because one determination has been incorrect, for example, then this can be accomplished in the statistical overview (see chapter 25, page 229).

Input range	2 to 20
Default value	3

[Save automat.]

Define the settings for automatically saving the determination (see chapter 16.5.6, page 180).

[Properties]

Display the properties of the current method and sign electronically (see chapter 16.5.4, page 178).

[Note]

Enter a note on the method (see chapter 16.5.5, page 179).

[Start/Stop options]

Edit functions which are carried out immediately after the method is started or stopped, see chapter Start options, page 171 and chapter Stop options, page 173.

[Sample data]

Specify method-specific settings for the sample data, e.g. designation of the sample identifications, limits for the sample size (see chapter 16.5.3, page 174).

16.5.1 Start options

Main dialog: **Edit parameters ► Method options ► Start/Stop options ► Start options**

In the dialog **Method options / Start options**, you can configure the functions which are carried out when a determination is started.

Selection	no Selection of the configured sample racks
Default value	no

Increase sample variable automatically

on | **off** (Default value: **on**)

The sample variable describes the current position of the sample on the rack of the Sample Processor. If this option is activated, then the sample variable is automatically increased by 1 at the end of each determination.

This option must be deactivated if the sample variable is to be altered specifically with a **SAMPLE** command.

Request rack position at start

on | **off** (Default value: **off**)

If this option is activated, then the rack position of the first sample to be processed is requested.

16.5.2 Stop options

Main dialog: **Edit parameters** ► **Method options** ► **Start/Stop options** ► **Stop options**

In the dialog **Method options / Stop options**, you can define the actions which are carried out when a method is canceled.

The method can be canceled as follows:

- Manual stop with the fixed key []
- Stop because of an error
- Stop by remote signal via the Control Remote Box

Switch off pumps

on | **off** (Default value: **on**)

If this option is activated, then all of the connected pumps will be switched off.

Switch off stirrers

on | **off** (Default value: **on**)

If this option is activated, then all of the connected stirrers will be switched off.

Set remote lines

Selection of the signal out of the templates or entering the required bit pattern. Templates are defined under **System** ► **Templates** ► **Output lines**.



NOTE

A line set active is not being reset automatically, not even at the end of the determination.

Entering the bit pattern:

- 0 = line inactive
- 1 = line active
- * = retain line status
- p = set pulse

The output lines are always numbered from right to left, i.e. with the signal *******1** line 0 is set. With a pulse, the length is set to 200 ms. If you wish to set pulses with other lengths, you have to define a corresponding template.

Entry	Bit pattern containing exactly 14 characters or a max. of 24 characters for the name of the template
Default value	*****
Selection	Selection of the templates defined

16.5.3 Sample data

Main dialog: **Edit parameters** ► **Method options** ► **Sample data**

You can define the following method-specific settings in the dialog

Method options / Sample data:

- Modify the designation for the sample identifications
- Define a fixed sample size
- Define limits for the sample size

If you define limits, these will be monitored at the start and at the end of the determination.

Below you will find an example demonstrating the effects on the main dialog of the settings in this dialog (see *"Modifying sample data"*, page 176).

Identification 1, Identification 2

on | off (Default value: **on**)

If this option is activated, then the input field for the sample identification in the main dialog will be displayed.

Designation

A designation of your own choosing can be defined for the input field for each method.

Entry	16 characters maximum
-------	-----------------------

Fixed sample size

on | off (Default value: **off**)

Activate this parameter if the same sample amount is always to be used for all determinations. If this is case, then you can define the sample size here. It will be displayed in the main dialog, but it can no longer be changed there.

Sample size

This parameter can only be edited when **Fixed sample size** is activated.

Input range	–999999999 to 999999999
Default value	1.0

Unit of sample size.

Selection	g mg µg mL µL pieces
Default value	g



NOTE

You can also use the parameter **Fixed sample size** if you would like to establish a default for the sample size.

Example: For a particular determination, you require 10 mL of sample. Now and again, however, you have too little sample available and you must perform the determination with 5 mL, for instance. In this case, proceed as follows:

1. Activate the parameter **Fixed sample size**.
2. Define the value and unit for the sample size, in the above example, for instance, **[10 mL]**.
3. Deactivate the parameter **Fixed sample size** again.
4. Save method.

⇒ Each time you load this method, 10 mL is entered as the sample size, but in the example mentioned you can enter the lesser quantity of 5 mL.

[Sample size limits]

This button is accessible only if **Fixed sample size** is deactivated.

Define the limit values for the sample size (see "Defining limit values for the sample size", page 177).

Modifying sample data

This instruction is for the purpose of clarifying the settings which are possible in this dialog. The following settings should be made:

- Change the designation for the input field of the first sample identification.
- Hide the input field for the second sample identification.
- Define a fixed sample size.

Proceed as follows:

1 Change the designation for the input field

- Tap on the input field **Identification 1**.
The text editor opens.
- Change the designation in **Batch** and confirm with **[OK]**.
- Deactivate the parameter **Identification 2**.

2 Define a fixed sample size

- Activate the parameter **Fixed sample size**.
- Enter the value **10 mL** as **Sample size**.

3 Displaying the main dialog

Tap on the fixed key [].

The data is saved and the adjusted main dialog appears.

New method

09:08:13

User



Batch



Sample size



System

Load method

Control/ Logout

Edit parameters

Results

Only Identification 1 is shown with the title "Batch." The fixed sample size "10 mL" is shown, but cannot be edited.

Defining limit values for the sample size

In the dialog **Sample data / Sample size limits**, you can define the limit values for the sample size. These values are not monitored during sample data input, but rather:

- when the determination is started.
- when the determination is finished.
- upon automatic sample data query via the command **REQUEST**.
- when the determination is being recalculated.

An entry is made in the message list automatically when these limits are infringed (see chapter 24.2, page 216).

Monitoring sample size limits

on | off (Default value: **off**)

If this parameter is activated, then the sample size will be monitored.

Lower limit

Input range	–999999999 to 999999999
Default value	0

Upper limit

Input range	–999999999 to 999999999
Default value	999999999

Display message

on | off (Default value: **off**)

If this parameter is activated, a message is displayed in case the limits are infringed. You can select whether you want to continue with the determination or cancel the run.

If this parameter is deactivated, then the message will be entered only in the message list of the determination.

Timeout

The message is displayed during this time. Afterwards, the determination is continued automatically.

Input range	0 to 999999 s
Default value	30 s
Selection	off

off

The determination will not be continued until after the message has been confirmed.

- **Reviewed on**
Date and time at which the method has been signed at level 1.
- **Reason**
Reason for the signing of the method.
- **Comment**
Comment for the signing of the method.

The following information is only shown if the method has been signed at level 2:

- **Released by**
User who signed the method at level 2.
- **Released on**
Date and time at which the method has been signed at level 2.
- **Reason**
Reason for the signing of the method.
- **Comment**
Comment for the signing of the method.

[Delete signature]

Delete all of the signatures for the method or the determination. This button is inactive if level 2 has not yet been signed (*see chapter 17.2, page 183*).

[Sign]

Sign the method or determination. This button is active only when working with activated login function and password (*see chapter 17.1, page 182*).

16.5.5 Note

Main dialog: **Edit parameters ► Method options ► Note**

In this dialog, you can enter a short text, e.g. important information about carrying out the determination.

[Display options]

Definition when the note is being displayed.

Automatically after loading the method

on | off (Default value: **off**)

If this parameter is activated, the note is displayed when loading the method. It can otherwise only be read in this dialog.

16.5.6 Saving a determination automatically

Main dialog: **Edit parameters** ► **Method options** ► **Save automat.**

In the dialog **Method options / Save automatically**, you can specify whether the determination is to be saved automatically and/or whether a PC/LIMS report is to be created.

Save determination automatically

on | off (Default value: **off**)

If this parameter is activated, then the determination will be saved automatically.

Memory

Memory location where the determination will be saved.

Selection	External memory 1 External memory 2
	Shared memory
Default value	External memory 1

Shared memory

The determination will be saved in a shared directory on the network. The shared directory is selected in the **Edit device / Shared memory** dialog (see chapter 11.3.3, page 86).

Group

Directory where the determination will be saved.

Entry	32 characters maximum
Selection	Selection of available directories empty
Default value	empty

File name

File name of the determination. The file name is always extended to include the date and the time of day (YYYYMMDD-hhmmss) in order to ensure that the file names are unique for all determinations.

Entry	max. 16 characters
Selection	Identification 1 Identification 2 Method
Default value	Identification 1

Identification 1

The first 16 characters of the text which was entered in the main dialog for Identification 1 + YYYYMMDD-hhmmss.

Identification 2

The first 16 characters of the text which was entered in the main dialog for Identification 2 + YYYYMMDD-hhmmss.

Method

The first 16 characters of the text, of the method name +
YYYYMMDD-hhmmss.

Write protection

on | off (Default value: **on**)

If this parameter is activated, then the file cannot be saved, deleted, or renamed. This is only an internal write protection feature and is independent of the write protection of the operating system on your computer. This write protection feature protects the saved determination data against accidental modifications or modifications by unauthorized persons.

Create PC/LIMS report

on | off (Default value: **off**)

If this parameter is activated, then a machine-readable report with all of the important data for the determination is created, a so-called PC/LIMS report.

The settings are implemented in the device manager (*see chapter 11.3.2, page 84*).

When you save the PC/LIMS report as a file, the file name will be generated automatically as follows: **PC_LIMS_Report**-first 16 characters of the value of the parameter **File name-Date-Time.txt**. If the value of the parameter **File name** contains one or more of the following characters, then these will be automatically converted in the file name into the character "_": / \ : * ? " < > |.

Reason

Selection of the reasons for the signature. The list with reasons is compiled in the user administration (*see chapter 7.3.6, page 48*).

Selection	Selection of configured reasons
-----------	--

Comment

Entry	24 characters maximum
-------	------------------------------

[Change password]

Change the current password.

[Sign]

Sign the method or determination with the data entered above.

17.2 Deleting electronic signatures

If a released method must be revised, e.g. in order to adapt to altered conditions, then the signatures must be deleted. The method can otherwise not be modified. In the case of released determinations, the signatures can also be deleted. Signatures, however, cannot be deleted until after the method/determination has been signed at Level 2.

Signatures are deleted automatically in the following cases:

- A method signed at Level 1 is being altered.
- A determination signed at Level 1 is being recalculated.

User

Identification of the user who deletes the signatures.

Password

Password of the user.

Reason

Selection of the reasons for the deletion of the signature. The list with reasons is compiled in the user administration (*see chapter 7.3.6, page 48*).

Selection	Selection of configured reasons
-----------	--

Comment


Entry	24 characters maximum
-------	------------------------------

[Delete signature]

Delete all the signatures for the method or the determination. The method is assigned the status **saved**.

18 Control

Main dialog: **Control**

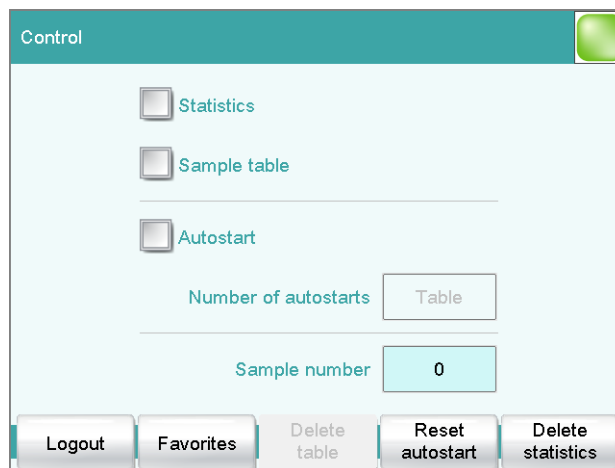

In the dialog **Control**, the settings for the execution of a single determination or of one sample series are defined.

If you work with the login function activated, you can log off the system in this dialog. The login dialog will then be shown immediately.

Statistics

on | off (Default value: **off**)

With this option, you can enable or disable the statistics calculations for individual determinations (see *"Statistics"*, page 186).



NOTE

The parameter **Statistics** in the method options will not be changed as a result. This setting here will be adjusted automatically to match the one in the method options.

Sample table

on | off (Default value: **off**)

If this option is activated, the sample data for a sample series can be entered in a table (see *chapter 21*, page 195).

Autostart

on | off (Default value: **off**)

If this option is activated, a new determination is started automatically at the end of a determination. This continues until the number specified has been reached (siehe **Number of autostarts**).



NOTE

If you deactivate this option while a determination is running, then that determination will be ended, although the next one will however also not be started. This can be used to cancel a sample series.

Number of autostarts

Number of automatic starts.

Input range	1 to 9999
Selection	Table
Default value	Table

Table

The number of automatic starts corresponds to the number of samples in the sample table.

Sample number

The sample number is primarily used for documentation. Each time a determination is started, it is increased by one. The sample number can, however, also be entered manually. Each time the Ti-Touch is switched on, the value is reset to zero.

Input range	0 to 9999
Default value	0

[Logout]

User logout (visible only when work is done with activated login function).

[Favorites]

Save methods/sample tables as a favorite (*see chapter 19, page 188*). Depending on whether you are working with or without login function, user-specific or common favorites can be created. On the main dialog, a button is created for each favorite. This way, a method or a sample table can be started with one single mouse click.

[Delete table]

Delete the sample table completely.

[Reset autostart]

Reset the value of the autostart counter to zero.


Continuing the sample series

19 Favorites

Main dialog: **Control ► Favorites**

Main dialog: **System ► System settings ► User admin. ► Edit ► Favorites**

This chapter describes how you can create favorites. Favorites are comparable to the favorites/bookmarks in your Internet browser. In the main dialog, a button is created for each favorite. By pressing a single key, you can trigger an action without having to navigate through different dialogs.

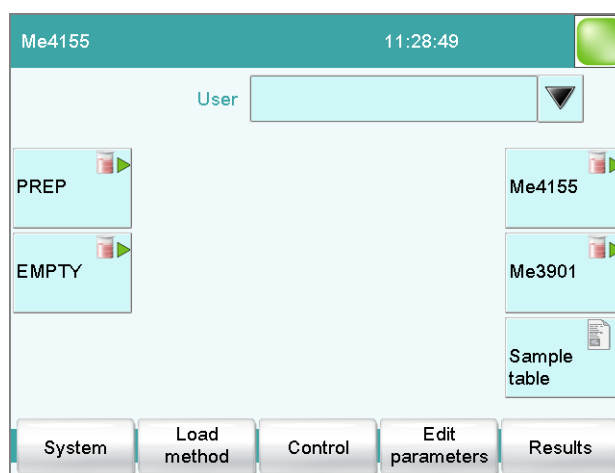

The following objects can be saved as favorites:

- Methods
- Sample tables

User-specific or common favorites can be created, depending on whether you are working with or without the login function.

- Operation with login function:
If you are working with the login function activated, **user-specific favorites** can be created. These can only be used by a certain user. User-specific favorites are created as follows:
 - in the user administration by a user with administrator rights (**System ► System settings ► User admin. ► Edit**).
 - in the dialog **Control** for the user who is logged in.
- Operation without login function:
If you are working without the login function, **common favorites** can be created. These favorites are available for all users. Common favorites are created in the **Control** dialog.

19.1 Creating favorites

Control / Favorites	
Name	Position
PREP	1
EMPTY	6
Me4155	5
Me3901	10
New Delete Edit	

The list of favorites gives the designation and position of the button on the main dialog of each favorite configured.

[New]

Add a new favorite to the list, see following chapter.

[Delete]

Delete the selected favorite from the list.

[Edit]

Edit the properties of the selected favorite, see following chapter.

19.1.1 Editing favorites

Favorites / Edit	
Position	5 ▼
Name	Me4155
Type	Method ▼
Memory	Internal memory ▼
Method / Sample table	Me4155 ▼
More options	

on | off (Default value: **off**)

Number of autostarts

Number of automatic starts.

Input range	1 to 9999
Default value	1
Selection	Table

The number of automatic starts corresponds to the number of samples in the sample table.

Sample identification. The sample identification can be used in calculations as the variable **C12**. If you have defined a sample identification list (see *chapter 15.1.1, page 151*) then the entries can be selected here.

Sample size. The value of the sample size can be used in calculations as the variable **C00**.



Me2115



01 REQUEST Data request

Identification 1

Sample size

10

mL

Continue

To accomplish this, you must insert and configure the command **REQUEST** in the method (see chapter 29.9.3, page 346).

If the **Hold sequence** parameter is activated, then the run will be paused and must be continued with **[Continue]** after the sample data has been entered. If the **Hold sequence** parameter is deactivated, then the determination will continue in the background until the measured data is required for further processing. This dialog will be displayed until the entering of the sample data is confirmed with **[Continue]**, even if the titration/measurement is already completed. This ensures that the sample data is available for calculations.

If a fixed sample size has been defined in the method (see *chapter 16.5.3, page 174*), then it will be displayed at the time of the entry but it will not be editable.

21 Sample table

21.1 General

Main dialog: **Sample table**

The sample table is a table in which the sample data for up to 999 samples can be entered. If you are processing larger-sized sample series, it makes sense to enter the sample data in this table.


Activating the sample table

Proceed as follows to enable the sample table:

1 Opening the control dialog

- In the main dialog, tap on **[Control]**.


The control dialog is displayed:

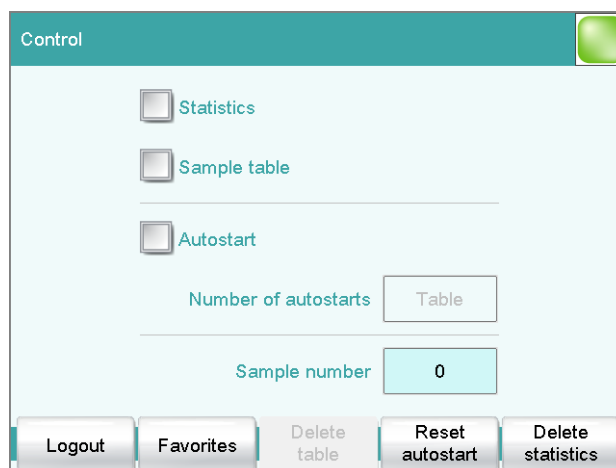

2 Activating the sample table

- Activate the **Sample table** option.
- Tap on the fixed key **[]**.

The main dialog is displayed:

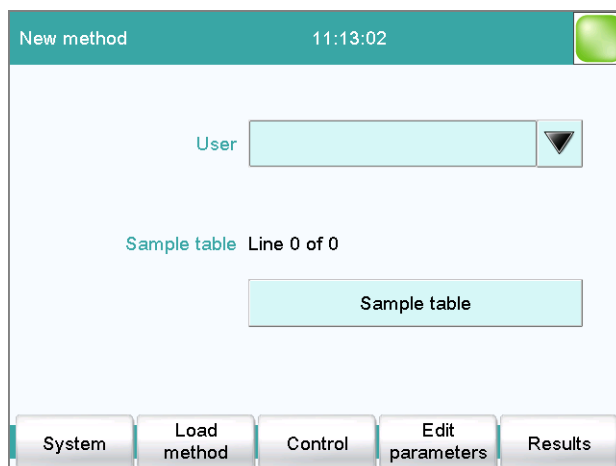

The display will show the number of determinations already carried out and the total number of sample lines containing data. The sample table is still empty in this example.

The sample table contains numbered lines; each line represents one determination. In the standard settings, the first sample identification and the sample size of each sample are displayed. In the properties, you define the columns to be displayed (*see chapter 21.3, page 200*).

Sample table		
No.	Identification 1	Sample size
1	#2370015	1.0 g
2	#2370015	1.0 g
3	#2370016	1.0 g
4	#2370016	1.0 g
5	...	

Load/Save

Properties

Insert line

Delete

Edit

[Load/ Save]

Load a saved sample table (see "Loading a sample table", page 198) or save the current sample table (see "Saving a sample table", page 197).

[Properties]

Edit the properties of the sample table (see chapter 21.3, page 200).

[Insert line]

Insert a new line above the line selected.

[Delete]

Delete the selected line from the sample table.

[Edit]

Edit the data of the selected line (see chapter 21.2, page 198).

Saving a sample table



NOTE

Sample tables can only be saved on external storage media or on a shared file system.

Proceed as follows to save a sample table:

1 Opening the dialog

- In the **Sample table** dialog, tap on the **[Load/ Save]** button.

The **Sample table / Load/Save** dialog is displayed.

If no external storage medium is plugged in or if no file system is shared, then **[Save]** is disabled.

2 Modify/apply a name

- Tap on **[Save]**.
A name is suggested for new sample tables. If the sample table has already been saved, the previous name will be displayed.

Applying the name:

- Tap on **[Save]**.

The sample table is saved.

Entering a new name:

- Tap on the input field **File name**.
The text editor opens.
- Enter a new file name (max. 32 characters) and confirm with **[OK]**.
- Tap on **[Save]**.

The sample table is saved.

Loading a sample table

Proceed as follows to load a sample table:

1 Display the list of saved sample tables

- In the **Sample table** dialog, tap on the **[Load/ Save]** button. The **Sample table / Load/Save** dialog is displayed. If no external storage medium is plugged in or if no file system is shared, then **[Load]** is disabled.
- Tap on **[Load]**. The selection of memory locations is displayed. This selection will be skipped if only one memory can be accessed.
- Select the memory location where the desired sample table is stored. The selection of file groups is displayed (*see chapter 12.1, page 124*). If only one group is available, then the list of the saved sample tables will be displayed directly.
- Select the group with the desired sample table.
- Tap on **[Show files]**.

The list with the saved sample table is opened.

2 Select the sample table

- Select the desired sample table.

3 Loading a sample table

- Tap on **[Load]**.

The sample table is now loaded.

21.2 Editing sample data

Sample table / Edit

Line number:

Method:

Identification 1:

Identification 2:

Sample size:

At the very top you will see the line number of the selected line. In this example, the sample data of the first line is displayed. You can scroll between individual data sets with the keys **[-]** and **[+]**.

Method

Method used for processing the sample.

If you use the sample assignment table (*see chapter 15.1.2, page 152*), then the method can no longer be defined here.

Selection	Selection of stored methods empty
Default value	empty

empty

The currently loaded method is used.

Identification 1

Sample identification. The sample identification can be used in calculations as the variable **CI1**. If you have defined a sample identification list (*see chapter 15.1.1, page 151*) then the entries can be selected here.

Entry	24 characters maximum
Default value	empty

Identification 2

Sample identification. The sample identification can be used in calculations as the variable **CI2**. If you have defined a sample identification list (*see chapter 15.1.1, page 151*) then the entries can be selected here.

Entry	24 characters maximum
Default value	empty

Sample size

Sample size. The value of the sample size can be used in calculations as the variable **C00**.

Input range	-99999999 to 999999999
Default value	1.0

Unit of sample size.

Entry	6 characters maximum
Selection	g mg µg mL µL pieces
Default value	g

Unit defined in method:

The unit defined in the method is ignored. The unit entered in the sample table is used.

Method

Only the assigned method is displayed for each sample.

Identification 1

The first sample identification and the sample size are displayed for each sample.

Identification 2

The second sample identification and the sample size are displayed for each sample.

Number of lines

Maximum number of lines which can be present in the sample table.

Input range	2 to 999
Default value	99

Lock editing

on | off (Default value: **off**)

If this option is activated, then the sample data can no longer be altered.

Display options

Configure the sample data to be displayed in the dialog **Sample table / Edit** (see *"Properties / Display options" dialog*, page 201).

Editing options

Specify the settings for the processing of the sample data (see *"Properties / Editing options" dialog*, page 201).

Identifications

Define the designations for the input fields of Identification 1 and Identification 2 (see *"Dialog "Properties / Identifications"*, page 202).

"Properties / Display options" dialog

Sample table: **Properties ► Display options**

In this dialog, you can configure which of the following sample data can be edited in the sample table:

- **Method**
- **Identification 1**
- **Identification 2**
- **Sample size**

"Properties / Editing options" dialog

Sample table: **Properties ► Editing options**

In this dialog, you can define the settings so that the input of the sample data becomes simpler and faster.



Properties / Editing options

End line with Sample size

Copy automatically into next line:

☒ Method ☒ Identification 2

☒ Identification 1 ☒ Sample size

End line with

Selection as to with which input the editing dialog for the next sample is automatically displayed.

Selection	Sample size Identification 1 Identification 2 manual
Default value	Sample size

manual

The editing dialog for the next sample can be displayed with the **[New sample]** button in the editing dialog.

If, for example, you have to enter similar data for each sample, then you can automatically copy the data which is identical for each sample into the next line:

- Method
- Identification 1
- Identification 2
- Sample size

Dialog "Properties / Identifications"

Sample table: **Properties** ► **Identifications**

In the dialog **Properties / Identifications**, you can alter the designations for the sample identifications in accordance with specific methods.

22 Determination run

This chapter provides a description of how to carry out a determination, i.e. how to process a single sample or a sample series. The current sample data is used for calculations.

22.1 Carrying out a single KF determination

When you carry out a determination, you can enter the sample data in various ways (see chapter 20, page 192). The following instructions describe how to enter the sample data in the main dialog.

To carry out a single determination, proceed as follows:

1 Loading the method

- See Chapter "Loading a method", Page 167.

2 Starting conditioning

Tap on the fixed key [].

Conditioning starts. **Conditioning...** is displayed until the endpoint is reached. The working medium is titrated to the endpoint. This is indicated by **Conditioning OK**. The status is kept stable.

3 Starting a determination

- If **Conditioning OK** is displayed, tap on the [] fixed key. Conditioning is stopped. The request for adding the sample will be displayed.
- Tap on the fixed key [].

4 Adding sample



NOTE

Calculate the amount of the sample in such a way that a titrant consumption of 10...90% of the cylinder volume will result.

- Add sample into the titration cell. Afterward, the request for the sample size is displayed.

5 Entering the sample size

- Tap on the input field **Sample size** in the dialog. The number editor opens.
- Enter the sample size and confirm with **[OK]**.

The KF titration is already executed in the background during the entry of the sample size.

After the confirmation of the sample size, the live curve of the titration is displayed.

Sequence of the determination

- After a determination is started, a check is made as to
 - whether all of the devices required are connected and ready for operation.
 - whether all titrants and sensors are present.
 - whether the monitoring intervals for titrants, sensors and common variables are still valid.
 - if a sample assignment table is activated (*see chapter 15.1.2, page 152*):
whether the entered sample identification is defined in the sample assignment table. If this is the case, then the assigned method will be loaded automatically and the actual determination will be started.

22.2 Performing a sample series

You can use the following functions with sample series:

- **Statistics for multiple determinations**

The statistics calculations are enabled in the method options (see chapter 16.5, page 170).

- **Sample table**

You can enter the sample data for an entire sample series in the sample table (see *chapter 21, page 195*).

Sequence of the sample series

- After a determination is started, a check is made as to
 - whether all of the devices required are connected and ready for operation.
 - whether all titrants and sensors are present.
 -
 - whether the monitoring intervals for titrants, sensors and common variables are still valid.
 -
 - if a sample assignment table is activated (*see chapter 15.1.2, page 152*):
 - whether the entered sample identification appears in the sample assignment table. If this is the case, then the assigned method will be loaded automatically and the actual determination will be started. If a different method has been entered in the sample table, then it will be ignored.
- At the end of the determination,
 - the line in the sample table is deleted.
 - the next determination is started (if you are working with the autostart function enabled (*see chapter 18, page 184*)).

22.3 Canceling determinations manually

You can cancel a determination at any time with the fixed key []. The command being carried out at that moment is canceled and only calculations and reports are executed or created.

If you are working with the autostart function activated (*see chapter 18, page 184*), the entire series will be canceled.



NOTE

If you do not wish to cancel the entire determination, but only the current command, use the function **[Cancel command]** in the live display.

The subsequent method commands are carried out normally.

23.2 Editing the sample table while a determination is running

You can insert new lines or delete existing ones or edit sample data while a determination is running.



NOTE

To ensure that no problems arise during the run and that the current data is always available for the calculation, we recommend that the input dialogs always be closed.

Editing the sample table

Proceed as follows to edit the sample table:

1 Display the main dialog

- Tap on the fixed key [].

The main dialog is displayed. The determination continues to run in the background.

2 Select the sample data

- Tap on **Sample table**.
The sample table is displayed. The first line always contains the sample data of the running determination.
- Select the desired line.
- Tap on the button **[Edit]**.

3 Edit the sample data

- Edit the sample data and apply with **[OK]**.



NOTE

In addition to the sample data, the method can also be modified, except in cases where the determination is running.


4 Display the live dialog

- Tap on the fixed key [].
The main dialog is displayed.
- Tap on **[Live display]**.

The live dialog is displayed once again.

23.3 Live display

Main dialog: **Live display**

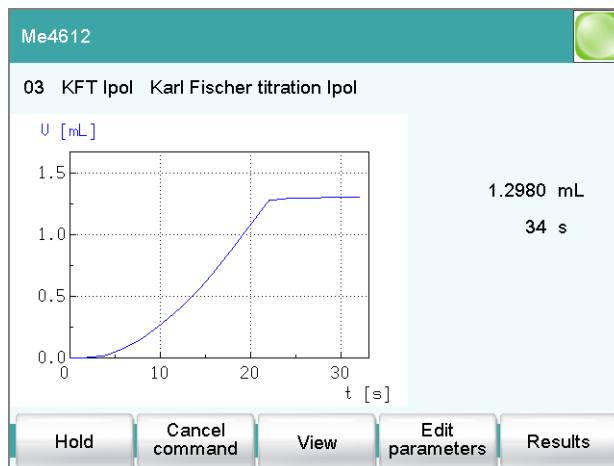


The current method and the system status are displayed in the title line. The currently running command is displayed directly underneath the title line.

The following functions are available to you while a determination is running:

- **[Hold]**
Pause the determination. The label changes to **[Continue]**. This can be used to continue the run.
- **[Cancel command]**
Cancel the current method command. The next command is started.
- **[View]** (only for titrations/measurements)
To modify the settings for the curve and measured value display, see below.
- **[Sample data]** (visible only if a titration or a measurement is not currently running)
Switch to the main dialog in order to edit the sample data.
- **[Edit parameters]**
Edit the method options and live parameters (*see chapter 23.4, page 210*).
- **[Results]**
Display the results dialog (*see chapter 24, page 212*).

Dialog "View"

Main dialog: **Live display** ► **View**

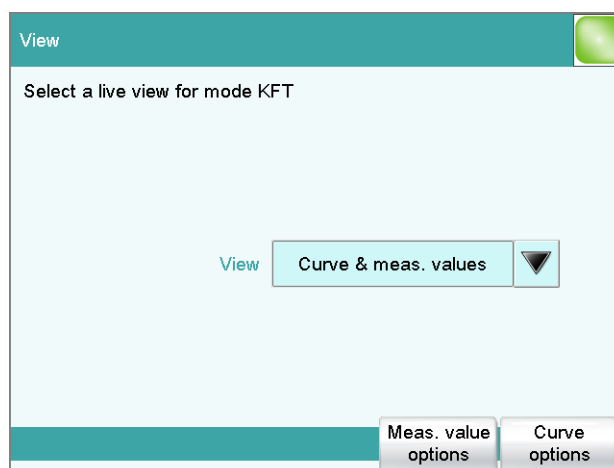

View

Type of curve display.

Selection	Curve Curve & meas. values Measured values Curve & sample data
Default value	Curve & meas. values

Curve

Only the curve is displayed.

Curve & meas. values

The measured values are displayed to the right next to the curve.

Measured values

Only the measured values are displayed.

Curve & sample data

The sample data is displayed to the right next to the curve.

[Meas. value options]

Define settings for the display of the measured values.

A maximum of three different values can be displayed at the same time. The selection depends on the titration/measuring mode. The settings are saved separately for each mode.

[Curve options]

Define the settings for the curve display (see "Curve options", page 226). The settings are saved separately for each titration/measuring mode and apply not only for the live display but also for the curve display in the results dialog.

2 Edit the method parameters

- Change the desired parameters accordingly.

3 Display the live dialog


- Tap on the fixed key [].
The main dialog is displayed.
- Tap on **[Live display]**.

The live dialog is displayed once again.

24 Results and more determination data

Main dialog: **Results**

Detailed specifications concerning the current determination are displayed in the **Results** dialog:

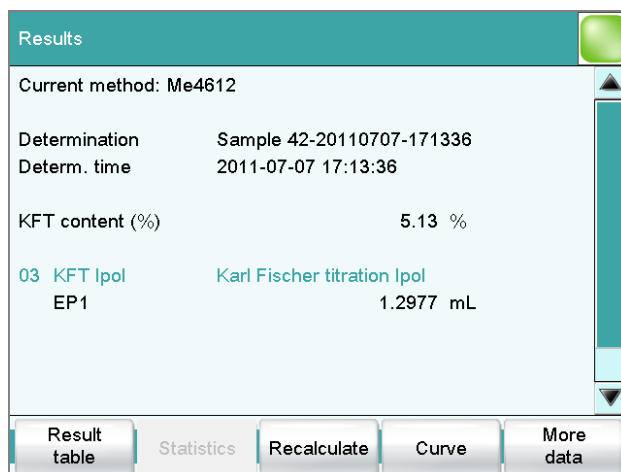

- **Current method**
Method with which the determination was carried out.
- **Determination**
Name of the determination. Displayed only if the determination has been saved.
- **Determination time**
Date and time at which the determination was carried out.
- **Results**
Calculated results. If monitoring of the result limits is activated, then the result will be highlighted in color:
 - Green lettering: The result is within the limits
 - Red lettering: The result is outside the limits
- **Commands**
List of the data-generating commands (titrations, measurements, etc.). The most important data is displayed for each command.

[Result table]

Open result table (*see chapter 26, page 234*). The result table is a table in which the results for several determinations can be saved.

[Statistics]

Display the statistical overview of a determination series (*see chapter 25, page 229*).

[Recalculate]

Recalculate the determination. The procedure will be executed immediately. Recalculation cannot be undone (*see chapter 24.8, page 227*).

[Sign]

Sign the current determination.

If the method contains a command **SIGN**, then this button will be displayed instead of **[Recalculate]**. The determination can be signed directly on level 1 (with a second command also on level 2).


[Curve]

Display the curve of the current determination (see chapter 24.7, page 225).

[More data]

Display detailed specifications for the determination (determination properties, measuring point list, endpoint details, etc.), save determination (*see chapter 24.1, page 213*).

24.1 More determination data

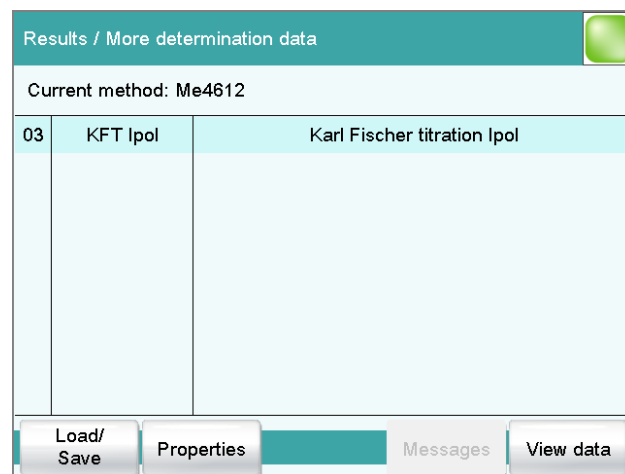Results dialog: **More data**

All of the commands for which you can display other data (titrations, measurements, evaluations, etc.) are listed in the **Results / More determination data** dialog.

[Load/ Save]

Load saved determinations(see *chapter 24.5, page 219*) and save current determination (see *chapter 24.6, page 224*).

[Properties]

Display the properties of the current determination (see *chapter 24.4*, *page 217*).



[Local com. var.]

(available only with loaded determinations)

Display the list of the common variables defined at the time of the determination (see *chapter 24.3, page 216*).

[Messages]

Display list with all of the messages which occurred during the determination (see chapter 24.2, page 216).

[\[View data\]](#)

Display all of the details for the selected command.

24.1.1 Details

Results dialog: **More data** ► **View data**

More determination data / View data		
03 KFT Ipol Karl Fischer titration Ipol		1M
Stop criterion	Regular stop	MTS
Titration duration	27.8 s	MCD
Start volume	0.0000 mL	MSV
End volume	1.2975 mL	MCV
Temp. measurement	manual	MTM
Initial measured value	636.7 mV	MIM
Initial temperature	25.0 °C	MIT
Start duration	1.2 s	MSD
Start measured value	636.7 mV	MSM
Start temperature	25.0 °C	MST
End measured value	247.3 mV	MCM
Endpoint details		Measuring point list

All of the variables generated by the determination are listed in the **More determination data / View data** dialog. These can be used in calculations.

[Titrant data]

(available only with loaded determinations)

Display titrant data which was valid at the time of the determination (see "Titrant data", page 215).

[Endpoint details]

Display details for each endpoint found (see "Endpoint details", page 215).

[Measuring point list]

Display measuring point list of titrations and measurements (see "*Measuring point list*", page 215).

Measuring point list

Results dialog: **More data ► View data ► Measuring point list**

Results dialog: **Curve ► Display curve ► Measuring point list**

View data / Measuring point list			
03 KFT Ipol Karl Fischer-Titration Ipol			1M
Time [s]	Volume [mL]	Meas.value [mV]	
0.0	0.00000	636.1	
2.0	0.00500	636.0	
4.0	0.02125	635.9	
6.0	0.07450	634.7	
8.0	0.15625	633.2	
10.0	0.26300	630.7	
12.0	0.39125	627.7	
			View

The individual measured values are listed in the **View data / Measuring point list** dialog. You can configure which columns are displayed (**[View]** button). These settings are saved separately for each mode (KFT, MEAS, etc.).

Column 1 / Column 2 / Column 3

Quantity displayed in the columns 1, 2 and 3. Both the selection of the quantities and the default setting depend on the mode.

Selection	Volume	Measured value	Time	Temperature	Volume drift	Signal drift	Index
-----------	--------	----------------	------	-------------	--------------	--------------	-------

Endpoint details

Results dialog: **More data ► View data ► Endpoint details**

Results dialog: **Curve ► Display curve ► Endpoint details**

Detailed information concerning each equivalence point or endpoint determined is displayed in the **View data / Endpoint details** dialog.

The volume, the measured value, the time and the temperature are specified for each equivalence point or endpoint determined.

Titrant data

Results dialog: **More data ► View data ► Titrant data**

The titrant data which was current at the time the determination was performed is displayed in the **View data / Titrant data** dialog.



NOTE

Recalculation:

The titer of the titrant can be modified retroactively. The values displayed here are used for the recalculation of a loaded determination. The values saved in the dialog **System / Titrants** are not changed as a result.

24.2 Messages

Results dialog: **More data** ► **Messages**

If messages occurred during a determination, they will be displayed in the **More determination data / Messages** dialog. The following data for the message is indicated under **[Details]**:

- **Message number**
Unique identification number.
- **Time**
Moment that the message occurred.
- **Message**
Message text.

24.3 Local common variables

Results dialog: **More data** ► **Local com. var.**

The list of common variables as it existed at the time of the determination is displayed in the **More determination data / Local common variables** dialog.



NOTE

Recalculation:

The common variables can be modified retroactively. The values displayed here are used for the recalculation of a loaded determination. The list of common variables in the **System / Common variables** dialog is not changed as a result.

24.4 Determination properties

Results dialog: **More data ► Properties**

Detailed information concerning the determination is displayed in the **More determination data / Properties** dialog:

- **User**
User who carried out the determination.
- **Sample number**
Number which was entered in the dialog **Control** at the start of the determination.
- **Start counter**
This value, just like the sample number, is increased by one each time a determination starts. In contrast to the sample number, this value cannot be modified by the user.
- **Duration determ.**
Duration of the entire determination.
- **Method status**
Current status of the method.
 - **new**
The method has been newly created and not yet saved.
 - **saved**
The current method has been saved.
 - **modified**
The current method has been modified but not saved yet (e.g. live modification).
 - **reviewed**
The current method has been signed at level 1.
 - **released**
The current method has been signed at level 2.
- **Determination**
Name of the determination. Displayed only with saved determinations.
- **End of determin.**
Way in which the determination was ended.
 - **Regular without errors**
The determination was automatically ended at the end of the method.
 - **Manual stop**
The determination has been canceled with the fixed key [].
 - **Error**
The determination has been canceled due to an error.

- **Status of deter.**
 - **original**
The determination has not been recalculated.
 - **recalculated**
The determination has been recalculated.
 - **loaded**
The determination has been loaded.
 - **saved**
The determination has been saved. The version number indicates how many times the file has been saved.
 - **reviewed**
The determination has been signed at level 1.
 - **released**
The determination has been signed at level 2.

The following data is displayed only if the option **Recalculating determination** under **Login options / Modification options** is activated:

- **Recalculated by**
User who recalculated the determination.
- **Recalculated on**
Date and time at which the determination was recalculated.
- **Reason**
Reason for the recalculation of the determination.
- **Comment**
Comment for the recalculation of the determination.

The following information is displayed only if the determination has been signed at level 1:

- **Reviewed by**
User who signed the determination at level 1.
- **Reviewed on**
Date and time at which the determination was signed at level 1.
- **Reason**
Reason for the signing of the determination.
- **Comment**
Comment for the signing of the determination.

The following information is displayed only if the determination has been signed at level 2:

- **Released by**
User who signed the determination at level 2.
- **Released on**
Date and time at which the determination was signed at level 2.
- **Reason**
Reason for the signing of the determination.

- **Comment**

Comment for the signing of the determination.

- **Used devices**

Hardware that was used for the determination.

- Ti-Touch with program version of the software and the serial number.
- Measuring inputs with ADC type und serial number.
- MSB connector with peripheral devices connected to it (dosing device, stirrer).

[Delete signature]

Delete all of the signatures for the method or the determination. This button is inactive if level 2 has not yet been signed (*see chapter 17.2, page 183*).

[Sign]

Sign the method or determination. This button is active only when working with activated login function and password (*see chapter 17.1, page 182*).


24.5 Loading a determination

Proceed as follows to load a determination:

1 Opening the determination list

- In the results dialog, tap on **[More data]**.
- Tap on **[Load/ Save]**.
- Tap on **[Load]**.
The selection of memory locations is displayed. This selection will be skipped if only one external storage medium is available.
- Select the memory location where the desired determination is saved.
The selection of file groups is displayed (*see chapter 12.1, page 124*). If only one group is available, then the determination list will be displayed directly.
- Select the group with the desired determination.
- Tap on **[Show files]**.

The determination list with the stored determination opens:

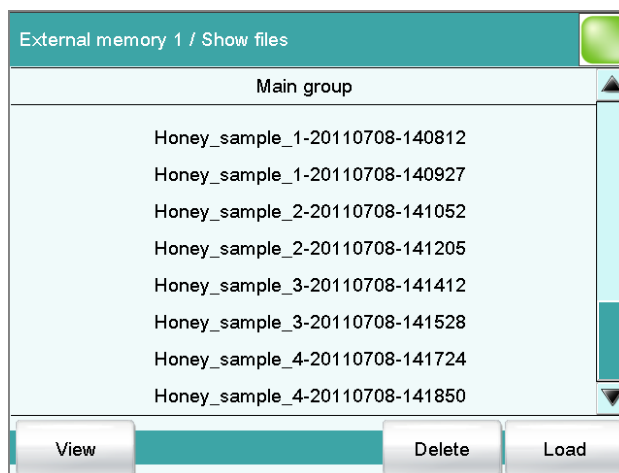

2 Selecting a determination

- Select the desired determination.

3 Loading the determination

- Tap on **[Load]**.

The determination is now loaded. The properties are displayed in the results dialog.



NOTE

When you load a determination, the method with which the determination was carried out will be loaded automatically.

For this reason, save any modifications which may have been made to the method currently loaded before you load a determination.

24.5.1 Determination list

Results dialog: **More data ► Load/ Save ► Load ► Show files ► View**

Show files / View			
Sorted by Determ. time			27
No.	Identification 1	Result 1	
20	Honey_sample_1	16.72 %	 
21	Honey_sample_1	16.89 %	
22	Honey_sample_2	16.62 %	
23	Honey_sample_2	16.70 %	
24	Honey_sample_3	16.67 %	
25	Honey_sample_3	16.84 %	
26	Honey_sample_4	16.87 %	

Previous 100

Next 100

Properties

Delete

Load

In the **Show files / View** dialog you can depict the individual determinations in greater detail. The sorting criterion and the number of saved determinations are displayed at the very top. The list contains numbered lines; each line represents one determination. In the default settings, the first sample identification and the first calculated result of each determination are displayed. In the properties, you define which data is displayed in the columns and the criteria according to which the determinations are sorted (see *"Properties"*, page 221). In addition, you can define search filters to aid you in finding the desired determinations (see *"Filter criteria"*, page 223).

[Previous 100]

Display the previous 100 determinations.

[Next 100]

Display the next 100 determinations.

[Properties]

Configure the display of the columns and define filter criteria.

[Delete]

Delete the selected file.

[Load]

Load the selected file.

Properties

Determination list: **Properties**



View / Properties

Column 1

Identification 1

Column 2

Result 1

Sort by

Determ. time

Sort sequence

descending

Filter

In the **View / Properties** dialog, you can define which data is displayed in the columns and the criteria according to which the determinations are to be sorted.

Column 1

Selection of what is to be displayed in the first column.

Selection	Identification 1 Identification 2 Determ. time Result name Method Sample size User
Default value	Identification 1

Column 2

Selection of what is to be displayed in the second column.

Selection	Result 1...9 The numbering of the results corresponds to the order in which the results were calculated in the determination run.
-----------	--

Sort by

Selection according to which criterion the list is to be sorted.

Selection	Identification 1 Identification 2 Determ.
Default value	time Result Method Sample size User Determ. time

Identification 1

Sort alphabetically according to the first sample identification.

Identification 2

Sort alphabetically according to the second sample identification.

Determ. time

Sort chronologically according to the determination time.



no filter

The list itself will not be filtered.

Identification 1

In the input field **is equal to**, select or enter the first sample identification according to which filtering is to be carried out.

Identification 2

In the input field **is equal to**, select or enter the second sample identification according to which filtering is to be carried out.

Date _____

In the input fields **Date from ... to**, specify the period of time according to which filtering is to be carried out.

Result

In the input fields **Value from ... to**, specify the value range according to which filtering is to be carried out.

Method

In the field **is equal to**, select or enter the method name according to which filtering is to be carried out.

Sample size

In the input fields **Value from ... to**, specify the value range according to which filtering is to be carried out.

User

In the field **is equal to**, select or enter the user according to whom filtering is to be carried out.



NOTE

The filter can be applied to empty fields as well, for example, if no sample identification has been entered. In this case, select only the desired filter criterion and leave the fields below empty.

24.6 Saving a determination

Results dialog: **More data** ► **Load/ Save** ► **Save**

Determinations can be saved as follows:


- Manually in this dialog.
The suggested default setting for file names is the first 16 characters of Identification 1, followed by the determination time in the format YYYYMMDD-hhmmss. Extending the file name with the determination time has the advantage that the file names of all determinations are unique.
- Automatically at the end of the determination (*see chapter 16.5.6, page 180*).

The following data is stored for each determination:

- Calculated results
- Measuring point lists, equivalence points, curves and variables which were generated during the determination
- Determination properties
- Statistics (if carried out)
- Method with which the determination was carried out
- Common variables (values and properties of all common variables at the time of the determination)
- Titrant data of the titrant used

24.7 Curves

Results dialog: **Curve**

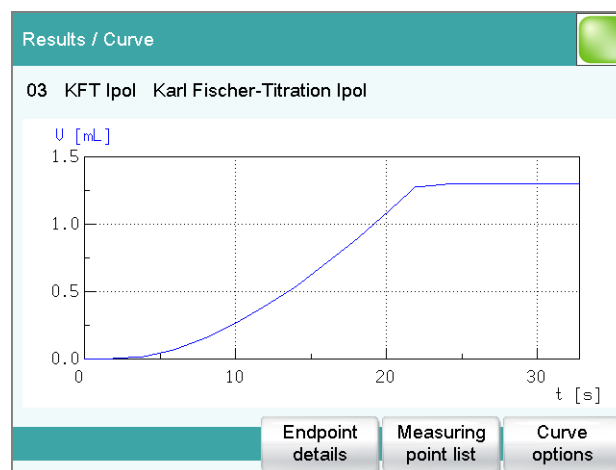

The curve of the current determination is displayed in the **Results / Curve** dialog. You can configure the curve display.

[Endpoint details]

Display details for each endpoint found (see *"Endpoint details"*, page 215).

[Measuring point list]


Display measuring point list of titrations and measurements (see *"Measuring point list"*, page 215).

[Curve options]

Define the settings for the curve display (see *"Curve options"*, page 226). The settings are saved separately for each titration/measuring mode and apply not only for the live display but also for the curve display in the results dialog.

Curve options

Results dialog: **Curve** ► **Curve options**

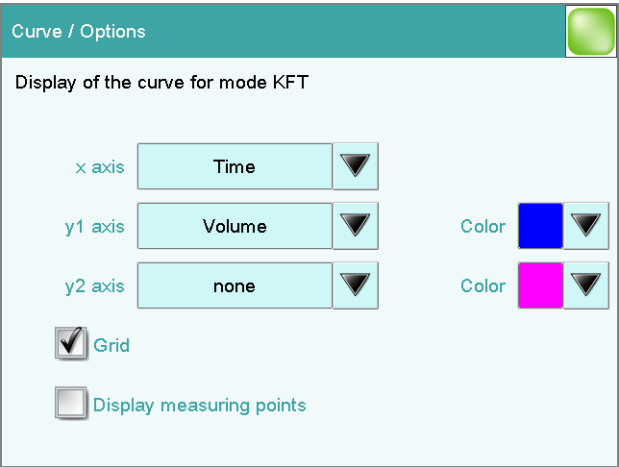

You can make the settings for the curve display in the **Curve / Options** dialog. These settings are saved separately for each mode (KFT, MEAS, etc.) and also apply for the live display. Both the selection of the quantities and the default setting depend on the mode.

x axis

Quantity which is shown on the x axis of the graph.

Selection	Volume Measured value Time Temperature Volume drift Signal drift
-----------	---

y1 axis

Quantity which is shown on the y1 axis of the graph.

Selection	Volume Measured value Time Temperature Volume drift Signal drift
-----------	---

y2 axis

Quantity which is shown on the y2 axis of the graph.

Selection	none Volume Measured value Time Temperature Volume drift Signal drift
Default value	none

Color

Selection of a color for the curves.

Grid

on | off (Default value: **on**)

If this parameter is activated, grid lines are shown.

Display measuring points

on | off (Default value: **off**)

If this parameter is activated, the curve is shown with the single measured values.

24.8 Recalculation and reevaluation

Results dialog: **Recalculate**

You can use the **[Recalculate]** function to recalculate and reevaluate determinations. This is necessary in the following cases:

- If you have modified the evaluation parameters in a way that changes the recognized endpoints or equivalence points.
- If you have modified the calculations.
- If you have modified the variables used in calculations, e.g. sample size, titer or common variables.



NOTE

Neither the recalculation nor the reevaluation can be undone.

We therefore recommend that the original determination be saved beforehand (*see chapter 24.6, page 224*). The best way of doing this is to use the **Save determination automatically** function (*see chapter 16.5.6, page 180*).

The following method commands are executed once again with the **[Recalculate]** function:

- EVAL commands:
All evaluations are carried out again. You can also insert new evaluations retroactively.
- CALC command:
All calculations are carried out again. You can also define new calculations retroactively.
If you subsequently modify the titrant in a method, then the variable TITER of the "new" titrant will be used. If a result is assigned to the variable TITER in a calculation, then the titer of the "new" titrant will also be overwritten.

If a determination has been recalculated, this will be documented in the determination properties (*see chapter 24.4, page 217*).

Recalculation/reevaluation of loaded determinations

With loaded determinations, it is assumed that they were carried out a long time ago and that the **titrant data** and **common variables** of the

system no longer correspond to those of the determination. This is why the most important titrant data of the titrants used and the list of **common variables** are stored with each determination. This data will be used if you recalculate a loaded determination. If you wish to recalculate the result with a corrected **common variable** or with a corrected titer, you must therefore modify the respective values in the **More determination data / Local common variables** or **More determination data / View data** dialog.

If a result is assigned to the variable TITER or to a **common variable**, then the titer or **common variable** which is currently valid in the system will not be overwritten until after a query has been confirmed.

If you modify the titrant retroactively in the method, this will not be taken into account.

The recalculated results of loaded determinations are not saved in the result table.

25 Statistics

Main dialog: **Results ► Statistics**


In the dialog **Results**, you can use **[Statistics]** to display the statistical overview of a determination overview. You can statistically evaluate a maximum of nine results calculated in a determination. A statistical series can contain a maximum of 20 determinations.

Activating the statistics function

Proceed as follows to activate the statistics function:

1 Editing the method options

- In the main dialog, tap on **[Edit parameters]**.
- Tap on **[Method options]**.

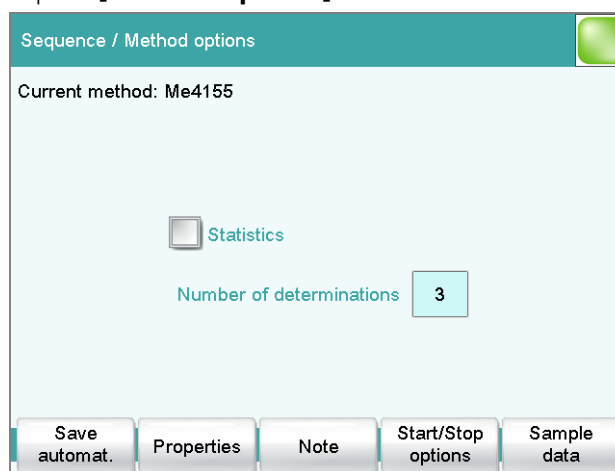

- Activate the **Statistics** option.
- Under **Number of determinations**, enter the desired number of determinations.

The results for which statistical calculations are performed are defined in the calculation command **CALC** (see "Variable for mean value", page 156).

25.1 Displaying details for a result

Results dialog: **Statistics ► Details**

Statistics / Details		
Result name: KFT Titer Mean value5.0080 mg/mL s abs0.02571 mg/mL s rel0.51 % SMN1 n=03		
No.	Sample size	Result
1	1.2 g	5.0170 mg/mL
2	1.7 g	4.9790 mg/mL
3	1.4 g	5.0280 mg/mL
Sample data Determ. on/off Result on/off		

In the dialog **Statistics / Details**, the following details are displayed:

- Result name
- Mean value
On the right-hand side of the dialog the number of single results from which the mean value was calculated is displayed. In this example, it is 3.
- Absolute standard deviation (**s abs**)
- Relative standard deviation (**s rel**)
- Result and sample size of each single determination
Results that could not be calculated are indicated with **invalid** and ignored.

[Sample data]

Display the sample data of the selected determination.

[Determ. on/off]

Remove all of the results of the selected determination from the statistics. All of the entries of this line are marked with an asterisk (*), the statistics will be recalculated automatically.

[Result on/off]

Remove the selected result from the statistics. The result is marked with an asterisk (*), the statistics will be recalculated automatically.



NOTE

If the mean value has been assigned to a common variable or to the variable TITER (see command **CALC**), then the determination must be recalculated manually (*see chapter 24.8, page 227*).

25.2 Deleting statistical data

In the dialog **Control**, you can delete all of the statistics data manually with the function **[Delete statistics]**. The statistics data is deleted automatically in the following cases:

- when all of the determinations of the determination series have been carried out and a new determination has been started afterwards.
 - when a new method is loaded (even if it is the same method)
- If the **sample table** or the **sample assignment table** is used, the statistics data is only deleted when another method is loaded.

25.3 Adding a determination to a determination series

You can use the function **[Add determ.]** to add an additional sample to a determination series, e.g. because a determination was faulty and had to be removed from the statistics.

Results / Statistics	
Current method: Me4155 3 of 3 (+2)	
Result name	Mean value
KFT Titer	5.0080 mg/mL

The number of determinations by which the set statistics counter has been increased is shown above the list in brackets. This modification is also displayed in the main dialog and in the statistics report.



The number of samples defined in the method options does not, however, change as a result.

[Delete]

Delete the selected line from the result table.

[Details]

Display details concerning the selected determination.


"Result table / Details" dialog

The following determination data is displayed in this dialog:

- **User**
User who carried out the determination.
- **Method**
Method with which the determination was carried out.
- **Determ. time**
Date and time at which the determination was carried out.
- Sample data (Identification 1 and 2 and sample size).
- All of the results of this determination saved in the result table.

26.1 Properties

Result table: **Properties**

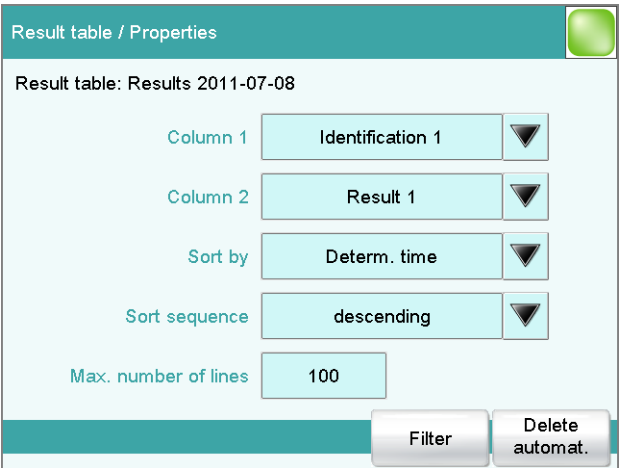

In the **Result table / Properties** dialog, you can select the data to be displayed with each determination as well as the sorting criterion for the list.

Column 1

Selection of what is to be displayed in the first column.

Selection	Identification 1 Identification 2 Determ. time Result name Method Sample size User
Default value	Identification 1

Column 2

Selection of what is to be displayed in the second column.

Selection	Result 1...9 The numbering of the results corresponds to the order in which the results were calculated in the determination run.
-----------	--

Sort by

Selection according to which criterion the list is to be sorted.

Selection	Identification 1 Identification 2 Determ.
	time Result Method Sample size User
Default value	Determ. time

Identification 1

Sort alphabetically according to the first sample identification.

Identification 2

Sort alphabetically according to the second sample identification.

Determ. time

Sort chronologically according to the determination time.

Result

Sort according to the numerical value of the result.

Method

Sort alphabetically according to the method name.

Sample size

Sort according to the numerical value of the sample size.

User

Sort alphabetically according to the user who carried out the determination.

Sort sequence

Selection	ascending descending
Default value	descending

Max. number of lines

If the result table contains the quantity of lines defined here, then no further results can be saved. If this is the case, save the result table and create a new one.

Input range	10 to 200
Default value	100

[Filter]

Specify the filter criteria for the list.

[Delete automat.]

Define the settings for the automatic deletion of the result table (see *"Deleting the result table automatically"*, page 238).

Filter criteria

Result table: **Properties** ► **Filter**

Properties / Filter

Filter

Date

is equal to

Date from

2011-07-01

to

2011-07-08

Value from

to

You can set filters in the **Properties / Filter** dialog. A filter is a rule defining what is to be displayed in the result table. If a filter has been set, this is indicated with a corresponding symbol in the **Results / Result table** dialog.

Filter

Define the criterion according to which the list is to be filtered. The only lines which will be displayed are those which fulfil the filter criterion.

Selection	no filter Identification 1 Identification 2 Date Result Method Sample size User
Default value	no filter

no filter

The list itself will not be filtered.

Identification 1

In the input field **is equal to**, select or enter the first sample identification according to which filtering is to be carried out.

Identification 2

In the input field **is equal to**, select or enter the second sample identification according to which filtering is to be carried out.

Date _____

In the input fields **Date from ... to**, specify the period of time according to which filtering is to be carried out.

Save result table

on | off (Default value: **on**)

If this option is activated, the result table will be deleted as soon as you have saved it.

26.2 Saving the result table



NOTE

When you save the result table, the complete determination will not be saved with each entry. If you wish to save the determination data for each determination, use the **Save determination automatically** function (see chapter 16.5.6, page 180).

26.3 Loading the result table

Proceed as follows to load a result table:

1 Opening the current result table

- In the results dialog, tap on **[Result table]**.
- Tap on **[Load/ Save]**.
- Tap on **[Load]**.
The selection of memory locations is displayed. This selection will be skipped if only one external storage medium is available.
- Select the memory location where the desired result table is saved.
The selection of file groups is displayed (see chapter 12.1, page 124). If only one group is available, then the list of the saved result tables will be displayed directly.
- Select the group with the desired result table.
- Tap on **[Show files]**.

The list with the saved result tables opens.

2 Selecting the result table

- Select the desired result table.

3 Loading the result table


- Tap on **[Load]**.

The result table is now loaded and will be displayed immediately.

More reports

Fixed key : **More reports**

The dialog **Print / More reports** displays a list of all of the available reports.

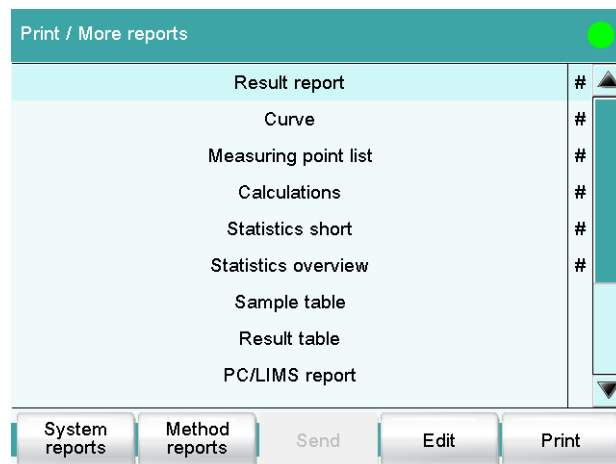

NOTE

If you are using a **CUSTOM Neos** printer, you can print only the reports marked with a #. The designation on the reports is only visible once the **CUSTOM NEOS** printer is selected (see screenshot above).

[System reports]

Opening the list of reports for the system settings.



NOTE

If you are using a **CUSTOM NEOS** printer, you cannot print system reports.

[Method reports]

Opening the list of reports for the current method.



NOTE

If you are using a **CUSTOM NEOS** printer, you cannot print method reports.

[Send]

Create a PC/LIMS report.

[Edit]

Edit the settings of the selected report.

[Print]

Print the selected report.

27.1 General report options

Fixed key []: **Report options**

The settings for the report printout can be adjusted in the dialog **Print / Report options** to indicate, for example, whether a report header, a signature line or a frame is to be printed.

Report header

Defining the output of the report header. The report header contains general information such as device type, serial number, program version and the printing date.



NOTE

You can also define your own report header, which is printed in addition to this report header (**System ► Templates ► Report header**).

Selection	off once on each page
Default value	on each page

off

The report header will not be printed.

once

The report header will only be printed on the first page.

on each page

The report header will be printed on every page.

Signature line

Output of a special line for date and signature. This line will be printed at the very bottom of every page.

Selection	off once on each page
Default value	off

off

No signature line will be printed.

once

The signature line will only be printed on the last page.

on each page

The signature line will be printed on every page.

Frame

on | off (Default value: **on**)

If this parameter is activated, a frame is printed as lateral marking.

27.2 Settings of the individual reports

Fixed key []: **More reports ► Report ► Edit**

Command REPORT: **Report ► Edit**


Settings can be edited for the following reports (for details, see online help):

- **Result report**
- **Curve**
Definition of the curve size, display of the individual measuring points and grid lines, etc.
- **Measuring point list**
Definition of the method command to which the measuring point list is to be printed.
- **Calculations**
Definition of the accuracy with which the variables used as well as the results calculated are printed.
- **Statistics short**
Definition of the time at which the report is to be printed (with each determination or only at the end of a sample series).
- **Statistics overview**
Definition of the time at which the report is to be printed (with each determination or only at the end of a sample series).
- **Sample table**
- **Result table**
Definition of the time at which the report is to be printed (with each determination, at the end of a sample series or at the end of the sample table).
- **Used devices**
- **Variables**
- **Monitoring**

Report	Contents
PC/LIMS report	Machine-readable report with all of the important data for a determination. This report can be saved as a TXT or UTF-8 file or sent to a terminal program or a LIMS via an RS-232 interface. The definition is carried out in the device manager (<i>see chapter 11.3.2, page 84</i>).
Sample data	
Current sample data	Sample data of the current determination, <i>context-sensitive only</i> from the main dialog.
Sample table	List of all determinations in the sample table with the respective sample data, as entered in the sample table.
Result table	
Result table	List of all determinations in the result table with results and with the determination data, as saved in the result table.
More reports ► Method reports	
Method sequence	Method properties and the list of all method commands.
Parameters full	Method properties and options, all method commands with all parameters. All parameters which no longer have their default settings will be printed in bold . All parameters which have been modified in comparison to the stored version of the method are indicated by an *.
Titration & measuring param.	Method properties; titration and measuring commands with all parameters. All parameters which no longer have their default settings will be printed in bold . All parameters which have been modified in comparison to the stored version of the method are indicated by an *.
Modified parameters	Method properties, all method parameters which have been modified in comparison to the stored version of the method.
Non default parameters	Method properties, all method parameters which no longer have any default settings.
More reports ► System reports	
System settings	
System settings	Settings for acoustic signals and accuracy of the measured value display.

Report	Contents
Common variable list	List of all the common variables defined in the system, together with their most important data (name, value, status).
All common variable properties	Properties of all common variables (name, value, validity, status).
Common variable properties	Properties of the selected common variable (name, value, validity, status), <i>context-sensitive only</i> from the editing dialog.
Local Common variable list	List with the common variables available at the time of the determination, together with their most important data (name, value, status), <i>context-sensitive only</i> from the dialog "More determination data / Local common variables".
Local common var. properties	Properties of the selected common variable (name, value, validity, status), <i>context-sensitive only</i> from the dialog "More determination data / Local common variables".
Templates	
Templates sample data	Sample identification list and sample assignment table.
Result template list	List with all of the user-generated result templates.
All result templates details	Details of the all of the user-generated result templates (calculation formula, result options, note).
Result template details	Details of the selected result template, <i>context-sensitive only</i> from the editing dialog.
Input/Output lines	List with all of the defined input and output lines at the remote interface (name, bit pattern).
Rack tables	
Sample rack list	List with all of the sample racks configured in the system, together with designation, number of positions and rack code.

28 Manual control

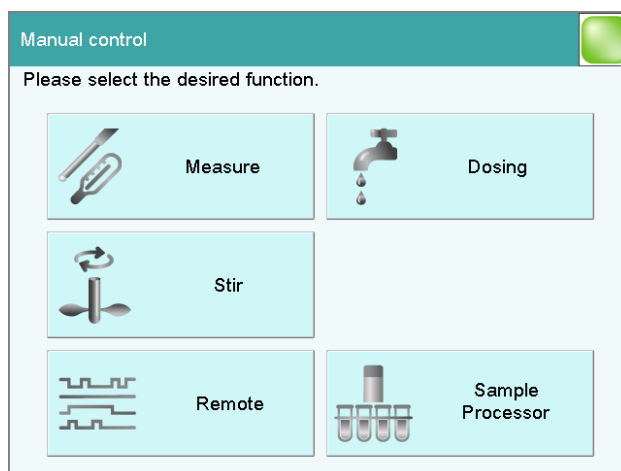

The following functions are available in the manual control:

[Measure]

Carrying out manual measurements (*see chapter 28.2, page 250*).

Measuring modes:

- **Ipol** (voltametric measurement with selectable polarization current)
- **Upol** (amperometric measurement with selectable polarization voltage)
- **T** (temperature measurement)

[Dosing]



NOTE

Only active if a Dosino or dosing unit is connected.

Manual dosing (*see chapter 28.3, page 252*).

The following dosing functions are available:

- Preparing the exchange unit or dosing unit
- Emptying the dosing unit
- Filling the dosing cylinder of the exchange unit/dosing unit
- Dosing a specified volume
- Dosing continuously

[Stir]

Controlling the stirrer manually (*see chapter 28.4, page 261*).

[Remote]

Scanning the input lines on the remote interface and setting the output lines (see chapter 28.5, page 263).

[Sample Processor]

Controlling the Sample Processor manually (see chapter 28.6, page 264).


28.1 Opening and closing the manual control

Opening the manual control

1 Opening the dialog window

- Tap on the fixed key [].

Manual control opens:

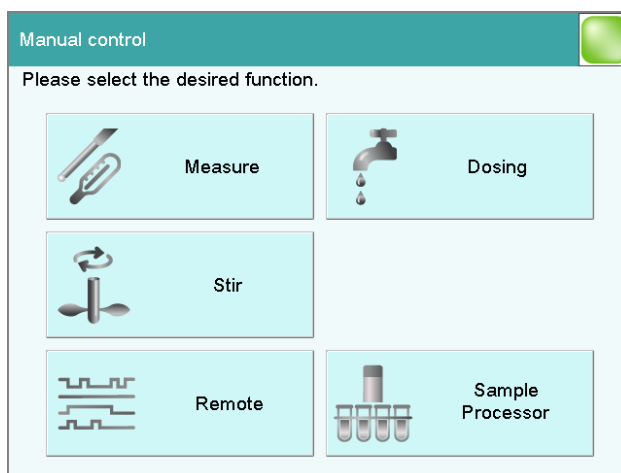

2 Select the function

- Tap on the button for the desired function.

If the function can be carried out by several devices, a dialog to select the device will appear first; otherwise the dialog for the selected function will appear immediately.

Closing the manual control

- 1 ▪ Tap on the fixed key [].



NOTE

Manual control can also be exited when a manual sequence has been started but has not yet ended. The fact that a sequence has been started in manual control can be recognized by the corresponding symbol in the title line (*see chapter 6.2.3, page 30*).

28.2 Measuring

Manual control ► Measure

Manual control / Measure

Device: Ti-Touch

Sensor: Metal electrode

Measuring input: 1

Measuring mode: Ipol

Temperature: 25 °C

Pol. current: - 5.0 + μ A

Info sensor Start

Manual measurements can be carried out with the function **[Measure]**.

Proceed as follows:

1 Selecting a sensor

- Select the desired sensor from the sensor list.
All of the sensors available in the sensor list are displayed. Sensors are defined under **System ► Sensors**.
- Select the measuring input at which the sensor is connected.

2 Select a measuring mode

- Select the desired measuring mode.
Only those measuring modes are displayed which are advisable for the selected sensor.


3 Enter a measuring temperature

- Enter the measuring temperature if no temperature sensor is connected. If a temperature sensor is connected, then the temperature will be measured automatically.

4 Starting the measurement

- Tap on **[Start]**.

The measured value display is opened:

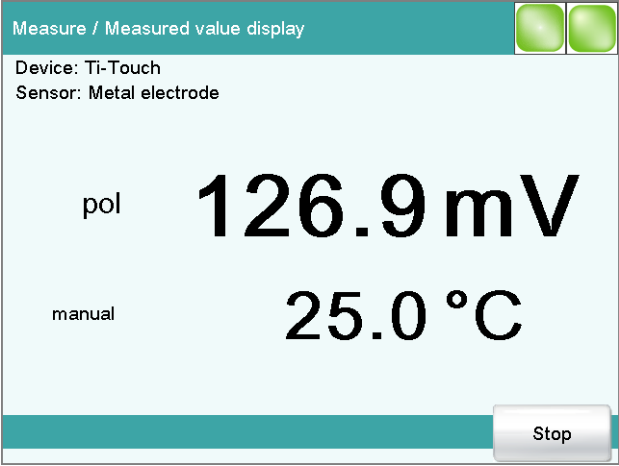

The ongoing measurement is canceled with **[Stop]**. The label switches once again to **[Start]**.

28.2.1 Parameter description

"Manual control / Measure" dialog

Sensor

Selection of the sensor from the sensor list. Sensors are defined under **System ► Sensors**.

Selection	Selection of configured sensors
-----------	---------------------------------

Measuring input

Selection of the measuring input the sensor is connected to. The selection depends on whether the control device has one or two measuring interfaces.

Selection	1 2
Default value	1

Measuring mode

Selection of the measuring mode. Only those measuring modes are displayed which are advisable for the selected sensor.

Selection	Ipol Upol T
-----------	-----------------

Ipol
voltametric measurement with selectable polarization current

Upol

amperometric measurement with selectable polarization voltage

T

temperature measurement

Temperature

Temperature entered manually. If a temperature sensor is connected, then the temperature will be measured continuously.

Input range	-20.0 to 150.0 °C
Default value	25.0 °C

[Info sensor]

Display information on the sensor.

28.3 Dosing


Manual control ► Dosing

The following dosing functions are available:

- Preparing the exchange or dosing unit (*see chapter 28.3.3, page 259*)
- Emptying the dosing unit (*see chapter 28.3.4, page 260*)
- Filling the dosing cylinder of the exchange unit/dosing unit (*see chapter 28.3.5, page 260*)
- Dosing a specified volume (*see chapter 28.3.2, page 256*)
- Dosing continuously (*see chapter 28.3.1, page 254*)

Selecting the dosing device

If several dosing devices are connected to a control instrument, then the dialog for selecting a dosing device will be opened.

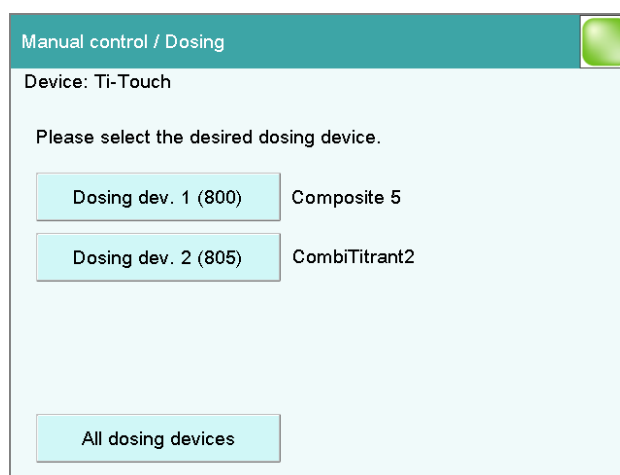

[Dosing device X]

Selection of the dosing device. The MSB connector and the type of the connected dosing device are displayed.

Various information is shown next to the button:

- **Titrant**
For dosing devices of the type 8XX with an intelligent exchange unit or dosing unit the titrant name is read off the data chip.
- **Cylinder volume**
For dosing devices of the type 8XX with an exchange unit without data chip only the cylinder volume is displayed.
- **---**
For exchange or dosing units which are not or not correctly attached.
- **Manual busy**
When a dosing function is being carried out.


[All dosing devices]

Carrying out the functions **[Prepare]**, **[Empty]** and **[Fill]** at the same time with several dosing devices of the control device.

Selecting the dosing function

If only one dosing device is connected, then the properties of this dosing device will be shown directly. This information varies according to the type of dosing device:

- **Type 8XX with integrated data chip:**
If an exchange or dosing unit with integrated data chip is attached, the data stored on the data chip is displayed.

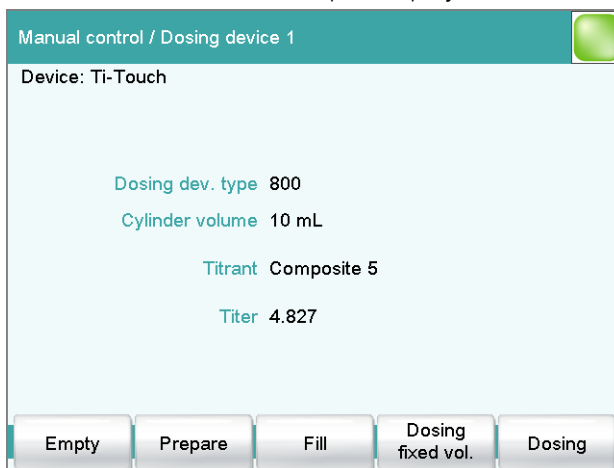

- Type 8XX without integrated data chip:
If an exchange unit without integrated data chip is attached, then the titrant can be selected from the titrant list. We recommend always selecting the titrant. This ensures that the parameters defined for the titrant for the preparing are used. Titrants are defined under **System ► Titrants**.

Manual control / Dosing device 2

Device: Ti-Touch

Dosing dev. type 805

Cylinder volume 10 mL

Titrant Titrant 5 ▼

Prepare
Fill
Dosing fixed vol.
Dosing

[Empty]

Empty the cylinder and the tubings of the dosing unit (see *chapter 28.3.4, page 260*).

[Prepare]

Rinse the cylinder and the tubings of the exchange / dosing unit (see chapter 28.3.3, page 259).

[Fill]

Fill the cylinder of the exchange / dosing unit (see chapter 28.3.5, page 260).

[Dosing fixed vol.]


Dose a specified volume (see chapter 28.3.2, page 256).

[Dosing]

Dose continuously (see chapter 28.3.1, page 254).

28.3.1 Continuous dosing

Manual control: **Dosing** ► **Dosing**

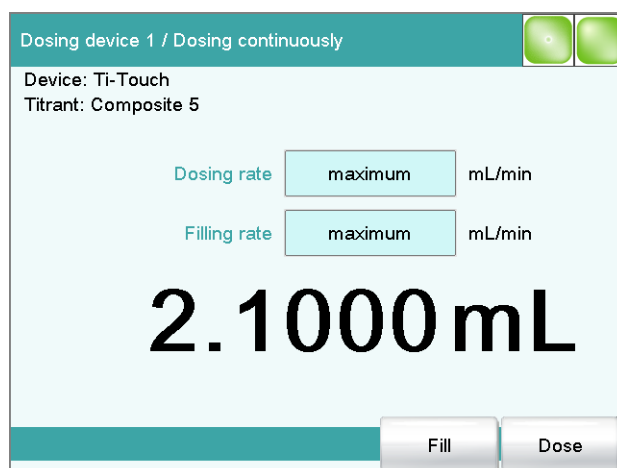

With the function **[Dosing]**, dosing is continuous for as long as you hold the button **[Dose]** pressed down.

Proceed as follows:

1 Configure the dosing function



NOTE

- The dosing and filling rates should be decreased for viscous liquids.
 - The maximum dosing and filling rates depend on the cylinder volume (*see chapter 32, page 356*).
- Enter the dosing rate.
 - Enter the filling rate.

2 Start dosing

- Tap on **[Dose]**.

The volume dosed is displayed. When the volume of one cylinder has been dosed, the dosing cylinder will be refilled automatically.

3 Fill the cylinder

- Tap on **[Fill]**.

The dosing cylinder is filled. The displayed volume value will be reset to 0.0000 mL.



NOTE

If you exit the dialog with the fixed key [, the dosing cylinder is filled automatically.

Parameter description

Dosing rate

Rate at which it is dosed.

Input range	0.01 to 166.00 mL/min
Selection	maximum dynamic
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (see chapter 32, page 356).

dynamic

This setting can only be selected when the dosing cylinder is filled. The dosing is being carried out faster and faster until the maximum dosing rate is reached (starting with 1 mL/min and then doubling the dosing rate every 1.5 s).

Filling rate

Rate at which the dosing cylinder is filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum filling rate depends on the cylinder volume (see chapter 32, page 356).

28.3.2 Dosing fixed volumes

Manual control: **Dosing ► Dosing fixed vol.**

Dosing device 1 / Dosing fixed volume

Device: Ti-Touch
Titrant: Solub.promoter: Formamid

Volume mL

Dosing rate mL/min

Filling rate mL/min

Fill automatic. ☒

You can dose a particular volume with the **[Dosing fixed vol.]** function.
Proceed as follows:

1 Configure the dosing function




NOTE

- The dosing and filling rates should be decreased for viscous liquids.
 - The maximum dosing and filling rates depend on the cylinder volume (*see chapter 32, page 356*).
- Enter the desired volume.
 - Enter the dosing rate.
 - Enter the filling rate.

2 Start dosing

- Tap on **[Start]**.

The volume display is shown. When the volume of one cylinder has been dosed, the dosing cylinder will be refilled automatically.

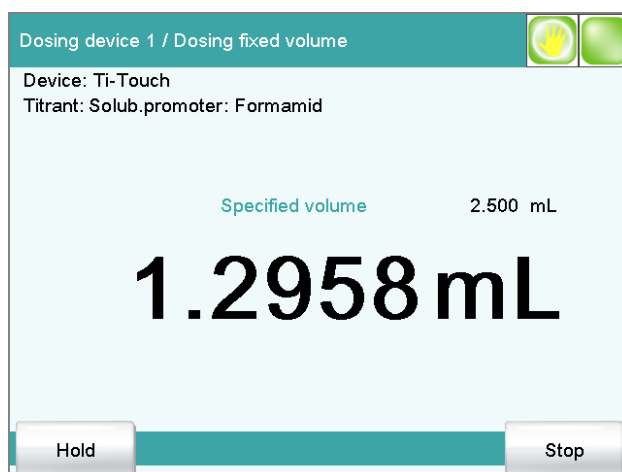

Continuous dosing is paused with **[Hold]**. The label changes to **[Continue]**. It can be used to resume the run.

3 Fill the cylinder

With the default settings (see "Fill automatic.", page 259) the dosing cylinder is filled automatically.

Otherwise:

- Tap on **[Fill]**.

The dosing cylinder is filled. The displayed volume value will be reset to 0.0000 mL.

Parameter description

Volume

Volume which is dosed.

Input range	0.00000 to 99999.9 mL
Default value	0.10000 mL

Dosing rate

The rate at which dosing takes place.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32, page 356*).

Filling rate

Rate at which the dosing cylinder is filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum filling rate depends on the cylinder volume (*see chapter 32, page 356*).

Fill automatic.

on | off (Default value: **on**)

If this parameter is activated, then the dosing cylinder is being automatically filled at the end of the dosing.

28.3.3 Preparing

Manual control: **Dosing ► Prepare**

The **Prepare** function is used to rinse and fill the cylinder and tubings of the exchange or dosing unit air bubble-free. You should carry out this function before the first determination or once a day.

If the titrant is selected (*see "Selecting the dosing function", page 253*), the parameters defined for the titrant for preparing/emptying and the tubing parameters will be used (dosing unit only). If the titrant is not selected, default parameters will be used (*see Chapter 32.1.2, page 356 and Chapter 32.2.2, page 357*).

The following two figures show the live display of an exchange unit and a dosing unit:

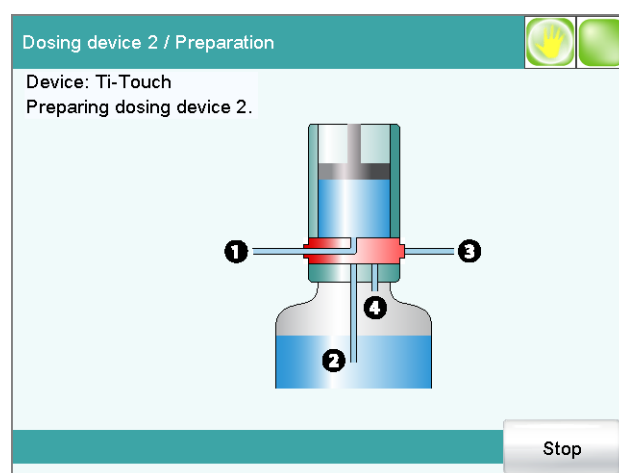


Figure 13 Live display "Preparing the dosing unit"

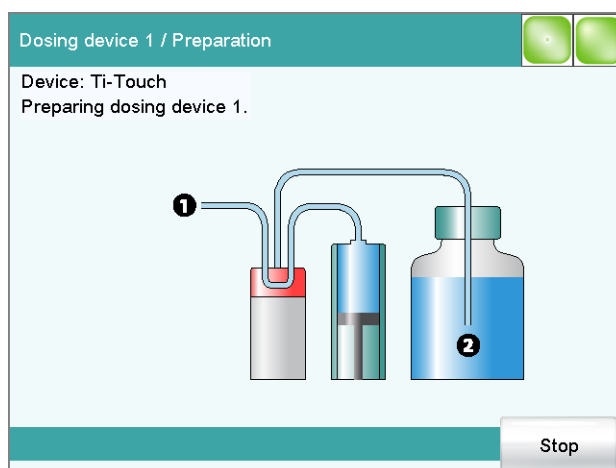


Figure 14 Live display "Preparing the exchange unit"

28.3.4 Emptying

Manual control: **Dosing** ► **Empty**

NOTE

The function **Empty** is possible only with dosing units.

The cylinders and the tubings of the dosing unit are emptied with this function.

If the titrant is selected (see *"Selecting the dosing function"*, page 253), the parameters defined for the titrant for preparing/emptying and the tubing parameters will be used. If the titrant is not selected, default parameters will be used (see *Chapter 32.2.2*, page 357 and *Chapter 32.1.2*, page 356).

28.3.5 Filling

Manual control: **Dosing** ► **Fill**

You can use the function **[Fill]** to fill the dosing cylinder of the exchange unit/dosing unit manually. The filling rate cannot be configured; the maximum filling rate will be used.

28.3.6 Replacing reagent



NOTE


You can easily change the reagent in a dosing unit without no contact with the chemicals using the functions **Prepare** and **Empty**.

Proceed as follows:

- 1 Execute the function **[Empty]**.
- 2 Attach the dosing unit to the bottle with a suitable solvent.
- 3 Execute the function **[Prepare]**.
The last traces of the reagent in the cylinder and tubings are removed with the solvent.
- 4 Execute the function **[Empty]** again.
- 5 Attach the dosing unit to the bottle with the new reagent.
- 6 Execute the function **[Prepare]**.
The dosing unit can now be used.

28.4 Stirring

Manual control ► Stir

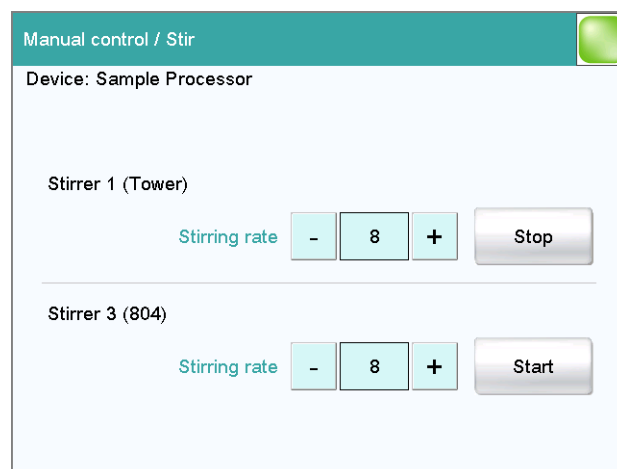


The MSB connector and type of stirrer is indicated for each connected stirrer.

You can control a connected stirrer manually with the function **[Stir]**.

28.5 Remote

Manual control ► Remote

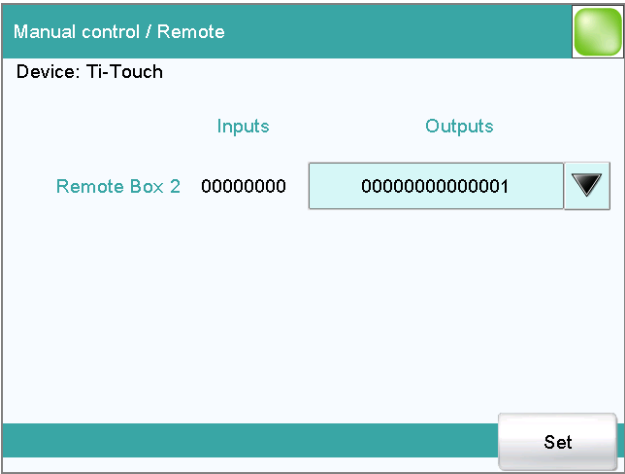

With the function **[Remote]**, you can define manual output signals to the remote interface of a connected Remote Box. It is not possible to define input signals; the current status of the input lines is, however, displayed.

Inputs

Current status of the input lines.

Outputs

Selection of the signal out of the templates or entering the required bit pattern. Templates are defined under **System ► Templates ► Output lines**.

Entering the bit pattern:

- 0 = line inactive
- 1 = line active
- * = retain line status
- p = set pulse

The output lines are always numbered from right to left, i.e. with the signal **000000000000001** line 0 is set. With a pulse, the length is set to 200 ms. If you wish to set pulses with other lengths, you have to define a corresponding template.

Entry	Bit pattern containing exactly 14 characters or a max. of 24 characters for the name of the template
Default value	000000000000001
Selection	Selection of the templates defined

**NOTE**

Note that the rack has to be set to a valid position for all lift movements! Otherwise a corresponding error message will be displayed.

Move to a valid rack position with **[Previous]** or **[Next]**. Alternatively, you can carry out a MOVE command.

[Home position]

Move the lift to the home position (corresponds to lift position 0 mm).

[Work position]

Move the lift to the work position of the current rack position or the external position.

The work position can directly be defined here (see *"Assigning lift positions"*, page 267) or in the device manager (see *"Lift positions Tower 1/2"*, page 103).

[↑]

Continuously move the lift upwards as long as the button is pressed down. The lift rate used is defined in the dialog **Sample Processor / Lift**.

[↓]

Continuously move the lift downwards as long as the button is pressed down. The lift rate used is defined in the dialog **Sample Processor / Lift**.

[Previous]

Move to the previous rack position.

[Next]

Move to the next rack position.

[Reset rack]

Initialize the rack. This function is identical with the command **RACK**. The following actions are carried out:

- The lift is moved upwards.
- The sample rack is rotated to the starting position.
- The rack code of the rack attached is being read out and the corresponding rack data is transferred to the Sample Processor.
- The sample variable is being reset to the value 1.

[Pump]

Switch on and off the pumps mounted or connected to the tower.

[Lift]

Move the lift to any or a specific lift position and assigning lift positions (see chapter 28.6.1, page 266).

[Rack position]

Move the rack to any sample position or to defined positions which are defined as special beakers (*see chapter 28.6.2, page 268*).

[External position]

This button is only enabled when a Swing Head is mounted to the Sample Processor.

Swing the robotic arm to an external position outside the sample rack and assign the swing angle and the lift position (*see chapter 28.6.3, page 269*).

28.6.1 Moving the lift

Manual control: **Sample Processor ► Lift**

Sample Processor / Lift

Tower 1

Current lift pos. 0 mm

Move to lift pos. mm

Lift rate mm/s

Assign lift pos. Start

In the dialog **Sample Processor / Lift**, you can move the lift to any lift position. You can also define frequently used positions (work position, rinse position, etc.) as specific lift positions and then to move to them with ease (see *"Assigning lift positions"*, page 267). This definition is possible in this dialog for the attached rack, and in the device manager for all of the racks in the list.

Current lift pos.

Current lift position.

Move to lift pos.

Selecting a predefined lift position or enter any lift position.

Input range	0 to 'maximum stroke path' mm A lift position of 0 mm corresponds to the "home position", i.e. the lift is located at the upper stop position. The maximum stroke path is defined in the properties of the tower (<i>see "Maximum stroke path", page 94</i>). If a higher value is entered, an error message will be displayed.
Selection	Work position Shift position Rinse position Special position
Default value	Work position

Lift rate

Rate at which the lift is moved in the manual control.

Input range	5 to 25 mm/s
Default value	25 mm/s

[Assign lift pos.]

Assign the current lift position to a preset lift position.

Assigning lift positions

Specific lift position are stored separately for every sample rack used as well as for tower 1 and tower 2. The following specific lift positions are available:

- **General rack positions**
For general rack positions you can define one work position, one shift position and one rinse position as well as one special position.
- **Special beaker positions**
A specific work position can be defined for each special beaker. The definitions of the general rack positions are used for the shift position, rinse position and special position of the affected tower.
- **External positions** (only with robotic arm and mounted Swing Head, *see Chapter 28.6.3, page 269*)
For each of the four possible external positions, a specific work position can be defined. Shift position and rinse position can only be defined for all four external positions commonly. A special position is not possible.

Proceed as follows to assign the current lift position to a specific lift position:

1 Move to lift position

- Enter the desired lift position in mm and tap on **[Start]**.

The lift moves to the desired lift position.

2 Assign lift position

- Tap on **[Assign lift pos.]**.

The **Lift / Assign lift position** dialog is displayed:

Lift / Assign lift position	
Current lift pos. 100 mm	
Position name	Position height
Work position	130 mm
Shift position	50 mm
Rinse position	100 mm
Special position	0 mm
Assign	

- Select the desired specific lift position and tap on **[Assign]**.

The current lift position is assigned to the specific lift position.

28.6.2 Moving to a rack position

Manual control: **Sample Processor** ► **Rack position**

Sample Processor / Rack position

Tower 1

Current rack pos. 1

Rack position ▼

Shift rate %/s

Shift direction ▼

Swing rate %/s

Start

In the dialog **Sample Processor / Rack position**, you can move to any position on the attached rack.

Current rack pos.

Current rack position. If the rack is not set to a valid position, e.g. after a rack reset, -- is displayed to indicate this.

Rack position

Desired rack position.

Input range	1 to Number of positions on the rack attached.
Default value	1
Selection	Special beaker 1...16

Special beaker 1...16

The selection depends on how many rack positions are defined as special beakers (see *"Editing special beakers"*, page 105).

Shift rate

Rate at which the sample rack is moved.

Input range	3 to 20 °/s
Default value	20 °/s

Shift direction

Direction in which the rack is moved.

Selection	auto + -
Default value	auto

auto

A shift direction, with which the shorter way has to be passed, is automatically selected.

+

Counterclockwise rotation.

—

Clockwise rotation.

Swing rate

Rate at which the robotic arm is swung in the manual control.

Input range	10 to 55 °/s
Default value	55 °/s

28.6.3 External positions

Manual control: **Sample Processor ► External position**



Sample Processor / External position

Tower 1

Lift position 100 mm Angle 105.0 °

°/s

Move to ext. pos. External position 1

NOTE

This dialog is accessible only if you have a Swing Head with robotic arm attached to the Sample Processor.

In the dialog **Sample Processor / External position**, you can swivel the robotic arm to an external position outside of the sample rack. As with the rack positions, here you can also define up to four frequently used external positions as specific positions. You can assign specific lift positions to these positions, analogous to the rack positions.

External positions are required if you have, for example, a measuring vessel or a titration cell mounted in addition to the rack.

Lift position

Current lift position.

Angle

Current angle position of the robotic arm.

[↑]

Only enabled, when the robotic arm has been moved to a specific external position with **[Start]**.

Continuously move the lift upwards as long as the button is pressed down. The lift rate used is defined in the dialog **Sample Processor / Lift**.

[↓]

Only enabled, when the robotic arm has been moved to a specific external position with **[Start]**.

Continuously move the lift downwards as long as the button is pressed down. The lift rate used is defined in the dialog **Sample Processor / Lift**.

[⇐]

Swing the robotic arm in clockwise direction as long as the button is pressed down.

[⇒]

Swing the robotic arm in counterclockwise direction as long as the button is pressed down.

Swing rate

Rate at which the robotic arm is swung in the manual control.

Input range	10 to 55 °/s
Default value	55 °/s

Move to ext. pos.

Selecting a predefined position or entering any swing angle.

Input range	(Offset) to (Offset + max. swing range) ° The offset is made up of a design-dependent angle (approx. 8...9°) together with the robotic arm offset from the robotic arm properties. The maximum swing range is also defined under the robotic arm properties (see "Properties – Robotic arm", page 97).
Selection	External position 1...4

[Assign Ext. pos.]

Assign the current angle position of the robotic arm to an external position as swing angle.

[Assign lift pos.]

Assign the current lift position to a preset lift position.

Assigning swing angles and lift positions

For each of the four possible external positions, a specific work position can be defined. Shift position and rinse position can only be defined for all four external positions commonly. A special position is not possible.

Proceed as follows to assign a swing angle to an external position and to assign the current lift position to a specific lift position:

1 Move to the external position

- Swing the robotic arm to the desired position with the arrow keys [⇐] or [⇒].

2 Assign the external position

- Tap on **[Assign Ext. pos.]**.

The **External position / Assign swing angle** dialog is displayed:

External position / Assign swing angle	
Current robotic arm angle 98.0 °	
External position	Angle
1	105.0 °
2	60.0 °
3	60.0 °
4	60.0 °

- Select the desired external position and tap on **[Assign]**.

The current angle position of the robotic arm is assigned to the external position.

3 Move to the external position

- Tap on [←].

The next higher-level dialog is displayed.

- Select **Move to ext. pos. = External position X** (X = 1...4) and tap on **[Start]**.

The robotic arm swings to the desired position.

4 Move to lift position

- Move the lift to the desired height with the arrow keys [**↑**] or [**↓**].

5 Assign lift position

- Tap on **[Assign lift pos.]**.

The **Lift / Assign lift position** dialog is displayed:

Lift / Assign lift position	
Current lift pos. 100 mm	
Position name	Position height
Work position	130 mm
Shift position	50 mm
Rinse position	100 mm
Special position	0 mm
Assign	

- Select the desired specific lift position.
- Tap on **[Assign]**.

The current lift position is assigned to the specific lift position.

28.7 885 Compact Oven SC

Manual control ► Sample Processor

Manual control / Sample Processor			
Tower 1 ▼			
Lift position	0 mm	Rack position	--
Home position	↑	Previous	Next
Work position	↓		
Reset rack		Heater Gas flow	

You can control a connected 885 Compact Oven SC manually with the **[Sample Processor]** function. The following functions are available:

- Initializing the rack
- Settings for the heating and the gas flow

Lift position

Current lift position.

Rack position

Current rack position. If the rack is not set to a valid position, e.g. after a rack reset, -- is displayed to indicate this.



NOTE

Note that the rack has to be set to a valid position for all lift movements! Otherwise a corresponding error message will be displayed.

Move to a valid rack position with **[Previous]** or **[Next]**. Alternatively, you can carry out a MOVE command.

[Home position]

Move the lift to the home position (corresponds to lift position 0 mm).

[Work position]

Move the lift to the work position of the current rack position.

The work position for sample beakers and conditioning beakers is defined in the device manager (*see chapter 11.7.2, page 110*).



Move the lift continuously upwards as long as the button is pressed down.
The lift moves in 6 mm steps.



Move the lift continuously downwards as long as the button is pressed down. The lift moves in 6 mm steps.

[\[Previous\]](#)

Move to the previous rack position.

[Next]

Move to the next rack position.

[Reset rack]

Initialize the rack. This function is identical with the command **RACK**. The following actions are carried out:

- The lift is moved upwards.
- The sample rack is rotated to the starting position.
- The sample variable is reset to the value 1.

[Heater Gas flow]

Display, control and settings of the current temperature and current gas flow.

valve

The integrated inlet valve is used for the gas flow. The inlet valve in the 885 Compact Oven SC is connected to compressed air or nitrogen.

Gas type

Gas type of the carrier gas used.

Selection	air nitrogen
Default value	air

Settings for the oven temperature

Temperature

Temperature to which the oven is heated.

Input range	50 to 250 °C (Increment: 1)
Default value	110 °C

29 Parameters



NOTE

For most of the numerical parameters, the result of a calculation can also be utilized in place of a number. A description of this can be found in the appendix (*see chapter 32.6, page 362*).

29.1 Volumetric Karl Fischer titrations (KFT)

29.1.1 Start conditions

The parameters that are carried out before the start of titration are defined under **[Start conditions]**.

Pause 1

Waiting time, e.g. for the electrode to settle down before a start volume is added.

Input range	0 to 999999 s
Default value	0 s

Start volume

Volume that is dosed prior to the start of the titration.

Input range	0.00000 to 9999.99 mL
Default value	0.00000 mL

Dosing rate

Rate at which the start volume is dosed.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (*see chapter 32, page 356*).

Pause 2

Waiting time, e.g. for the electrode to settle down after the start or a reaction time after the dosing of a start volume.

Input range	0 to 999999 s
Default value	0 s

29.1.2 Control parameters

The control parameters for the endpoint are defined under **[Control parameters]**.

Endpoint at

Measured value for the endpoint.

The values that are actually measurable may vary from the values specified here, depending on the hardware. The binding measurement value ranges can be found in the .

Measuring mode Ipol:

Input range	–2,000.0 to 2,000.0 mV
Default value	250.0 mV
Selection	off

Measuring mode Upol:

Input range	–200.00 to 200.00 µA
Default value	25.00 µA
Selection	off

Titration rate

Three predefined sets of parameters can be selected for the titration rate.

Selection	slow optimal fast user
Default value	optimal

slow

For samples with a low water content or samples which release their moisture only slowly.

optimal

For all standard titrations. The parameters have been optimized for the most frequent applications.

fast

For uncritical samples with high water content.

user

The individual titration parameters can be modified.

The settings of the individual titration rates are listed in *table 9, page 280*.

User-defined parameters

These parameters are only accessible when **Titration rate = user**.

Dynamics

This parameter defines the control range before the specified endpoint. Individual volume steps are dosed in the control range, the dosing is finely controlled. The closer the endpoint, the slower the dosing until the vol-

ume increment defined under **Min. increment** has been reached. The larger the control range, the slower the titration. Outside the control range, dosing is carried out continuously, and the dosing rate is defined under **Max. rate**.

The values that are actually measurable may vary from the values specified here, depending on the hardware. The binding measurement value ranges can be found in the .

Measuring mode I_{pol} :

Input range	0.1 to 2,000.0 mV
Default value	100.0 mV
Selection	off

Measuring mode Upol:

Input range	0.01 to 200.00 μA
Default value	10.00 μA
Selection	off

Max. rate

Rate at which dosing is carried out outside of the control range.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing rate depends on the cylinder volume (see chapter 32, page 356).

Min. increment

Minimum volume increment that is dosed at the beginning of the titration and in the control range at the end of the titration. This parameter has a decisive influence on the titration rate and thus also on the accuracy. The smaller the selected minimum increment, the slower the titration.

Input range	0.10 to 99.90 µL
Selection	minimum
Default value	minimum

Table 9 Default values of the predefined titration rates for KFT

	Titration rate		
	slow	optimal	fast
Dynamics			
– I _{pol}	300.0 mV	100.0 mV	30.0 mV
– U _{pol}	40.00 µA	10.00 µA	5.00 µA

	Titration rate		
	slow	optimal	fast
Max. rate	1.00 mL/min	maximum	maximum
Min. increment	minimum (= cylinder volume/10,000)	minimum (= cylinder volume/10,000)	5.00 µL

Stop criterion

Stop criterion

The titration is stopped when the endpoint has been reached and this stop criterion has been fulfilled. If no stop criterion has been selected then the titration will not be stopped. The stop conditions (*see chapter 29.1.4, page 283*) always lead to a stop, even if the stop criterion has not been reached.



NOTE

In older instructions the delay time was usually defined as the stop criterion. The same delay time may however result in different stopping time points – because of different smallest possible volume increments (depending on the cylinder volume). In contrast, if the stop drift is used, the titration will always be stopped at the same curve slope dV/dt .

Selection	drift time rel. drift off
Default value	drift

drift

The titration is stopped when the stop drift has been reached.

time

The titration is stopped if the endpoint has been exceeded during a certain time period (**Delay time**).

rel. drift

The titration is stopped when the sum of the drift at the start of the titration and the relative stop drift has been reached.

off

The titration will not be stopped until the stop conditions have been fulfilled.

Stop drift

This parameter can only be edited with **Stop criterion = drift**.

The titration is stopped when the endpoint and the stop drift have been reached.

Input range	0 to 999999 s
Default value	0 s

Temperature

Temperature entered manually. The temperature is being continuously measured when a temperature sensor is connected and when **Temp. measurement = automatic** or **continuous** is defined (see sensor dialog of the command).

Input range	–20.0 to 150.0 °C
Default value	25.0 °C

Time interval MP

Time interval for entering a measuring point in the measuring point list. The measuring point list is limited to 1000 measuring points.

Input range	0.1 to 999999.0 s
Default value	2.0 s

29.1.4 Stop conditions

The conditions for canceling the titration are defined under **[Stop conditions]**, if this does not occur automatically. This could be the case when the endpoint set is not reached or if the stop criterion (see "Stop criterion", page 281) is not fulfilled.

Stop volume

The titration is stopped when the specified volume has been dosed since the start of the titration. This volume should be adjusted to the size of the titration vessel in order to prevent the contents from running over.

Input range	0.00000 to 9999.99 mL
Default value	100,000 mL
Selection	off

Stop time

The titration is stopped when the specified time has elapsed following the termination of the start conditions.

Input range	1 to 999999 s
Selection	off
Default value	off

Filling rate

Rate at which the dosing cylinder is filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum filling rate depends on the cylinder volume (*see chapter 32, page 356*).

29.1.5 Conditioning

The conditions required for conditioning are defined under **[Conditioning]**.

Conditioning

on | off (Default value: **on**)

If this parameter is activated, then the first time the titration is started the working medium will be titrated to the endpoint with the specified control parameters. The status is kept stable. The actual method run does not begin until [] has been pressed once more.

Start drift

Conditioning OK will be displayed as soon as this drift has been reached and the titration can be started.

Input range	1 to 999 µL/min
Default value	20 µL/min

Drift correction

The endpoint volume can be corrected by taking the drift value into account. For this, the drift is multiplied with the drift correction time and this value is then subtracted from the endpoint volume. The drift correction time is the time interval between the end of conditioning and the end of the determination.

Selection	auto manual off
Default value	off

auto

The value of the current drift is automatically applied at the start of the titration.

manual

If the drift is known throughout a longer period of time, this can be entered manually.

off

No drift correction takes place.

Drift value

This parameter can only be edited with **Drift correction = manual**.

Drift for manual drift correction.

Input range	0.0 to 99.9 µL/min
Default value	0.0 µL/min

Measured value display

on | off (Default value: **off**)

If this parameter is activated, the currently measured value is displayed during the conditioning.

Conditioning options [Cond. options]

Cond. stop volume

Maximum permissible volume that can be dosed during conditioning. Conditioning is stopped when the specified volume is dosed. If conditioning is continued by pressing [] once again, then the titrant volume that has already been dosed will not be taken into account; i.e. the dosing starts again at zero. The stop volume should be adjusted to the size of the titration cell in order to prevent any overflow.

Input range	0.00000 to 9999.99 mL
Default value	20.0000 mL
Selection	off

Cond. stop time

Maximum permissible time over which conditioning may take place. Conditioning is stopped when the specified time has elapsed.

Input range	1 to 999999 s
Selection	off
Default value	off

Delay reconditioning

on | off (Default value: **off**)

If this parameter is activated, the conditioning will not start before all the method commands have been processed. Otherwise, the conditioning will immediately start after the titration command.



NOTE

This parameter must be activated if commands follow which must not run simultaneously with the conditioning.

Example: Emptying the titration cell followed by adding new working medium.

determination is started there is a check whether the sensor is contained in the sensor list.

Selection	Selection of configured sensors
-----------	---------------------------------

I(pol)

The polarization current is the current that is applied to a polarizable electrode during voltametric measurement. This parameter is available only with I(pol) determinations.

The values that are actually adjustable may vary from the values specified here. The binding values can be found in the manual, in the "Technical specifications" chapter.

Input range	–125.0 to 125.0 μA (Increment: 2.5)
Default value	50.0 μA
Selection	–1.0 μA 1.0 μA

U(pol)

The polarization voltage is the voltage applied to the polarizable electrode during an amperometric measurement. This parameter is available only with U(pol) determinations.

The values that are actually adjustable may vary from the values specified here. The binding values can be found in the manual, in the *Technical specifications* chapter.

Input range	–1,250 to 1,250 mV (Increment: 25)
Default value	400 mV

Electrode check

on | off (Default value: **off**)

For the following electrodes, an electrode check can be carried out:

- Metal electrodes

A check is made that the electrode is properly connected and that no short-circuit is present. The electrode check is carried out when this command is started.

Temp. measurement

Type of temperature measurement.

Selection	continuous automatic off
Default value	automatic

continuous

A temperature sensor must be connected. The temperature is measured continuously.

automatic

If a temperature sensor is connected then the temperature will be measured continuously. Otherwise, the temperature entered manually will be used (see dialog of the titration and measuring parameters).

off

The temperature will not be measured. The temperature entered manually is used (see dialog of the titration and measuring parameters).

29.1.9 Dosing device

The parameters for the dosing device are edited under **[Dosing device]**.

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

Titrant

Selection of the titrant from the titrant list. We recommend always selecting the titrant. This ensures that the correct data (titer, etc.) is always used for the calculation. Titrants are defined under **System ► Titrants**. You can also enter a name which is not contained in the titrant list. When a determination is started there is a check whether the solution is contained in the list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether the type of dosing drive matches. For buret units without an integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, the validity of the titer, the working life of the titrant and the time interval for the GLP test of the buret unit are checked at the start of the determination.

Selection	Selection of configured titrants not defined
Default value	not defined

not defined

No check takes place.

29.1.10 Stirrer

The parameters for the stirrer are edited under **[Stirrer]**.

Stirrer

Selection of the MSB connector the stirrer is connected to. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4 off
Default value	1

off

No stirrer will be used.

Stirring rate

Setting the stirring rate. It can be set in steps of –15 to +15. The default setting **8** corresponds to 1000 rpm. The formula for calculating the rotational speed is specified in *chapter 32.3, page 357*. The optimum stirring rate can be tested in the manual control.

The algebraic sign of the stirring rate changes the direction in which the stirring is done. When the stirrer is viewed from above, this means:

- "+": counterclockwise rotation
- "–": clockwise rotation

Input range	–15 to 15
Default value	8

Switch off automatically

on | off (Default value: **on**)

If this parameter is activated, the stirrer will be switched off automatically at the end of the titration, measurement, etc.

29.2 Measurements (MEAS)

29.2.1 Measuring parameters

Under **[Measuring parameters]**, the parameters influencing the run of the entire measurement are defined.

Signal drift

Maximum permissible drift for the measured value acceptance, i.e. maximum change of the measured value per minute. For this value, a suitable maximum waiting time is automatically calculated except you already have modified this waiting time.



NOTE

A constant measured value is often only reached after a certain time, as mixing and the reaction itself require a certain time. The response time of an electrode can also increase with time, i.e., reaching a constant measured value takes longer and longer. Drift-controlled measured value acceptance is particularly advisable in such cases, as the measured values are only accepted when equilibrium has almost been reached.

Measuring mode I_{pol} :

Input range	0.1 to 999.0 mV/min
Default value	10.0 mV/min
Selection	off

off

Measured value acceptance will take place after the maximum waiting time has elapsed.

Measuring mode Upol:

Input range	0.01 to 99.90 $\mu\text{A}/\text{min}$
Default value	10.00 $\mu\text{A}/\text{min}$
Selection	off

off

Measured value acceptance will take place after the maximum waiting time has elapsed.

Measuring mode T:

Input range	0.1 to 999.0 °C/min
Default value	0.5 °C/min
Selection	off

off

Measured value acceptance will take place after the maximum waiting time has elapsed.

Min. waiting time

The measured value is not accepted until the minimum waiting time has elapsed, even if the signal drift has already been reached. The minimum waiting time is only important for drift-controlled measurements.

Input range	0 to 999999 s
Default value	0 s

Max. waiting time

If the signal drift has been switched off or has not yet been reached, then the measured value will be accepted when the maximum waiting time has elapsed. As default value, a suitable waiting time for the signal drift is automatically calculated.

Input range	0 to 999999 s
	The default value depends on the measuring mode.
Selection	off
off	The measurement is continued endlessly.

Stop meas. value

The measurement is canceled when the specified measured value has been reached since the start of the measurement.

The values that are actually measurable may vary from the values specified here, depending on the hardware. The binding measurement value ranges can be found in the .

Measuring mode I_{pol}:

Input range	–2000.0 to 2000.0 mV
Selection	off
Default value	off

Measuring mode U_{pol}:

Input range	–200.00 to 200.00 µA
Selection	off
Default value	off

Measuring mode T:

Input range	–20.0 to 150.0 °C
Selection	off
Default value	off

Temperature

This parameter is not available with the **MEAS T** command (temperature measurement).

Temperature entered manually. The temperature is being continuously measured when a temperature sensor is connected and when **Temp. measurement = automatic** or **continuous** is defined (see sensor dialog of the command).

Input range	–20.0 to 150.0 °C
Default value	25.0 °C

Time interval MP

Time interval for entering a measuring point in the measuring point list.
The measuring point list is limited to 1000 measuring points.

Input range	0.1 to 999999.0 s
Default value	2.0 s

29.2.2 Control instrument

The control instrument the determination is carried out with is selected under **[Control device]**. Control instruments are defined under **System ► Device manager**.



NOTE

This button will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control instrument.

Control device

Selection of the control instrument from the list of devices. Only those devices are displayed which are able to carry out the command.

Selection	Selection of configured control instruments
-----------	---

29.2.3 Sensor

The parameters for the sensor are edited under **[Sensor]**.

The parameters are valid for all measuring modes except for **Temp. measurement**. This parameter is not available with the **MEAS T** command (temperature measurement).

Measuring input

Selection of the measuring input the sensor is connected to. The selection is not dependent on whether the control device has one or two measuring interfaces.

Selection	1 2
Default value	1

Sensor

Selection of the sensor from the sensor list. The selection depends on the measuring mode. Sensors are defined under **System ► Sensors**. You can also enter a sensor name which is not contained in the sensor list. When a determination is started there is a check whether the sensor is contained in the sensor list.

Selection	Selection of configured sensors
-----------	---------------------------------

I(pol)

The polarization current is the current that is applied to a polarizable electrode during voltametric measurement. This parameter is available only with I(pol) determinations.

The values that are actually adjustable may vary from the values specified here. The binding values can be found in the manual, in the *Technical specifications* chapter.

Input range	–125.0 to 125.0 μA (Increment: 2.5)
Default value	5.0 μA
Selection	–1.0 μA 1.0 μA

U(pol)

The polarization voltage is the voltage applied to the polarizable electrode during an amperometric measurement. This parameter is available only with U(pol) determinations.

The values that are actually adjustable may vary from the values specified here. The binding values can be found in the manual, in the *Technical specifications* chapter.

Input range	–1,250 to 1,250 mV (Increment: 25)
Default value	400 mV

Electrode check

on | off (Default value: **off**)

For the following electrodes, an electrode check can be carried out:

- Metal electrodes

A check is made that the electrode is properly connected and that no short-circuit is present. The electrode check is carried out when this command is started.

Temp. measurement

Type of temperature measurement.

Selection	continuous automatic off
Default value	automatic

continuous

A temperature sensor must be connected. The temperature is measured continuously.

automatic

If a temperature sensor is connected then the temperature will be measured continuously. Otherwise, the temperature entered manually will be used (see dialog of the titration and measuring parameters).

off

The temperature will not be measured. The temperature entered manually is used (see dialog of the titration and measuring parameters).

29.2.4 Stirrer

The parameters for the stirrer are edited under **[Stirrer]**.

Stirrer

Selection of the MSB connector the stirrer is connected to. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4 off
Default value	1

off

No stirrer will be used.

Stirring rate

Setting the stirring rate. It can be set in steps of -15 to $+15$. The default setting **8** corresponds to 1000 rpm. The formula for calculating the rotational speed is specified in *chapter 32.3, page 357*. The optimum stirring rate can be tested in the manual control.

The algebraic sign of the stirring rate changes the direction in which the stirring is done. When the stirrer is viewed from above, this means:

- "+": counterclockwise rotation
- "-": clockwise rotation

Input range	-15 to 15
Default value	8

Switch off automatically

on | off (Default value: **on**)

If this parameter is activated, the stirrer will be switched off automatically at the end of the titration, measurement, etc.

29.3 Evaluations (EVAL)

Various **additional evaluations** (EVAL commands) can be carried out for titrations and measurements. The evaluation always refers to the last titration or measurement prior to the EVAL command. Only those evaluation commands which are available for the last titration or measurement prior to the EVAL command can be inserted in the command list. If a titration or measuring command is deleted before the EVAL command, then it will be shown in red in the list of commands, as the reference is missing.

EVAL commands can also be inserted in the method run at a later time and the determination can be reevaluated (**[Recalculate]** function in the results dialog).

Additional evaluations for KFT titrations

The following additional evaluations are possible:

- **EVAL FIX-EP** (fixed endpoints)
Measured quantities:
 - Measured value
 - Time
 - Volume
- **EVAL MIN/MAX** (minimum/maximum)
- **EVAL RATE** (average dosing rate)

Additional evaluations for measurements

The following additional evaluations are possible:

- **EVAL FIX-EP** (fixed endpoints)
Measured quantities:
 - Measured value
 - Time
- **EVAL MIN/MAX** (minimum/maximum)
- **EVAL BREAK** (break points)

29.3.1 Fixed endpoint evaluation (EVAL FIX-EP)

For the fixed endpoint evaluation, the associated values are interpolated from the measuring point list for one quantity (measured value, volume, etc.). Up to nine fixed endpoints can be evaluated with one command.

Fixed quantity

Selection of the quantity to which the associated value is interpolated from the measuring point list.

Selection	Measured value Time Volume
Default value	Measured value

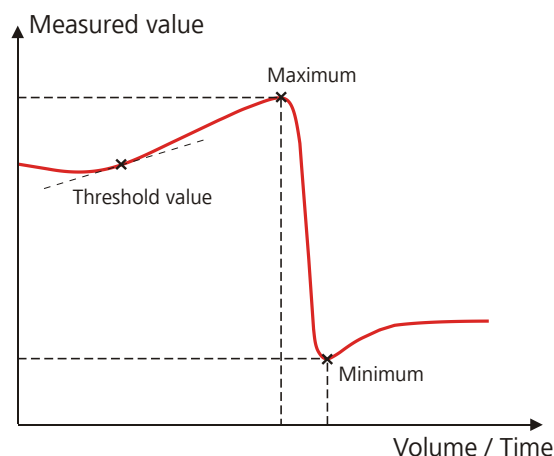


Figure 15 Evaluation of minimum and maximum

Evaluation

Selection whether the minimum or the maximum of the titration curve is being evaluated. With a command either the minimum or the maximum can be evaluated. If both values are needed, a second command has to be defined in the method.

Selection	Minimum Maximum
Default value	Maximum

Threshold value

The evaluation of the minimum or maximum begins as soon as the slope of the curve exceeds the set threshold value. Use a lower threshold value if the minimum or maximum is not found.

Measuring mode Ipol:

Input range	1.0 to 2,000.0
Default value	25.0
	With titrations the unit is mV/mL, with measurements mV/s.

Measuring mode Upol:

Input range	0.5 to 10.0
Default value	5.0
	With titrations the unit is $\mu\text{A/mL}$, with measurements $\mu\text{A/s}$.

Measuring mode T:

Input range	0.1 to 20.0 °C/s
Default value	1.0 °C/s

29.3.3 Rate evaluation (EVAL RATE)

The rate evaluation provides the median dosing rate in a defined range (so-called window). The median dosing rate is determined by means of linear regression through a minimum of three measuring points. If the dosing is carried out with only one dosing device and if more than one cylinder volume is dosed, then the median dosing rate will be reduced as a result of the filling time.

When the command is opened for the first time, a window over the entire time range is already defined. In contrary to the potentiometric evaluation of a DET/MET titration the individual windows may overlap. A maximum of 9 windows can be defined.

[New]

Define a new window.

[Delete]

Delete the selected window.

[Edit]

Edit the settings of the selected window.

Lower limit

Time for the lower limit.

Input range	0 to 999999 s
Default value	0 s

Upper limit

Time for the upper limit.

Input range	0 to 999999 s
Default value	999999 s

29.4 Calculations

The following calculation commands are available:

CALC	Defining calculations of a determination.
CALC LIVE	Defining the calculation whose result will be displayed in the live display (so-called live result) during a titration or measurement.

29.4.1 Calculations (CALC)

Calculations are defined with the command **CALC**. A method may contain a maximum of nine calculation commands. A maximum of nine calculations can be defined per command. A series of variables (raw data from the determination, previously calculated results, etc.) is available for the calculations.

The command contains the following options, among others:

- Save result as titer
- Save result as common variable
- Save result in result table
- Define result limits

Sequence / Edit command	
02 CALC	Calculation
Result	Result name
R1	Content
New Delete Edit	

The result name is specified in the list for each calculation.

[New]

Define a new calculation or load an existing template, see following chapter.

[Delete]

Delete the selected calculation from the list.

[Edit]

Edit the data of the selected calculation (*see chapter 29.4.1.2, page 302*).

29.4.1.1 Loading a result template


The most common calculations are already available, the so-called result templates, in order to facilitate the definition of a new calculation. You also have the option, however, of defining templates of your own (see dialog **System / Templates**).

Proceed as follows to load an existing template:

1 Select a result template

- In the **Sequence / Edit command** dialog, tap on the **[New]** button.

The list of Metrohm result templates is displayed:

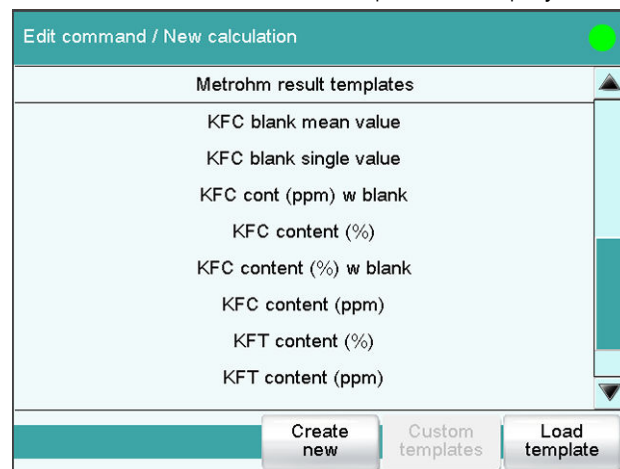


- Select the desired Metrohm result template and tap on **[Load template]**.

or

Tap on **[Custom templates]** and select a template you generated yourself.

If a statistic variable SMN# (# = 1...9) is defined in the template which is already used in another calculation, then the following message will be displayed:

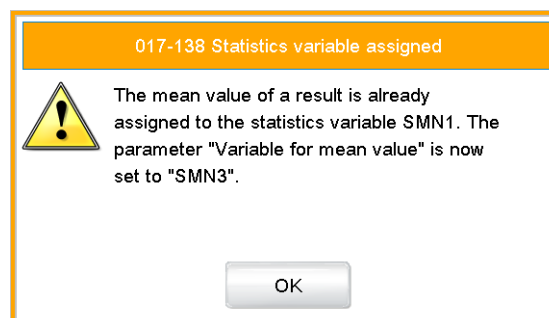

2

- Tap on **[OK]**.

The Note dialog is displayed. This note is defined in the dialog **Edit result template / Note for wizard**.

New calculation / Load template

Enter the following variable:

F1= Water content of the standard in mg/g

Please consider the information
in the note under edit result.

Cancel

Back

Next

3

- Tap on **[Next]**.

The list of wildcards is displayed:

New calculation / Load template

F1=10.02

F2=

F3=

F4=

F5=

F6=

F7=

F8=

F9=

Cancel

Back

Next

- Enter the numerical values of the wildcards.
- Tap on **[Next]**.

The editing dialog of the new calculation is displayed:

Edit command / Edit calculation

03 CALC

Calculation

Result name

KFT Titer

Calc. formula R1

C00/EP1*10.02

Decimal places

4

▼

Result unit

mg/mL

▼

Note

Result variable

Result limits

Result options

The calculation can now be edited further, e.g. define result limits, see following chapter.

29.4.1.2 Editing a calculation

Result name

The result name is the text which will be shown in the result display and in the report. The default name corresponds to the result variable.

Entry	max. 24 characters
Default value	R1...R9

Calc. formula R1...R9

Shows the calculation formula. A special editor is opened for the definition (see chapter 29.4.3, page 308). The designation **R1...R9** corresponds to the result variable. Under **[Result variable]** this can later be modified.

Entry	max. 100 characters
Default value	empty

Decimal places

Number of decimal places used to display the result.

Input range	0 to 5
Default value	2

Result unit

The result unit is displayed and saved together with the result.

Entry	max. 10 characters
Selection	% mol/L mmol/L g/L mg/L mg/mL mg/100 g ppm g mg µg mL µL mg/ piece °C µg/min mL/min µL/min
Default value	%

[Note]

Entering a note on the calculation.

[Result variable]

Modifying the result variable.

[Result limits]

Defining the limits for the monitoring of the result.

[Result options]

Defining additional settings for the calculation.

Dialog "Edit calculation / Note"

In this dialog you can enter a short text, e.g. to describe the variables used.

Dialog "Edit calculation / Result variable"

When a new calculation is created, the result variable will be automatically issued. In this dialog, however, it can be modified at any time.

Result variable

To each calculation belongs an unambiguous result variable. With this result variable you can use this result in additional calculations.

Selection	R1 R2 R3 R4 R5 R6 R7 R8 R9
	The selection contains only those result variables which are not issued in this calculation command yet.

Dialog "Edit calculation / Result limits"

For each result, limit values can be defined. These result limits are monitored when the calculation is carried out. If monitoring is activated, the results are shown as follows in the result display:

- Green, if the result falls within the limit values
- Red, if the result lies outside the limit values

Monitoring result limits

on | off (Default value: **off**)

Enabling and disabling the monitoring function.

Lower limit

When the monitored quantity falls below this value, the action defined in the following is triggered.

Input range	0 to 9999999999
Default value	0

Upper limit

When this value is exceeded, the action defined in the following is triggered.

Input range	0 to 9999999999
Default value	9999999999

Action

Selection of the action, which is carried out when the result limits are exceeded.

Selection	Display message Document message Cancel determination
Default value	Display message

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the result limits have been exceeded.

Cancel determination

The determination is stopped.

Dialog "Edit calculation / Result options"

In the dialog **Edit calculation / Result options**, settings for how to process the calculated result are defined.

Variable for mean value

If the statistic calculations have been activated (see method options), the mean value of the single results will be saved as variable SMN1 to SMN9. As default value, always the first free variable is being displayed.

Selection	SMN1 ... SMN9	off
-----------	----------------------	------------

off

For the result not statistic calculations are carried out.

Save as titer

The result can be saved as titer for the selected titrant.



NOTE

The titer is assigned to the titrant defined in the last titration command before the calculation. Please note that the CALC command containing the titer assignment is inserted after the titration command with which the titer is determined.

Selection	Single value Mean value off
Default value	off

Single value

The result of the current determination is saved as titer.

Mean value

If statistic calculations are carried out for the result, the current mean value of the determination series is saved.

off

The result is not saved as titer.

Save as common variable

on | off (Default value: **off**)

The calculated result can be saved as a method-independent variable, a so-called common variable. The result is then also available in other methods for calculations. In all cases, only the single value will be saved, even if the statistics function is enabled.

Variable

Selection of the common variable to which the result is being assigned.

Selection	CV01 ... CV25
-----------	----------------------

Display result

on | off (Default value: **on**)

If you deactivate this parameter, the result is neither displayed in the result dialog nor printed in the result report. This can be advisable for intermediate results.

Save result in result table

on | off (Default value: **off**)

The calculated result can be saved in the result table. This may be advisable if e.g. the results of all determinations carried out on a particular day are to be displayed clearly. A maximum of nine results from a determination can be saved in the result table.

Precision

Setting, with which accuracy the result is used in additional calculations.

Selection	Round Truncate Full precision
Default value	Round

Decimal places

Number of decimal places used to display the result.

Input range	0 to 5
Default value	2

Result unit

The result unit is displayed and saved together with the result.

Entry	max. 10 characters
Selection	% mol/L mmol/L g/L mg/L mg/mL mg/100 g ppm g mg µg mL µL mg/ piece °C µg/min mL/min µL/min
Default value	%

[Note]

Entering a note on the calculation.

Dialog "Edit calculation / Note"

In this dialog you can enter a short text, e.g. to describe the variables used.

Displaying the live result

Proceed as follows to display the live result in the live display:

1 Starting a determination

- Tap on [].

The determination starts and the live display is shown:

2 Define the display options

- Tap on the **[View]** button.
- Tap on the **[Meas. value options]** button.
A maximum of three measured values can be shown in the live display.
- Select the option **Live result** at one of the three parameters.

3 Display the live display

- Tap twice on [].

The current result is displayed,

If the method contains more than one data-generating command (titration, measurement etc.) before the calculation command, then you will have to indicate the command identification in front of the variable:

Variable	Description
#M.	Titration or measuring command
#C.	Calculation command
#E.	Evaluation command

"#" stands for a sequential number that you must enter manually. The number for each command identification will be counted separately and does not necessarily correspond to the number of the command line.

Example: The variable **1M.TITER** corresponds to the titer of the titrant selected in the first titration command.



CAUTION

The command identifications are not adjusted automatically when additional commands are subsequently inserted into the method. Check the formula for correctness in such cases.

[Variables]

Opening the list of additional variables (*see chapter 29.4.3.1, page 309*).

[Math. functions]

Opening the list with mathematical functions (*see chapter 29.4.3.2, page 310*).

29.4.3.1 Variables

All of the variables which can be used for calculations are divided up according to subject in the dialog **Calculation formula / Variables**:

- **System variable**
List of the system-specific variables. These variables describe the current status of the system.
- **Result and statistics variables**
List of the variables which are required for calculating results and statistics for the current determination.
- **Common variables**
List of common variables. All of the common variables are listed, no matter whether they possess a value or not.

You will find a detailed description of all of the variables in *chapter 29.4.3.3, page 310*.


Variable	Description
DD	Duration of the entire determination
Titration	
TITER	Titer of the titrant selected in the titration command
CONC	Concentration of the titrant selected in the titration command
Titration, measurements, calibrations	
EP#	Volume/quantity of endpoint EP# (# = 1 - 9)
EC#	Charge at endpoint EP# (# = 1 - 9)
EM#	Measured value of endpoint EP# (# = 1...9)
EF#	ERC of endpoint EP# (# = 1 - 9)
ET#	Temperature at endpoint EP# (# = 1...9)
ED#	Time at endpoint EP# (# = 1...9)
ESI#	Recognition of endpoint EP# (# = 1 - 9) Endpoint found = 1, no endpoint = 0
FMN	Mean gas flow during MEAS T/F ON to MEAS T/F OFF
TMN	Mean temperature during MEAS T/F ON to MEAS T/F OFF
TMI	Minimum temperature during MEAS T/F ON to MEAS T/F OFF
TMA	Maximum temperature during MEAS T/F ON to MEAS T/F OFF
MIM	Initial measured value, i.e. measured value prior to the processing of the start conditions
MIT	Initial temperature, i.e. temperature prior to the processing of the start conditions
MSA	Volume for start volume
MSP	Volume for start measured value
MSS	Volume for start slope
MSV	Volume for all the start conditions
MSD	Duration of start conditions
MSM	Start measured value, i.e. measured value after the processing of the start conditions
MST	Start temperature, i.e. temperature after the processing of the start conditions
MCD	Total duration of the titration, measurement or calibration
MTM	Type of temperature measurement (Pt1000, NTC or manually) Format = Text
MDD	Duration of effective dosing, i.e. without filling times, pauses
MDC	Drift for drift correction
DDC	Time for drift correction
MCV	End volume, i.e. total dosed volume at the end of the titration

Variable	Description
BP#	Volume of break point BP# (# = 1 - 9)
BM#	Measured value of break point BP# (# = 1 - 9)
BT#	Temperature at break point BP# (# = 1 - 9)
BD#	Time at break point BP# (# = 1 - 9)
RD#	Mean rate in time slot # (# = 1 - 9)
RDS#	Standard deviation for RD# (# = 1 - 9)
RDC#	Correlation coefficient for RD# (# = 1 - 9)
RM	Mean rate for whole range
RMS	Standard deviation for RM
RMC	Correlation coefficient for RM
Results, statistics	
R#	Result (# = 1–9)
SMN#	Mean value of result R# (# = 1–9)
SSA#	Absolute standard deviation of SMN# (# = 1 - 9)
SSR#	Relative standard deviation of SMN# in % (# = 1 - 9)
SNR#	Number of results from which the mean value SMN# was calculated (# = 1 - 9)
SSD	Number of determinations that are carried out for the statistics calculations
SNT	Statistics status (statistics activated = 1, statistics deactivated = 0)
Common variables	
CV#	Common variable (# = 1...25)
System variables	
%RN	Sample number, i.e. number that is increased by one at each start of a determination
%SC	Start counter, see dialog More determination data / Properties
%AS	Autostart status (autostart activated = 1, autostart deactivated = 0)
%AC	Autostart counter, i.e. number of autostarts already carried out
%AD	Autostart setpoint counter, i.e. number of autostarts to be carried out
%SS	Status of the sample table (sample table activated = 1, sample table deactivated = 0)
%SL	Number of the current sample table line
%SE	Last sample in the sample table processed (yes = 1, no = 0)

29.5 Reports (REPORT)

The reports to be printed out are defined with the command **REPORT**. A maximum of ten reports can be defined for each command. The following reports can be inserted, among others:

- Result report
- Calculation report
- Parameter report
- System reports (system settings, sensor list, GLP data etc.)
- Method reports (method run, method parameters etc.)

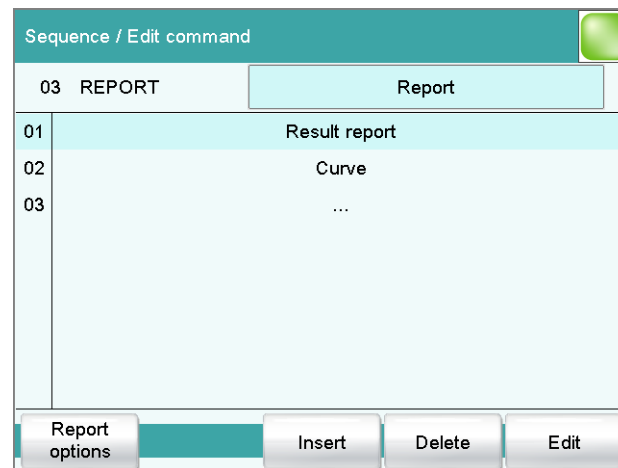

[Report options]

Defining the general settings for the report output.

[Insert]

Add a new report to the list.

[Delete]

Delete the selected report from the list.

[Edit]

Edit the settings of the selected report.

29.5.1 General report options

The settings for the report printout can be adjusted in the dialog **Edit command / Report options** to indicate, for example, whether a report header, a signature line or a frame is to be printed.

Report header

Defining the output of the report header. The report header contains general information such as device type, serial number, program version and the printing date.



Selection	off once on each page
Default value	on each page

The report header will not be printed.

The report header will only be printed on the first page.

The report header will be printed on every page.

Output of a special line for date and signature. This line will be printed at the very bottom of every page.

No signature line will be printed.

The signature line will only be printed on the last page.

The signature line will be printed on every page.

on | off (Default value: **on**)

If this parameter is activated, a frame is printed as lateral marking.

Fixed key [📄]: **More reports** ► **Report** ► **Edit**Command REPORT: **Report** ► **Edit**

Settings can be edited for the following reports (for details, see online help):

- **Result report**
- **Curve**
Definition of the curve size, display of the individual measuring points and grid lines, etc.

- **Measuring point list**
Definition of the method command to which the measuring point list is to be printed.
- **Calculations**
Definition of the accuracy with which the variables used as well as the results calculated are printed.
- **Statistics short**
Definition of the time at which the report is to be printed (with each determination or only at the end of a sample series).
- **Statistics overview**
Definition of the time at which the report is to be printed (with each determination or only at the end of a sample series).
- **Sample table**
- **Result table**
Definition of the time at which the report is to be printed (with each determination, at the end of a sample series or at the end of the sample table).
- **Used devices**
- **Variables**
- **Monitoring**

29.5.3 List of reports

The following reports can be inserted in the command **REPORT**:



NOTE

If you are using a **CUSTOM NEOS** printer, you can print only the reports that are marked with a **#** on the display.

Report	Contents
Result report	# Report with determination properties, sample data, calculated results, etc. If there are multiple determinations, the statistics will also be printed out.
Curve	# Curve report.
Measuring point list	# Measuring point list report.
Calculations	# Details concerning the calculations carried out (parameter settings in accordance with the CALC command and the calculated results).
Used devices	The devices used for the determination, as displayed in the More determination data / Properties dialog.
Variables	All determination variables, as displayed in the results dialog.
Monitoring	Details for the monitored quantities.
Statistics	

Report	Contents
Statistics short	# Summary of the statistics calculations. The number of determinations, the mean value, and the absolute and relative standard deviations are printed out for each result.
Statistics overview	# Detailed statistical overview. The sample data and all individual results are printed out for each determination. The number of determinations, the mean value, and the absolute and relative standard deviations are printed out for each result.
Sample data	
Sample table	List of all determinations in the sample table with the respective sample data, as entered in the sample table.
Result table	
Result table	List of all determinations in the result table with results and with the determination data, as saved in the result table.
More reports ► Method reports	
Method sequence	Method properties and the list of all method commands.
Parameters full	Method properties and options, all method commands with all parameters. All parameters which no longer have their default settings will be printed in bold . All parameters which have been modified in comparison to the stored version of the method are indicated by an *.
Titration & measuring param.	Method properties; titration and measuring commands with all parameters. All parameters which no longer have their default settings will be printed in bold . All parameters which have been modified in comparison to the stored version of the method are indicated by an *.
Modified parameters	Method properties, all method parameters which have been modified in comparison to the stored version of the method.
Non default parameters	Method properties, all method parameters which no longer have any default settings.
More reports ► System reports	
System settings	
System settings	Settings for acoustic signals and accuracy of the measured value display.
Dialog options	Settings for routine dialog and expert dialog.
Titrants	
Titrant list	List of all the titrants configured in the system.
All titrant data short	The most important titrant data of all titrants (name, concentration, titer, last titer determination).
All titrant data full	All of the titrant data of all of the titrants (working life, titer options, data concerning the exchange unit/dosing unit).

Report	Contents
Sensors	
Sensor list	List of all the sensors configured in the system.
All sensor data short	The most important sensor data for all sensors.
All sensor data full	All of the sensor data for all of the sensors (name, working life, etc.).
Device manager	
Device list	List of all devices configured in the system.
All device properties	Properties of all the devices configured in the system.
GLP manager	
GLP data	All data stored in the GLP manager.
Common variables	
Common variable list	List of all the common variables defined in the system, together with their most important data (name, value, status).
All common variable properties	Properties of all common variables (name, value, validity, status).
Templates	
Templates sample data	Sample identification list and sample assignment table.
Result template list	List with all of the user-generated result templates.
All result templates details	Details of the all of the user-generated result templates (calculation formula, result options, note).
Input/Output lines	List with all of the defined input and output lines at the remote interface (name, bit pattern).
Custom calibration buffers	Temperature tables for all defined custom calibration buffers.
Rack tables	
Sample rack list	List with all of the sample racks configured in the system, together with designation, number of positions and rack code.
Miscellaneous	
Form feed	If this entry is inserted between two reports, then these will each be printed out on a separate page.

29.6 Dosing and Liquid Handling

The following dosing commands are available:

Preparing the exchange unit or dosing unit (PREP)	Rinsing the cylinder and the tubings of the exchange unit / dosing unit.
Emptying the dosing unit (EMPTY)	Emptying the cylinder and the tubings of the dosing unit.
Dosing a fixed volume (ADD)	Dosing a specified volume.
Liquid Handling (LQH)	Carrying out complex dosing tasks with a Dosino.

29.6.1 Preparing an exchange or dosing unit (PREP)

The **PREP** command is used to rinse and fill air bubble-free the cylinder and tubings of the exchange or dosing unit. You should carry out this function before the first determination or once a day.

Control device

This parameter will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control device. Control devices are defined under **System ► Device manager**.

Selection of the control device from the list of devices.

Selection	Selection of configured control devices
-----------	---

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

Titrant

Selection of the titrant from the titrant list. We recommend always selecting the titrant. This ensures that the parameters defined for the titrant for the preparing as well as the tubing parameters (dosing unit only) are used. Titrants are defined under **System ► Titrants**. You can also enter a name which is not contained in the titrant list. When a determination is started there is a check whether the solution is contained in the list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether

the type of dosing drive matches. For buret units without integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, the validity of the titer, the working life of the titrant and the time interval for the GLP test of the buret unit are checked at the start of the determination.

Selection	Selection of configured titrants not defined
Default value	not defined

not defined

No check takes place. The command is carried out with default parameters (see *Chapter 32.1.2, page 356* and *Chapter 32.2.2, page 357*).

Cylinder volume to

This parameter can only be edited with **Titrant = not defined**.

Selection of the port via which the titrant is ejected. This setting is only relevant for dosing units, in case of preparing an exchange unit, this setting is ignored.

Selection	Port 1 Port 2 Port 3 Port 4
Default value	Port 1

29.6.2 Emptying a dosing unit (EMPTY)

The **EMPTY** command is used to empty the cylinder and tubings of the dosing unit.

Control device

This parameter will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control device. Control devices are defined under **System ► Device manager**.

Selection of the control device from the list of devices.

Selection	Selection of configured control devices
-----------	---

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

Titrant

Selection of the titrant from the titrant list. We recommend always selecting the titrant. This ensures that the parameters defined for the titrant for the preparing as well as the tubing parameters (dosing unit only) are used.

Titration are defined under **System ► Titrants**. You can also enter a name which is not contained in the titrant list. When a determination is started there is a check whether the solution is contained in the list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether the type of dosing drive matches. For buret units without integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, the validity of the titer, the working life of the titrant and the time interval for the GLP test of the buret unit are checked at the start of the determination.

Selection	Selection of configured titrants not defined
Default value	not defined

not defined

No check takes place. The command is carried out with default parameters (see Chapter 32.1.2, page 356 and Chapter 32.2.2, page 357).

Air inlet

This parameter can only be edited with **Titrant = not defined**.

Selection of the port via which air is aspirated. This setting is only relevant for dosing units, in case of preparing an exchange unit, this setting is ignored.

Selection	Port 1 Port 2 Port 3 Port 4
Default value	Port 4

29.6.3 Dosing a specified volume (ADD)

You can dose a specified volume with the command **ADD**.

29.6.3.1 Dosing parameters

The parameters for the dosing are defined under **[Dosing parameters]**.

Volume

Volume which is dosed.

Input range	0.00000 to 99999.9 mL
Default value	10.0000 mL

Dosing rate

The rate at which dosing takes place.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

determination is started there is a check whether the solution is contained in the list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether the type of dosing drive matches. For buret units without an integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, the validity of the titer, the working life of the titrant and the time interval for the GLP test of the buret unit are checked at the start of the determination.

Selection	Selection of configured titrants not defined
Default value	not defined

not defined

No check takes place.

Tandem dosing

A second dosing device can be defined under **[Tandem dosing]** in order to enable uninterrupted dosing. Dosing is carried out with a combination of two dosing devices so that the second dosing device is dosing while the first one is being filled and vice versa.

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4 off
Default value	off

Titrant

Selection of the titrant from the titrant list. We recommend always selecting the titrant. Titrants are defined under **System ► Titrants**. You can also enter a name which is not contained in the titrant list. When a determination is started there is a check whether the solution is contained in list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether the type of dosing drive matches. For buret units without integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, only the time interval for the GLP test of the buret unit is checked at the start of the determination.



NOTE

The concentration, validity of the titer and the working life of the titrant are ignored.

Selection	Selection of configured titrants not defined
Default value	not defined

not defined

No check takes place.

Filling rate

Rate at which the dosing cylinder is filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum filling rate depends on the cylinder volume (*see chapter 32, page 356*).

However, in order to ensure uninterrupted dosing, the following points must be taken into account:

- Use the highest possible filling rate in order to keep the filling times as short as possible. The filling rate should be decreased for viscous liquids.
- When you use two buret units with different dosing cylinders, the filling rate for the larger sized cylinder must be at the minimum:

$$v_{2,Fill} \geq v_{1,Fill} \cdot \frac{V_{Cyl.2}}{V_{Cyl.1}}$$

 $v_{2,Fill}$ = filling rate in mL/min for the larger sized cylinder
$$v_{1,Fill} = \text{filling rate in mL/min for the smaller cylinder}$$

$V_{Cyl,2}$ = cylinder volume in mL of the buret unit of the second dosing device

$V_{Cyl.1}$ = cylinder volume in mL of the buret unit of the first dosing device

Example:

Dosing device 1: volume = 20 mL, filling rate = 50 mL/min

Dosing device 2: volume = 50 mL

$$\text{Filling rate } 2 \geq 50 \text{ mL/min} \cdot 50 \text{ mL} / 20 \text{ mL} \geq 125 \text{ mL/min}$$

- The dosing rate must not exceed 75% of the filling rate of the smaller cylinder. These values are listed in the following table, valid at maximum filling rate:

Table 10 Maximum dosing rate for different dosing cylinders

Cylinder volume	maximum dosing rate	
	Exchange unit	Dosing unit
1 mL	2.25 mL/min	–
2 mL	–	5.00 mL/min
5 mL	11.25 mL/min	12.50 mL/min
10 mL	22.50 mL/min	25.00 mL/min
20 mL	45.00 mL/min	50.00 mL/min
50 mL	112.50 mL/min	124.50 mL/min

29.6.3.4 Stirrer

The parameters for the stirrer are edited under **[Stirrer]**.

Stirrer

Selection of the MSB connector the stirrer is connected to. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4 off
Default value	1

off

No stirrer will be used.

Stirring rate

Setting the stirring rate. It can be set in steps of –15 to +15. The default setting **8** corresponds to 1000 rpm. The formula for calculating the rotational speed is specified in *chapter 32.3, page 357*. The optimum stirring rate can be tested in the manual control.

The algebraic sign of the stirring rate changes the direction in which the stirring is done. When the stirrer is viewed from above, this means:

- "+": counterclockwise rotation
- "–": clockwise rotation

Input range	–15 to 15
Default value	8

Change port

The stopcock is moved to the given port, but no piston movement takes place.

Compensate

Because of the fact that the dosing units are interchangeable, the coupling of the Dosino piston rod (spindle) exhibits a low mechanical tolerance that can be noticed when the piston changes its direction of movement. This tolerance can be compensated with this function. A short piston movement is first made in the same direction as the previous movement, which is then followed by a piston movement in the reverse direction.

End volume

The whole cylinder content is ejected. The piston is lowered to the maximum volume mark i.e. until 10'000 pulses have been dosed. This command should be used for pipetting functions for emptying the cylinder.

In-/Outlet

Selection of the port via which the Liquid Handling command is carried out.

Selection	Port 1	Port 2	Port 3	Port 4
Port 1	Default value, when Function = Dose, Aspirate and End volume .			
Port 2	Default value, when Function = Fill, Exchange position and Change port .			
Port 4	Default value, when Function = Eject and Compensate .			

Volume

This parameter can only be edited with **Function = Dose** and **Aspirate**.
Volume which is dosed or aspirated.

Function = Dose:

Input range	0.00000 to 99999.9 mL
Default value	1.00000 mL

Function = Aspirate:

Input range	0.00000 to 50.0000 mL
Default value	1.00000 mL

Flow rate

This parameter can only be edited with **Function = Dose, Fill, Aspirate, Eject, Exchange position** and **End volume**.

Rate at which it is dosed or filled.

Input range	0.01 to 166.00 mL/min
Selection	maximum
Default value	maximum

maximum

The maximum dosing and filling rates depend on the cylinder volume (see chapter 32, page 356).

29.6.4.2 Control instrument

The control instrument the determination is carried out with is selected under **[Control device]**. Control instruments are defined under **System ► Device manager**.



NOTE

This button will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control instrument.

Control device

Selection of the control instrument from the list of devices. Only those devices are displayed which are able to carry out the command.

Selection	Selection of configured control instruments
-----------	---

29.6.4.3 Dosing device

The parameters for the dosing device are edited under **[Dosing device]**.

Dosing device

Selection of the MSB connector to which the dosing device is connected. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

Titrant

Selection of the titrant from the titrant list. We recommend always selecting the titrant. This ensures that the correct data (titer, etc.) is always used for the calculation. Titrants are defined under **System ► Titrants**. You can also enter a name which is not contained in the titrant list. When a determination is started there is a check whether the solution is contained in the list.

For buret units with integrated data chip, a check is made in the method run to verify whether the correct titrant has been attached and whether

the type of dosing drive matches. For buret units without an integrated data chip, the cylinder volume and the type of dosing drive are checked. For the selected titrant, the validity of the titer, the working life of the titrant and the time interval for the GLP test of the buret unit are checked at the start of the determination.

Selection	Selection of configured titrants not defined
Default value	not defined
not defined	
No check takes place.	

29.7 Communication

The following communications commands are available:

Scanning remote lines (SCAN)	Defining remote signals which are awaited before the next method command is started.
Setting remote lines (CTRL)	Defining remote signals.
Receiving RS-232 commands (SCAN RS)	Defining RS-232 commands which are awaited before the next method command is started.
Sending RS-232 commands (CTRL RS)	Defining the RS-232 commands which are sent.

29.7.1 Scanning remote lines (SCAN)

The command **SCAN** can be used to define input signals at the remote interface which are awaited before the next method command is started.

Control device

This parameter will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control device. Control devices are defined under **System ► Device manager**.

Selection of the control device from the list of devices.

Selection	Selection of configured control devices
-----------	--

Remote Box

Selection of the MSB connector the Remote Box is connected to. The selection does not depend on the number of MSB connectors of the control instrument. The remote signal defined is requested on this Remote Box.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

29.7.2 Setting remote lines (CTRL)

The command **CTRL** can be used to define output signals at the remote interface.

Control device

This parameter will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control device. Control devices are defined under **System ► Device manager**.

Selection of the control device from the list of devices.

Selection	Selection of configured control devices
-----------	---

Remote Box

Selection of the MSB connector the Remote Box is connected to. The selection does not depend on the number of MSB connectors of the control instrument. The remote signal defined is sent via this Remote Box.

Selection	1 2 3 4
Default value	1

Output signal

Selection of the signal out of the templates or entering the required bit pattern. Templates are defined under **System ► Templates ► Output lines**.



NOTE

A line set active is not being reset automatically, not even at the end of the determination.

Entering a bit pattern:

- 0 = line inactive
- 1 = line active
- * = retain line status
- p = set pulse

The output lines are always numbered from right to left, i.e. with the signal *******1** line 0 is set. With a pulse, the length is set to 200 ms. If it is required to set pulses with other lengths, you have to define a corresponding template.



NOTE

We recommend masking lines that are of no interest or for which no defined condition can be predicted with an asterisk (*).

Entry	Bit patterns containing exactly 14 characters or a max. of 24 characters for the name of the template
Default value	*****
Selection	Selection of the templates defined

29.7.3 Scanning the RS-232 interface (SCAN RS)

The command **SCAN RS** can be used to define RS-232 commands which are awaited before the next method command is started.

Serial port

Selection of the serial interface the peripheral device is connected to. The RS-232 command defined is requested on this interface.

Selection	COM1 COM2 COM3 COM4 COM5 COM6 COM7 COM8
Default value	COM1

Character string

Input of the RS-232 command as character string. All characters of the ASCII code page can be used. Control characters (Esc, FF etc.) have to be entered as a three-place, decimal ASCII code, starting with a slash. Each command is automatically terminated with the ASCII characters **CR** and **LF**. The * can be used for one or more characters.

Entry	ASCII string with a max. of 24 characters
Default value	!*R

Timeout

When this time interval has expired without recognizing the RS-232 command, the action defined in the following is triggered.

Input range	0 to 999 s
Default value	0 s

Action

Selection of the action which is carried out when the time interval has expired.

Selection	Display message Document message Cancel determination
Default value	Display message For all three options it is documented in the determination data (see dialog More determination data / Messages), that the time interval has been expired.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.

Document message

In the determination data it will be documented, that the time interval has been expired.

Cancel determination

The determination is stopped.

29.7.4 Defining RS-232 commands (CONTROL RS)

The command **CONTROL RS** can be used to define the RS-232 commands which are sent.

Serial port

Selection of the serial interface the peripheral device is connected to. The RS-232 command defined is sent via this interface.

Selection	COM1 COM2 COM3 COM4 COM5 COM6 COM7 COM8
Default value	COM1

Character string

Input of the RS-232 command as character string. All characters of the ASCII code page can be used. Control characters (Esc, FF etc.) have to be entered as a three-place, decimal ASCII code, starting with a slash. Each command is automatically terminated with the ASCII characters **CR** and **LF**.

Entry	ASCII string with a max. of 24 characters
Default value	&M;\$G

29.8 Automation

29.8.1 Rotating sample rack (MOVE)

The command **MOVE** is used to move to a rack position or an external position.

Tower

Selection of the tower, with which the command is carried out. Both Tower 1 and Tower 2 can always be selected, even if your Sample Processor has only one tower.

Selection	1 2
Default value	1

1

View from the front, the right tower.

2

View from the front, the left tower.

Destination

Selection of the required target position.

Selection	Sample Ext. position Special beaker Rack position Current sample + Current sample - Next position Prev. position Calibration pos. rotate + rotate - swing + swing -
Default value	Sample

Sample

Rack position, which is defined by the sample variable. The sample variable is automatically increased by 1 after every method run or it can selectively be modified with the command **SAMPLE**.

Ext. position

One of the four external positions which can be moved to with the Swing Head. Enter the number of the external position (**1...4**) in the field to the right.

Special beaker

Special beakers, which are defined in the rack table, can be directly moved to. Enter the number of the special beaker (**1...16**) in field to the right.

Rack position

Any rack position. Enter the number of the rack position (**1...999**) in the field to the right.

Current sample +

Starting from the current sample (defined by the sample variable) the rack is moved forward by the number of rack positions (**1...999**) entered in the field to the right.

Current sample -

Starting from the current sample (defined by the sample variable) the rack is moved backward by the number of rack positions (**1...999**) entered in the field to the right.

Next position

Starting from the current rack position, the rack is moved forward by one position.

Prev. position

Starting from the current rack position, the rack is moved backward by one position.

Calibration pos.

For automatic calibrations with a Sample Processor (relevant only if the Ti-Touch supports calibrations).

rotate +

Moving the rack forward by a certain increment. The rotation increment is defined in the properties of the tower.

rotate -

Moving the rack backward by a certain increment. The rotation increment is defined in the properties of the tower.

swing +

Swinging the robotic arm outward by a certain increment (towards higher angular degrees). The swing increment is defined in the properties of the Swing Head.

swing -

Swinging the robotic arm towards the center of the rack by a certain increment (towards lower angular degrees). The swing increment is defined in the properties of the Swing Head.

Beaker test action

This parameter can only be edited with **Destination = Sample, Next position** or **Prev. position**.

Selection of the action which is carried out when the beaker sensor does not detect any vessel at the sample position moved to. In addition, enable the beaker sensor in the rack table.

Selection	Rotate rack Display message
Default value	Display message



Rotate rack

The rack is moved to the next vessel on a sample position. **Exception:** If the sample variable is newly defined in the method before this command (command **SAMPLE**), then the option **Display message** is automatically valid.

Display message

A message is displayed. You can select whether you want to continue with the determination or cancel the run.



NOTE

This parameter is ignored with the 885 Compact Oven SC.

Options

Shift rate

Rate at which the sample rack is moved.

Input range	3 to 20 °/s
Default value	20 °/s



NOTE

This parameter is ignored with the 885 Compact Oven SC.

Shift direction

Direction in which the rack is moved.

Selection	auto + -
Default value	auto

auto

A shift direction, with which the shorter way has to be passed, is automatically selected.

+

Counterclockwise rotation.

—

Clockwise rotation.

Swing rate

Rate at which the robotic arm is swung.

Input range	10 to 55 °/s
Default value	55 °/s

29.8.2 Moving the lift (LIFT)

The lift of a Sample Processor is moved with the command **LIFT**. The movement can, however, only be carried out if the Sample Processor is located at a valid rack position. This will not be the case, for example, after a Rack Reset.

Tower

Selection of the tower, with which the command is carried out. Both Tower 1 and Tower 2 can always be selected, even if your Sample Processor has only one tower.

Selection	1 2
Default value	1

1

View from the front, the right tower.

2

View from the front, the left tower.

Lift position

Selecting a predefined lift position or enter any lift position. The predefined lift positions are defined in the device properties of the Sample Processor. They can, however, also be assigned in the manual control.

Input range	0 to 235 mm
Selection	Work position Shift position Rinse position Special position Home position
Default value	Work position

Home position

The lift is being located at the upper stop position, i.e. at the position 0 mm.

For the 885 Compact Oven SC, the position 0 mm is approached when **Shift position**, **Rinse position** and **Special position** is selected.

Lift rate

Speed at which the lift is moved.

Input range	5 to 25 mm/s
Default value	20 mm/s



NOTE

This parameter is ignored with the 885 Compact Oven SC.

29.8.4 Resetting the rack (RACK)

The following actions are carried out with the command **RACK**:

- Rack, lift and robotic arm are being reset.
- The rack code of the rack attached is being read out and the corresponding rack data is transferred to the Sample Processor.
- The sample variable is reset to the value 1.

Check rack

on | off (Default value: **off**)

Activate this parameter when the rack attached is to be checked. For this you, however, have to additionally select the rack in the dialog **Method options / Start options** at the option **Check rack**.



NOTE

This option is not available with the 885 Compact Oven SC. Therefore, this parameter must be set to **no**.

29.8.5 Defining the sample variable (SAMPLE)

The sample variable describes the current position of the sample on the rack of the Sample Processor. It is automatically increased by 1 after the end of a method run. In the following cases, the sample variable is automatically reset to the value 1:

- if the autostart counter is reset.
- if the sample rack is reset.

The sample variable can be altered specifically with the command **SAMPLE**.

Sample variable

Modifying the sample variable.

Selection	= + -
Default value	+

$$=$$

The sample variable corresponds to the number entered in the field **Value**.

+

The sample variable corresponds to the current value plus the number entered in the field **Value**.


—

The sample variable corresponds to the current value minus the number entered in the field **Value**.

29.8.6.1 Loading / saving a subsequence

[Load/ Save] is used to load stored subsequences or to save new subsequences. This means that they are available for all methods and only have to be created once.

Loading a subsequence

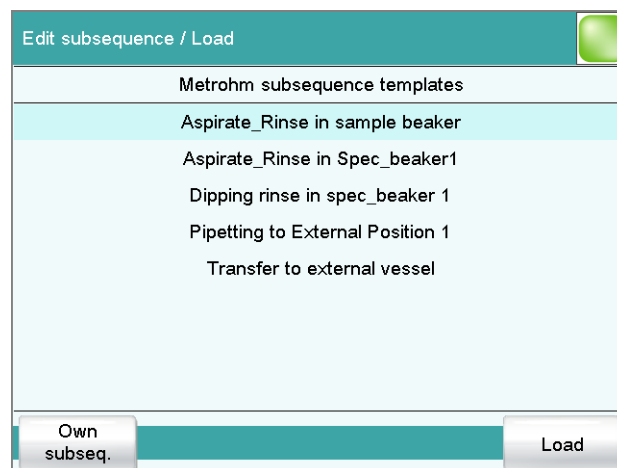

[Own subseq.]

Opening the list of own subsequences.

[Load]

Loading the selected subsequence.

Saving a subsequence

Custom-made subsequences are saved in the internal memory.



NOTE

In order to be able to use the subsequences also for other titration systems, you have to create a backup. You can only restore the subsequences from this backup on another system.

File name

File name of the subsequence.

Entry	max. 32 characters
-------	---------------------------

[Save]

Saving the subsequence to the internal memory.

[Display options]

Definition when the note is being displayed.

Automatically after loading the subsequence

on | off (Default value: **off**)

If this parameter is activated, the note is displayed when loading the subsequence. It can otherwise only be read in this dialog.

29.8.7 Controlling the heating (HEATER)

The heating on the 885 Compact Oven SC is controlled with the **HEATER** command.

Status

Switch the heater on or off.

Selection	on off
Default value	on

on

The heater is switched on. The command is displayed in the method run until the specified temperature is reached. The heater runs in the background until it is explicitly switched off. The current temperature can be seen in the Manual control.

off

The heater is switched off.

Temperature

Temperature to which the oven is heated.

Input range	50 to 250 °C (Increment: 1)
Default value	110 °C

29.8.8 Controlling the gas flow (FLOW)

The gas flow is controlled on the 885 Compact Oven SC with the **FLOW** command.

Status

Switch the gas flow on or off.

Selection	on off
Default value	on

on

The gas flow is switched on. The gas flow runs until it is switched off. The current gas flow can be seen in the Manual control.

off

The gas flow is switched off.

Flow rate

Flow rate for the gas flow.

Input range	10 to 150 mL/min (Increment: 1)
Default value	50 mL/min

Gas supply

Selection of the gas supply.

Selection	pump valve
Default value	pump

pump

The integrated air pump is used for the gas flow.

valve

The integrated inlet valve is used for the gas flow. The inlet valve in the 885 Compact Oven SC is connected to compressed air or nitrogen.

Gas type

Gas type of the carrier gas used.

Selection	air nitrogen
Default value	air

29.8.9 Starting the measurement for temperature and gas flow (MEAS T/F ON)

Data acquisition and evaluation of the data of the 885 Compact Oven SC is started with the **MEAS T/F ON** command.

The following values are output as raw data and are available as variables:

- Mean gas flow (evaluation starts after 10 s).
- Mean temperature
- Minimum temperature
- Maximum temperature

After a **MEAS T/F ON** command, there must also be a **MEAS T/F OFF** command present in the method run. If one of the two commands is missing or the order is not followed, the command is depicted in red. No parameters can be edited for this command.

29.8.10 Exiting the measurement for temperature and gas flow (MEAS T/F OFF)

The data acquisition evaluation of the data of the 885 Compact Oven SC is finished with the **MEAS T/F OFF** command. There must be a **MEAS T/F ON** command present before a **MEAS T/F OFF** command in the method run. The command will be depicted in red if this is not fulfilled. No parameters can be edited for this command.

29.9 Miscellaneous commands

29.9.1 Controlling a stirrer (STIR)

A stirrer is controlled with the command **STIR**.

Control device

This parameter will only be shown when, in addition to the Ti-Touch, a Sample Processor has been configured as control device. Control devices are defined under **System ► Device manager**.

Selection of the control device from the list of devices.

Selection	Selection of configured control devices
-----------	---

Stirrer

Selection of the MSB connector the stirrer is connected to. The selection does not depend on the number of MSB connectors of the control instrument.

Selection	1 2 3 4
Default value	1

Status/Duration

Switching the stirrer on or off.

Input range	0 to 999 s (Increment: 1) During this time the stirrer is running.
Selection	on off
Default value	on

on

The stirrer is switched on. The stirrer runs until it is explicitly switched off. If the determination is canceled with the key [], the stirrer is also switched off.

off

The stirrer is switched off.

Stirring rate

Setting the stirring rate. It can be set in steps of -15 to $+15$. The default setting **8** corresponds to 1000 rpm. The formula for calculating the rotational speed is specified in *chapter 32.3, page 357*. The optimum stirring rate can be tested in the manual control.

The algebraic sign of the stirring rate changes the direction in which the stirring is done. When the stirrer is viewed from above, this means:

- "+": counterclockwise rotation
- "-": clockwise rotation

Input range	-15 to 15
Default value	8

29.9.2 Pausing the method run (WAIT)

The method run can be paused with the command **WAIT**.

Hold sequence

on | off (Default value: **off**)

If this parameter is activated, then the method run is paused until it is continued manually.

Waiting time

This parameter can only be edited when **Hold sequence** is deactivated.

The method run is automatically continued after this waiting time.

Input range	0 to 99999 s
Default value	30 s

Message

on | off (Default value: **off**)

If this parameter is activated, the text message defined in the following is displayed during the waiting time.

Message text

This parameter can only be edited when **Message** is activated.

Text which is displayed during the waiting time.

Entry	max. 28 characters
Default value	empty

29.9.3 Scan data (REQUEST)

The following data can be scanned in the method run with the command **REQUEST**:

- Sample data
 - Sample identification 1
 - Sample identification 2
 - Sample size (value and unit)
- Common variable

You can select thereby whether the sequence should be paused or resumed in the background.

Sample ident.

Selection of the sample identification that is queried in the method run.

Selection	off Identification 1 Identification 2 ID1 & ID2
Default value	off

Sample size

on | off (Default value: **on**)

If this parameter is activated, then the value for the sample size will be requested.

Sample size unit

on | off (Default value: **off**)

If this parameter is activated, then the unit for the sample size will be requested.

Common variable

Selection of the common variable which is queried in the method run.

Selection	CV01...CV25 off
Default value	off

Hold sequence

on | off (Default value: **on**)

If this parameter is activated, then the run will be paused during the request. If this parameter is deactivated, then the method continues in the background until the following measurement or titration has been completed.

29.9.4 Defining an acoustic signal (BEEP)

The **BEEP** command can be used to produce an acoustic signal.

Duration

Approximate playing time in seconds of the permanently stored melody.

Input range	1 to 9
Default value	1

29.9.5 Signing a determination (SIGN)

A determination can be signed with the **SIGN** command. The method run is automatically paused as soon as this command is reached. At the same time a switch is made to the result view. The run will not be continued until the determination has been signed or the **SIGN** command has been canceled. You can specify whether a routine user is allowed to cancel the **SIGN** command in the setting **System ► System settings ► Dialog options ► Routine dialog**. No parameters can be edited for this com-

mand. The parameters for signing determinations are described in *chapter 17.1, page 182*.



NOTE

Methods with this command can only be started if you are working with login and password protection.

The command should always be inserted at the end directly before the **REPORT** command or at the end of the method sequence if a PC/LIMS report is to be sent automatically.

29.9.6 Canceling the method run (END)

The method run will cancel as soon as the **END** command is reached. This is a good idea if you wish to test only the first part of a method. No parameters can be edited for this command.

30.1 System initialization



We recommend creating a backup of the system at regular intervals in order to avoid data losses.

Initialization is started. After the initialization has been completed, the 915 KF Ti-Touch is started automatically.

31 Troubleshooting

31.1 Editing methods

Problem	Cause	Remedy
The method command cannot be inserted.	No control instrument which supports the method command is configured in the device manager.	Add a control instrument which supports the desired method command to the list of devices (see chapter 11, page 80).
The additional evaluations (EVAL commands) cannot be inserted.	A titration or measuring command to which the evaluation can be applied is missing.	1. Insert a titration or measuring command. 2. Insert the EVAL command directly after it.
A MEAS T/F ON command or MEAS T/F OFF command is depicted in red in the command list.	The corresponding MEAS T/F ON command or MEAS T/F OFF command is missing or not inserted in the correct order in the method run. To be able to execute a method, both MEAS T/F commands are needed in the method run.	Insert the missing MEAS T/F command and observe the logical order of the commands.
An EVAL command is depicted in red in the command list.	The associated titration or measuring command has been deleted or replaced by a different one.	1. Delete the EVAL command. 2. Insert a titration or measuring command. 3. Insert the EVAL command directly after it.

31.2 Sample table

Problem	Cause	Remedy
The method cannot be entered in the sample table.	The option Use sample assignment table is activated. The method assigned to the sample identification will be loaded automatically at the start of the determination.	Deactivate the option (see chapter 15.1, page 150).

31.3 Results/Statistics

Problem	Cause	Remedy
No statistics are carried out for a result.	<i>Statistics is not activated.</i>	In the method options, activate the option Statistics (see chapter 16.5, page 170).
	<i>No statistics variable has been assigned to the calculation.</i>	In the result options (command CALC), define a variable for the mean value (see "Variable for mean value", page 156).
	<i>Statistics is not activated.</i>	In the control dialog, activate the Statistics option (see chapter 18, page 184).
A result has been removed from the statistics, the mean value has not been reassigned to the common variable or to the TITER variable.	<i>Assignment is not carried out automatically with retroactive modifications.</i>	Recalculate the determination manually (see chapter 24, page 212).
The result is not displayed in the result table.	<i>The column display in the result table is incorrectly configured.</i>	Modify the settings for the display accordingly (see chapter 26.1, page 235).
	<i>The parameter Save result in result table is not activated.</i>	In the result options (command CALC), activate the parameter (see "Save result in result table", page 157).

31.4 Printing

Problem	Cause	Remedy
Tapping on the fixed key [] has no effect.	<i>The fixed key is disabled.</i>	In the Dialog options / Fixed keys dialog, activate the Print option (<i>see chapter 7.2, page 35</i>).
	<i>A determination is running.</i>	Wait until the determination is finished.

31.5 Manual control

Problem	Cause	Remedy
Tapping on the fixed key [] has no effect.	<i>The fixed key is disabled.</i>	In the Dialog options / Fixed keys dialog, activate the Manual control option (see chapter 7.2, page 35).
A button in the manual control is disabled.	<i>The required hardware is either not connected or incorrectly connected.</i>	1. Switch off the 915 KF Ti-Touch. 2. Connect the hardware correctly. 3. Switch the 915 KF Ti-Touch back on.
	<i>The required hardware is being used in the ongoing determination.</i>	Wait until the determination is finished.

31.6 Titration stand with pump

Problem	Cause	Remedy
Solvent is added at the same time during aspiration.	<i>The adsorber tube is blocked.</i>	▪ Enlarge the hole of the drying tube cover to 2 mm at least. ▪ Refill the adsorber tube, pack not too tightly, perhaps remove some of the molecular sieve.
The pump conveys no liquid.	<i>The add/aspirate key is not pressed tightly enough.</i>	▪ Usually the bottle attachments are not screwed tightly enough to the reagent and waste bottles. ▪ Check all connections for leak-tightness.
The solvent flows on into the titration vessel after the addition, without the add key being pressed.	<i>The solvent bottle is not positioned correctly.</i>	Position the solvent bottle in such a way that the liquid level in the bottle is lower than that in the titration vessel.

31.7 Miscellaneous

Problem	Cause	Remedy
The 915 KF Ti-Touch cannot be switched on.	<i>The Ti-Touch is connected to an MSB socket on the control instrument.</i>	Plug in the Ti-Touch at the Controller connector (see manual for control instrument).

31.8 Volumetric Karl Fischer Titration

Problem	Cause	Remedy
The drift is very high during conditioning.	<i>The titration cell is leaking.</i>	- Check the seals and the septum. Replace if necessary. - Replace the molecular sieve.
The drift becomes greater after each titration.	<i>The sample releases water very slowly.</i>	- Adjust the method. - Add solubility promoter. - Increase the temperature (possibly using a KF oven). - See technical literature.
	<i>A side reaction is taking place.</i>	- Use special reagents. - Adjust the method (increase/decrease the temperature, external extraction). - See technical literature.
	<i>The pH value is no longer in the optimum range.</i>	Add buffer, see technical literature.
The titration will not be finished.	<i>The titration cell is leaking.</i>	- Check the seals and the septum. Replace if necessary. - Replace the molecular sieve.
	<i>The minimum increment is too low.</i>	Define Titration rate = user and increase the minimum volume increment (Min. increment) (see chapter 29.1.2, page 279).
	<i>The stop criterion is unsuitable.</i>	Adjust the control parameters (see chapter 29.1.2, page 279): - Increase the stop drift. - Select a short delay time.
	<i>See also: The drift becomes greater after each titration.</i>	
The sample is over-titrated.	<i>The increments at the end of the titration are too high.</i>	Define Titration rate = user and reduce the dosing rate (Max. rate) (see chapter 29.1.2, page 279). The following experiment provides a clue for the optimum dosing rate: During conditioning, display the drift and add sample without starting the titration. Select a value below the highest drift as dosing rate.

Problem	Cause	Remedy
		Stir faster.
	<i>The amount of methanol in the working medium is too low.</i>	- Replace the working medium. - Reduce the amount of solubility promoter, if working with solvent mixtures, see technical literature.
	<i>The electrode may be covered.</i>	Wipe off the electrode with ethanol or a suitable solvent.
The solution becomes darker after each titration.		Replace the working medium.
	<i>The electrode may be covered.</i>	Wipe off the electrode with ethanol or a suitable solvent.
	<i>The electrode has a short circuit.</i>	- Check the Pt wires. - Activate the electrode check.
The endpoint is reached too quickly.	<i>The dosing rate outside the control range is too high.</i>	Define Titration rate = user and reduce the dosing rate (Max. rate) (see chapter 29.1.2, page 279).
The titration times with volumetric titration are constantly longer.	<i>The buffer capacity of the solvent can be exhausted for two-component reagents.</i>	Replace the working medium.

32.2 Exchange unit

32.2.1 Maximum dosing and filling rate

The maximum dosing rate and maximum filling rate for the exchange unit depend on the cylinder volume:

Cylinder volume	Maximum rate
1 mL	3.00 mL/min
5 mL	15.00 mL/min
10 mL	30.00 mL/min
20 mL	60.00 mL/min
50 mL	150.00 mL/min

Independent of the cylinder volume, values ranging from 0.01 to 166.00 mL/min can always be entered. When the function is carried out the rate will be, if necessary, decreased automatically to the highest possible value.

32.2.2 Default parameters for preparing (PREP)

The **PREP** command is used to rinse and fill the cylinder and tubings of the exchange unit air bubble-free. You should carry out this function before the first determination or once per day.

If no titrant is selected in the command, preparing will be carried out with the following default parameters:

- The entire cylinder volume is dosed twice at the maximum dosing rate.
- The cylinder is filled at the maximum filling rate.

32.3 Stirring rate

The stirring rate can be adjusted in steps from -15 to +15.

The approximate rotational speed for the internal magnetic stirrer (depends on the product version) can be calculated with the following formula:

$$\text{Rotational speed/min (r/min)} = 125 \cdot \text{Stirring rate}$$

Example:

Configured stirring rate: 8

Rotational speed in revolutions per minutes = $125 \cdot 8 = 1,000$

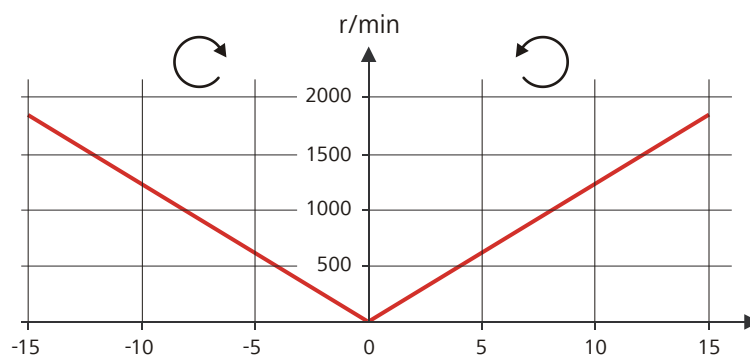


Figure 16 Rotational speed depending on the stirring rate

The information on the separately connectable 802 propeller stirrer can be found in the "802 Stirrer" manual.

32.4 Balance

The sample size and the associated unit can be sent from a connected balance. The sample size is transmitted as a number with up to ten characters (including algebraic sign and decimal point).

Sample size and unit are sent as a single character string. They are separated by a space character. The string is terminated with the ASCII characters **CR** and **LF**.

If the balance sends a negative sample size (e.g. when you are reweighing a sample), then the algebraic sign is adopted. The algebraic sign is, however, ignored for the calculations.



NOTE

With some balances, the sample identification and the method can be sent in addition to the sample size.

Make sure that the balance does not send the sample size until the end.

Mettler AX

For the Mettler AX balance, the fields that contain the sample identification or the method must be designated as follows:

- Designation for the field with the method name: **METHOD**
- Designation for the field with sample identification 1: **ID1**
- Designation for the field with sample identification 2: **ID2**

32.5 Remote interface

The 6.2148.010 Remote Box allows devices to be controlled which cannot be connected directly to the MSB interface of the Ti-Touch.

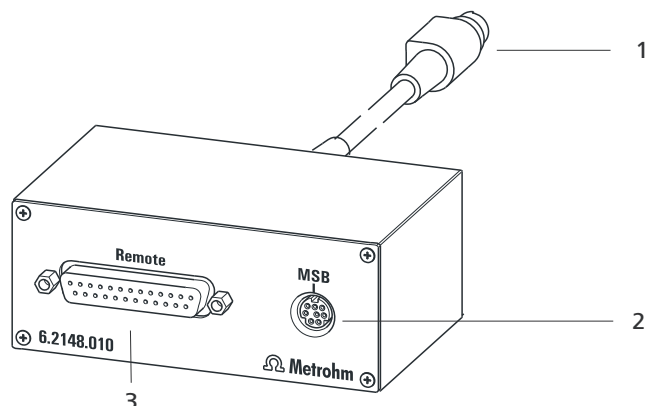


Figure 17 Connectors of the Remote Box

1 Cable

For connecting to the Ti-Touch.

2 MSB connector

Metrohm Serial Bus. For connecting external dosing devices or stirrers.

3 Remote connector

For connecting instruments with a remote interface.

32.5.1 Pin assignment of the remote interface


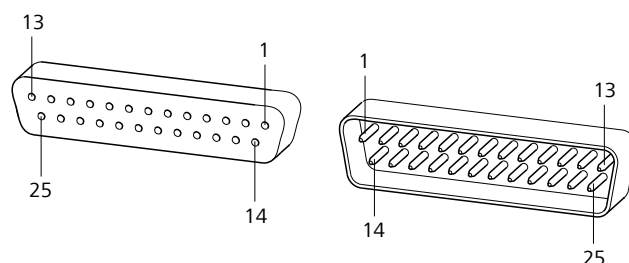


Figure 18 Pin assignment of remote socket and plug

The above figure of the pin assignment of a Metrohm remote interface applies not only for the Remote Box, but also for all Metrohm devices with 25-pin D-Sub remote connector.

Inputs

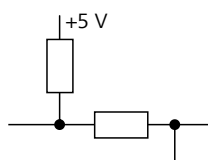

approx. 50 k Ω Pull-up

$t_p > 20$ ms

active = low, inactive = high

Assignment	Pin No.	Function
Output 13	20	
0 volts / GND	14	
+5 volts	15	
0 volts / GND	25	

Table 12 Explanation of the individual functions

Function	Explanation
Start	The current method is started at the time of activation. $t_{\text{pulse}} > 100 \text{ ms}$
Stop	The current method is canceled (Stop) at the time of activation. $t_{\text{pulse}} > 100 \text{ ms}$
Quit	The current command in the determination run will be canceled at the time of activation. $t_{\text{pulse}} > 100 \text{ ms}$
Ready	The instrument is ready to receive a start signal.
Conditioning OK	The line is set when Conditioning with KFT titration is at OK. The line remains set until the determination is started with [START] .
Determination	The instrument performs a data-generating determination.
EOD	End of Determination. Pulse ($t_{\text{pulse}} = 200 \text{ ms}$) after a determination.
Error	The line is set for error message display.
Warning	The line is set for warning message display.

32.6 Result variables as parameter setting

For most of the numerical parameters, the result of a calculation can also be utilized in place of a number. This result must be calculated in a preceding calculation command (command **CALC**). A typical application is the use of a relative start volume.

The following example is used to explain how this method is used to define a sample size-dependent start volume for a titration. Proceed as follows:

1 Create a new calculation

- Insert the calculation command **CALC** in front of the titration command.
- Enter the calculation formula, e.g. $R1 = C00*3$.

01

CALC

Calculation

Result name

Relative start volume

Calc. formula R1

C00*3

Decimal places

3

▼

Result unit

mL

▼

Note

Result variable

Result limits

Result options

2 Define a result variable as parameter setting

- In the titration command, enter the variable **R1** as the start volume.

The selection of the result variables is displayed in the number editor.



NOTE

If the method before the command with the result variables applied contains more than one calculation commands (possibly with result variables bearing the same name), then the result variable of the preceding calculation command will always be used.

Please note that the result must be within the input range of the parameter, as otherwise the determination run will be stopped when the parameter is accessed.


32.7 Using AuditTrailViewer

When you work in accordance with the *FDA Guidance 21 CFR Part 11*, logs with detailed user activities (Audit Trail) are important. The *AuditTrailViewer* supports you in your efforts to compile, filter and export logs of this type as TXT files.

First you must install the *AuditTrailViewer* on your PC. Afterwards, you must create a backup of your Ti-Touch (see chapter 12.3, page 130). The Audit Trail is generated automatically when you create the backup. The file name of the Audit Trail is **log.madt**. In the following you will find detailed information about how you must proceed.

32.7.1 Installing AuditTrailViewer

You will find the file for installing the *AuditTrailViewer* on the USB flash drive provided.

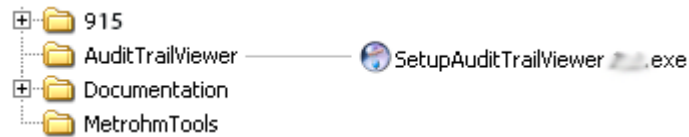


Proceed as follows:

- 1 Double-click on the file **SetupAuditTrailViewer 1.2.3.exe**.
The dialog **Welcome to the AuditTrailViewer Setup Wizard** is displayed.
- 2 Click on **[Next >]**.
The **License Agreement** dialog is displayed.
- 3 Click on **[I Agree]** to accept the contractual conditions.
The **Choose Install Location** dialog is displayed. A memory location for the *AuditTrailViewer* is suggested in this dialog.
- 4 If you accept the suggested memory location, click on **[Next >]**.
or
Select a different memory location and click on **[Next >]**.
The **Confirm Installation** dialog is displayed.
- 5 Click on **[Install]** to install the *AuditTrailViewer*.
As soon as the installation is completed, the dialog **Installation Complete** is displayed.
- 6 Click on **[Next >]**.
The dialog **Completing the AuditTrailViewer Setup Wizard** is displayed.
- 7 If the *AuditTrailViewer* should not be started immediately, deactivate the check box **Run AuditTrailViewer**.
- 8 Click on **[Finish]**.
If the preceding instruction step has been skipped, then the *AuditTrailViewer* will start.

32.7.2 Opening the Audit Trail

Create a backup of the data and of the settings of your Ti-Touch (see chapter 12.3, page 130).

1 Open the backup on your PC.

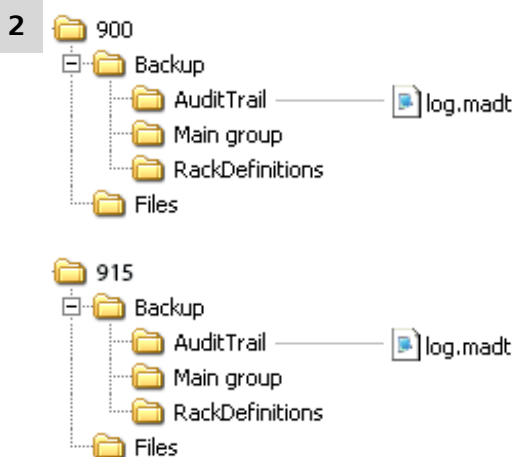

Open the file **log.madt** with the *AuditTrailViewer*.

Column	Content										
Symbol	Classification of the entry:  Actions that are neither relevant to safety nor alter the determination data.  Actions, e.g. changes to the login options and recalculation of determinations.  Errors that occur, e.g. entry of an incorrect password.										
Date	Precise time of the event.										
User	User who triggered the action.										
Category	Category to which the entry belongs.										
Action	Designation of the action.										
Details	Details of the action.										
<table> <tr> <td colspan="2">Menu bar</td></tr> <tr> <td> Print</td><td> Print Audit Trail. An Audit Trail can be printed by means of the  symbol. </td></tr> <tr> <td> Update</td><td> Update view. The list is updated automatically only at the time the dialog window is opened. </td></tr> <tr> <td> Quick filter</td><td>Filter the Audit Trail.</td></tr> <tr> <td> Show all</td><td>Display all entries once again.</td></tr> </table>		Menu bar		 Print	Print Audit Trail. An Audit Trail can be printed by means of the  symbol.	 Update	Update view. The list is updated automatically only at the time the dialog window is opened.	 Quick filter	Filter the Audit Trail.	 Show all	Display all entries once again.
Menu bar											
 Print	Print Audit Trail. An Audit Trail can be printed by means of the  symbol.										
 Update	Update view. The list is updated automatically only at the time the dialog window is opened.										
 Quick filter	Filter the Audit Trail.										
 Show all	Display all entries once again.										

32.7.4 Filtering the Audit Trail

The entries in the Audit Trail can be filtered. You can set the following filter criteria:

- **Date**
- **User**
- **Category**
- **Action**

Proceed as follows:

- 1 Click in a cell that contains the desired filter criterion (e.g. **Category** = **Method**).

- **Touch screen test**
Test the function of the touch-sensitive screen (*see chapter 32.8.5, page 371*).
- **Software update**
Update the software of Ti-Touch and of the control instruments (*see chapter 32.8.6, page 372*).
- **Service**
Special functions for the service technician (*see chapter 32.8.7, page 375*).

32.8.1 LCD test

System ► Diagnosis ► LCD test

You can use the LCD test to check the display for faulty pixels. To do this, various test images are displayed one after the other.



NOTE

The [↔] fixed key is used to display the preceding test image; the test can be stopped at any time with the [🏠] fixed key.

Proceed as follows:

- 1 ▪ Tap on **[LCD test]**.
The entire display turns white.
- 2 ▪ Check the display for pixel errors and other irregularities.
 - Continue the test with [▶].
 Other colors and patterns will be displayed one after the other.
- 3 ▪ Repeat step 2 until the **System / Diagnosis** dialog is displayed again.
- 4 ▪ Contact your Metrohm representative in the event of faults or irregularities.

32.8.2 Formatting an external storage medium

System ► Diagnosis ► Format storage medium

With this function, you have the option of formatting an external storage medium directly on the 915 KF Ti-Touch (quick formatting).

The following file systems can be selected:

- **FAT**

- FAT32
- ExFAT



CAUTION

If you format the storage medium with the **ExFAT** file system, it is possible that it will no longer be recognized by a PC. Read the pertinent information from Microsoft Support at <http://support.microsoft.com>.

Do not use this file system unless you need to be able to save more than 999 files in a single group.

32.8.3 Removing an external storage medium

System ► Diagnosis ► Remove storage medium

As long as the stored data is not being accessed (reading/saving data), you can plug in and disconnect the storage medium without difficulty at any time. The **[Remove storage medium]** function offers additional protection. This ensures that the storage medium cannot be removed while any data is being transferred.

32.8.4 Adjusting the touch screen

System ► Diagnosis ► Touch adjustment

With time the touch-sensitive screen may not register touches at their precise position. You can readjust the screen in this dialog.

Proceed as follows:

- 1 In the main dialog, tap on **[System]**.

- 2** Tap on **[Diagnosis]**.

- 3** Tap on **[Touch adjustment]**.

A message in English and a crosshair will appear in the center of the screen.

- 4 Touch a stylus (a special pen-shaped instrument for operating devices with touch screens) to the center of the crosshair and maintain this position until the cross hair moves to a new position on the screen.

- 5 Repeat this process as long as the crosshair continues to be repositioned on the screen.

A message in English will appear.

6 Tap at any location on the screen to confirm the adjustment.

7 Press [] to finish the adjustment.

32.8.5 Testing the touch screen

System ► Diagnosis ► Touch screen test

You can use this test to check the function of the touch-sensitive monitor screen.

A calibration is not possible with this test. If the Touch screen does not work correctly anymore, it will need to be replaced. For that please contact your Metrohm representative.



NOTE

You can use the [] fixed key to stop the test at any time.

Proceed as follows:


1 Tap on **[Touch screen test]**.

A red square is displayed at the upper left.

2 Tap on the square as close to the center as possible.

The actual detected point of contact will be displayed by a crosshair, and the deviation from the center will also be displayed on the lower edge of the image.

The next red square will also be displayed.

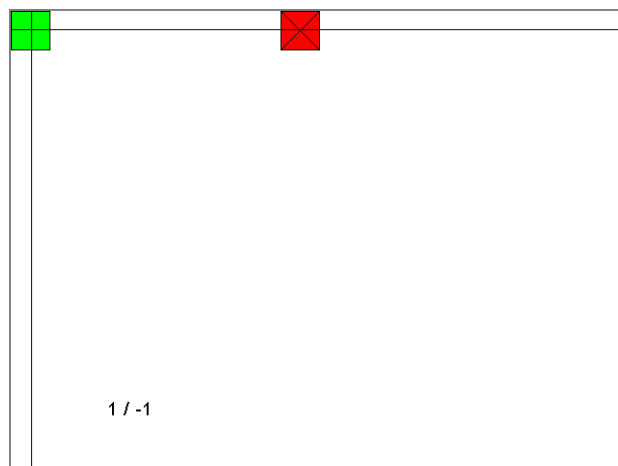

Program files

The files are instrument-specific. The file names are structured as follows:

- Control software for 915 KF Ti-Touch:
5XXXyyyy.BIN
5XXXyyyy.MBIN
 - XXX = device type (i.e. "915" for the 915 KF Ti-Touch)
 - yyyy = program version
- Firmware of the control instrument:
5XXXyyyy.BIN
 - XXX = device type (e.g. 814 for the 814 USB Sample Processor)
 - yyyy = program version

Language files

Language files may have the following content:

- one or more additional dialog languages
- the online help for one or more additional dialog languages
- supplements for existing dialog languages/online help documents

They can be recognized by means of the two-digit language code in the file name. The file name has the following structure:

- **5XXXyZZZML.BIN**
5XXXyZZZML.MBIN
 - XXX = device type (i.e. "915" for the 915 KF Ti-Touch)
 - y = version number of the language package
 - ZZZ = program version

32.8.6.1 Carrying out a software update



CAUTION

Make sure that the power supply is guaranteed during the entire update process. Otherwise there is the chance that the instrument will no longer be able to be switched on and that it will need to be sent in for repair.



NOTE

Make sure that no USB/RS-232 adapter is connected!

Proceed as follows:


1 Copying files to an external storage medium



NOTE

In the case of the program versions and language files for the Ti-Touch, it is imperative that both the BIN file and the associated MBIN file be copied.

- Copy the files to the "SwUpdates" directory.
 - Copy also the "MetrohmTools" folder with the "MUpdate-Tool.exe" file to the external storage medium.
- Make sure that the file is stored at the top level of the external storage medium. This also applies to the "915" folder.

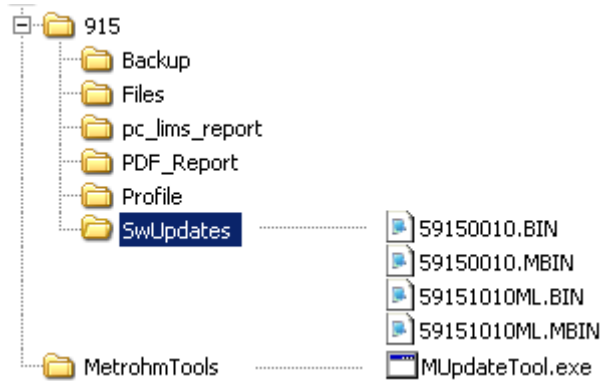

For the program versions 5.915.0043 to 5.915.0045, an additional folder that is on the same level as Metrohm-Tools with the name SRVPRF900 containing the file St900_SrPrf.cfg is required

If this structure is not maintained, then the files will not be found, because the directories are accessed directly during the update process.

- Plug in the external storage medium at the Ti-Touch.

2 Opening the update dialog

- Under **System ► Diagnosis**, tap on the button **Software update**.

3 Selecting the instrument

- Tap on the **Control device** list box and select the instrument to be updated.

4 Selecting file

- Tap on the **Binary file** list box.
The selection list with the BIN files saved in the "SwUpdates" directory is opened.
- Select the required file.
- Tap on **[Select]**.
- Update of the Ti-Touch: Continue with Step 5.
Update of the Sample Processor: Continue with Step 6.

5 Starting the update (915 KF Ti-Touch)

- Tap on the **[Start]** button (NOT on the [] fixed key).
- The message **023-102 Program update** is displayed.
- Confirm the message with **[Yes]**.
The update process is started, it runs automatically. The device is automatically switched off and back on during this process, possibly several times. No user intervention is required.
- **After successfully updating the software, remove the external storage medium (USB flash drive) from the instrument.**

6 Starting the update (control instruments)

- Tap on the **[Start]** button (NOT on the [] fixed key).
If **Result = No errors** is displayed, then the update has been successful.
- **Remove the external storage medium (USB flash drive) from the instrument.**
- Switch the 915 KF Ti-Touch off and then back on again.

32.8.7 Service

System ► Diagnosis ► Service

Access to the service functions is password-protected and only accessible to Metrohm service technicians. We recommend that you have maintenance performed on the system regularly by the Metrohm Service Department. You can automatically monitor the time interval after which the next service visit is due (*see chapter 13.5.1, page 142*).

33.2 Measuring inputs

The measuring cycle is 100 ms for all measuring modes.

33.2.1 Polarizer

One measuring input (**Pol.**) for polarizable electrodes.

Measuring mode
I_{pol} Determination with adjustable polarization current.

Polarization current –120 to +120 μA (increment: 1 μA)
–125 to –121 μA / +121 to +125 μA : non-guaranteed values,
dependent on reference voltage +2.5 V

Measuring range –1200 to +1200 mV

Resolution 0.1 mV

Measuring accuracy ± 0.2 mV
(± 1 digit, without sensor error, under reference conditions)

Measuring mode
U_{pol} Determination with adjustable polarization voltage.

Polarization voltage –1200 to +1200 mV (increment: 10 mV)
–1250 to –1210 mV / +1210 to +1250 mV: non-guaranteed values,
dependent on reference voltage +2.5 V

Measuring range –120 to +120 μA

Resolution 0.01 μA

Measuring accuracy –

33.2.2 Temperature

A measuring input (**Temp.**) for temperature sensors of the Pt1000 or NTC type with automatic temperature compensation.

R (25 °C) and B value can be configured for NTC sensors.

Measuring range

Pt1000 –150 to +250 °C

NTC –5 to +250 °C

(For an NTC sensor with R (25 °C) = 30,000 Ω and B (25/50) = 4,100 K)

Resolution

Pt1000 0.1 °C

NTC 0.1 °C

<i>Pt1000</i>	±0.2 °C (Applies for measuring range –20 to +150 °C)
<i>NTC</i>	±0.6 °C (Applies for measuring range +10 to +40 °C)

33.6 Energy supply

Device

Inlet

Nominal voltage 24 V DC

Power consumption (typical) 10 W

External power supply unit 6.2164.010

Inlet

Nominal voltage range 100–240 VAC ($\pm 10\%$)

Nominal frequency range 50–60 Hz ($\pm 3\%$)

Max. current 1.5 A

Outlet

Nominal voltage 24 V DC

Max. current 2.7 A

Power output 65 W

33.7 Declarations of Conformity

Up-to-date information on the *EU conformity* and the *requirements for electromagnetic compatibility* for your product can be found on the Internet under <https://www.metrohm.com/en-us/products-overview/>.

Using the article number, you can download the *Declaration of conformity for ...* under the documents of the product.

Base Steel, stainless, coated

33.11 Storage capacity

Random access memory 128 MB

Internal memory 1 GB of which 900 MB available.

Required memory 3–6 KB
per method

800 Dosino	15
801 Stirrer	16
803 Ti Stand	16
804 Ti Stand	16
805 Dosimat	15

Acoustic signals	51
ADD	321
Audit Trail	49
Archive	368
Export	368
Filter	367
Open	365
AuditTrailViewer	
Installing	363
Automatic system test	133
Automation	334
Axial distance	95

Backup	130
Interval	142
Balance	18, 116, 358
Barcode reader	
Connect	19
Beaker radius	
Sample vessel	102
Special beaker	106
Beaker sensor	
Sample vessel	103
Special beaker	107
BEEP	347

CALC	299
CALC LIVE	306
Calculation	299
Add	299
Delete	299
Edit	299
Load template	300
Note	303
Result options	304
Cell	
KFT	286
Command	
ADD	321
BEEP	347
CALC	299

CALC LIVE	306
CONTROL RS	333
CTRL	331
EMPTY	320
END	348
EVAL FIX-EP	295
EVAL MIN/MAX	296
EVAL RATE	298
FLOW	343
HEATER	343
Insert	169
KFT	277
LIFT	337
LQH	326
MEAS	289
MEAS T/F OFF	344
MEAS T/F ON	344
MOVE	334
PREP	319
PUMP	338
RACK	339
REPORT	314
REQUEST	346
SAMPLE	339
SCAN	329
SCAN RS	332
SIGN	347
STIR	345
SUBSEQ	340
WAIT	346
Command list	168
Common variable	
Delete	144
Edit	144
Properties	146
Validity	147
Communication	329
Conditioning	
KFT	284
Connect	
Balance	18
Barcode reader	19
Dosing device	15
MSB devices	14
PC keyboard	19
Power supply unit	11
Printer	17
Remote Box	16
Sample Processor	20

Stirrer or titration stand	16
USB devices	17
USB hub	17
Control	184
Control instrument	88
Parameters	
.....	286, 292, 322, 328
Control parameters	
KFT	279
Control Remote Box	83
CONTROL RS	333
CTRL	331
Curve options	226
Curves	225
Custom logo	163

Delete	
Signature	183
Statistics data	186
Determination	
Cancel	205
Carry out	203
Data	212
Define filters	223
Delete signatures	183
Endpoint details	215
List of saved determinations	220
Load	219
Measuring point list	215
Message	216
Properties	217
Recalculate	227
Reevaluate	227
Save	224
Save automatically	180
Send PC/LIMS report	180
Sign	182
Titrant data	215
Determination run	203
Device	
885 Compact Oven SC	108
Add	81
Balance	116
Configure	81
Delete	80
Edit	80
New	80

- Printer 112
- Sample Processor 92
- Ti-Touch titrator 88
- USB/RS-232 adapter 118
- Device manager 80
- 885 Compact Oven SC 108
- Diagnosis 368
- Dialog language 33
- Dialog options 35
- Dialog settings
- System-specific 35
- Directory structure 129
- Dosing 252
- Continuous 254
- Fixed volume 256
- Dosing curve 225
- Dosing device
- Connect 15
- Parameters 288, 322, 328
- Dosing parameters
- ADD 321
- Dosing port 58
- Dosing rate
- Maximum 356, 357
- Dosing unit
- Data 57
- GLP test 67
- Preparing parameters 58
- Replace reagent 260
- Tubing parameters 60
- E**
- E-mail 83
- Electrode 76
- Parameters 286, 292
- Electrostatic charge 7
- EMPTY 320
- Default parameters 356
- Emptying 260
- Default parameters 356
- END 348
- Endpoint details 215
- Ethernet connector 24
- EVAL FIX-EP 295
- EVAL MIN/MAX 296
- EVAL RATE 298
- Evaluation 295
- Fixed endpoint 295
- Maximum 296
- Minimum 296
- Rate 298
- Exchange unit
- Data 63
- GLP test 67
- Preparing parameters 65
- Tubing parameters 66
- Expiry date
- Sensor 78
- Titrant 56
- External position
- Swing angle 99
- F**
- Favorite 188
- Add 189
- Delete 189
- Edit 189
- FDA Guidance 21 CFR Part 11
- 38, 49
- File
- Copy 126
- Properties 127
- Rename 126
- File manager 124
- Filling 260
- Filling rate
- Maximum 356, 357
- Filter
- Determination 223
- Result table 237
- Fixed endpoint
- Evaluation 295
- Fixed key 9
- FLOW 343
- Format storage medium 369
- Formula editor 308
- G**
- GLP manager 132
- GLP test "Measurement" .. 134
- GLP test "Titration" 134
- System validation 138
- Test tool 133
- GLP test
- Dosing unit 67
- Exchange unit 67
- H**
- HEATER 343
- I**
- ID profile
- Create 43
- Identification profile
- Create 43
- Initialization 349
- Input line
- Template 158
- Installation
- Set up instrument 11
- Instrument
- Switch off 28
- Switch on 27
- Ti-Touch 82
- Intervention limit
- Titer 71
- K**
- Key functions 120
- Keyboard
- Connect 19
- KFT 1, 26, 277
- L**
- LCD test 369
- LIFT 337
- Lift rate 267
- Limit values
- Sample size 177
- Liquid Handling parameters 326
- Live display 208
- Live modifications 206
- Live parameters 210
- Live result
- Note 307
- Load
- Determination 219
- Sample table 198
- Load result template 300
- Login
- Options 44
- Logout automatically 45
- LQH 326
- M**
- Magnetic stirrer
- Connect 16
- Manual control 248
- 885 Compact Oven SC 273
- Continuous dosing 254
- Dosing 252
- Dosing fixed volume 256
- Emptying 260
- Filling 260
- Measuring 250
- Preparing 259
- Remote lines 263
- Sample Processor 264
- Stirring 261
- Troubleshooting 353
- Max. swing angle 97
- Maximum evaluation 296
- Maximum stroke path 94

O	
Operation	
General	28
Output line	
Template	160
Oven module	
Properties	111

Power supply unit	
Connect	11
PREP	319
Default parameters ...	356, 357
Parameters for the dosing unit	58

Rack	
Copy	100
Delete	100
Edit	100
Load	100
RACK	339
Rack adjustment	107
Rack name	93
Rack offset	103
Rate evaluation	298
Reagent	72
Add	72
Delete	72
Edit	72
Replace	260
Reagent monitoring	73
Reagent replacement	74
Reason	48
Modification option	48
Recalculation	227
Reevaluation	227
Remote	
Input	359
Interface	359
Output	360
Remote Box	
Connect	16

- Pin assignment 359
- Remote lines
- Manual control 263
- Report 240
- Print 240
- REPORT 314
- Report header 163
- REQUEST 346
- Resolution
- Measured value display 51
- Result 212
- Calculation formula 302
- Limits 303
- Save as a common variable 148
- Template 154
- Result table 234
- Define filters 237
- Delete automatically 238
- Save 239
- Result variable 303
- As parameter setting 362
- Robotic arm 97
- Max. swing angle 97
- Rotation offset 98
- Swing direction 98
- Swing increment 97
- Swing offset 97
- Swing radius 98
- Swing rate 269, 271, 336
- Rotation increment 95
- Rotation offset 98
- Routine dialog 37
- Standard configuration 37
- S**
- Safety 5
- Safety instructions 6
- SAMPLE 339
- Sample assignment table 150
- Sample data 192
- Modify live 206
- Settings 174
- Sample identification
- Modify designation 174
- Modify the designation 202
- Sample identification list 150
- Sample Processor 264
- Connect 20
- Sample rack 100
- Sample series
- Perform 204
- Sample size
- Fixed 174
- Limit value 177
- Limits 174
- Send from a balance 358
- Sample table 195
- Save as favorite 188
- Troubleshooting 351
- Save
- Determination 224
- Result table 239
- Sample table 197
- SCAN 329
- SCAN RS 332
- Sensor 76
- Add 76
- Delete 76
- Edit 76
- Expiry date 78
- Parameters 286, 292
- Working life 78
- Sensor list 76
- Serial number
- Dosing cylinder 58, 64
- Dosing unit 57
- Exchange unit 64
- Service 6, 375
- Set the date 34
- Set the time 34
- Set up instrument 11
- Sign
- Determination 182
- Method 182
- SIGN 347
- Signature
- Delete 183
- Software update 372
- Special beakers
- Edit 105
- Spooler 116
- Start conditions
- KFT 277
- Start options 171
- Statistics 186, 229
- Delete 186
- Troubleshooting 352
- STIR 345
- Stirrer
- Connect 16
- Parameters 288, 294, 325
- Stirring 261
- Stirring rate 357
- Stop conditions
- KFT 283
- Stop criterion
- KFT 281
- Stop options 173
- SUBSEQ 340
- Subsequence 340
- Load 341
- Note 342
- Options 342
- Save 341
- Supply voltage 7
- Swing direction 98
- Swing Head
- Properties 95
- Swing increment 97
- Swing offset 97
- Swing radius 98
- Swing rate 269, 271, 336
- Switch off 28
- Switch on 27
- System initialization 349
- System settings 33
- System test 133
- System validation 138
- T**
- Tandem dosing 323
- Template
- Input line 158
- Output line 160
- Report header 163
- Result template 154
- Sample assignment table .. 150
- Sample identification list ... 150
- Templates 150
- Text input 31
- Ti-Touch 82
- Time zone 34
- Titer
- History 70
- Intervention limit 71
- Properties 69
- Warning limit 71
- Titrant 52
- Add 53
- Delete 53
- Edit 53
- Expiry date 56
- Working life 56
- Titrant data 215
- Titrant list 52
- Titration curve 225
- Titration mode
- KFT 1
- Titration parameters
- KFT 282

Index

Titration stand	16
Connect	16
Troubleshooting	353
Touch screen	28
Test	371
Tower	
Properties	94, 110
Tubing parameters	
Dosing unit	60
Exchange unit	66
U	
Uninterrupted dosing	323
Update	372

USB	
Connect devices	17
USB hub	
Connect	17
USB/RS-232 adapter	118
User profile	
Create	43
V	
Validation interval	140
Validity	
Common variable	147
Titer	69

Valve disk	
Shift direction	63
W	
WAIT	346
Warning limit	
Titer	71
Water determination according to Karl Fischer	26
Working life	
Sensor	78
Titrant	56
Write protection	
External memory	128