



Application Note AN-PAN-1037

Online measurement of the acid number (AN) in oils with thermometric titration (ASTM D8045)

Metrohm has partnered with industry leaders to develop an alternative standard for the measurement of acid number (AN) in crude oil and petroleum products to overcome shortcomings in ASTM D664. This standard method (ASTM D8045) describes the use of thermometric catalytic titration for the analysis. Results agree closely with those from ASTM D664, but the thermometric catalytic titration method is far superior in terms of reproducibility and speed of analysis, with complete determinations in one minute. This alternative standard method also decreases the

amount of solvent used, reducing waste disposal costs. Comparison studies show very close data correlation between ASTM D8045 and traditional potentiometric AN titration methods, making it possible to seamlessly implement this method into a refinery that relies on historic data.

This Process Application Note presents a method to regularly monitor AN online in crude oil according to ASTM D8045. Online monitoring helps operators avoid corrosion issues in refinery processes.

However, instrumental methods that use glass membrane pH electrodes suffer from the difficulty of working in a water-free environment which results in

Metrohm thermometric titration uses a maintenance-free temperature sensor that does not require calibration nor rehydration and is resistant against fouling and matrix effects. The procedure requires minimal sample preparation and less solvent than traditional methods, which saves on waste disposal costs. Results agree closely with those from the potentiometric titrimetric procedure according to ASTM D664, but the thermometric catalytic titration method is far superior in terms of reproducibility and speed of analysis, with determinations completed within **one minute**. The thermometric endpoint titration (TET) method utilizes the same titrant (0.1 mol/L KOH in isopropanol) as ASTM D664. Online process analyzers from Metrohm Process Analytics can continuously sample and monitor AN in a refinery (**Figure 1**) with this fast, robust thermometric technique.

The 2060 TI Process Analyzer (Figure 2) can monitor acidity of crude oil according to ASTM D8045 testing procedures. By monitoring the acidity of crude oil and associated products, billions of dollars are saved annually by avoiding unexpected shutdowns and preserving expensive treatment chemicals.



Figure 2. 2060 TI Process Analyzer, suitable for ASTM D8045.

APPLICATION

ASTM D8045 describes the nonaqueous thermometric catalytic titration of weakly acidic species in crude oil. A defined sample amount is introduced via the sample loop and dissolved in 30–35 mL of a 3:1 xylene/isopropanol mixture. Paraformaldehyde is added as a catalytic indicator

before the titration is performed with 0.1 mol/L KOH in isopropanol. The endpoint, indicated by a temperature change with the robust Metrohm Thermoprobe, is identified by the second derivative (Figure 3). Benzoic acid is used as a standard in this method.

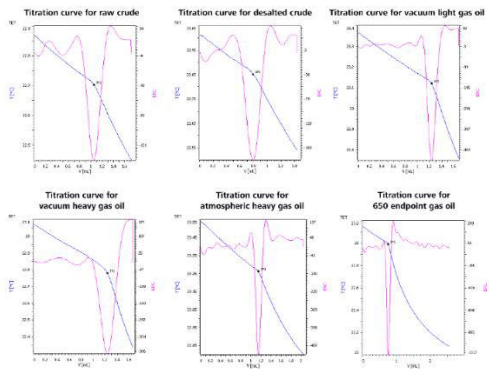


Figure 3. Acidity in crude oils and petroleum products by Metrohm thermometric catalytic titration, according to ASTM D8045.

Table 1. Typical range for AN in petroleum.

Parameter	Range [mg KOH/g]
AN	0.1–16

If the sample is not in liquid form, preconditioning at 65 °C is permitted to decrease sample viscosity. In thermometric titration, enthalpy change rather than potential of the reaction is monitored. Catalytically enhanced titrations using paraformaldehyde as catalyst are based on the endothermic hydrolysis of

paraformaldehyde in the presence of an excess of hydroxide ions. It is recommended to use the paraformaldehyde specified in the given ASTM method, as not every type is suited for the catalysis of this reaction.

CONCLUSION

The Metrohm Process Analytics 2060 TI Process Analyzer can reliably measure the acid number in crude oil and petroleum products according to ASTM

D8045. Additionally, it offers automated analysis results for different parts of a refinery process and helps to safeguard plant operations.

RELATED APPLICATION NOTES

[AN-PAN-1014 Automated online determination of salt in crude oil according to ASTM D3230](#)

[AN-PAN-1026 Analysis of hydrogen sulfide and](#)

[mercaptans derived from crude oil – Online determination according to ASTM D3227 and UOP163](#)

BENEFITS FOR ONLINE ANALYSIS IN PROCESS

- **More savings per measurement**, making results more cost-effective.
- **Increased product throughput**, reproducibility, production rates, and profitability.
- **Guarantee compliance** with government standards.
- **Protection of company assets** with built-in alarms at specified warning limits to prevent corrosion.



CONTACT

Metrohm ANZ
56 Buffalo Road
2111 Gladesville

info@metrohm.com.au

CONFIGURATION



2060 Process Analyzer

The 2060 Process Analyzer is an online wet chemistry analyzer that is suitable for countless applications. This process analyzer offers a new modularity concept consisting of a central platform, which is called a «basic cabinet».

The basic cabinet consists of two parts. The upper part contains a touch screen and an industrial PC. The lower part contains the flexible wet part where the hardware for the actual analysis is housed. If the basic wet part capacity is not sufficient enough to solve an analytical challenge, then the basic cabinet can be expanded to up to four additional wet part cabinets to ensure enough space to solve even the most challenging applications. The additional cabinets can be configured in such a way that each wet part cabinet can be combined with a reagent cabinet with integrated (non-contact) level detection to increase analyzer uptime.

The 2060 process analyzer offers different wet chem techniques: titration, Karl Fischer titration, photometry, direct measurement and standard additions methods.

To meet all project requirements (or to meet all your needs) sample preconditioning systems can be provided to guarantee a robust analytical solution. We can provide any sample preconditioning system, such as cooling or heating, pressure reduction and degassing, filtration, and many more.