



## Application Note AN-C-183

# Potassium in potassium bicarbonate effervescent tablets for oral solution

## Method validation according to the U.S. Pharmacopoeia

Potassium bicarbonate effervescent tablets are used to prevent low levels of potassium in the blood [1]. Pharmaceutical manufacturers and laboratories are obliged to use the monographs from the United States Pharmacopeia (USP) and the National Formulary (NF) to test drugs and formulations. The USP has embarked on a global initiative to modernize many of their existing monographs. As an alternative to atomic absorption spectroscopy (AAS), ion chromatography (IC) with non-suppressed conductivity detection has been approved by the USP

as a validated method to quantify the potassium content in potassium bicarbonate effervescent tablets for oral solution [2].

The Metrosep C 6 - 150/4.0 (L76) column provides the required separation of potassium. All acceptance criteria from the USP monograph «Potassium Bicarbonate Effervescent Tablets for Oral Solution» are fulfilled. The present IC method has been validated according to USP General Chapter <621> Chromatography [3].


## SAMPLE AND SAMPLE PREPARATION

Potassium bicarbonate effervescent tablets for oral solution (Effer-K 25 mEq, potassium 978 mg, unflavored) from two different lots were evaluated in this application study.

Sample stock solutions with nominally 4890 mg/L potassium were prepared from at least 20 finely powdered tablets. A small amount (47.2 mg) of the resulting powder was transferred to a 2000 mL volumetric flask. 200 mL of ultrapure water was

added, and the flask was swirled until effervescence ceased. Then the solution was diluted to volume with ultrapure water and mixed well.

Sample solutions with nominally 15.0 mg/L potassium were prepared by transferring 1.533 mL of the sample stock solution to a 500 mL volumetric flask, then diluted to volume with ultrapure water and mixed well.

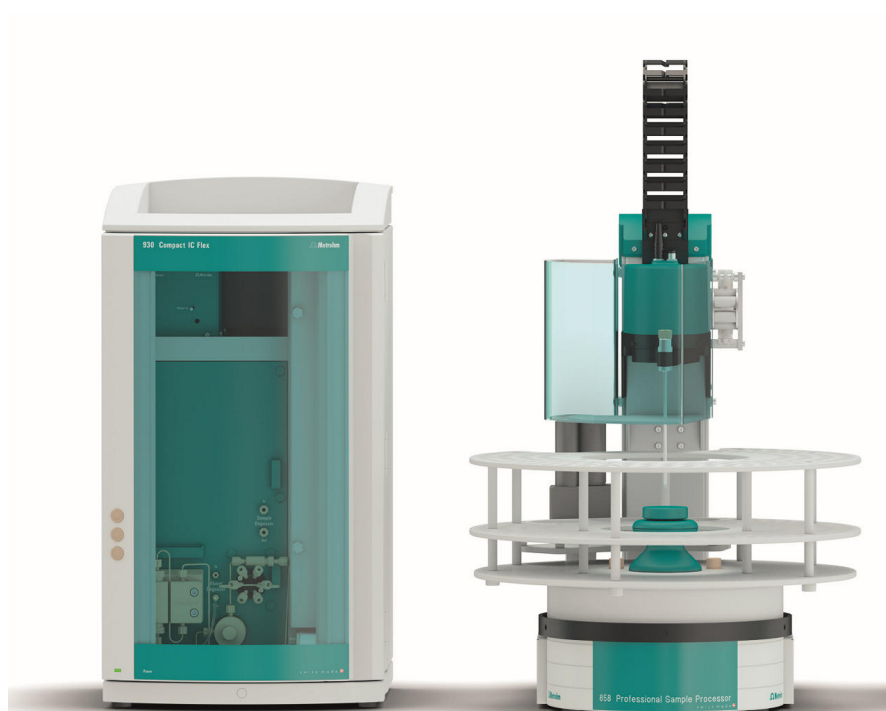

**Figure 1.** Instrumental setup including a 930 Compact IC Flex Oven/Deg and an 858 Professional Sample Processor.

## EXPERIMENTAL

The working standard solution of 15 µg/mL potassium was prepared from a certified 1000 µg/mL potassium USP reference standard.

Samples and standard solutions were injected directly into the IC using an 858 Professional Sample Processor (**Figure 1**). Potassium was separated from all

other cations using a Metrosep C 6 - 150/4.0 column (L76) and the signal was recorded with a conductivity detector.

The calibration was performed using a single standard at 15 µg/mL that was injected six times. The samples were analyzed in duplicate.

**Table 1.** IC method parameters as per USP monograph «Potassium Bicarbonate Effervescent Tablets for Oral Solution» [2].

| Column with L76 packing | Metrosep C 6 - 150/4.0 |
|-------------------------|------------------------|
| Eluent                  | 4 mmol/L nitric acid   |
| Flow rate               | 0.9 mL/min             |
| Temperature             | 30 °C                  |
| Injection volume        | 20 µL                  |
| Detection               | Direct conductivity    |

## RESULTS

The IC assay for potassium content was validated according to USP Monograph «Potassium Bicarbonate Effervescent Tablets for Oral Solution» [2]. The accuracy of the potassium determination was calculated as 100% (Figure 2).

All acceptance criteria were fulfilled, e.g., asymmetry (tailing factors) for the potassium peak was <2, or the relative standard deviation (%RSD) of standard solutions was <0.5% (n=6) (Table 2).

**Table 2.** Required acceptance criteria as per USP monograph «Potassium Bicarbonate Effervescent Tablets for Oral Solution» [2].

| Parameter                     | Actual | USP requirement | Status |
|-------------------------------|--------|-----------------|--------|
| % RSD standard solution (n=6) | 0.05   | NMT 0.5         | Pass   |
| Tailing factor                | 1.5    | NMT 2.0         | Pass   |
| Resolution                    | 3.89   | NLT 3.0         | Pass   |
| Solution stability            | 0.08%  | NMT 1.0%        | Pass   |
| Assay recovery                | 100.0% | 90–110%         | Pass   |
| % RSD assay (n=6)             | 0.15%  | NMT 1.0%        | Pass   |



**Figure 2.** Chromatogram of 15.0 µg/mL potassium in sample solution (100.0% recovery of the nominal concentration).

## CONCLUSION

The presented IC method for potassium in potassium bicarbonate effervescent tablets for oral solution with the Metrosep C 6 column (packing material L76) is officially included into the USP. Robustness and

reliability of the method was demonstrated following the guidelines of the USP General Chapter <621> [3]. The presented setup is suitable to quantify potassium according to USP requirements.

## REFERENCES

1. Kardalas, E.; Paschou, S.S.; Anagnostis, P.; et al. Hypokalemia: a clinical update. *Endocrine Connections* **2018**, 7 (4), R135–R146.  
<https://www.ncbi.nlm.nih.gov/pmc/articles/PMC5881435/>
2. U.S. Pharmacopeia. USP-NF Potassium Bicarbonate Effervescent Tablets for Oral Solution. *Monograph*.  
[https://doi.usp.org/USPNF/USPNF\\_M67194\\_02\\_01.html](https://doi.usp.org/USPNF/USPNF_M67194_02_01.html)
3. <621> Chromatography.  
[https://doi.org/10.31003/USPNF\\_M99380\\_01\\_01](https://doi.org/10.31003/USPNF_M99380_01_01)

Internal references: AW IC IN6-1888-062018

## CONTACT

Metrohm Vietnam Co., LTD  
Park IX Building 08 Phan  
Dinh Giot Street, Ward Tan  
Son Hoa  
700000 Ho Chi Minh City

info@metrohm.vn

## CONFIGURATION



### 930 Compact IC Flex Oven/Deg

The 930 Compact IC Flex Oven/Deg is the intelligent Compact IC instrument with **column oven**, **without suppression** and with built-in **degasser**. The instrument can be used with any separation and detection methods.


Typical areas of application:

- Anion and cation determinations without suppression with conductivity detection
- Simple applications with UV/VIS or amperometric detection



### 858 Professional Sample Processor – Pump

The 858 Professional Sample Processor – Pump processes samples from 500 µL to 500 mL. The sample transfer takes place either with the installed bidirectional two-channel peristaltic pump or with an 800 Dosino.

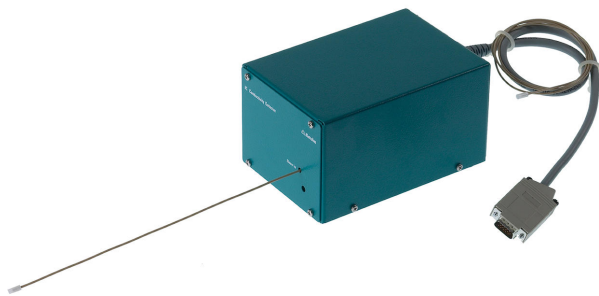

### IC Conductivity Detector

Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.