



Application Note AN-RA-008

# Easy detection of enzymes with the electrochemical-SERS effect

Activation of enhanced features of silver electrodes helps obtain characteristic Raman bands of biological molecules

Raman spectroscopy is one of the most promising chemical analysis techniques. This is due to its inherent fingerprint properties which allow the identification of different species present in a studied system.

Although low sensitivity has limited its use as a detection method, the surface-enhanced Raman scattering (SERS) effect has improved its effectivity for analytical use. Enhancement of the Raman signal has

allowed the development of many sensing applications. In particular, the energy provided by the 638 nm laser ensures a balance between risking damage to the sample and the generation of fluorescence, making this laser popular for most biological applications.

In this Application Note, aldehyde dehydrogenase and cytochrome c are analyzed by Raman spectroelectrochemistry as a proof of concept.

## INSTRUMENTATION AND SOFTWARE

Measurements were performed using a SPELEC RAMAN 638 instrument (638 nm laser), a Raman probe corresponding to the laser wavelength, and

spectroelectrochemical cells for screen-printed electrodes (**Figure 1a**) as well as for conventional electrodes (**Figure 1b**).

a)



b)



**Figure 1.** SPELEC RAMAN 638 instrument and Raman probe used in combination with a Raman spectroelectrochemical cell for (a) screen-printed and (b) conventional electrodes.

Silver screen-printed electrodes (Ag SPEs, C013) were used as SERS substrate because of their electrochemical activation properties. In addition, a conventional silver working electrode (6.09395.044) was also used as SERS substrate in combination with steel (6.0343.110) and Ag/AgCl (6.0728.120) counter and reference electrodes, respectively.

The SPELEC RAMAN 638 was controlled with DropView SPELEC, a dedicated software that provides spectroelectrochemical information and includes tools to perform adequate treatment and analysis of the collected data. All hardware and software used for this study is compiled in **Table 1**.

**Table 1.** Hardware and software equipment overview.

| Equipment                                                     | Article number  |
|---------------------------------------------------------------|-----------------|
| Instrument                                                    | SPELECRAMAN638  |
| Probe                                                         | RAMANPROBE638   |
| Raman spectroelectrochemical cell for SPEs                    | RAMANCELL       |
| Raman spectroelectrochemical cell for conventional electrodes | RAMANCELL-C     |
| Silver SPE                                                    | C013            |
| Connection cable for SPEs                                     | CAST            |
| Silver working electrode                                      | 6.09395.044     |
| Steel counter electrode                                       | 6.0343.110      |
| Ag/AgCl reference electrode                                   | 6.0728.120      |
| Connection cable for conventional electrodes                  | CABSTAT         |
| Software                                                      | DropView SPELEC |

## EC-SERS EFFECT OF SCREEN-PRINTED ELECTRODES: DETECTION OF ALDEHYDE DEHYDROGENASE

Raman spectroelectrochemistry was employed for the detection of aldehyde dehydrogenase (ALDH) in aqueous solution. Electrochemical SERS (EC-SERS) protocol requires two steps in a single experiment: the electrochemical activation of SERS features of Ag SPEs, and then the spectroscopic detection of the sample.


Electrochemical activation consists of scanning the potential from +0.50 V to produce the initial oxidation of the silver surface, followed by its reduction up to -0.60 V to generate nanostructures with the SERS effect. The spectroelectrochemical experiment is performed in presence of the analyte to be detected (in this case ALDH), but chloride ions are also present in solution to favor the generation of silver SERS

substrate [1].


In order to evaluate the evolution of Raman bands during the entire experiment, and therefore to optimize the optical detection, SPELEC RAMAN 638 performs measurements in operando mode. This way, spectra are recorded continuously during the whole experiment and not only at certain potentials.

**Figure 2a** displays the cyclic voltammogram of ALDH, and the characteristic Raman spectrum of ALDH is shown in **Figure 2b**. The results were obtained for 1 mg/mL ALDH in 0.1 mol/L KCl aqueous solution. Although the acquisition of spectra was continuously performed, the Raman spectrum in **Figure 2b** was recorded at -0.50 V as this potential provides the highest Raman intensity.

a)

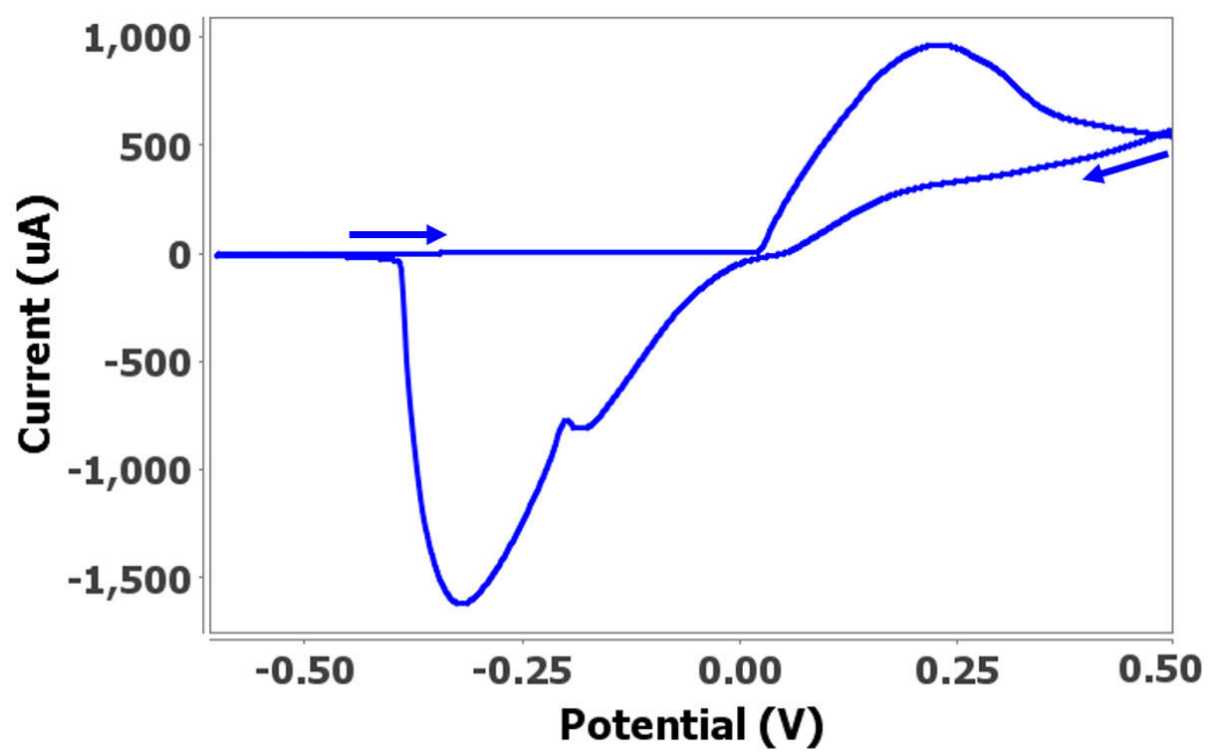

b)

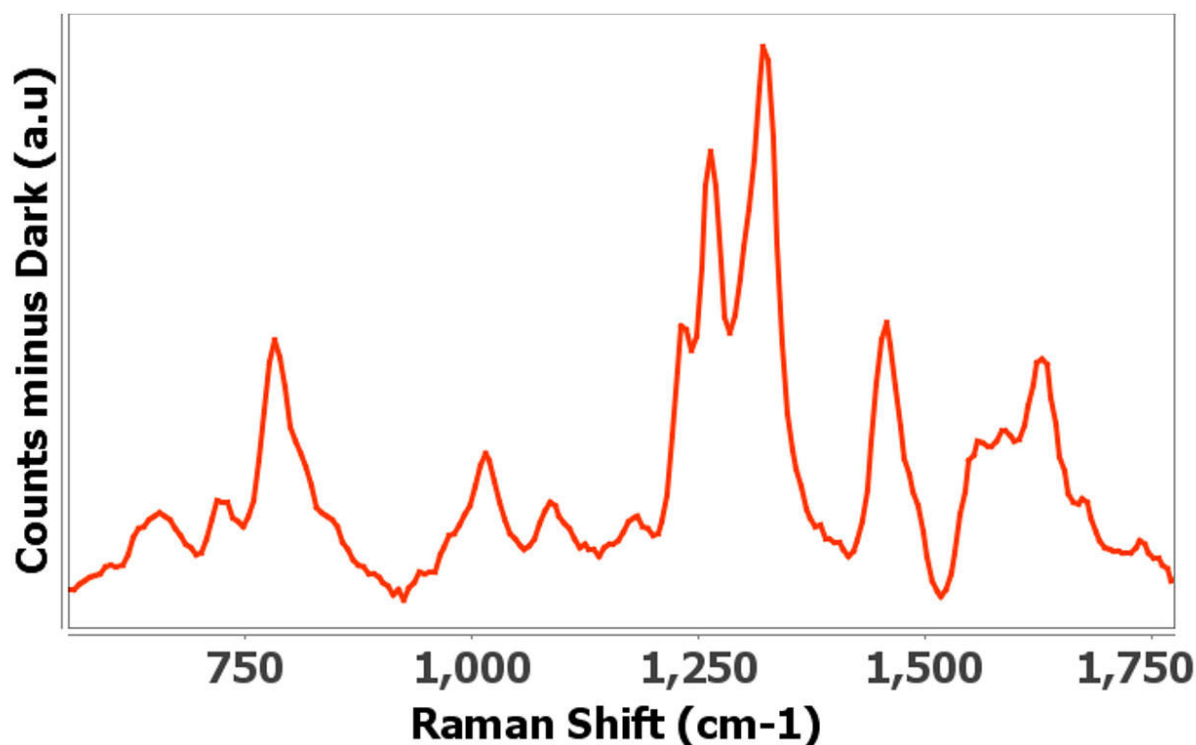

**Figure 2.** (a) Cyclic voltammogram and (b) Raman spectrum of 1 mg/mL aldehyde dehydrogenase in 0.1 mol/L potassium chloride aqueous solution using C013 electrodes. Integration time was 2000 ms.

The combination of the proposed EC-SERS procedure for the easy activation of Ag SPEs with a Raman spectroelectrochemical instrument (638 nm laser)


offers a fast and interesting alternative for the characterization of ALDH in solution not previously reported in the literature.

## EC-SERS EFFECT OF CONVENTIONAL ELECTRODES: DETECTION OF CYTOCHROME C

Conventional electrodes were also used for the EC-SERS detection of different enzymes such as cytochrome c. The protocol for this type of electrode follows the same steps as for SPEs. An initial oxidation of the surface is required followed by the subsequent reduction to generate silver nanostructures with SERS effect. The potential window was adjusted according to the electrodes used in this experiment. The best

results were obtained by scanning the potential of 0.1 mg/mL cytochrome c in 0.1 mol/L KCl aqueous solution from +0.80 V to -0.80 V. The Raman spectrum with the highest intensity (**Figure 3**) was recorded at -0.70 V.

Assignments of vibrational modes of each Raman band of cytochrome c are listed in **Table 2**.

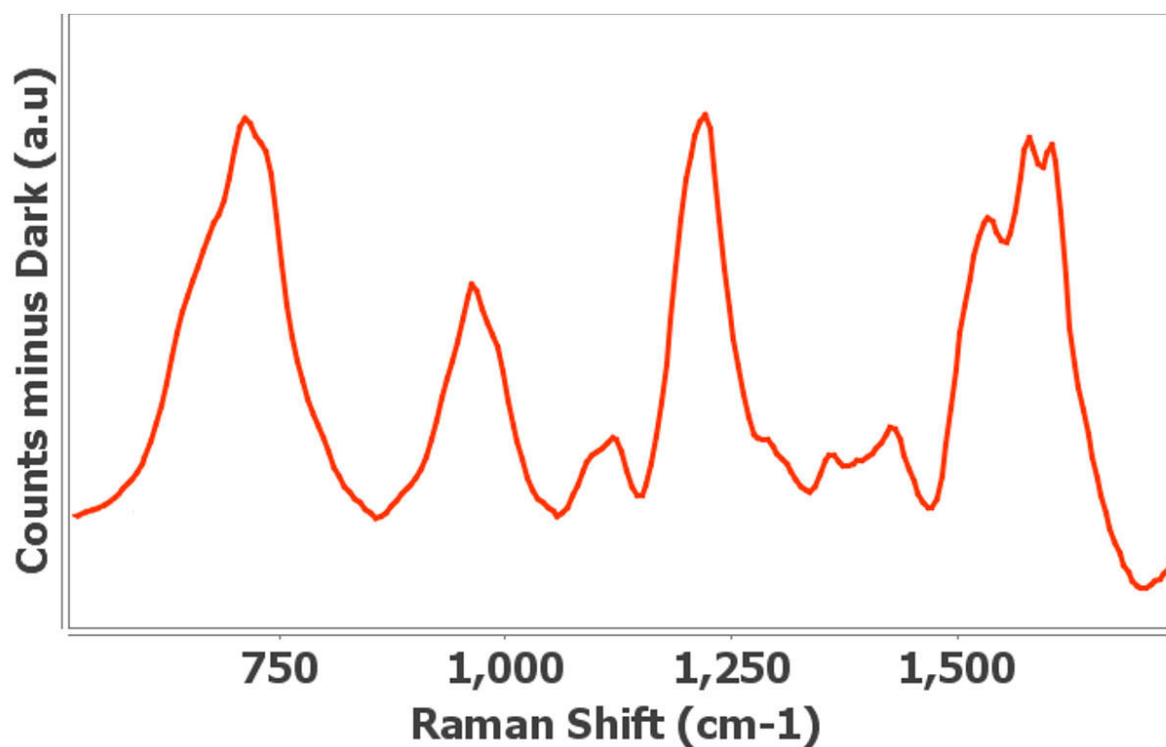

**Figure 3.** Raman spectrum obtained of 0.1 mg/mL cytochrome c in 0.1 mol/L potassium chloride aqueous solution using a conventional silver electrode. Integration time was 2000 ms.

**Table 2.** Vibrational assignment of cytochrome c [2,3].

| Cytochrome c SERS bands (cm <sup>-1</sup> ) | Assignment                                                      |
|---------------------------------------------|-----------------------------------------------------------------|
| 713                                         | Heme breathing                                                  |
| 969                                         | Asymmetric pyrrole deformation                                  |
| 1123                                        | C <sub>β</sub> -C <sub>1</sub>                                  |
| 1220                                        | Asymmetric pyrrole half-ring                                    |
| 1358                                        | Symmetric pyrrole half-ring                                     |
| 1426                                        | Pyrrole quarter-ring                                            |
| 1528                                        | C <sub>α</sub> C <sub>m'</sub> , C <sub>α</sub> N               |
| 1578                                        | C <sub>β</sub> C <sub>β'</sub> , C <sub>α</sub> C <sub>m'</sub> |
| 1604                                        | C <sub>α</sub> C <sub>m'</sub> , C <sub>α</sub> C <sub>β</sub>  |

Considering that cytochrome c exists as interconvertible reduced and oxidized forms, the oxidation state of the Fe ion can be determined according to the position of their characteristic Raman bands. Reduced cytochrome c shows one Raman band centered at  $1604\text{ cm}^{-1}$ , while the oxidized form displays an upshifted band centered at  $1636\text{ cm}^{-1}$ .

## CONCLUSION

Raman spectroelectrochemistry is a multi-response technique that provides outstanding results when studying a wide variety of systems. In the particular case of the 638 nm laser, the energy provided at this wavelength is suitable for biological applications.

In this work, the proposed electrochemical protocols for the activation of SERS properties of different silver electrodes (screen-printed as well as conventional)

According to the Raman spectrum (**Figure 3**) and the vibrational assignments (**Table 2**), the reduced form of the analyte is detected during the EC-SERS experiment. This experiment demonstrates the possibility to detect cytochrome c as well as characterize its oxidation state.

offers an easy and fast procedure to enhance the Raman intensity. This enhancement of Raman intensity allows the detection of different molecules present in solution. Fingerprint Raman bands of ALDH, not previously reported, are defined in this work. In the study of cytochrome c, the characteristic Raman bands are defined in addition to characterization of the redox state of the Fe ion.

## REFERENCES

1. Martín-Yerga, D.; Pérez-Junquera, A.; González-García, M. B.; et al. Quantitative Raman Spectroelectrochemistry Using Silver Screen-Printed Electrodes. *Electrochimica Acta* **2018**, *264*, 183–190.  
<https://doi.org/10.1016/j.electacta.2018.01.060>.
2. Brazhe, N. A.; Evlyukhin, A. B.; Goodilin, E. A.; et al. Probing Cytochrome c in Living Mitochondria with Surface-Enhanced Raman Spectroscopy. *Sci Rep* **2015**, *5* (1), 13793.  
<https://doi.org/10.1038/srep13793>.
3. Hu, S.; Morris, I. K.; Singh, J. P.; Complete Assignment of Cytochrome c Resonance Raman Spectra via Enzymic Reconstitution with Isotopically Labeled Hemes. *J. Am. Chem. Soc.* **1993**, *115* (26), 12446–12458.  
<https://doi.org/10.1021/ja00079a028>.

## RELATED APPLICATION NOTES

**AN-RA-006** New strategies for obtaining the SERS effect in organic solvents

**AN-SEC-001** Spectroelectrochemistry: an autovalidated analytical technique – Confirm results

via two different routes in a single experiment

**AN-SEC-002** Gathering information from spectroelectrochemical experiments – Calculation of electrochemical parameters from data

## CONTACT

Metrohm France  
13, avenue du Québec - CS  
90038  
91978 VILLEBON  
COURTABOEUF CEDEX

[info@metrohm.fr](mailto:info@metrohm.fr)

## CONFIGURATION



### Logiciel DropView SPELEC

DropView SPELEC est un logiciel de spectroélectrochimie qui pilote l'appareil SPELEC en offrant une synchronisation parfaite des mesures optiques et électrochimiques et des outils avancés pour le traitement de données.