



Application Note AN-V-179

# Iron in boiler feed water

Straightforward, sensitive, and cost-efficient determination using adsorptive stripping voltammetry (DHN method)

The iron concentration in boiler feed water must be monitored to ensure reliable and safe operation of the water-steam circuit in thermal power plants. Iron ions, even in low concentrations, are an indicator for corrosion and therefore can signal potential safety issues. Various guidelines set limits for the maximum iron content in boiler feed water.

The concentration of total iron in boiler feed water can be determined with high sensitivity using adsorptive stripping voltammetry (AdSV) using 2,3-

dihydroxynaphthalene (DHN) as complexing agent. Concentrations of total iron in water samples can be determined down to approximately 0.1 µg/L with this method.

The AdSV method is simple to perform, specific, and free of interferences. It is a viable, less sophisticated alternative to atomic absorption spectroscopy (AAS) or inductively coupled plasma (ICP) requiring only a moderate investment in hardware and low running costs.

# SAMPLE

Boiler feed water

# EXPERIMENTAL

The water sample, DHN solution, and the buffer solution are pipetted into the measuring vessel. The determination of total iron is carried out with an 884 Professional VA using the parameters specified in Table 1. The concentration is determined by two additions of an iron standard addition solution.



Figure 1. 884 Professional VA.

Table 1. Parameters for adsorptive stripping voltammetric (AdSV) analysis of Fe in boiler feed water

| Parameter            | Setting                 |
|----------------------|-------------------------|
| Working electrode    | HMDE                    |
| Mode                 | DP – Differential Pulse |
| Deposition potential | -0.1 V                  |
| Deposition time      | 30 s                    |
| Start potential      | -0.2 V                  |
| End potential        | -1.2 V                  |
| Peak potential Fe    | -0.7 V                  |

# ELECTRODES

- Working electrode: Multi-Mode Electrode pro with standard glass capillaries
- Reference electrode: Ag/AgCl/KCl (3 mol/L) reference electrode with electrolyte vessel.  
Bridge electrolyte: KCl (3 mol/L)
- Auxiliary electrode: Platinum rod electrode

RESULTS

The determination of total iron in boiler feed water samples can be carried out in a simple and straightforward manner with adsorptive stripping voltammetry (AdSV).

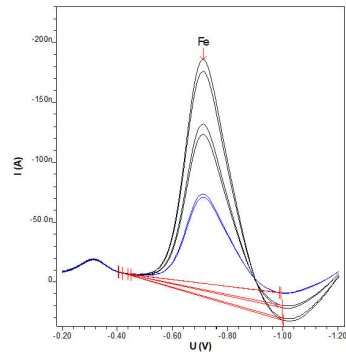


Figure 2. Determination of total Fe in boiler feed water.

Table 2. Results of total Fe determination with the 884 Professional VA

| Sample            | Total Fe [µg/L] |
|-------------------|-----------------|
| Boiler feed water | 3.0             |

Internal reference: AW VA CH4-0513-072012

CONTACT

CÔNG TY TNHH METROHM  
VIỆT NAM  
Tòa nhà Park IX, số 08  
Đường Phan Đình Giót,  
Phường Tân Sơn Hòa  
null Thành phố Hồ Chí Minh

info@metrohm.vn

## CONFIGURATION



### 884 Professional VA semiautomated für Multi-Mode Electrode (MME) with 2 Dosinos

884 Professional VA semiautomated for Multi-Mode Electrode (MME) is a convenient high-end routine analyzer for trace determinations with voltammetry and polarography with the Multi-Mode Electrode pro or the scTRACE Gold. The proven Metrohm electrode methods in combination with a high-performance potentiostat/galvanostat and the extremely flexible **viva** software open up new perspectives for the determination of heavy metals. The potentiostat with a certified calibrator readjusts itself automatically before each measurement, thus guaranteeing maximum precision.

Determinations with rotating disc electrodes can also be performed with the instrument, e.g. determinations of organic additives in electroplating baths with "Cyclic Voltammetric Stripping" (CVS), "Cyclic Pulse Voltammetric Stripping" (CPVS), and chronopotentiometry (CP). The replaceable measuring head enables rapid changes between the various applications with different electrodes.

2x 800 Dosinos (supplied) permit the automatic addition of auxiliary solutions during the determination, e.g., electrolyte, buffer or standard solutions.

The **viva** software is required for control, data collection, and evaluation.

The 884 Professional VA semiautomated for Multi-Mode Electrode (MME) is supplied with extensive accessories and a measuring head for the Multi-Mode Electrode pro. Electrode set and **viva** license need to be ordered separately.



### VA electrode equipment with Multi-Mode Electrode pro for Professional VA instruments

Complete electrode set for polarographic and voltammetric determinations. Includes Multi-Mode Electrode pro, reference electrode, platinum auxiliary electrode, measuring vessel, stirrer, electrolyte solution and additional accessories for setting up and operating the Multi-Mode Electrode.

### viva 2.1 Full CD: 1 license

Computer program for control and data management of Professional VA/ CVS instruments. The software permits checks, data acquisition, evaluation, and monitoring, as well as report generation.

Graphic user interface for routine operations, extensive database programs with reevaluation, graphic method editor with numerous templates, system configuration, very flexible user administration, extensive data export functions, individually configurable report generator.

Includes the measuring modes "Differential Pulse" (DP), "Square Wave" (SQW), "Cyclic Voltammetric Stripping" (CVS), "Cyclic Pulse Voltammetric Stripping" (CPVS), and Chronopotentiometry (CP). The following calibration techniques are supported: standard addition, external calibration, DT (Dilution Titration), LAT (Linear Approximation Technique), MLAT (Modified Linear Approximation Technique), RC (Response Curve).

Dialog languages: German, English, Chinese.

The following instruments are supported:

884 Professional VA, 894 Professional CVS

1 license for a maximum of 4 Professional VA/ CVS instruments on a single PC.