



Application Note AN-C-198

醋酸钙胶囊中的醋酸钙测定

Method validation according to the U.S. Pharmacopeia (USP)

Calcium acetate functions as a phosphate binder in the gastrointestinal tract, helping to lower high phosphate levels in individuals with kidney disease who are receiving dialysis treatment [1,2]. To meet the stringent quality standards for pharmaceutical products, manufacturers and laboratories must use validated methods from the United States Pharmacopeia – National Formulary (USP-NF). Previously, such methods included titration or liquid chromatography (LC) with UV detection. As part of their

modernization efforts, the USP has updated the calcium monograph to include ion chromatographic (IC) analysis, which is more straightforward and sensitive than previous methods. For the calcium acetate assay, the USP specifies ion chromatography using a cation-exchange column with L76 column material and non-suppressed conductivity detection to quantify the amount of calcium ions in calcium acetate capsules [3].

The current IC method employs a Metrosep C 6 -

150/4.0 column (L76) to separate calcium from other ions in calcium acetate capsules. This method has been validated according to USP General Chapters <621> Chromatography [4] and <1225> Validation of Compendial

SAMPLE AND SAMPLE PREPARATION

The standard solution nominally contains 0.08 mg/mL of USP Calcium Acetate Reference Standard (Cat# 1086334) in water. It is prepared by accurately weighing 80.0 mg of USP Calcium Acetate Reference Standard and transferring it into a clean 1000 mL volumetric flask. It is dissolved and made up to the mark with ultrapure water (UPW).

For sample stock solutions nominally containing 6.7 mg/mL of calcium acetate, an appropriate portion of the contents of at least 20 capsules are transferred into a 2000 mL volumetric flask.

EXPERIMENTAL

Samples are directly injected into the IC with a 919 IC Autosampler plus (Figure 1).

Cations are separated using a Metrosep C 6 - 150/4.0 column (L76) and detected with non-suppressed conductivity (Table 1). The run time was 40 minutes which complies with the USP requirements of 1.5 times the retention time of the calcium peak (here: 24 minutes).

A one-point calibration with the standard solution of 0.08 mg/mL of USP Calcium Acetate was used for quantification. Samples were evaluated as triplicates. Repeatability studies are done with 6-fold injections.

Procedures [5]. All acceptance criteria for the calcium acetate assay in the USP monograph «Calcium Acetate Capsules» are met.

UPW is added to about 40% of the final volume of the volumetric flask, and the solution is then sonicated for 30 minutes with intermittent shaking. Afterwards, the solution is made up to the mark with UPW and filtered through 0.2 µm filter paper.

For sample solutions nominally containing 0.08 mg/mL of calcium acetate, 5.97 mL of sample stock solution are transferred into a clean 500 mL volumetric flask. This is diluted and made up to the mark with UPW. All solutions are sonicated for 5 minutes before injection.



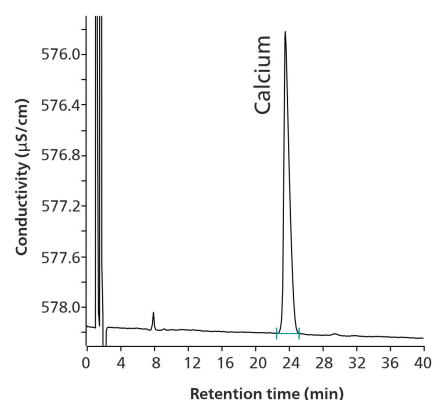
Figure 1. Instrumental setup including a 930 Compact IC Flex Oven/Deg and a 919 IC Autosampler plus.

Table 1. Parameters for the IC method as per USP Monograph «Calcium Acetate Capsules» [3].

Column with L76 packing	Metrosep C 6 - 150/4.0
Eluent	0.75 mmol/L dipicolinic acid + 1.7 mmol/L nitric acid
Flow rate	0.9 mL/min
Temperature	35 ° C
Injection volume	10 µL
Detection	Direct conductivity

RESULTS

IC method parameters followed the requirements of the USP assay for calcium acetate (**Table 2**). Chromatograms did not show any interferences or contaminations, and the calcium peak eluted after 24 minutes (**Figure 2**). All results for validation testing were within the specified USP requirements. Column efficiency was exceptional with >5900 theoretical plates. Relative standard deviation for a 6-fold standard injection was 0.4% (USP requirement <2.0%). The measured amount of calcium acetate in the capsules correlated well with the labeled amount of calcium acetate, e.g., the calculated percentage was 102.6% of the labeled content (90.0110.0% is acceptable) (**Table 2**). Thus, the IC method was suitable to determine calcium in calcium acetate capsules.

**Figure 2.** Chromatogram of a calcium peak from a calcium acetate capsule sample, containing 0.082 mg/mL calcium acetate (102% recovery).**Table 2.** Exemplary results and USP requirements from the IC method validation for calcium acetate in calcium acetate capsules as per USP [3].

Parameter	Result	USP requirement
Theoretical plates	5909	NLT 1000
RSD % (n = 6)	0.406%	NMT 2.0%
Percentage of labeled amount	102.6%	90.0110.0%

CONCLUSION

According to the USP monograph for calcium acetate capsules [3], the assay for calcium acetate involves determining calcium content using ion chromatography (IC) on a separation column with L76 packing material (here: Metrosep C 6). The validation results met all

requirements of the monograph and adhered to the guidelines specified in USP General Chapters <621> Chromatography and <1225> Validation of Compendial Procedures [4,5]. The described IC method is appropriate for quantifying calcium in calcium acetate capsules.

REFERENCES

1. Biruete, A.; Gallant, K. M. H.; Lindemann, S. R.; et al. Phosphate Binders and Nonphosphate Effects in the Gastrointestinal Tract. *Journal of Renal Nutrition* **2020**, 30 (1), 4–10. DOI:10.1053/j.jrn.2019.01.004
2. Sekar, A.; Kaur, T.; Nally, J. V.; et al. Phosphorus Binders: The New and the Old, and How to Choose. *CCJM* **2018**, 85 (8), 629–638. DOI:10.3949/ccjm.85a.17054
3. U.S. Pharmacopeia. USP-NF Calcium Acetate Capsules. *Monograph*. DOI:10.31003/USPNF_M11403_08_01
4. <621> Chromatography, General Chapter; U.S. Pharmacopeia/National Formulary: Rockville, MD.
5. <1225> Validation of Compendial Procedures; General Chapter; U.S. Pharmacopeia/National Formulary: Rockville, MD. DOI:10.31003/USPNF_M99945_04_01

CONTACT

瑞士万通中国
北京市海淀区上地东路1号
院1号楼7层702
100085 北京

marketing@metrohm.com.cn

CONFIGURATION



Metrosep C 6 - 150/4.0

高容量的 C 6 材料使分离柱 Metrosep C 6 - 150/4.0 成为在浓度差异很高时的正常保留时间内分离标准阳离子的解决方法。可用此柱测定铵含量很高的饮用水。



930 Compact IC Flex Oven/Deg

930 Compact IC Flex Oven/Deg 是智能型带柱加热炉、无抑制的 Compact IC 仪器,并且内置脱气装置。该仪器可使用各种分离和检测方法。


典型的应用范围:

- 阴离子和阳离子测定,无抑制的电导检测
- 使用 UV-VIS 或安培检测的简单应用



919 IC Autosampler plus

919 IC Autosampler plus 满足实验室内处理中等样品量的要求。使用该仪器,可实现万通产品的各种离子色谱分析自动化。

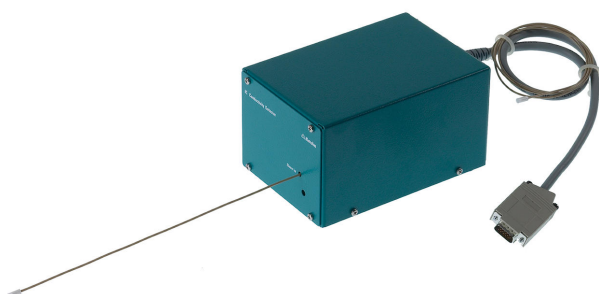

IC Conductivity Detector

用于智能型离子色谱仪的紧凑智能型高性能电导检测器。不凡的温度稳定性,受保护的检测器端子板内的整个信号处理过程以及新一代的 DSP(数字式信号处理)均能确保测量的准确性。归功于动态工作范围,无需进行范围变更(也不是自动进行)。