



Application Note AN-PAN-1063

用拉曼光谱对硼酸盐和硫酸盐溶液进行在线分析

Boron is a semimetal found in the form of borax (sodium tetraborate) and other oxides in nature [1]. Boric acid (H_3BO_3) is derived from borax and is used in several industrial applications such as glass manufacturing, electronics, detergents, food preservatives, and more [2].

Boric acid can be produced from borax through various processes. Sulfuric acid is mainly used during borax synthesis since it is considered to have minimal environmental impact.

In 2021, the boric acid market was estimated at US\$ 706.52 million and is expected to reach US\$


1,169.89 million by 2030 [3]. As the market grows, so does the need for a more cost-efficient and environmentally friendly production process.

This Process Application Note shows the excellent performance of the PTRam, the single-channel Raman instrument for process development from Metrohm Process Analytics, while measuring boric acid and sodium sulfate solutions inline at low concentrations (<100 mg/L).

INTRODUCTION

There are several ways to manufacture boric acid from borax. Some of these processes use strong acids like nitric acid or chloric acid, however, production costs using these chemicals are high

due to equipment wear. Out of these acids, sulfuric acid is primarily used (**Reaction 1**) since it is considered to have the smallest environmental footprint.

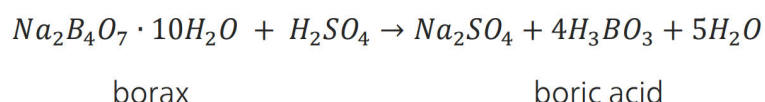

Reaction 1. Boric acid can be produced from the reaction between borax and sulfuric acid.

Maximizing production efficiency and reducing costs in a boric acid refinery is possible by monitoring and controlling the sodium sulfate (Na_2SO_4) chemistry in the crystallization stage

(**Figure 1**). If the reagent concentrations are outside of the set limits, chemical dosage is not controlled, waste increases, and production costs are higher.

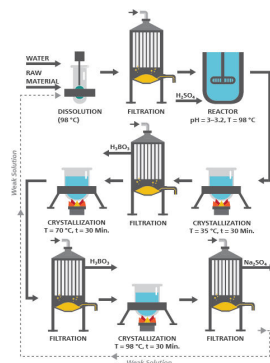


Figure 1. Schematic illustration of boric acid production from borax (adapted from [4]).

Gravimetric analysis is traditionally used for quantifying desired constituents of mixed chemicals or solutions by weight after separation [5]. This conventional method can monitor H_3BO_3 [6] and Na_2SO_4 concentrations

[7]. However, practical challenges arise through laborious sample preparation methods and manual data analysis. Additionally, gravimetric analysis does not provide users with real-time process information.

For optimal boric acid production, multiple parameters must be monitored in a safer, more efficient, and faster manner. This is possible via inline process analysis with reagent-free spectroscopy (e.g., Raman).

Metrohm Process Analytics offers the PTRam Process Analyzer (Figure 2) which enables direct

comparison of «real-time» spectral data from the process to a reference method (e.g., titration). This allows operators to create a simple, yet indispensable calibration model used to produce quantitative results during the boric acid manufacturing process.



Figure 2. The PTRam Analyzer is suitable for quantitative inline process analysis.

APPLICATION

Laser used: 785 nm. Samples of H_3BO_3 and Na_2SO_4 single salt and mixed salt solutions were taken for this study. Only a minimal number of reference measurements were required for calibration and modeling.

One of the many functions of Raman spectroscopy is its use for material identification. Most materials can be identified by their Raman signature, as they exhibit sharp, distinctive peaks serving as a molecular fingerprint. The spectrum

contains information not only regarding sample composition, but also about the concentrations of its constituents, which are directly proportional to the intensity of the spectrum. Because of the spectral differences, Raman analyzers are known to be capable of identification and confirmation of chemical substances used in a variety of industries (e.g., semiconductors, food, pharmaceuticals, etc.).

RESULTS

In this application study, the boric acid salt and sodium sulfate salt exhibit intense Raman bands at 880 cm^{-1} and 993 cm^{-1} , respectively (**Figure 3**). The limits of detection (LODs) for H_3BO_3 and Na_2SO_4 solutions are 15 mg/L (15 mg/L BO_3^{3-})

and 10 mg/L (7 mg/L SO_4^{2-}) (**Table 1**). This clearly demonstrates the ability of inline Raman spectroscopy for accurate quantitative analysis at low analyte concentrations.

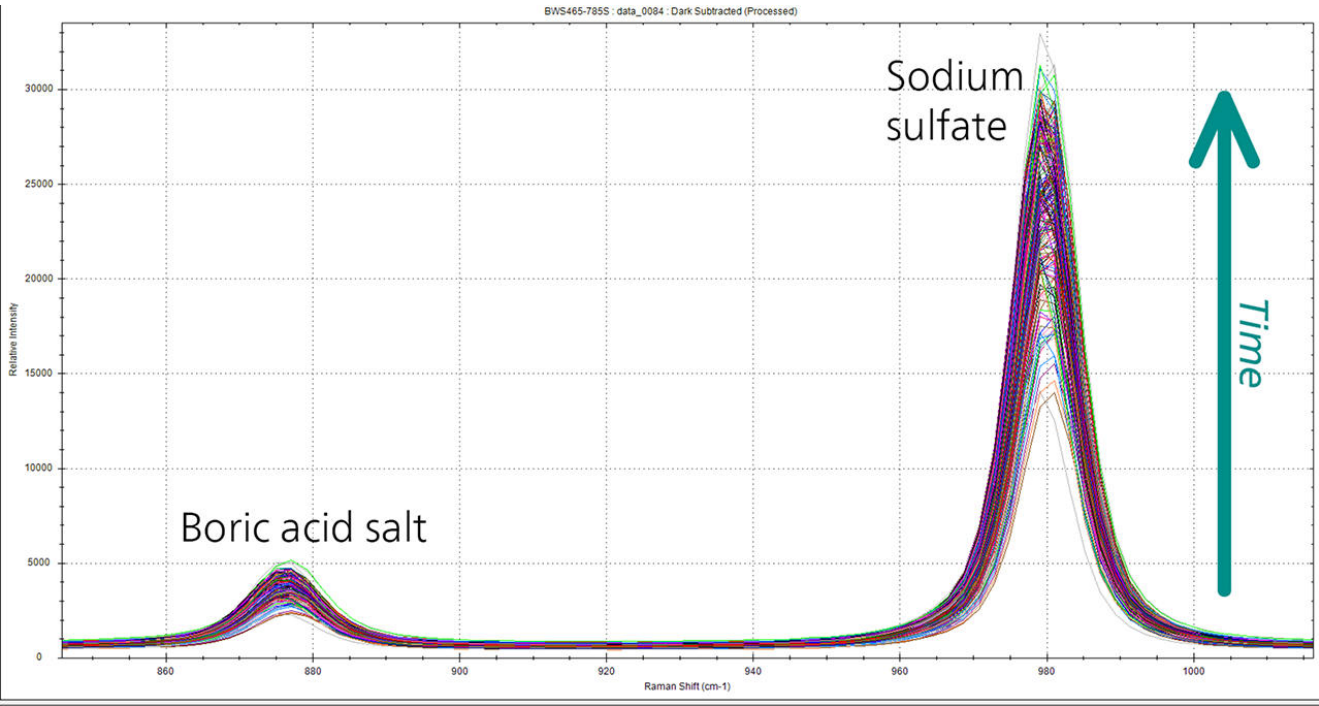


Figure 3. Raman spectra from the reaction stage. Boric acid salt (left) and sodium sulfate (right) exhibit clearly defined peaks.

Table 1. LOD (limit of detection in mg/L) of H_3BO_3 and Na_2SO_4 solutions with inline Raman spectroscopy.

† Standard Error of Prediction without bias correction; ‡ H_3BO_3 and Na_2SO_4 , mg/L; * BO_3^{3-} and SO_4^{2-} , mg/L. st: salt weight. io: ion weight.

Parameters	Factors	Concentration range (mg/L)	SEP [†]	LOD _{st} [‡]	LOD _{io} [*]
H_3BO_3 solution	2	0–80	4.6	15.2	14.5
Na_2SO_4 solution	2	0–80	3.1	10.2	6.9
H_3BO_3 in 1 g/L Na_2SO_4	2	0–80	10.1	33.3	31.6
Na_2SO_4 in 5 g/L H_3BO_3	2	0–80	3.5	1.6	7.8

CONCLUSION

In conclusion, the use of Raman spectroscopy for inline analysis of borate and sulfate solutions presents significant advantages in the context of boric acid production. As the boric acid market expands, the demand for efficient and eco-friendly production intensifies. Inline analysis, facilitated by reagent-free spectroscopy techniques such as Raman, enables real-time

monitoring of process parameters, as demonstrated by the PTRam Process Analyzer. By identifying distinctive molecular signatures and offering accurate quantitative results, Raman spectroscopy presents a robust approach for optimizing boric acid production, addressing challenges associated with traditional methods.

BENEFITS FOR RAMAN IN PROCESS

- Increased product throughput, reproducibility, production rates, and profitability.
- Gain insight into chemical reactions occurring in the manufacturing process.



REFERENCES

1. *Boron | Properties, Uses, & Facts | Britannica.*
<https://www.britannica.com/science/boron-chemical-element> (accessed 2023-08-21).
2. *Boron Mining: Sources And Major Producers | Borates Today.*
<https://borates.today/boron-mining-sources-and-major-producers/> (accessed 2023-08-21).
3. *Boric Acid Market Size to Worth Around US\$ 1,169.89 Million by 2030.*
<https://www.precedenceresearch.com/boric-acid-market> (accessed 2022-07-22).
4. Mergen, A.; Demirhan, M. H.; Bilen, M. Processing of Boric Acid from Borax by a Wet Chemical Method. *Adv. Powder Technol.* **2003**, 14 (3), 279–293.
<https://doi.org/10.1163/15685520360685947>.
5. *Gravimetric analysis | Definition, Steps, Types, & Facts | Britannica.*
<https://www.britannica.com/science/gravimetric-analysis> (accessed 2023-08-21).
6. Childs, M. P. Quantification of Boric Acid Concentration and Losses Due to Vaporization in the PASTA Facility., Texas A&M University, 2016.
7. Sodium Sulphate for Industrial Use — Determination of Sulphates Content — Calculation Method and Barium Sulphate Gravimetric Method, July 1975.

CONTACT

瑞士万通中国
北京市海淀区上地东路1号
院1号楼7层702
100085 北京

marketing@metrohm.com.cn

CONFIGURATION



PTRam Analyzer

PTRam Analyzer为一款 785 nm 进程开发拉曼分析仪,其设计旨在用于分析实验室和试点工厂环境中的产品和进程开发。这是一个高性能、准确、结实且可靠的拉曼系统,其具有自校正和自动性能校验功能,来确保每次测量的有效性。

这一单样本通道系统包括了一个实验室光纤探头,以及用户可更换的框架。该 PTRam 为 19“ 可固定样品盘。其 PTRam 可通过 Vision 软件进行操作,并可以与 2060 Human Interface 连接。