



Application Note AN-RS-026

婴幼儿食品中二苯基胺的痕量检测

Protecting consumer safety with Misa

Diphenylamine (DPA) is used as a dye fixative and antioxidant in industrial applications and as a produce preservative in agricultural operations. Despite its low toxicity to humans and the environment, there is controversy about its application on fruits and vegetables post-harvest to ensure long-term storage. Food safety advocates are concerned that daily ingestion of DPA, particularly in foods meant for babies, could have negative effects on children's health. To mitigate potentially toxic effects of DPA, both

the US and EU stipulate a maximum residue limit (MRL) of 5 µg/g for whole pears and a stringent MRL of 10 ng/g for all processed baby foods.

Misa (Metrohm Instant SERS Analyzer) provides a user-friendly and cost-effective alternative to traditional analytical methods used for detecting DPA in foods, such as GC-MS and GC-NPD. With Misa, trace level detection of DPA in infant food is achieved quickly and easily in a format that is easily adapted to diverse sample matrices.

INTRODUCTION

This application note describes a simulated test procedure for detecting DPA in pureed pears marketed as food for infants. The assay is based

on the acquisition of SERS-specific spectra for diphenylamine in chloroform extracts using Misa and gold nanoparticles (Au NPs).

REFERENCE MATERIAL AND LIBRARY CREATION

To establish a reference spectrum, a pure DPA standard (100 µg/g in ethanol) was analyzed using Au NPs. The unique SERS spectrum shown

in **Figure 1** can be used to create a library entry for DPA.

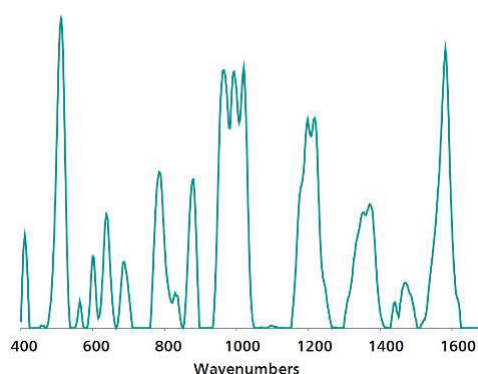


Figure 1. Standard gold nanoparticle SERS reference spectrum of diphenylamine.

EXPERIMENT

Serial dilutions of a stock solution of DPA in ethanol were spiked into purchased pureed pears sold as baby food to provide a range of test samples: 100, 50, 25, 10, 5, 2.5, and 1 µg/g. Samples were mixed with a vortex, then 0.5 mL chloroform was added, and samples were mixed again for 1 minute. After allowing phase separation for 10 minutes, 100 µL of the bottom chloroform layer was siphoned with a pipette and transferred to a glass vial for evaporative heating. The extract residue was suspended in 450 µL of Au NPs and 50 µL of 0.5 mol/L NaCl. Each vial was vigorously shaken and immediately placed into the vial attachment on Misa for measurement.




Table 1. Experimental parameters

Instrument		Acquisition	
Firmware	0.9.33	Laser Power	5
Software	Misa Cal V1.0.15	Int. Time	5 s
Misa Vial Attachment	6.07505.040	Averages	10
ID Kit - Au NP	6.07506.440	Raster	ON

RESULTS

Overlaid of baseline-corrected spectra acquired for each tested concentration shows high-

resolution DPA detection down to 1 µg/g (Figure 2).

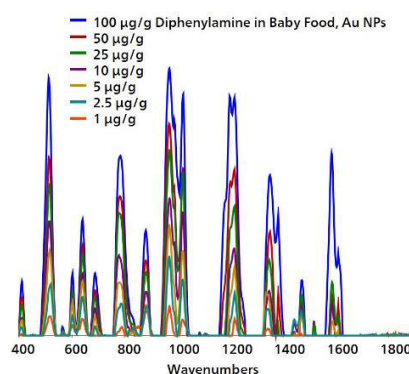**Figure 2.** Gold NP SERS spectral profile for diphenylamine in pureed pears, demonstrating detection down to 1 µg/g.

FIELD TEST PROTOCOL

Detection of diphenylamine in the field

Add baby food to glass vial until approximately 1/4 full. Using a pipette, add chloroform to the same vial until halfway full. Shake mixture vigorously for 1 minute, and then rest for 10 minutes to permit separation. Using a pipette, carefully remove half of the *bottom layer*, avoiding the top layer. Transfer this solution into

a *clean glass vial* and heat until dry on a hotplate (5 minutes). Using clean pipettes, add Au NPs to vial until halfway full, then add 2 drops of NaCl solution. Cap vial and shake to mix, then allow to rest for 5 minutes before inserting into the vial attachment on Misa for testing.

Table 2. Requirements for field test protocol

ID Kit - Au NP	6.07506.440
includes:	Gold nanoparticles (Au NP)
	Scoop
	Disposable pipettes
	2 mL glass vials
Reagents	
Chloroform	
NaCl solution	3 g NaCl in 100 mL water
Test settings	Use ID Kit OP on MISA

CONCLUSION

With successful, sensitive detection of diphenylamine in baby food, Misa is demonstrated to be a reliable QC tool for monitoring food safety, especially safeguarding the health of children. Misa is a portable and fully contained system, ideal for on-site testing and uniquely tailored for operation in a wide

range of different environments.

Additionally, operators in postharvest produce treatment facilities who are at risk of high-level occupational exposure, for instance, would benefit greatly from a mobile monitoring system for such chemical hazards.

CONTACT

瑞士万通中国
北京市海淀区上地东路1号
院1号楼7层702
100085 北京

marketing@metrohm.com.cn

CONFIGURATION



MISA Advanced

Metrohm Instant SERS Analyzer (MISA) 是一款高性能、便携式分析系统,可快速检测/鉴定非法物质、食品添加剂和微量食品污染物。MISA 的特点是配备了 Metrohm 的轨道光栅扫描 (ORS) 技术设备的光谱仪。其空间需求小和并且电池寿命有所延长,适用于现场测试或移动实验室应用。MISA 提供各种 1 级激光附件,可灵活选择取样。该分析仪可通过 Bluetooth 或 USB 连接运行。

MISA Advanced 套件是一个完整套件,其作用是让用户能用 Metrohms 纳米颗粒溶液和 P-SERS 试剂条进行分析。

MISA Advanced 套件包含了一个 MISA 小管附件、一个 P-SERS-附件、一个 ASTM 校正标准件、一个 USB 迷你电缆、一个 USB 供电单元和用于运行 MISA 仪器的 MISA Cal 软件。还随供了一个用来安全保管仪器和附件的坚固保护箱。



ID – Au NP

ID 套件 - Au NP 包含了 Mira/Misa 用户使用胶体金溶液进行 SERS 分析所需的组件。该套件包含了一个一次性抹刀、一个移液管、样品小瓶和一个含金胶体的瓶子。