



Application Note AN-RS-032

拉曼和 SERS 法鉴定阿片类药物复方制剂

Two chemical analyses from one tablet

Handheld Raman spectrometers are valued for their ability to provide onsite material identification in seconds. A simple point-and-shoot analysis of bulk materials with the handheld MIRA DS results in chemical identification accompanied by both statistical support and relevant color-coded warnings. In the case of combination pharmaceuticals, a single tablet contains more than one active ingredient in different proportions. MIRA DS is uniquely capable of identifying multiple

compounds in such tablets by using Raman to identify the major component and SERS (surface-enhanced Raman spectroscopy) for the minor component.

This application describes quick, dual analysis of a prescription medication containing acetaminophen and hydrocodone. The application is easily extrapolated to the study of street drugs. With MIRA DS, forensic analysis of tablets reaches a whole new level of accuracy.

INTRODUCTION

Norco® is a combination pharmaceutical containing 5 mg of hydrocodone, an opioid pain medication, and 325 mg of acetaminophen, a less potent analgesic that enhances the effects of hydrocodone. As with all opioids, Norco® has the potential for deadly side effects and may be habit-forming. Such medications require restricted access for children and care to avoid potential recreational use. Raman spectroscopy

may be used to identify tablets for the purpose of controlling recreational use of prescription drugs. A Raman scan of Norco® would simply report the bulk material: acetaminophen. Identification of hydrocodone component requires trace analysis with SERS. MIRA DS offers both Raman and SERS analyses with its signature compact size, ease-of-use, speed, and sampling flexibility.

SAMPLING WITH RAMAN

The beauty of sampling with MIRA DS is its simplicity. A Long Working Distance (LWD) Smart Tip (8 mm focal length) was used to penetrate beyond the tablet coating. Once the LWD was attached to MIRA DS, it automatically activated Smart Acquire routines that prompted the instrument to optimize acquisition parameters, take a sample, and match to

internal libraries—all in **less than a minute**.

Figure 1 contains Library Match information as it is seen in Mira Cal DS. The red Norco® sample spectrum matches very nicely to the library spectrum of acetaminophen, with a Hit Quality Index (HQL, a measure of spectral correlation) of 0.96.

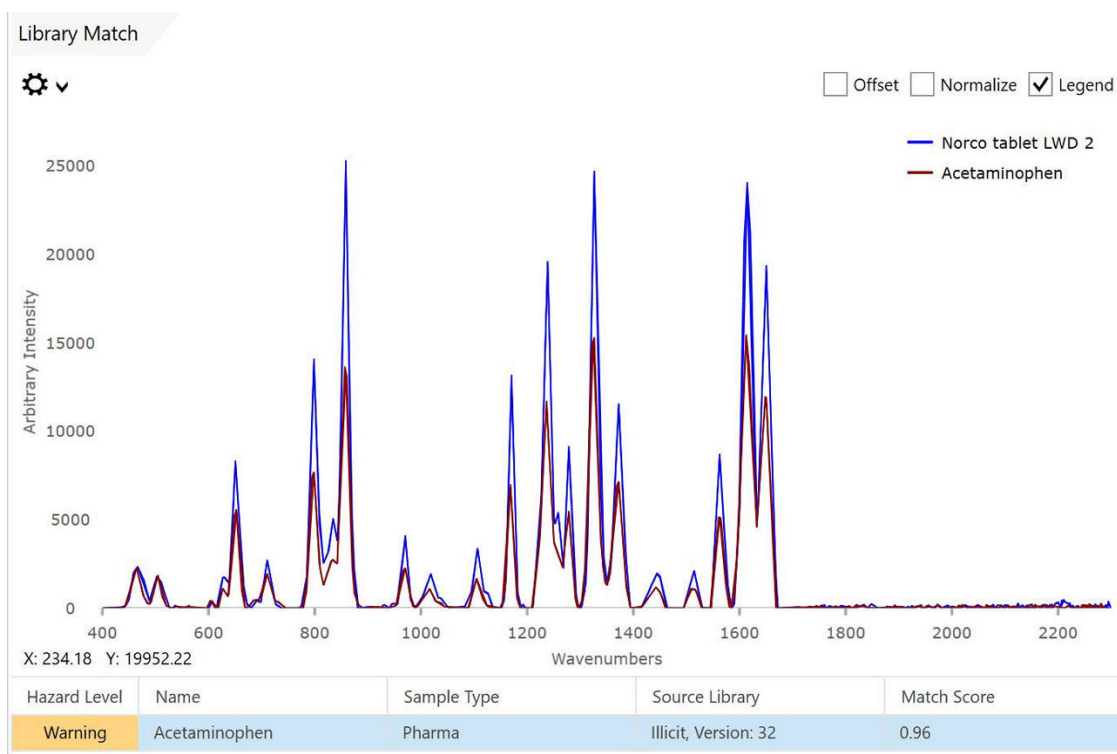


Figure 1. Mira Cal DS spectrum of Norco® tablet with acetaminophen library match supported by high HQI = 0.96.

SAMPLING WITH SERS

Crushed Norco® was dissolved in methanol to prepare a 1 mg/mL solution. Then, 100 µL each of the Norco® solution and 50 mmol/L NaCl solution were added to a glass vial with 800 µL of proprietary SERS gold colloid solution (available from Metrohm Raman). The vial was capped, shaken, then inserted into the Vial

Holder Smart Tip on MIRA DS. Selection of the SERS OP, laser activation, and sample acquisition was completed in **seconds**.

The presence of hydrocodone in Norco® is clearly indicated by peak agreement between sample and SERS library spectra. (**Figure 2**).

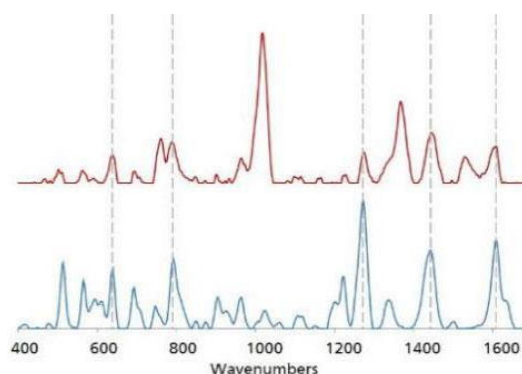


Figure 2. The gold SERS spectrum of Norco® (red), is overlaid with a gold SERS library spectrum of hydrocodone (blue). Dotted lines indicate peak agreement between the two spectra.

RESULTS AND DISCUSSION

This application is an excellent illustration of two important functions of the MIRA DS handheld Raman system. A good way to characterize Raman sampling is, «If you can see it, Raman can identify it.» Raman also excels at through-container sampling, meaning that information about a material can be collected with no sample contact. This may not be crucial for a sample like Norco®, but during encounters with synthetic opioids, which can be deadly in microdoses, it is essential.

In contrast, SERS is used primarily as a method for detection of trace materials or contaminants in a larger sample. SERS can also be used to study materials in a sample that would



otherwise fluoresce.

In this application, both the bulk acetaminophen and trace hydrocodone were determined in a single tablet of Norco® using MIRA DS.

CONCLUSION

Dual analysis with Raman and SERS can be an effective tool in the fight against abuse of prescription medications, which is responsible for nearly a third of annual opioid-related deaths. It is also an excellent solution for law enforcement professionals who encounter unknown street drugs containing potentially

illegal or even lethal substances.

With MIRA DS, one compact instrument can be used—even by non-technicians in the field—to quickly and easily perform both bulk Raman identification and SERS trace analysis with great accuracy.

CONTACT

瑞士万通中国
北京市海淀区上地东路1号
院1号楼7层702
100085 北京

marketing@metrohm.com.cn

CONFIGURATION



MIRA DS Advanced

Metrohm Instant Raman Analyzer (MIRA) DS 是一款坚固的高性能手持式拉曼光谱仪,用于快速、无损地确定非法物质,例如毒品、炸药、起始原料和有害物质。尽管仪器体积小,但 Mira DS 非常坚固耐用,提供有高效光谱仪,并配备了我们独特的栅格环绕扫描 (ORS) 技术。

该高级套件包括校正标准件、供瓶装或袋装分析或者直接分析用的通用附件以及直角附件,非常适合在表面和/或袋中运行样品。3B 类操作。MIRA DS 支持瑞士万通手持式拉曼功能库。



(LWD)

Mira M-3/P/DS 高级套件包含一个适用于 7.6 mm 焦距长距离的附加透镜 (LWD)。3B 类操作。



ID – Au NP

ID 套件 - Au NP 包含了 Mira/Misa 用户使用胶体金溶液进行 SERS 分析所需的组件。该套件包含了一个一次性抹刀、一个移液管、样品小瓶和一个含金胶体的瓶子。