



907 Titrando

2.907.0020

High-end titrator for potentiometric and volumetric Karl Fischer titration with two measuring interfaces and Dosino dosing units.

- up to four dosing device systems of the 800 Dosino type
- dynamic (DET), monotonic (MET) and endpoint titration (SET), enzymatic and pH-STAT titrations (STAT), Karl Fischer titration (KFT)
- "iTrode" intelligent electrodes
- Measurement with ion-selective electrodes (MEAS CONC)
- Dosing functions with monitoring, liquid handling
- four MSB connectors for additional stirrers or dosing device systems
- USB connector
- For use with OMNIS Software, *tiamo* software, or Touch Control
- Complies with GMP/GLP and FDA regulations such as 21 CFR Part 11, if required

Below, the accessories are grouped into Scope of delivery and Optional accessories.

Please keep this printout at hand for ordering replacement material.

These lists may be subject to change.

Scope of delivery 2.907.0020

Qt.	Order no.	Description	
1 PCS	1.907.0020	907 Titrande	
		High-end titrator for use with intelligent electrodes – iTrodes – with	
		• up to four dosing systems of the 800 Dosino type • dynamic (DET), monotonic (MET) and endpoint titration (SET), enzymatic and pH-STAT titrations (STAT), Karl Fischer titration (KFT) • Measurement with ion-selective electrodes (MEAS CONC) • Dosing functions with monitoring (DOS), Liquid Handling • four MSB connectors • two galvanically isolated measuring interfaces • USB connector	
2 PCS	6.2043.005	Holding clip for bottles	
		Holding spring for reagent bottles in exchange units.	

Optional accessories

Order no.	Description
6.2148.020	RS-232/USB box Extends the Titrando, USB Sample Processors, Dosing Interface, Titrosampler by two RS-232 interfaces for connecting e.g. balances, PC, LIMS. Requires 6.2151.030 cable. 
2.800.0010	800 Dosino The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm). 
2.801.0040	801 Stirrer with stand Magnetic stirrer including base plate, support rod and electrode holder for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimat and 780/781 pH meters as well as the 856 and 867 measuring modules. With permanently attached cable for MSB (Metrohm Serial Bus). 
2.802.0040	802 Stirrer for 804 Ti Stand Rod stirrer including 6.1909.010 stirring propeller. 

2.804.0040

804 Ti Stand with stand

Titration stand and controller for 802 Rod Stirrer. The 804 Ti Stand together with the optional 802 Rod Stirrer provides an alternative to the magnetic stirrer. Ti Stand with base plate, support rod and electrode holder.



2.805.0010

805 Dosimat

Dosing device for the Titrando and Sample Processors with read/write hardware for intelligent Exchange Units. With permanently attached cable. Without Exchange Unit.



2.846.0010

846 Dosing Interface

USB-capable control unit for connecting a maximum of four 800 Dosinos or 805 Dosimats for dosing and Liquid Handling tasks. A Touch Control or the connection to a PC with OMNIS Software, tiamo™, MagIC Net, viva, or 797 VA Computrace is required for operation.



6.2104.020

Electrode cable / 1 m / F

For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



6.2104.600

Electrode cable for plug in head U/plug F, 2x2 mm B, 1m

For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2148.010

Remote Box MSB

Additional remote interface for the connection of devices that can be controlled via remote lines. With permanently attached cable.



6.3026.110

Exchange Unit 1 mL

Exchange Unit with integrated data chip with 1 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.



6.3026.150

Exchange Unit 5 mL

Exchange Unit with integrated data chip with 5 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.



6.3026.210

Exchange Unit 10 mL

Exchange unit with integrated data chip with 10 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



6.3026.220

Exchange Unit 20 mL

Exchange unit with integrated data chip with 20 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



6.3026.250

Exchange unit 50 mL

Exchange unit with integrated data chip with 50 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion tip and standard amber glass reagent bottle.



6.3032.120

807 Dosing Unit 2 mL

807 Dosing Unit with integrated data chip with 2 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.150

807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.210

807 Dosing Unit 10 mL

807 Dosing Unit with integrated data chip with 10 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.220

807 Dosing Unit 20 mL

807 Dosing Unit with integrated data chip with 20 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.




6.3032.250

807 Dosing Unit 50 mL

807 Dosing Unit with integrated data chip with 50 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



Operating unit for the Titrandos, USB Sample Processors, 856 Conductivity Module, 867 pH Module and 846 Dosing Interface. Touch-sensitive, high-resolution color display, simple and intuitive operation, thanks to Favorites for direct method access. With integrated Ethernet interface for direct connection to the Internet and USB interface for connecting USB printers or a USB memory stick.

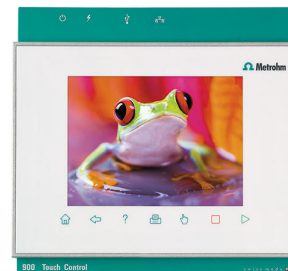

Dialog languages: German, English, Chinese, French, Spanish, Portuguese, Russian, Korean, Polish and Italian.

The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).



Mini USB adapter cable for connecting one 854 iConnect to an 867 pH Module, 901 to 907 Titrandos and 916 Ti-Touch, 25 cm. The adapter cable remains plugged into the instrument as strain relief and the 854 iConnect is only plugged into and unplugged from the adapter cable.



For connecting plug B (4 mm) to 2 mm socket.



6.2103.140

Adapter black 2 mm plug / B socket 4 mm

For connecting plug B (4 mm) to 2 mm socket.

