



949 pH Meter without electrode

2.949.0210

Discover the 949 pH Tabletop Meter – Your compact pH measuring instrument for quick and reliable laboratory measurements.

- Complete pH measuring instrument with color display
- Precise pH, mV and ORP measurements with integrated temperature display
- Simple calibration with Metrohm, USA, and NIST buffers up to 3 points or by 2-point calibration with in-house buffers
- Clearly structured display of the calibrated values using symbols for simple read-off
- 3 stability criteria for measured value recognition for reliable results
- GLP conformity of the measured values with embedded date and time
- Overview of calibration data with monitoring function
- RS232 port enables connection of a printer.

Included in the package are the 949 pH Meter, the electrode holder, and the power supply unit for immediate use in the laboratory.

Below, the accessories are grouped into Scope of delivery and Optional accessories.

Please keep this printout at hand for ordering replacement material.

These lists may be subject to change.

Scope of delivery 2.949.0210

Qt.	Order no.	Description
1 PCS	1.949.0210	949 pH Meter with electrode holder and power supply unit



1 PCS	6.2001.140	Electrode holder with stand base
		Electrode holder for 949 pH Meter



1 PCS	6.2166.200	Power supply unit for 949 pH Meter
		Tabletop power supply unit, 100–240 V, 50–60 Hz, 5 V DC, for 949 pH Meter
		Efficiency Level VI



Optional accessories

Order no.	Description	
6.2307.230	Buffer solutions pH 4, 7 and 9 Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL	
6.2104.020	Electrode cable / 1 m / F For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).	
6.2104.600	Electrode cable for plug in head U/plug F, 2x2 mm B, 1m For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).	
6.1110.100	Pt1000 temperature sensor (installation length 12.5 cm) Pt1000 temperature sensor (class B) made of glass. This PT1000 temperature sensor is also available under the article number 6.1110.110 with an installation length of 17.8 cm.	

6.2104.140

Electrode cable / 1 m / 2 x 2 mm

For connecting temperature sensors with Metrohm plug-in head G with 780, 781, 826, 827 pH/Ion meters and all Titrandos



6.0228.000

Solitrode with Pt1000 (fixed cable 1.2 m)

Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m). This electrode is suitable for routine pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.



Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

6.00226.600

Spearhead electrode with Pt1000

Maintenance-free combined pH electrode (gel electrolyte) for piercing measurements of all types (e.g., with cheese, meat, dough) with integrated Pt1000 temperature sensor. The electrode is stored in saturated potassium chloride solution $c(\text{KCl}) = \text{sat.}$ (6.2308.000) and is not suitable for low-ion solutions. The ageing indicator gives early indication of when the electrode needs to be replaced.



6.0256.100

Flat membrane electrode

This combined pH electrode is well suited for:

- pH measurements on surfaces, such as paper, textiles, leather
- pH measurement/titration in small sample volumes (min. immersion depth = 1 mm)
- pH measurement/titration of aqueous suspension

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.



Compact impact printer with USB and RS232 interface for

- 900 Touch Control (+)
- 915 KF Ti-Touch (+)
- 916 Ti-Touch (+)
- 917 Coulometer (+)
- 91X Meter (+)
- Eco Dosimat/Titrator (+)
- 877/848 Titrino plus (additional 6.2151.100 required)
- 865/876 Dosimat plus (additional 6.2151.100 required)
- 862 Compact Titrosampler (additional 6.2151.100 required)
- 870 KF Titrino plus (additional 6.2151.100 required)
- 899 Coulometer (additional 6.2151.100 required)



(+) Cable included in scope of delivery

6.0257.600

Aquatrode Plus with Pt1000

Combined pH electrode with integrated Pt1000 temperature sensor for pH measurements/titrations in ion-deficient aqueous media (e.g., drinking water, process water). This electrode shows a very short response time in these samples.

The fixed ground-joint diaphragm is insensitive to contamination. When $c(\text{KCl}) = 3 \text{ mol/L}$ is used as bridge electrolyte, storage in storage solution is recommended. The bridge electrolyte can be easily replaced with a chloride-free electrolyte (e.g., potassium nitrate $c(\text{KNO}_3) = 1 \text{ mol/L}$ (6.2310.010)), storage in the used bridge electrolyte.



6.0239.100

Viscotrode

Combined pH electrode for pH measurement / titration of:

- viscous samples (e.g., nutrient solutions, emulsions)
- protein- or sulfide-containing media

The flexible ground-joint diaphragm is insensitive to contamination and easy to clean.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.



6.0224.100

Biotrode

Combined pH electrode for measurements in very small sample volumes (>50 µL) and biological samples.

Idrolyte (6.2308.040) is used as reference electrolyte and storage solution.



6.2307.100

Buffer solution pH 4 (500 mL)

Ready-to-use buffer solution pH = 4 (25 °C), colored, with closure seal



6.2307.110

Buffer solution pH 7 (500 mL)

Ready-to-use buffer solution pH = 7 (25 °C), colored, with closure seal



6.2307.120

Buffer solution pH 9 (500 mL)

Ready-to-use buffer solution pH = 9 (25 °C), colorless, with closure seal



Combined pH electrode with integrated Pt1000 temperature sensor for pH measurements/titrations. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements




The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m, diameter banana plug 2 mm). This electrode is particularly suitable

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

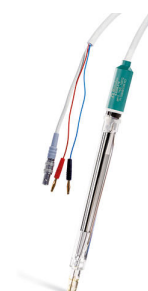

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

Combined platinum ring electrode with a ceramic pin diaphragm. This electrode is well suited for redox titrations when the pH value varies, e.g.:

- oxygen content according to Winkler
- determination of hydrogen peroxide with KMnO_4
- diazotization titrations



$c(\text{KCl}) = 3 \text{ mol/L}$ is used as reference electrolyte and storage solution.

6.2306.020 Redox standard 250 mV (250 mL)

Redox standard (with Ag/AgCl/c(KCl)=3 mol/L reference electrode gives
 $U = 250 \pm 5 \text{ mV}$ (20° C); can also be used as buffer pH=7)



6.2307.200 Buffer solution pH 4.00 (30 x 30 mL)

Buffer solution pH 4.00 (25 °C) in single use sachets, colorless



6.2307.210 Buffer solution pH 7.00 (30 x 30 mL)

Buffer solution pH 7.00 (25 °C) in single use sachets, colorless



6.2307.220 Buffer solution pH 9.00 (30 x 30 mL)

Buffer solution pH 9.00 (25 °C) in single use sachets, colorless

