


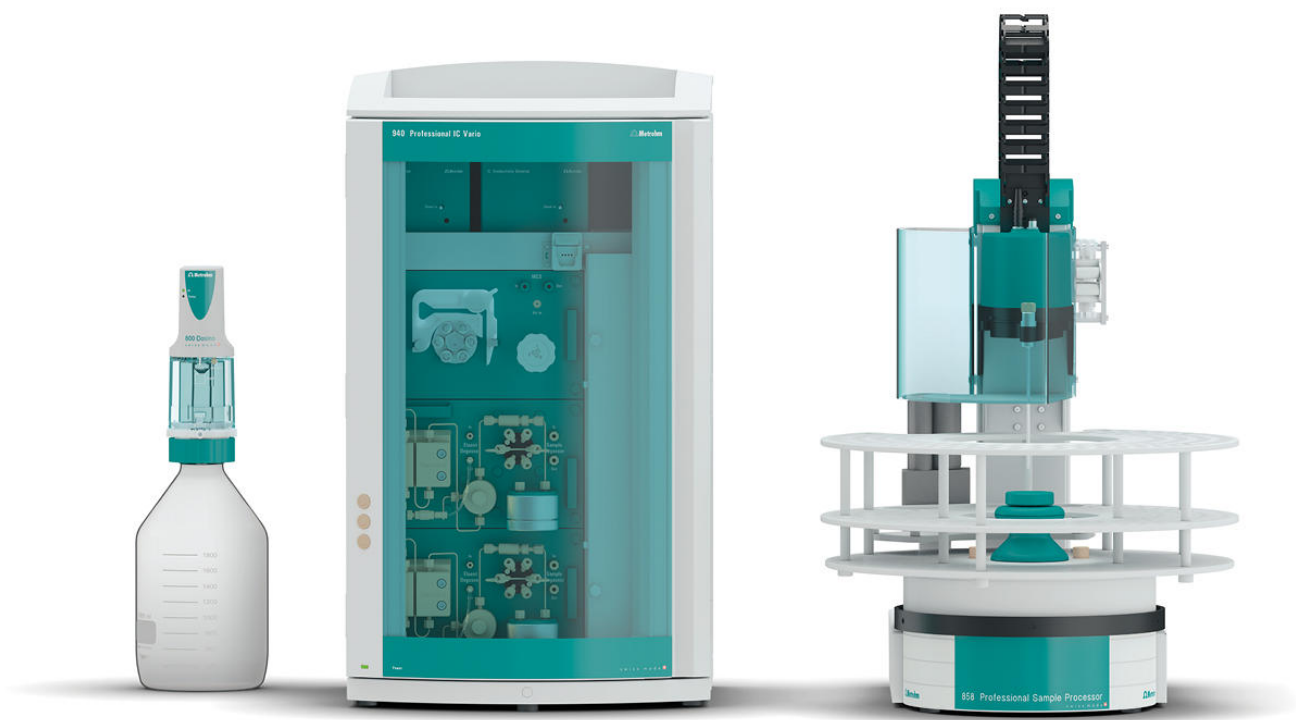

"ProfIC Vario 6 AnCat" – Professional IC Vario system with Inline Dilution and Inline Ultrafiltration

ProfIC-Vario-06-AnCat

The Professional IC Vario system with Metrohm Inline Dilution, Metrohm Inline Ultrafiltration, sequential suppression, and conductivity detection enables fully automatic determination of anions and cations that are present in very high concentrations and that have matrices containing particles.

Typical areas of application:

- Rinse, process, and waste water
- Extractions and digestion solutions
- Food samples

Schematic representation


Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery ProFI C-Vario-06-AnCat

Qt.	Order no.	Description
1 PCS	2.940.2500	940 Professional IC Vario TWO/SeS/PP
The 940 Professional IC Vario TWO/SeS/PP is the intelligent two-channel IC instrument with sequential suppression (one channel) and a peristaltic pump for suppressor regeneration. The instrument can be used with any separation and detection methods. Typical areas of application: • Standard instrument for parallel anion and cation determinations • Trace analysis for anions and cations • Online monitoring for anions and cations		
2 PCS	2.850.9010	IC Conductivity Detector
Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.		
1 PCS	2.858.0020	858 Professional Sample Processor – Pump
The 858 Professional Sample Processor – Pump processes samples from 500 µL to 500 mL. The sample transfer takes place either with the installed bidirectional two-channel peristaltic pump or with an 800 Dosino.		
1 PCS	2.741.0010	741 Magnetic Stirrer
Magnetic stirrer for sample changers.		



The magnetic stirrer is located beneath the sample rack.

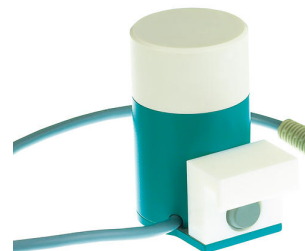

1 PCS

2.800.0010

800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



1 PCS

6.5330.120

IC equipment: Inline dilution

Accessory set for Inline Dilution. For use with the 858 Professional Sample Processor, 800 Dosino and 741 Magnetic Stirrer.




1 PCS

6.05330.110

IC equipment: Inline Ultrafiltration 2 - push mode

Accessory set for Inline Ultrafiltration 2 - push mode. For use with the 858 Professional Sample Processor/ 919 IC Autosampler plus.

Recommended set for the combination with Inline Dilution (MIDT).

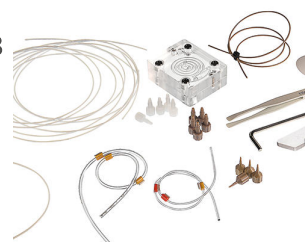

Optional accessories

Order no.	Description	
6.2041.760	Sample rack 54 x 11 mL + 1 x 300 mL Sample rack with small diameter for 54 x 11 mL sample vessels in two rows and additionally 1 x 300 mL bottle.	
6.2041.440	Sample rack, 148 x 11 mL and 3 x 300 mL (PVC) Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, suitable for 148 sample vessels and 3 rinsing bottles. The following sample vessels and rinsing bottles can be used: 6.2743.050, 6.1608.080, 6.1627.000, 6.2743.040. Plastic: Polyvinyl chloride (PVC)	
6.2743.050	Sample vessel 11 mL 2000 pieces. For IC Sample Processors and VA Autosamplers. The IC sample vessels are made of polypropylene (PP) plastic for health-care applications. Production and packaging both take place under sterile and clean room conditions (Clean Room Class 7, in accordance with DIN ISO 14644). This guarantees minimal degrees of contamination. The sample vessels are subjected to random checks for leachable anions and cations.	
6.2743.070	Stopper with perforation 2000 pieces. For IC Sample Processors	
6.1050.420	Metrosep C 4 - 150/4.0	

The Metrosep C 4 - 150/4.0 is the universal standard column in cation analysis for accomplishing high separating efficiency in a short time. The Metrosep C 4 - 150/4.0 is the ideal separation column for the analysis of alkaline and earth alkaline metals in aqueous media.



6.1050.500 Metrosep C 4 Guard/4.0

The Metrosep C 4 Guard/4.0 contains the C-4 column material and is used to protect all Metrosep cation columns that have a substrate based on silica gel. Particles and contamination are reliably retained, considerably prolonging the service life of the analytical separation column. The economical price is an additional plus. The Metrosep C 4 Guard/4.0 also works based on the "On Column Guard System" and is screwed directly onto the separation column with virtually no dead volume.



6.2842.000 MSM-HC Rotor A

Suppressor rotor for all IC instruments with MSM-HC (Metrohm Suppressor Module with high capacity)



6.2844.000 MSM-LC Rotor A

Suppressor rotor for use with 2 mm columns.



6.2832.000 MSM Rotor A

Suppressor rotor for all IC instruments with MSM (Metrohm Suppressor Module)



6.2842.020 Adapter sleeve, MSM

Not only the MSM Rotor A but also the MSM-LC Rotor A can be installed in the rotor housing with the aid of the adapter sleeve.



2.941.0010 941 Eluent Production Module

The 941 Eluent Production Module enables the automatic production of an eluent. It enables continuous working without manual intervention and guarantees stable retention times. It can be combined with all Metrohm IC instruments and is monitored by MagIC Net.



6.5330.090 IC equipment: Additional eluent on the Eluent Production Module

Accessory set for Eluent Production Module. Contains all parts necessary for expansion by one eluent. The level sensor must be ordered in accordance with the application.



2.800.0010 800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



6.2769.110 EMPTY sensor for level sensor for 2 L bottles

This sensor, mounted in the 6.1626.000 bottle attachment, results in the level sensor EMPTY for 2 L bottles. Included are two sensor rods, one aspiration tubing and two holders for fixing the tubing to one of the sensor rods.



6.3032.150 807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.2744.080 M6 thread / UNF 10/32 adapter

With M6 outer thread. For IC instruments.



6.6059.422 MagIC Net 4.2 Professional: 1 license

Professional chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Auto samplers, 800 Dosino, 771 Compact Interface, etc. The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses.

Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.



MagIC Net Professional complies fully with FDA Regulation 21 CFR Part 11 as well as GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

1 license

The installation and documents are supplied on a USB flash drive.

6.6059.423

MagIC Net 4.2 Multi: 3 licenses

Client/server-compatible chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Autosamplers, 800 Dosino, 771 Compact Interface, etc.

The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses. Data and method are centrally stored on a Server.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Multi complies fully with FDA Regulation 21 CFR Part 11 as well as with GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

Performance for client/server systems with servers outside the in-house company network has been improved.

Client/server version with 3 licenses
