


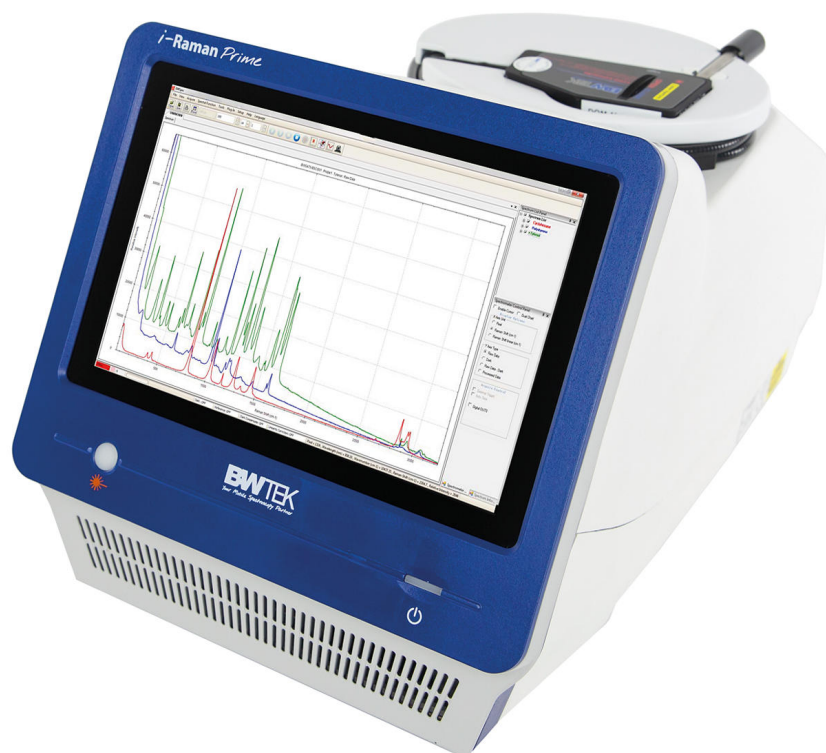

i-Raman Prime 785H Portable Raman Spectrometer


BWT-840000674

The i-Raman[®] Prime 785H is a low-noise, high-throughput, and fully integrated Raman system with an embedded tablet computer and a fiber-optic sampling probe. Using a high-quantum-efficiency CCD array detector with TE deep cooling (-25 °C) and high dynamic range, this portable Raman spectrometer delivers research-grade Raman analysis capabilities, including real-time quantitation and identification. The high throughput gives Raman spectra with excellent signal-to-noise ratio, making it possible to measure rapid processes and to measure even the slightest Raman signals, detecting subtle sample differences.

In addition to its portable construction, the i-Raman Prime 785H features the unique combination of wide spectral coverage and high resolution, thus enabling measurements from 150 cm⁻¹ to 2,800 cm⁻¹. The i-Raman Prime can be battery-operated for easy portability, providing research-grade Raman analysis capabilities for high-precision qualitative and quantitative work wherever needed. The system is optimized for use with our STRaman[®] technology for analyses through non-transparent packaging.

Below, the accessories are grouped into Scope of delivery and Optional accessories.
 Please keep this printout at hand for ordering replacement material.
 These lists may be subject to change.

Scope of delivery BWT-840000674

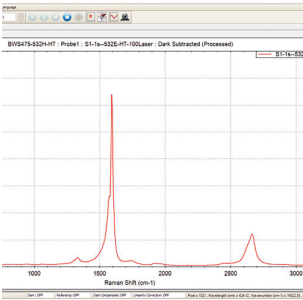Qt.	Order no.	Description
1 PCS	607510140	High Throughput Lab-grade Raman Probe (785 nm) High-throughput lab-grade fiber optic Raman probe assembly with trigger button. 1.5 meter fiber length. Compatible with i-Raman NxG 785S and 785H models. 
1 PCS	607510040	Polystyrene Validation Cap Validation cap with integrated polystyrene standard for lab-grade Raman probes with shaft diameter of 9.5 mm. 
1 PCS	607510130	Laser Safety Goggles, 785 nm wavelength Laser Safety Goggles, 785 nm Wavelength. 
1 PCS	BWT-810000044	BWSpec Software BWSpec[®] is B&W Tek's general spectroscopy software for instrument control and data acquisition, including real-time peak analysis and trending. BWSpec is the operating software included with the purchase of all B&W Tek portable Raman systems and spectrometer products. It is designed with features for broad range applications, performing complex measurements and calculations at the click of a button. It supports 

multiple data formats and provides the option to optimize measurement parameters, such as integration time and laser output power control. In addition to data acquisition and data processing, it also offers automatic dark removal, spectral smoothing, baseline correction, as well as peak monitoring and trending.

1 PCS

BWT-810000721 Carrying Case for i-Raman Pro/Prime Series

Carrying Case for i-Raman Pro and Prime series. Includes a soft, black carrying case with wheels and custom-cut foam for packing and safely carrying the i-Raman Pro, i-Raman Prime, STRam, or QTRam.

CCRM



1 PCS




BWT-840000296 Stylus/pen

Stylus/pen combo for touch screens.


NR-STY



Optional accessories

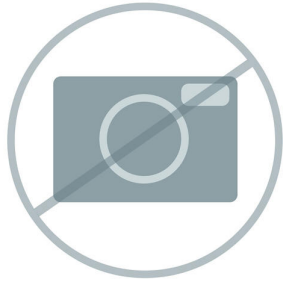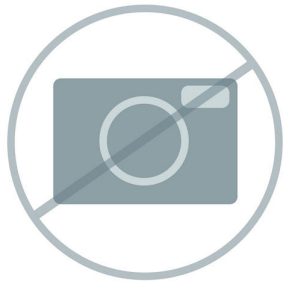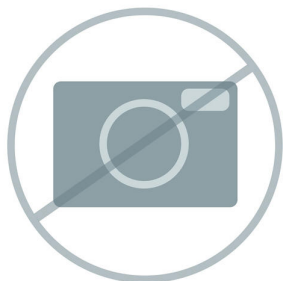Order no.	Description	
BWT-840001202	Hastelloy immersion shaft for 785nm Raman probe RIS100-HS-785-08: Hastelloy immersion shaft for BAC102/BAC100B Raman probes with excitation wavelength 785nm. Features a 0.5" (12.7 mm) O.D. x 8" (203.2 mm) L Hastelloy C-276 shaft and a gold- sealed sapphire ball lens. Working distance 0.4 mm in air, 0.6 mm in water. Shaft operating temperature is up to 250 degrees C, pressure up to 4000 psi. Shaft must be purchased with B&W Tek Raman probe. RIS100-HS-785-08	
BWT-810000719	Power cord for EU power for use with portable Raman systems Power cord for EU power for use with portable Raman systems. BWPC-EU	
BWT-810000468	Connection cleaner Cleaner for SMA905 and FC/PC connections.	
BWT-810000469	Cleaner for ferrule end surfaces FC/PC and FC/APC receptacle cleaner recommended for i-Raman Pro /Prime, STRam, QTRam. Do not use on SMA. FRC	
BWT-810000623	Microscope head tripod adaptor	

BAC151x-TRI tripod adapter for the BAC151x sampling system with Raman video microscope.

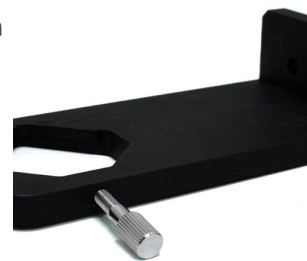


BWT-840000117 Raman Flow Cell (BAC160-SS)

Raman flow-through cell made of stainless steel, for inline Raman measurements of liquid samples by interfacing with standard laboratory probe (9.5 mm diameter; probe is not included in scope of delivery).




BWT-840000126 BWID Standard Raman ID Software

Data acquisition, material identification and reporting software for B&W Tek portable Raman spectrometers to rapidly identify or verify materials based on user-built libraries, proprietary B&W Tek and third-party libraries. Includes system performance tests. One user license provided. Free upgrades included within one year of purchase date.

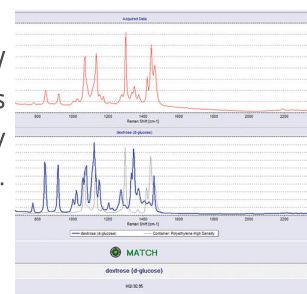


BWT-840000127 BWID Pharma Raman ID Software

Data acquisition and material identification software for B&W Tek portable Raman spectrometers to rapidly identify or verify materials based on user-built libraries, proprietary B&W Tek and third-party libraries. Meets all the requirements of FDA 21 CFR Part 11, including system access security, audit trail of data activities, and system performance tests. One user license provided. Free upgrades included within one year of purchase date.

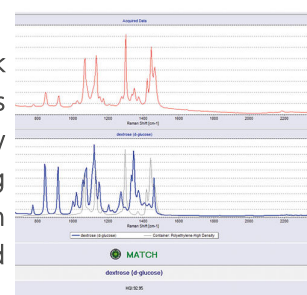


BWT-840000130 Immersion shaft for portable Raman instruments

Optional immersion shaft for BAC100/BAC102 Raman probes with an excitation wavelength of 532 nm or 785 nm. Equipped with a 76.2 mm long body made of 316L stainless steel with an outer diameter of 12.0 mm and a silica window that is sealed with an O-ring made of perfluoro rubber. Operating distance 5.0 mm in air. Can be used at -55 °C to 200 °C. The seal is liquid-tight at low pressure (1 bar).

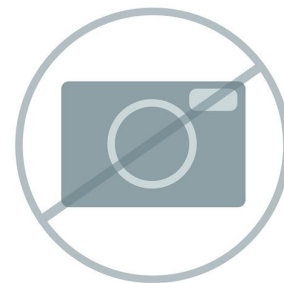


BWT-840000131 Immersion shaft for portable Raman instruments

Optional immersion shaft for BAC100/BAC102 Raman probes with an excitation wavelength of 532 nm or 785 nm. Equipped with a 76.2 mm long body made of 316L stainless steel with an outer diameter of 12.0 mm and a sapphire-glass window that is sealed with an O-ring made of perfluoro rubber. Operating distance 5.0 mm in air. Can be used at -55 °C to 200 °C. The seal is liquid-tight at low pressure (1 bar).

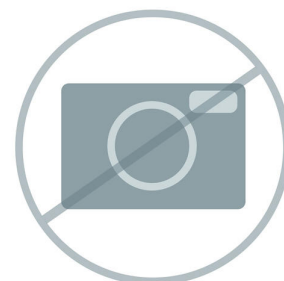

BWT-840000280 GemID Software

Data acquisition and identification software for the analysis of gemstones and minerals with B&W Tek portable Raman spectrometers. Includes Raman Spectral Library and images of more than 450 gemstones from Gem Expert and Japan Germany Gemmological Laboratory. One user license provided. Free upgrades included within one year of purchase date.

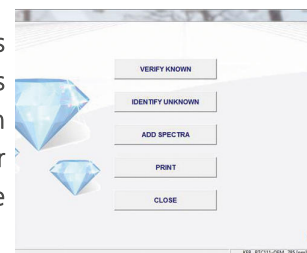

BWT-840000295 Vials for liquids, 15 mm in diameter


Borosilicate glass vials (diameter 15 mm), package of 6, with pouch.



BWT-840000322 Video Microscope Objective, 10x

Microscope objective, infinity-corrected, 10x, working distance (mm) = 16, focal length (mm) = 20, numerical aperture (NA) = 0.3.


RML110A

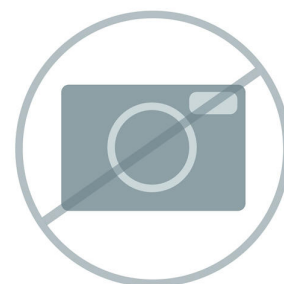

BWT-840000323 Video Microscope Objective, 20x (RML120A)

Microscope objective, infinity-corrected, 20x, working distance (mm) = 12, focal length (mm) = 10, numerical aperture (NA) = 0.4.


RML120A

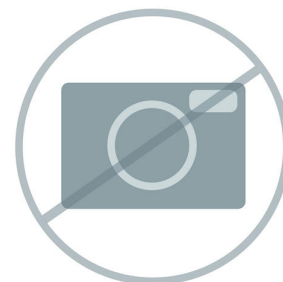

BWT-840000325 Video Microscope Objective, 50x (RML150A)

Microscope objective, infinity-corrected, 50x, working distance (mm) = 9.15, focal length (mm) = 4, numerical aperture (NA) = 0.55.


RML150A

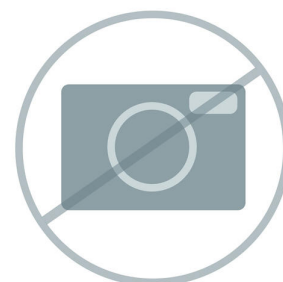

BWT-840000328 Video Microscope Objective, 100x

Microscope objective, infinity-corrected, 100x, working distance (mm) = 3.2, numerical aperture (NA) = 0.8, Focal length (mm) = 2.

RML1100A

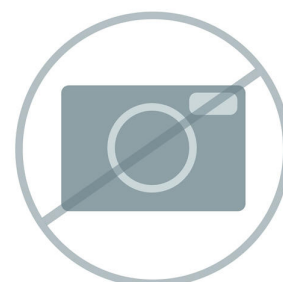

BWT-840000332 Tablet Holder

Tablet holder compatible with all models of BAC100/BAC102 lab-grade fiber-optic Raman probe with 9.5 mm shaft diameter.



BWT-840000361 Handheld Laser Power Meter

Laser Power Meter, Pen Type, Calibrated. Range: 5 micro-Watts to 1 Watt. Used in OQ/PQ and annual certification laser power measurement.



BW840000463 BAC151x probe adapter

BAC151x Probe Adaptor to connect B&W Tek Lab grade Raman probes to BAC151x Video microscope sampling system. Remove the standard shaft and distance regulator, to use the adaptor and connect the probe head to the video microscope. 3 adaptors per pack.



BWT-840000679 Raman Zoom Lens (20-60 mm)

Zoom lens for our BAC100/BAC102 fiber-optic Raman probes. Working distance is adjustable from 20 mm to 60 mm. Suitable for a Raman excitation range from 500 nm to 850 nm.



BWT-840000711 Raman Zoom Lens (20-60 mm)

Zoom lens for our BAC100/BAC102 fiber-optic Raman probes. Working distance is adjustable from 20 mm to 60 mm. Suitable for a Raman excitation range from 785 nm to 1064 nm.



BWT-840000680 Raman Zoom Lens (60-600 mm)

Long-working-distance zoom lens for BAC100/BAC102 fiber-optic Raman probes. Working distance is adjustable from 60 mm to 600 mm. Suitable for a Raman excitation range from 500 nm to 850 nm.



BWT-840000712 Raman Zoom Lens (60-600 mm)

Long-working-distance zoom lens for BAC100/BAC102 fiber-optic Raman probes. Working distance is adjustable from 60 mm to 600 mm. Suitable for a Raman excitation range from 1000 nm to 1100 nm.



BWT-840000681 Raman Stand-off Lens (0.6-6 m)

A telescope Stand-off lens for BAC100/BAC102 fiber-optic Raman probes. Wide zoom range for sample measurement ranges from 0.6 meters to 6 meters. Suitable for a Raman excitation range from 500 nm to 850 nm.

RTS202-VIS-NIR



BWT-840000843 Raman Gemstone Spectral Library

Raman Gemstone Spectral Library for BWID and portable Raman instruments with > 450 spectra from Gem Expert and from Japan Germany Gemmological Laboratory. Library supported in BWID Version 2.02 or higher.



BWT-840000847 Raman Explosives Spectral Library

Raman Explosives Spectral Library for portable Raman instruments including more than 180 spectra. Supported in BWID Version 2.02 and higher.

Note: This item may be subject to Export control to some territories.



BWT-840000910 ST Probe

Lab-grade Raman fiber-optic probe module for analyses through packaging, with hand trigger, 105 μ m diameter, 0.22 NA fiber for 785 nm excitation, 300 μ m diameter, 0.22 NA for Raman pickup, optical

density >6 ; up to 150 cm^{-1} . 1.5 m fiber cable with FC/PC connector on the spectrometer end and FC/PC connector on the laser end. Adhesive sealed quartz window. Not for immersive use. Accessories included: Focus adapter, surface regulator, calibration cap.



BWT-840000961 Raman Video Micro-Sampling System (785 nm)

Video microscope sampling system for use with B&W Tek's lab and industrial Raman probes. Includes a 20x objective at a working distance of 16 mm. Offers manual rough and fine adjustment on the X, Y, and Z-axes, coaxial LED illuminator for target alignment, video camera for sample observation, and is compatible with standard microscope objectives. Probe not included, available separately. 785 nm configuration.



BAC151C-785

BWT-840000963 Raman Video Micro-Sampling System (532/785 nm dual)

Video microscope sampling system for use with B&W Tek's lab and industrial Raman probes. Includes a 20x objective at a working distance of 16 mm. Offers manual rough and fine adjustment on the X, Y and Z-axes, coaxial LED illuminator for target alignment, video camera for sample observation, and is compatible with standard microscope objectives. Probe not included, available separately. Dual 532/785 nm configuration.



BAC151C-532/785

BWT-840000964 Raman Video Micro-Sampling Head (785 nm)

Video microscope sampling system head for use with B&W Tek's lab and industrial Raman probes. With coaxial LED illuminator for target alignment and video camera for sample observation. Compatible with standard microscope objectives. Probe not included, available separately. Objective lens not included, available separately. 785 nm configuration.




BWT-840000966 Raman Video Micro-Sampling Head (532/785 nm dual)

Video microscope sampling system head for use with B&W Tek's lab and industrial Raman probes. With coaxial LED illuminator for target alignment and video camera for sample observation. Compatible with standard microscope objectives. Probe not included, available separately. Objective lens not included, available separately. Dual 532/785 nm configuration.



BWT-840001123 Bluetooth 1D barcode scanner

Bluetooth 1D barcode scanner for use with BWID software operating a portable Raman instrument

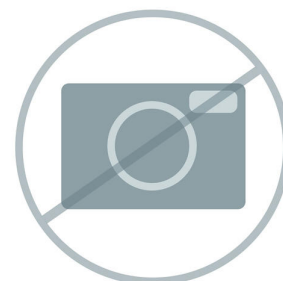

BWT-810000682 External rechargeable battery pack for portable Raman instruments

External Lithium-Ion battery that provides constant 12 VDC. (Compatible with the i-Raman-Plus, i-Raman Pro and i-Raman EX models).

EBP106



BWT-840001137 Vision 4.1

Vision is data acquisition and method development software solution for spectroscopic analysis and control of B&W Tek portable Raman instruments, Metrohm XDS laboratory and process NIRS instruments. A user-friendly graphical analysis interface enables a simple application of chemometric algorithms for creating identification, qualification, and quantification methods and running them in real-time. With Vision you can store, manage, reprocess and exchange data.



BWT-840001136 Vision Pharma 4.1

Vision Pharma is data acquisition and method development software solution for spectroscopic analysis and control of B&W Tek portable Raman instruments, Metrohm XDS laboratory and process NIRS instruments. A user-friendly graphical analysis interface enables a simple application of chemometric algorithms for creating identification, qualification, and quantification methods and running them in real-time.

With Vision you can store, manage, reprocess and exchange data. Vision Pharma is complies with the US FDA's 21 CFR Pt 11 regulations for electronic signatures and records.



BWT-840000934 Polystyrene Reference for QTRam

Polystyrene reference for QTRam.

QTR-PCC



BWT-840000940 QTRam Optics and Sample Holder Set – 2 mm

Set of 2 mm excitation aperture, collection aperture, large round sample holder and large square sample holder for QT-Sampler.



BWT-840000941 QTRam Optics and Sample Holder Set – 4 mm

Set of 4 mm excitation aperture, collection aperture, large round sample holder and large square sample holder for QT-Sampler.



BWT-840000942 QTRam Optics and Sample Holder Set – 6 mm

Set of 6 mm excitation aperture, collection aperture, large round sample holder and large square sample holder for QT-Sampler.



BWT-840000943 QTRam Optics and Sample Holder Set – 8 mm


Set of 8 mm excitation aperture, collection aperture, large round sample holder and large square sample holder for QT-Sampler.



BWT-840001177 Aperture Removal Tool

Aperture removal tool for easy interchanging different apertures for optimal Raman spot size for transmission Raman with QT-Sampler.


TOOL-APERTURE

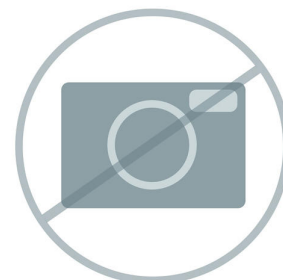

BWT-840001119 Focus adaptor

STRam 785 nm focus adaptor.


RST-FA-785

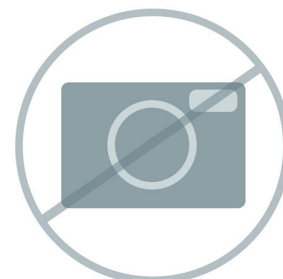

BWT-840001120 Light shield accessory for STRam


Light shield accessory for STRam.

RST-LS

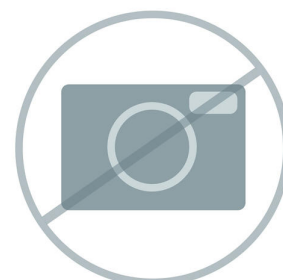

BWT-840001122 Surface regulator for ST probe

Surface regulator for ST probe.

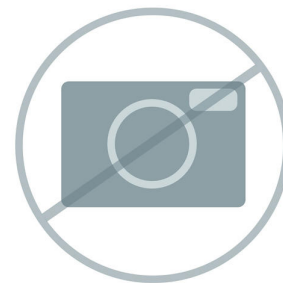

BWT-840001283 QT-Sampler

The QT-Sampler is an accessory for transmission Raman measurements, enabling content uniformity testing of pharmaceutical tablets. Compatible with i-Raman Prime 785S and H models. Apertures and sample holders are sold separately.



607506480 SERS Discovery Kit

Introductory kit for SERS analysis containing gold and silver P-SERS strips and gold and silver colloids.



BWT-840001265 BWID Comprehensive Raman Library

BWID Metrohm Comprehensive RamanLibrary - More than 27,700 spectral entries of chemicals, minerals, organic chemicals, solvents, polymers, dyes and pigments, TICs, TIMs, inorganics and organometallics, personal care products, and illicit materials.

PTL-MCRL



BWT-840001305 BWID MCRL - Polymers and Polymer Additives

Metrohm Comprehensive Raman Library (MCRL) Polymers and Polymer Additives for BWID.




BWT-840001309 BWID MCRL - Dyes, Pigments, and Stains

Metrohm Comprehensive Raman Library (MCRL) Dyes, Pigments, and Stains for BWID.



BWT-840001254 BWID Illicit and General Chemicals Library

BWID library of illicit materials includes drugs, TICs, TIMs, precursors, and general chemicals.

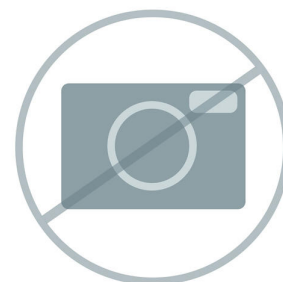

BWT-840001255 BWID CWA Library

Raman Chemical Warfare Agent (CWA) Library for BWID.



BWT-840001264 BWID Comprehensive USP Raman Library

Metrohm Comprehensive USP Raman Library for BWID containing ~ 900 traceable certified items.

