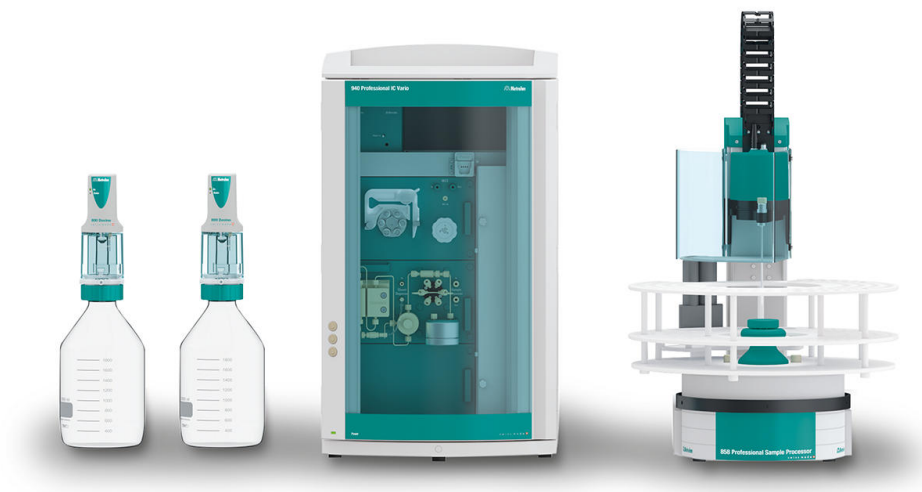




"ProfIC Vario 9 Anion" – Professional IC Vario system with Metrohm's intelligent preconcentration technique and Inline Matrix Elimination (MiPCT-ME)

ProfIC-Vario-09-Anion

The Professional IC Vario system with Metrohm's intelligent inline preconcentration technique, Inline Matrix Elimination, sequential suppression, and conductivity detection enables ultratrace analysis down to the lowest ppt range (ng/L) of anions or cations in complex matrices.

Typical areas of application:

- Ultratrace analysis in samples from nuclear power plants and conventional power plants
- Trace and ultratrace analysis in extraction agents and organics
- Trace analysis in snow and ice

Schematic representation

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery ProflC-Vario-09-Anion

Qt.	Order no.	Description
1 PCS	2.940.1500	940 Professional IC Vario ONE/SeS/PP
The 940 Professional IC Vario ONE/SeS/PP is the intelligent IC instrument with sequential suppression and a peristaltic pump for suppressor regeneration. The instrument can be used with any separation and detection methods. Typical areas of application: • Anion or cation determinations with sequential suppression and conductivity detection • Trace analysis for anions or cations • Online monitoring for anions or cations		
1 PCS	2.850.9010	IC Conductivity Detector
Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.		
1 PCS	2.858.0010	858 Professional Sample Processor
The 858 Professional Sample Processor processes samples from 500 µL to 500 mL. The sample transfer takes place either by means of a peristaltic pump on the 850 Professional IC system or with an 800 Dosino.		



2 PCS

2.800.0010

800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



1 PCS

6.5330.160

IC equipment: MiPCT-ME

Accessory set for assembling a Dosino for Metrohm intelligent Preconcentration Technique with Matrix Elimination (MiPCT-ME).



Optional accessories

Order no.	Description	
6.9920.191	Sample rack, 35 x 60 mL and 3 x 250 mL (PVC) Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, suitable for 35 sample vessels and 3 rinsing bottles. The following sample vessels and rinsing bottles can be used: 6.1608.100, 60 mL and 250 mL Nalgene FEP vessels. Plastic: Polyvinyl chloride (PVC)	
6.1608.100	PE bottle 50 mL Bottle for Extension Module Liquid Handling.	
6.2820.000	Al foil 1000 circular sheets of 80 mm diameter. Thickness: 0.010 mm. Used together with 6.2037.030 Foil holder for 6.1432.210 Sample beaker. For 674 KF Sample Changer	
6.2842.000	MSM-HC Rotor A Suppressor rotor for all IC instruments with MSM-HC (Metrohm Suppressor Module with high capacity)	

6.2844.000 MSM-LC Rotor A
Suppressor rotor for use with 2 mm columns.



6.2832.000 MSM Rotor A
Suppressor rotor for all IC instruments with MSM (Metrohm Suppressor Module)



6.2842.020 Adapter sleeve, MSM

Not only the MSM Rotor A but also the MSM-LC Rotor A can be installed in the rotor housing with the aid of the adapter sleeve.



2.941.0010 941 Eluent Production Module

The 941 Eluent Production Module enables the automatic production of an eluent. It enables continuous working without manual intervention and guarantees stable retention times. It can be combined with all Metrohm IC instruments and is monitored by MagIC Net.



Suppressor rotor for all IC instruments with MSM-HC drive (Metrohm Suppressor Module with high capacity) for cation suppression.



Professional chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Auto samplers, 800 Dosino, 771 Compact Interface, etc. The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Professional complies fully with FDA Regulation 21 CFR Part 11 as well as GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

1 license

The installation and documents are supplied on a USB flash drive.

Client/server-compatible chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Autosamplers, 800 Dosino, 771 Compact Interface, etc.

The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses. Data and method are centrally stored on a Server.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Multi complies fully with FDA Regulation 21 CFR Part 11 as well as with GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

Performance for client/server systems with servers outside the in-house company network has been improved.

Client/server version with 3 licenses
