



## 815 Robotic USB Sample Processor XL (1T/0P)

2.815.0030

Robotic USB Sample Processor XL with one workstation for the automatic processing of routine samples in series with large quantities and for complex sample preparation or parallel runs. Up to two pumps (membrane or peristaltic) and three dosing devices for Liquid Handling tasks can be connected.

Because of the multitude of application variants, rack, stirrer, titration head, robotic arm, Swing Head and sample vessels must be tailored to the application and ordered separately.

The control is "stand alone" using Touch Control. The following software products can be selected for the PC control: tiamo™ titration software, MagIC Net chromatography software, viva voltammetry software, or OMNIS.

Phạm vi giao hàng 2.815.0030

| Qt. | Order no. | Description |
|-----|-----------|-------------|
|-----|-----------|-------------|

1 PCS

1.815.0030

815 Robotic USB Sample Processor XL (1T/0P)

Robotic USB Sample Processor XL with one workstation for the automatic processing of routine samples in series with large quantities and for complex sample preparation or parallel runs. Up to two pumps (membrane or peristaltic) and three dosing devices for Liquid Handling tasks can be connected.

Because of the multitude of application variants, rack, stirrer, titration head, robotic arm, Swing Head and sample vessels must be tailored to the application and ordered separately.



2 PCS

6.1236.020

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm and O-ring.



1 PCS

6.1805.110

FEP tubing / M6 / 80 cm

With light and kink protection



1 PCS

6.2151.000

Cable USB A – mini-DIN 8-pin

Controller cable



1 PCS

6.2621.030

Hexagon key 4 mm



---

1 PCS

6.2621.070

5 mm hex key for IC Sample Processors



---

1 PCS

6.2621.140

Hexagon key 2.5 mm



---

3 PCS


6.2709.070

Guiding sleeve

For fixation of buret and dosing tips in titration heads and electrode holders



Safety shield for Sample Processors with 786 Swing Head

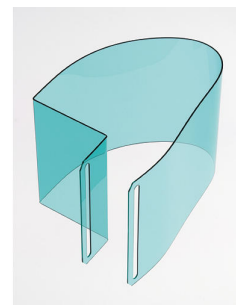

## Phụ kiện tùy chọn

| Order no.  | Description                                                                                                                                                                                                                                                                                     |                                                                                       |
|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------|
| 2.772.0120 | 772 Pump Unit – aspirate<br><br>Peristaltic pump for direct connection to Sample Processors. The "Aspiration" version contains extensive accessories for the aspiration of titrated-out solutions.                                                                                              |    |
| 2.772.0130 | 772 Pump Unit – rinse<br><br>Peristaltic pump for direct connection to Sample Processors. The "Rinse" version contains extensive accessories for rinsing.                                                                                                                                       |    |
| 2.786.0010 | 786 Swing Head with transfer head left<br><br>Additional equipment for Robotic Sample Processors XL for transferring / pipetting from small sample vessels into larger titration vessels on the rack or into an external titration cell; this can be located to the left of the work station.   |  |
| 2.786.0020 | 786 Swing Head with transfer head right<br><br>Additional equipment for Robotic Sample Processors XL for transferring / pipetting from small sample vessels into larger titration vessels on the rack or into an external titration cell; this can be located to the right of the work station. |  |

2.786.0030

### 786 Swing Head with titration swing arm, left- or right-swinging

Additional equipment for Robotic Sample Processors XL for direct titration in sample vessels of 75 mL and more.



2.786.0040

### 786 Swing Head

Additional equipment for Robotic Sample Processors XL without titration or transfer heads.



2.801.0010

### 801 Stirrer

Magnetic stirrer without stand for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimats and 780/781 pH Meters. With permanently attached cable for MSB (Metrohm Serial Bus).



2.802.0020

### 802 Stirrer for sample changer

Rod stirrer for sample changer and Sample Processor.

With 104 mm 6.1909.020 stirring propeller and permanently installed cable.



2.802.0040

#### 802 Stirrer for 804 Ti Stand

Rod stirrer including 6.1909.010 stirring propeller.



2.804.0040

#### 804 Ti Stand with stand

Titration stand and controller for 802 Rod Stirrer. The 804 Ti Stand together with the optional 802 Rod Stirrer provides an alternative to the magnetic stirrer. Ti Stand with base plate, support rod and electrode holder.



2.823.0020

#### 823 Membrane Pump Unit «aspire»

Membrane pump for direct connection to Sample Processors. The «Aspiration» version contains numerous accessories for aspirating off solutions.



2.823.0030

#### 823 Membrane Pump Unit «rinse»

Membrane pump for direct connection to Sample Processors. The «Rinsing» version contains numerous rinsing accessories.



2.843.0050

**843 Pump Station (membrane) - rinse/aspirate for Sample Processors**

The 843 Pump Station (membrane) has two built-in membrane pumps. These can be controlled directly from the Sample Processor via pump commands. The Rinse / Aspirate version is provided with all the accessories needed for automatically emptying the titration beaker and rinsing the titration equipment. The accessories are intended for use with the 814 and 815 Sample Processors.



2.843.0150

**843 Pump Station (peristaltic) - rinse/aspirate for Sample Processors**

The 843 Pump Station (peristaltic) has two built-in peristaltic pumps. These can be controlled directly from the Sample Processors via pump commands. The rinse/aspirate version is provided with all the accessories needed for automatically emptying the titration beaker and rinsing the titration equipment. The accessories are intended for use with the 814 and 815 Sample Processors.



2.849.0010

**849 Level Control**

Additional equipment for Sample Processors for monitoring the filling level of rinsing or waste canisters via Remote. Prevents pumps from running dry and/or canisters from overflowing and is suitable for use with aqueous solutions, solvents and suspensions.



6.1458.010

**Titration head, 6 x SGJ 14 and 3 x SGJ 9 openings.**

With 6xSGJ 14 and 3x SGJ 9 openings.



6.1458.020 Titration head / 4 x M10 connectors

With 4 M10 screw nipples, including stoppers.



---

6.1458.030 KF titration head / 4 x M10 connectors

With 4 M10 screw nipples, including stoppers.



---

6.1458.040 Titration Head, 3 x SGJ 14

Titration head for 120 mL beakers. With 3 openings SGJ 14 and 4 openings 6.4 mm, 2 dosing tips included. Suitable for rinsing and aspiration equipment.



---

6.1458.070 Titration head / 2 x SGJ 14

Titration head for 75 mL beakers. With 2 openings SGJ 14, 2 openings 6.4 mm, and 3 pre-mounted dosing tips. Suitable for rinsing and aspiration equipment.



6.1462.030      Robotic arm with transfer head for 786 Swing Head, left swinging  
Transfer robotic arm for liquid handling purposes



---

6.1462.040      Robotic arm with transfer head for 786 Swing Head, right swinging  
Transfer robotic arm for liquid handling purposes



---

6.1462.050      Robotic arm with titration head for 786 Swing Head, left or right swinging  
Titration head for Swing Head



---

6.1462.060      Robotic arm with holder for titration head, left swinging  
Titration head holder to 786 Swing Head



6.1462.070      Robotic arm with holder for titration head, right swinging  
Titration head holder to 786 Swing Head



---

6.1462.090      Robotic arm with Luer-Lock adapter for 786 Swing Head, right swinging  
Transfer robotic arm for liquid handling



---

6.1462.150      Robotic arm with titration head and beaker sensor for 786 Swing Head,  
left or right-swinging  
Titration head with beaker sensor for 786 Swing Head, 2 x M10, 1 x SGJ  
14.5, 2 x buret tip



---

6.1462.160      Robotic arm with holder for titration head and sensor for 786 Swing  
Head, left swinging  
Titration head holder with beaker sensor for 786 Swing Head



6.1462.170      Robotic arm with holder for titration head and sensor for 786 Swing Head, right swinging

Titration head holder with beaker sensor for 786 Swing Head



6.1462.260      Robotic arm with holder for titration head, left swinging, external

Titration head holder for 786 Swing Head with possibility of swinging to external positions



6.2001.070      Support stand for Sample Processors

Support stand for Sample Processors, for external titrations and dilutions (IC)



6.2041.350      Sample rack, 48 x 75 mL (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 48 sample beakers. The following sample beakers can be used: 6.1432.210, 6.1459.400.

Plastic: Polyvinyl chloride (PVC)



6.2041.400 Sample rack, 126 x 11 mL and 2 x 250 mL (PVC)

Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic" for 126 sample vessels and 2 sample beakers. The following sample vessels and sample beakers can be used: 6.2743.050, 6.1432.320, 6.1432.323, 6.1453.220, 6.1453.250, 6.2743.047, 6.2743.040.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.410 Sample rack, 141 x 11 mL and 1 x 300 mL (PVC)

Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic" for 141 sample vessels and 1 sample beaker. The following sample vessels and rinsing bottles can be used: 6.2743.050, 6.1608.080, 6.1627.000, 6.2743.040.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.430 Sample rack, 127 x 11 mL and 2 x 300 mL (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 127 sample vessels and 2 rinsing bottles. The following sample vessels and rinsing bottles can be used: 6.2743.050, 6.1608.080, 6.1627.000, 6.2743.040.

Plastic: Polyvinyl chloride (PVC)




---

6.2041.440 Sample rack, 148 x 11 mL and 3 x 300 mL (PVC)

Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, suitable for 148 sample vessels and 3 rinsing bottles. The following sample vessels and rinsing bottles can be used: 6.2743.050, 6.1608.080, 6.1627.000, 6.2743.040.


Plastic: Polyvinyl chloride (PVC)

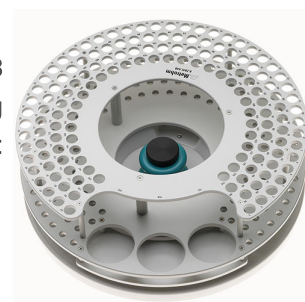

6.2041.450 Sample rack, 56 x 11 mL and 56 x 50 mL (PVC)

Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, suitable for 112 sample vessels. The following sample vessels can be used: 6.2743.050, 6.2747.010, 6.2743.040, 50 mL Sarstedt vials.

Plastic: Polyvinyl chloride (PVC)

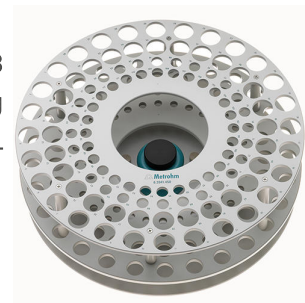

---

6.2041.710 Sample rack, 160 x 6 mL, (PVC)

Sample rack for 815 Robotic Sample Processor, 855 Robotic Titrosampler "Basic" suitable for 160 sample vessels. The following sample vessels can be used: 6.2419.000.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.800 Sample rack, 100 x 75 mL, (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic" suitable for 100 sample beakers. The following sample beakers can be used: 6.1432.210, 6.1459.400.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.810 Sample rack, 34 x 150 mL, (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic" suitable for 34 sample beakers. The following sample beakers can be used: 150 mL beaker (tall).

Plastic: Polyvinyl chloride (PVC)



6.2041.820      Sample rack, 28 x 250 mL, (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic" suitable for 28 sample beakers. The following sample beakers can be used: 6.14332.320, 6.1432.323, 6.1453.220, 6.1453.250.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.830      Sample rack, 28 x 200 mL, (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 28 sample beakers. The following sample beaker can be used: 6.1459.310.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.840      Sample rack, 59 x 120 mL, (PVC)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 59 sample beakers. The following sample beakers can be used: 6.1459.300, 6.1459.307.

Plastic: Polyvinyl chloride (PVC)



---

6.2041.850      Sample rack, 59 x 120 mL, (PP)

Sample rack for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 59 sample beakers. The following sample beakers can be used: 6.1459.300, 6.1459.307.

Plastic: Polyvinyl chloride (PVC)



6.2041.900 Sample rack insert, 54 x 75 mL, for water bath, (Al)

Sample rack insert for 815 Robotic USB Sample Processor XL, 855 Robotic Titrosampler "Basic", suitable for 54 sample beakers. The following sample beakers can be used: 6.1432.210, 6.1459.400.

Material: Aluminum (Al)



---

6.2058.020 Head extension 30 mm for Sample Processors

Head extension 30 mm for 789, 815 Sample Processors.



---

6.2058.030 Head extension 10 mm for Sample Processors

Head extension for use with the waterbath rack.



---

6.2058.110 Tower extension for Sample Processors

Height extension 130 mm for Sample Processors. Tower extension for COD applications.



6.2141.300

### Remote cable 843 Pump Station to Sample Processor

Remote cable for the direct connection of the 843 Pump Station to the pump connectors of the Sample Processors



6.2148.010

### Remote Box MSB

Additional remote interface for the connection of devices that can be controlled via remote lines. With permanently attached cable.



6.2711.070

### Drip pan for 789/815/855

Drip pan for 789, 815, 855 Robotic Sample Processors XL



6.2840.000

### Water bath for Robotic Sample Processors and Robotic Titrosampler

Water bath rack for Robotic Sample Processors and Titrosampler.



6.2841.000 External Polytron rinsing station for Robotic Sample Processors


The external rinsing station can be fitted on the side next to the Robotic Sample Processor tower, making it possible for external cleaning of the Polytron to be carried out during additional actions on the sample rack.



---

6.2841.010 External rinsing station for robotic arm with 6.1458.040 titration head

The external rinsing station can be fitted on the side next to the Robotic Sample Processor tower, making it possible for external cleaning of the titration equipment to be carried out during additional actions on the sample rack.

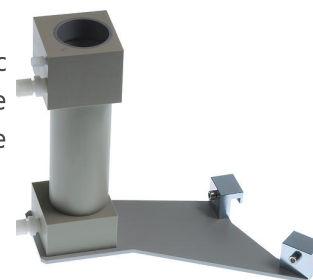

---

6.5619.000 Transfer equipment, 0.5 - 2 mL, for 789/815

Complete equipment for the transfer of solutions from 0.5 to 2.0 mL on the Robotic Sample Processor.



---

6.06000.000 OMNIS Software storage device

License-based PC program for controlling complex OMNIS analysis systems.

- Sample-oriented operating procedures with numerous templates
- Graphical method editor with numerous templates
- Clearly organized and professional sample lists
- Management of instruments, work systems, sensors, solutions
- Dialog languages: German, English, Chinese



6.06003.010 OMNIS Stand-Alone license

Enables stand-alone operation of the OMNIS software on a Windows™ computer.

Features:

- The license already includes one OMNIS instrument license.
- Must be activated via the Metrohm licensing portal.



---

6.2041.780 Sample rack 65 x 11 mL + 25 x 50 mL

Multi-row sample rack for 65 samples in 11 mL sample vessels (6.2743.050) and 25 samples in 50 mL sample vessels (6.02743.030)



---

6.2041.770 Sample rack 211 x 11 mL

Multi-row sample rack for 211 samples in 11 mL sample vessels (6.2743.050)

