



# 852 Titrando with Touch Control and generator electrode without diaphragm

2.852.0160

Coulometer including generator electrode without diaphragm, complete volumetric titration cell, and Touch Control.

Coulometry is the ideal method for water content determination in liquids, solids, and gases when it comes to water content determination in the trace range (10 µg to 10 mg absolute water). In addition, coulometry is an absolute method and thus no titer determination is necessary.

Recommended coulometry measuring range: 10 µg - 200 mg absolute water

The 852 Titrando controls not only coulometric, but also volumetric Karl Fischer titration.


For use with OMNIS Software, *tiamo* software, or Touch Control unit. Compliance with GMP/GLP and FDA regulations such as 21 CFR Part 11, if required.

Phạm vi giao hàng 2.852.0160

| Qt. | Order no. | Description |
|-----|-----------|-------------|
|-----|-----------|-------------|



Operating unit for the Titrandos, USB Sample Processors, 856 Conductivity Module, 867 pH Module and 846 Dosing Interface. Touch-sensitive, high-resolution color display, simple and intuitive operation, thanks to Favorites for direct method access. With integrated Ethernet interface for direct connection to the Internet and USB interface for connecting USB printers or a USB memory stick.

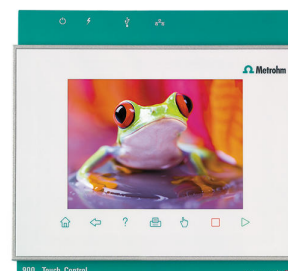

Dialog languages: German, English, Chinese, French, Spanish, Portuguese, Russian, Korean, Polish and Italian.

Indicator electrode for Karl Fischer titration which, thanks to its two fused platinum rings, is more robust and easier to clean than the conventional platinum pin electrode (6.0338.100).



Indicator electrode used for coulometric Karl Fischer titration.



1 PCS

6.0345.100

Generator electrode without diaphragm

Generator electrode without diaphragm for coulometric Karl Fischer titrations, SGJ 29/22. This generator electrode is well suited for most applications where the sample is highly soluble in alcohol.



1 PCS

6.1244.040

Set of seals

Set of seals for 6.1414.030. The set is comprised of one seal and various O-rings.



1 PCS

6.1403.030

KF adsorber tube for coulometer cell

For KF coulometers.



1 PCS

6.1403.040

KF Adsorber tube

Used with KF equipment. Including cover and O-ring.



1 PCS

6.1414.030

KF titration vessel lid with 5 openings

For volumetric KF titrations.



---

1 PCS

6.1415.220

Titration vessel / 20-90 mL



---

1 PCS

6.1415.250

Titration vessel / 50-150 mL



---

1 PCS

6.1437.000

SGJ stopper SGJ 14



1 PCS

6.1446.060

Stopper / B-14/15 / M10



---

2 PCS

6.1448.010

Septum 12 mm, 5 pieces

Set of 5 items.



---

2 PCS

6.1448.020

Septum 16 mm, 5 pieces

Set of 5 pcs.



---

1 PCS

6.1464.320

KF titration vessel / 80-250 mL / coulometric



2 PCS

6.1903.020

Stirring bar / 16 mm

Stirring Bar with magnetic core, PTFE coated, length 16 mm



---

3 PCS

6.1903.030

Stirring bar / 25 mm

Stirring bar with magnetic core, PTFE covering, length 25 mm.



---

1 PCS

6.2001.090

Base plate on both sides

Double stand for attaching 2 x 801, 803, 804 to a Titrando



---

2 PCS

6.2013.010

Clamping ring

For support rods with a diameter of 10 mm.



2 PCS

6.2016.030

Support rod / 250 mm



---

2 PCS

6.2043.005

Holding clip for bottles

Holding spring for reagent bottles in exchange units.



---

1 PCS

6.2047.020

Titration vessel holder

Titration vessel holder for coulometric cells. Used with 756, 831 KF Coulometer



---

2 PCS

6.2104.020

Electrode cable / 1 m / F

For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



1 PCS

6.2104.120

Electrode cable / 1 m / H

With plug H. For connecting generator electrode - KF Coulometer.



1 PCS

6.2412.000

Weighing boats

Weighing boat made of glass for the 6.1414.030 titration vessel lid or the 6.1455.31X, 6.1464.32X, 6.1465.320 titration vessel. For KF water determinations. Including protective tube.



1 PCS

6.2701.040

GL 18 screw cap with hole

Made of PBT. With a hole for the 6.1448.020 septum. For the 6.5405.000 cell without diaphragm and the 6.1464.32X and 6.1465.320 titration vessels.



3 PCS

6.2713.000

PTFE sleeve SGJ14

For fat-free ground-joint connection with the standard ground-joint SGJ 14.





1 PCS

6.2713.010

PTFE sleeve SGJ29

Sleeve for fat-free ground-joint connection with the standard ground-joint SGJ 29.



---

1 PCS

6.2713.020

PTFE sleeve SGJ19

For fat-free ground-joint connection with the standard ground-joint SGJ 19.



---

2 PCS

6.2730.010

Screw nipple

Made of PP. For 6.1448.010 Septum or for mounting 6.1403.040 Drying tube or 6.0338.100 Electrode in 6.1414.030 and 6.1414.040 KF Titration vessel lids. Used for KF titrations.



---

1 PCS

6.2730.020

Septum stopper

With E.301.0041 O-ring. For 6.1448.010 Septum and 6.1414.030 Titration vessel lid. Used for KF applications.



4 PCS

6.2730.030

Stopper

Complete with nipple and O-ring. For 6.1414.030 titration vessel lids. For KF instruments, VA measurement stands, VA Computrace and IC Dilution.



---

1 PCS

6.2738.000

Funnel

For the KF coulometer.



---

2 PCS

6.2816.030

Needle with Luer connector

Needle for filling in mercury into the Multi-Mode Electrode and for sample addition in Karl Fischer titration.



---

2 PCS

6.2816.090

Syringe 5 mL

5 mL syringe made of PP, with Luer connector.





## Phụ kiện tùy chọn

| Order no.  | Description                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                             |
|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 2.860.0010 | <p>860 KF Thermoprep</p> <p>The 860 KF Thermoprep is designed for thermal sample preparation in Karl Fischer titration.</p> <p>Numerous substances cannot be analysed by direct Karl Fischer titration as they are not soluble, react with the Karl Fischer reagent or release their water only very slowly or not until heated to high temperatures. Using tightly sealed sample vials, the sample is inserted in the oven. The samples can be subsequently analysed by any volumetric or coulometric KF titrator.</p> |
| 6.2148.020 | <p>RS-232/USB box</p> <p>Extends the Titrando, USB Sample Processors, Dosing Interface, Titrosampler by two RS-232 interfaces for connecting e.g. balances, PC, LIMS. Requires 6.2151.030 cable.</p>                                                                                                                                                                                                                                                                                                                    |
| 6.2103.140 | <p>Adapter black 2 mm plug / B socket 4 mm</p> <p>For connecting plug B (4 mm) to 2 mm socket.</p>                                                                                                                                                                                                                                                                                                                                                                                                                      |
| 6.2103.130 | <p>Adapter red 2 mm plug / 4 mm socket</p> <p>For connecting plug B (4 mm) to 2 mm socket.</p>                                                                                                                                                                                                                                                                                                                                                                                                                          |



6.2811.000

## Molecular sieve

Molecular sieve. Bottle containing 250 g. Pore size: 0.3 nm. Without moisture indicator. For Rancimats and Karl Fischer instruments.



---

6.2811.010

## Molecular sieve (US)

Molecular sieve. Bottle containing 250 g. Pore size: 0.3 nm. Without moisture indicator. For Rancimats and Karl Fischer instruments.

This product is intended for all users who are US residents. The material is rated hazardous in the US in accordance with OSHA Hazard Communication Standard (29 CFR 1910.1200).

