



870 KF Titrino plus and 860 KF Thermoprep

2.870.2010

This titration system is comprised of an 870 KF Titrino plus and an 860 KF Thermoprep. This means that samples unsuitable for direct Karl Fischer titration can be analyzed. With the aid of tightly sealed sample vials, the samples are moved into the oven of the 860 KF Thermoprep. The sample is heated up, the water it contains is vaporized and channeled into the titration cell. The water content of the sample is then determined with the 870 KF Titrino plus.

Phạm vi giao hàng 2.870.2010

Qt.	Order no.	Description
-----	-----------	-------------

The inexpensive 870 KF Titrino plus is a KF titrator for volumetric water content determination. With it, water contents from a few ppm to 100% can be determined reliably and precisely in solid, liquid and gaseous samples.

With its user interface, which is tailored to routine users, the 870 KF Titrino plus is so simple to operate that only brief orientation periods are required. It is the ideal titrator for routine determinations, also because of its robustness.

Titration vessel and electrode are included in the scope of delivery.



Compact titration stand with magnetic stirrer and built-in pump for solvent addition and aspiration of the titration vessel contents. Suitable for Karl Fischer titrations with Titrandos. MSB connection. Including comprehensive accessories, solvent bottle and waste bottle.



The 860 KF Thermoprep is designed for thermal sample preparation in Karl Fischer titration.

Numerous substances cannot be analysed by direct Karl Fischer titration as they are not soluble, react with the Karl Fischer reagent or release their water only very slowly or not until heated to high temperatures. Using tightly sealed sample vials, the sample is inserted in the oven. The samples can be subsequently analysed by any volumetric or coulometric KF titrator.



1 PCS

6.0338.100

Double Pt-wire electrode for volumetry

Indicator electrode used for volumetric Karl Fischer titration.



1 PCS

6.1244.040

Set of seals

Set of seals for 6.1414.030. The set is comprised of one seal and various O-rings.



1 PCS

6.1403.040

KF Adsorber tube

Used with KF equipment. Including cover and O-ring.



1 PCS

6.1414.030

KF titration vessel lid with 5 openings

For volumetric KF titrations.



1 PCS

6.1415.220

Titration vessel / 20-90 mL



1 PCS

6.1415.250

Titration vessel / 50-150 mL



2 PCS

6.1446.040

Threaded stopper / M6



1 PCS

6.1446.060

Stopper / B-14/15 / M10



1 PCS

6.1446.090

Stopper / B-14/15 / M8



1 PCS

6.1446.170

Heating tube stopper

for coulometric KF vessels in connection with heatable tubing connection



2 PCS

6.1448.010

Septum 12 mm, 5 pieces

Set of 5 items.




1 PCS

6.1448.057

Aluminum septum caps / 100 pieces

Including Al sealing.

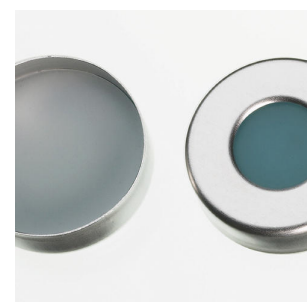

1 PCS

6.1543.110

Dosing tip / M8 thread

Transfer tip with M8 thread



1 PCS

6.1543.120

Aspiration tip / M8 thread

Aspiration tip if 6.1805.200 siphoning tubing is used



2 PCS

6.1602.105

Bottle attachment / GL 45 / green

For direct dosing from reagent bottles with GL 45 thread.



2 PCS

6.1602.145

Drying flask cap

Screw cap with GL 45 thread for the 6.1608.050 drying flask.



1 PCS

6.1608.023

Amber glass bottle / 1000 mL / GL 45

For Exchange Units. Bottle for auxiliary solutions.



1 PCS

6.1608.030

Round glass bottle / 1000 mL / GL 45



2 PCS

6.1608.050

Drying bottle / 100 mL / GL 45



1 PCS

6.1609.010

Adsorber tube with tubing nipple

Lid with tubing nipple.



1 PCS

6.1623.000

Bottle cap overflow protection

Bottle cap overflow protection for 1 waste bottle, prevents the waste bottle from overflowing when an aspiration pump is used



1 PCS

6.1801.120

PVC tubing / 4 mm / 6 mm / 2 m



1 PCS

6.1805.010

FEP tubing / M6 / 13 cm

With light and kink protection.



2 PCS

6.1805.080

FEP tubing / M6 / 25 cm

With light and kink protection.



1 PCS 6.1805.100 FEP tubing / M6 / 40 cm
With light and kink protection



2 PCS 6.1805.200 PTFE tubing / M8 / 0.5 m



1 PCS 6.1808.040 Thread adapter M6 outer / M8 inner
Outer thread M6, inner thread M8.



2 PCS 6.1808.050 Tubing adapter olive / M8 outer
1 M8 outer thread and 1 tubing olive. E.g. for thermostat jacket of
exchange units and stability measuring instruments.



1 PCS

6.1819.030

Tube

For M8 threaded holes. For 6.1602.1XX bottle attachments and 6.1608.0XX bottles



1 PCS

6.1819.050

Tube

For M8 threaded holes of 6.1602.10X Bottle attachments with 6.1608.030 Aspiration bottle. For 703 and 803 Ti Stand, 681 Pump Unit and 679 Rancimat



2 PCS

6.1821.040

Filter tube

Filter tube for 6.1608.050 Drying bottle. For Rancimats and Karl Fischer ovens.



1 PCS

6.1830.030

Heating tubing

Heating tubing with M6 thread for 860 KF Thermoprep, 885 Compact Oven Sample Processor and 874 USB Oven Sample Processor.



2 PCS

6.1903.020

Stirring bar / 16 mm

Stirring Bar with magnetic core, PTFE coated, length 16 mm



2 PCS

6.1903.030

Stirring bar / 25 mm

Stirring bar with magnetic core, PTFE covering, length 25 mm.




1 PCS

6.2001.060

Support stand without support rod

To attach an 801, 804, 803 Stirrer to a Titrando, Titrino plus, Dosimat plus or 856 and 867.

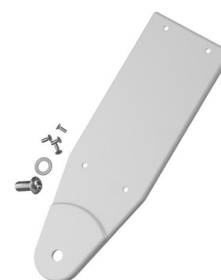

1 PCS

6.2013.010

Clamping ring

For support rods with a diameter of 10 mm.



1 PCS

6.2016.070

Support rod / 400 mm



1 PCS

6.2023.020

Clip for SGJ 14/15

Clip for SGJ 14/16



1 PCS

6.2049.040

Needle holder with Luer lock

Needle holder (58 mm) with Luer lock for the Karl Fischer ovens 860, 874 and 885.



1 PCS

6.2060.010

Cable holder

Used with 860 KF Thermoprep.



1 PCS

6.2063.010

Sample holder for 6.2419.007 Sample vials (6 mL)



1 PCS

6.2104.020

Electrode cable / 1 m / F

For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



1 PCS

6.2147.000

Numerical USB keypad

Numerical USB keypad for Titrino plus, Dosimat plus, 862 Compact Titrosampler, 869 Compact Sample Changer, 863 Compact Autosampler and 899 Coulometer. With two integrated USB connectors.



1 PCS

6.2151.100

Adapter USB MINI (OTG) - USB A

For connecting USB instruments.



1 PCS

6.2412.000

Weighing boats

Weighing boat made of glass for the 6.1414.030 titration vessel lid or the 6.1455.31X, 6.1464.32X, 6.1465.320 titration vessel. For KF water determinations. Including protective tube.



1 PCS

6.2419.007

Sample glass, 6 mL, 100 pieces

Sample glass for up to 6 mL sample, suitable for 6.1448.057 septum, 100 pieces



1 PCS

6.2621.070

5 mm hex key for IC Sample Processors



1 PCS

6.2621.110

Crimping tongs

Crimping tongs for sealing 6 mL sample vials (6.2419.000, 6.2419.007).



1 PCS

6.2621.120

Hexagon key 1.5 mm



1 PCS

6.2621.130

Hexagon key 2 mm

2 mm.



1 PCS

6.2627.000

Oven insert

Used with 860 KF Thermoprep.



1 PCS

6.2723.310

Protective cover for Titrino plus

Protective cover for display and keypad for all instruments of the Titrino plus family.

Suitable for: 848 Titrino plus, 862 Compact Titrosampler, 863 Compact Autosampler, 865 Dosimat plus, 869 Compact Sample Changer, 870 KF Titrino plus, 876 Dosimat plus, 877 Titrino plus, 899 Coulometer.



1 PCS

6.2724.010

Dust filter

Dust filter for Rancimats and Karl Fischer ovens.



2 PCS

6.2730.010

Screw nipple

Made of PP. For 6.1448.010 Septum or for mounting 6.1403.040 Drying tube or 6.0338.100 Electrode in 6.1414.030 and 6.1414.040 KF Titration vessel lids. Used for KF titrations.



1 PCS

6.2730.020

Septum stopper

With E.301.0041 O-ring. For 6.1448.010 Septum and 6.1414.030 Titration vessel lid. Used for KF applications.



3 PCS

6.2730.030

Stopper

Complete with nipple and O-ring. For 6.1414.030 titration vessel lids. For KF instruments, VA measurement stands, VA Computrace and IC Dilution.



1 PCS

6.2739.000

Wrench

For tightening connectors



1 PCS

6.2739.010

Wrench

For exchange units.



1 PCS

6.2811.000

Molecular sieve

Molecular sieve. Bottle containing 250 g. Pore size: 0.3 nm. Without moisture indicator. For Rancimats and Karl Fischer instruments.



1 PCS

6.2816.070

Septum piercing needle

Used with Karl Fischer ovens.



1 PCS

6.2816.080

Outlet needle

Used with KF Thermoprep and Oven Sample Processor.



1 PCS

6.3026.210

Exchange Unit 10 mL

Exchange unit with integrated data chip with 10 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



1 PCS

6.2048.030

Guidance rod

Guidance rod for 860 KF Thermoprep.



Phụ kiện tùy chọn

Order no.	Description
2.142.0100	Custom Q3X thermal printer Compact printer with USB interface for • 900 Touch Control • 915 KF Ti-Touch • 916 Ti-Touch • 917 Coulometer • 877 / 848 Titrino plus • 865 / 876 Dosimat plus • 91X Meter (cable 6.2151.140) • Eco Dosimat / Titrator • 862 Compact Titrosampler • 870 KF Titrino plus • 899 Coulometer Paper width 60 mm (40 characters). Including 6.2151.120 USB cable.

