


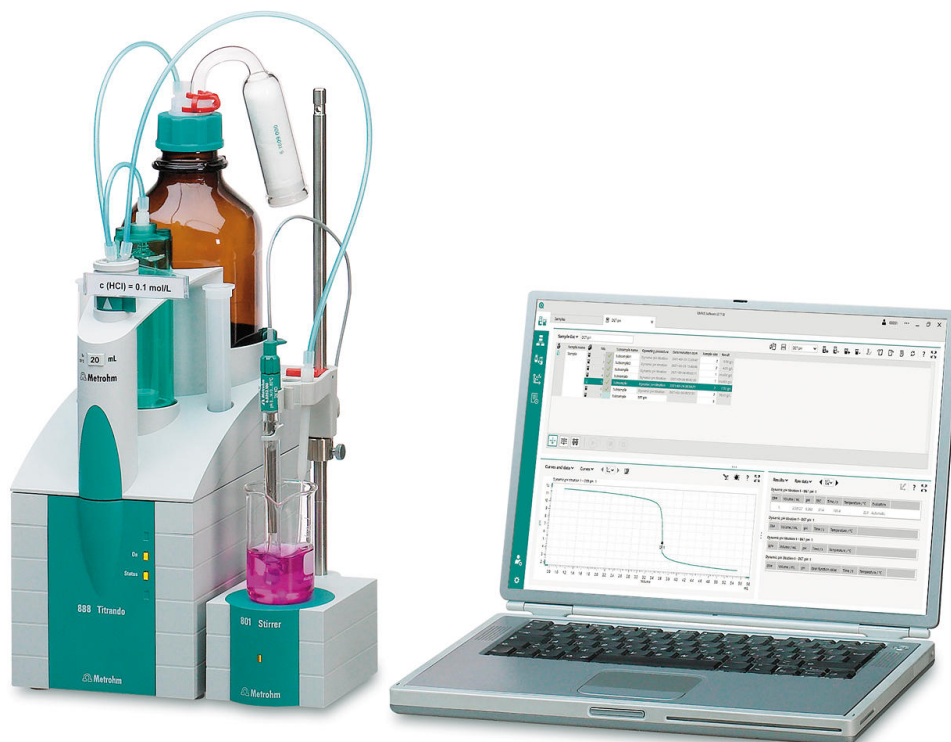

888 Titrando with OMNIS SW

2.888.0310

High-end titrator with built-in buret drive

- Includes OMNIS standalone license, 801 Stirrer magnetic stirrer, 20 mL exchange unit and combined "Ecotrode plus" pH glass electrode
- Dynamic (DET), monotone (MET), end point titration (SET), Karl-Fischer Titration (KFT)
- Four MSB connectors
- Optional second buret for titration
- A galvanically isolated measuring interface (also for intelligent electrodes "iTrodes")
- USB connector

Phạm vi giao hàng 2.888.0310

Qt.	Order no.	Description
-----	-----------	-------------

1 PCS

1.801.0010

801 Stirrer


Magnetic stirrer without stand for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimats and 780/781 pH Meters. With permanently attached cable for MSB (Metrohm Serial Bus).



1 PCS

1.888.0010

Titrando

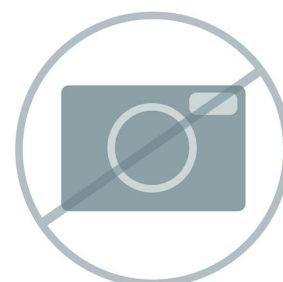

1 PCS

6.0262.100

Ecotrode Plus

Combined pH electrode for aqueous acid/base titration. The fixed ground-joint diaphragm is insensitive to contamination. Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.



1 PCS

6.1903.020

Stirring bar / 16 mm

Stirring Bar with magnetic core, PTFE coated, length 16 mm



1 PCS

6.1903.030

Stirring bar / 25 mm

Stirring bar with magnetic core, PTFE covering, length 25 mm.




1 PCS

6.2001.060

Support stand without support rod

To attach an 801, 804, 803 Stirrer to a Titrando, Titrino plus, Dosimat plus or 856 and 867.

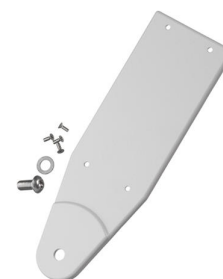

1 PCS

6.2013.010

Clamping ring

For support rods with a diameter of 10 mm.



1 PCS

6.2016.070

Support rod / 400 mm



1 PCS

6.2021.020

Electrode holder

Electrode holder for 4 electrodes and 2 buret tips



1 PCS

6.2103.130

Adapter red 2 mm plug / 4 mm socket

For connecting plug B (4 mm) to 2 mm socket.



1 PCS

6.2103.140

Adapter black 2 mm plug / B socket 4 mm

For connecting plug B (4 mm) to 2 mm socket.



1 PCS

6.2104.020

Electrode cable / 1 m / F

For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



1 PCS

6.2151.000

Cable USB A – mini-DIN 8-pin

Controller cable



1 PCS

6.2621.070

5 mm hex key for IC Sample Processors



1 PCS

6.2621.130

Hexagon key 2 mm

2 mm.



1 PCS

6.2739.010

Wrench

For exchange units.



Exchange unit with integrated data chip with 20 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



Enables stand-alone operation of the OMNIS software on a Windows™ computer.

Features:

- The license already includes one OMNIS instrument license.
- Must be activated via the Metrohm licensing portal.



License-based PC program for controlling complex OMNIS analysis systems.

- Sample-oriented operating procedures with numerous templates
- Graphical method editor with numerous templates
- Clearly organized and professional sample lists
- Management of instruments, work systems, sensors, solutions
- Dialog languages: German, English, Chinese



Phụ kiện tùy chọn

Order no.	Description	
6.5609.000	KF equipment Karl Fischer equipment for volumetry, complete set of accessories for Karl Fischer water determination	
6.3026.110	Exchange Unit 1 mL Exchange Unit with integrated data chip with 1 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.	
6.3026.150	Exchange Unit 5 mL Exchange Unit with integrated data chip with 5 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.	
6.3026.210	Exchange Unit 10 mL Exchange unit with integrated data chip with 10 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.	

6.3026.250

Exchange unit 50 mL

Exchange unit with integrated data chip with 50 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion tip and standard amber glass reagent bottle.



2.800.0020

800 Dosino

Drive with write/read hardware for intelligent dosing units. With permanently attached cable (length 0.65 m).



2.805.0010

805 Dosimat

Dosing device for the Titrando and Sample Processors with read/write hardware for intelligent Exchange Units. With permanently attached cable. Without Exchange Unit.



2.804.0040

804 Ti Stand with stand

Titration stand and controller for 802 Rod Stirrer. The 804 Ti Stand together with the optional 802 Rod Stirrer provides an alternative to the magnetic stirrer. Ti Stand with base plate, support rod and electrode holder.



2.854.0010

854 iConnect

The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).



6.3032.120

807 Dosing Unit 2 mL

807 Dosing Unit with integrated data chip with 2 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.150

807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.210

807 Dosing Unit 10 mL

807 Dosing Unit with integrated data chip with 10 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.220

807 Dosing Unit 20 mL

807 Dosing Unit with integrated data chip with 20 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.250

807 Dosing Unit 50 mL

807 Dosing Unit with integrated data chip with 50 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



2.803.0010

803 TI Stand with stirrer and pump

Compact titration stand with magnetic stirrer and built-in pump for solvent addition and aspiration of the titration vessel contents. Suitable for Karl Fischer titrations with Titrandos. MSB connection. Including comprehensive accessories, solvent bottle and waste bottle.

