



914 pH/Conductometer with accessories case

2.914.0120

Portable two-channel pH/conductivity measuring instrument for parallel measurement of pH/mV /conductivity/TDS/salinity and temperature. You will be optimally equipped for measurements outdoors with this battery-operated measuring instrument.

- Parallel measurement of pH value and conductivity
- Analog pH measuring input for Metrohm standard pH electrodes
- Analog conductivity measuring input for the 4-conductor conductivity measuring cells from Metrohm
- Robust, water-tight, and dust-tight housing (IP67) for tough outdoor and laboratory use
- LCD color display with background illumination making results easy to read
- USB interface for simple data export to PC or printer
- Large internal memory (10,000 data sets)
- Pin-protected User mode and Expert mode, prevents unwanted parameter changes
- GLP-compliant printout and data export with User ID and timestamp
- Robust transport case for routine measurements in the field

Phạm vi giao hàng 2.914.0120

Qt.	Order no.	Description
1 PCS	1.914.0020	914 pH/Conductometer
Two-channel pH/conductivity meter for routine use in the laboratory and on the road – version for connection of conventional pH electrodes. pH (or mV) and conductivity (or TDS, salinity) and temperatures can be measured in parallel and output to a large colour display with the 914 pH /Conductometer. Important information such as charge state, user, IDs can be clearly seen at a glance. A PIN-protected expert mode protects against unintentional changes of different parameters. The meter is furnished with an accumulator for mobile use that can be charged practically anywhere. It naturally also satisfies the requirements of IP67. A stand plate allows the mobile meter to be easily converted into a laboratory meter and vice versa. Very large measured value memory (10,000 data sets) and USB interface (GLP-compliant printout or data export with optional management of the data in tiBase) offer professional data handling.		
		
1 PCS	6.0228.030	Solitrade with Pt1000 (IP67, fixed cable 1.2 m)
Combined pH electrode with integrated temperature sensor (Pt1000) and fixed cable conforming to IP67 (1.2 m). This electrode is suitable for an introduction to GLP-compliant pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane. Additionally, this electrode offers a waterproof plug in accordance with IP67 for mobile use with Metrohm pH meters. Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.		
		

1 PCS

6.0917.080

Conductivity measuring cell $c = 0.5 \text{ cm}^{-1}$ with Pt1000
(fixed cable)

4-wire conductivity measuring cell with cell constant $c = 0.5 \text{ cm}^{-1}$ (guide value), with integrated Pt1000 temperature sensor and fixed cable for connecting to 912/914 Meters.

Thanks to the robust/break-proof plastic shaft made of PEEK, this sensor is mechanically very resistant. The sensor is suitable for measurements of medium conductivities ($15 \mu\text{S/cm}$ to 250 mS/cm), e.g., in:

- drinking water
- surface water
- wastewater



2 PCS

6.1446.000

SGJ stopper / B-14/(15)



2 PCS

6.1613.010

Bottle / 25 mL



1 PCS

6.1613.020

Bottle / 25 mL / label pH 7 (empty)



1 PCS

6.1613.030

Bottle / 25 mL / label pH 4 (empty)



1 PCS

6.1614.000

Wash bottle / 250 mL



1 PCS

6.2008.060

Holder for electrode storage vessel

Practical holder for fastening the electrode storage vessels to the 912, 913 or 914 pH/Conductometers.



1 PCS

6.2050.010

Carrying strap for 912/913/914

Carrying strap for 912/913/914 meters



1 PCS

6.2151.100

Adapter USB MINI (OTG) - USB A

For connecting USB instruments.



1 PCS

6.2151.110

Metrohm USB Mini B cable (OTG) - USB A, 1.8 m

For connecting USB instruments.



1 PCS

6.2166.100

USB power supply unit 5.25 V / 1.53 A

USB power supply unit for 912 / 913 / 914

Efficiency Level VI



1 PCS

6.2307.230

Buffer solutions pH 4, 7 and 9

Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL



1 PCS

6.2308.050

Electrolyte 3 mol/L KCl (50 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)

1 PCS

6.2324.110

Conductivity standard 100 $\mu\text{S/cm}$, 5 x 30 mL

Conductivity standard for calibration of conductivity measuring cells with cell constant = 0.1/cm.



1 PCS

6.2716.060

Case for 912 / 913 / 914

Case for the 912 / 913 / 914 pH/Conductometers.



1 PCS

6.2717.000

PP beaker, 100 mL

Beaker made of PP, 100 mL.



Phụ kiện tùy chọn

Order no.	Description
6.2001.130	Stand plate for 912/913/914 Stand plate for converting a mobile 912/913/914 pH/Conductometer into a laboratory meter. 
6.2166.500	12 V USB adapter for 912 / 913 / 914 pH/Conductometer 12 V USB adapter for 912 / 913 / 914 pH/Conductometer. 
6.0918.040	Conductivity measuring cell $c = 0.1 \text{ cm}^{-1}$ with Pt1000 (fixed cable) Conductivity measuring cell made of stainless steel with cell constant $c = 0.1 \text{ cm}^{-1}$ (guide value), with integrated Pt1000 temperature sensor and fixed cable for connecting to 912/914 Meters. This sensor is suitable for measurements of low conductivities ($0 \text{ }\mu\text{S/cm}$ to $300 \text{ }\mu\text{S/cm}$) in, e.g. deion. water. 
6.2325.000	pHit kit Maintenance kit for pH electrodes The kit contains: • 50 mL cleaning solution • 50 mL 3M KCl solution • 50 mL storage solution • 2 Storage vessels • Instructions for use 

6.0919.140

Conductivity measuring cell $c = 1.6 \text{ cm}^{-1}$ with Pt1000 (fixed cable)

3-ring conductivity measuring cell with cell constant $c = 1.6 \text{ cm}^{-1}$, with integrated Pt1000 temperature sensor and fixed cable for connecting to 912/914 Meters.

This sensor is suitable for measurements of high conductivities (0.1 to 1000 mS/cm), e.g., in:

- sea water
- flush water
- physiological solutions



6.2151.140

Y cable USB A St - USB B St - Mini B St

Y cable for connection of a USB printer to the pH/Conductometers 912 / 913 / 914. This cable allows a printer and the power supply unit to be connected to the measuring instrument at the same time.



2.142.0100

Custom Q3X thermal printer

Compact printer with USB interface for

- 900 Touch Control
- 915 KF Ti-Touch
- 916 Ti-Touch
- 917 Coulometer
- 877 / 848 Titrino plus
- 865 / 876 Dosimat plus
- 91X Meter (cable 6.2151.140)
- Eco Dosimat / Titrator
- 862 Compact Titrosampler
- 870 KF Titrino plus
- 899 Coulometer



Paper width 60 mm (40 characters). Including 6.2151.120 USB cable.


Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m, diameter banana plug 2 mm). This electrode is particularly suitable

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

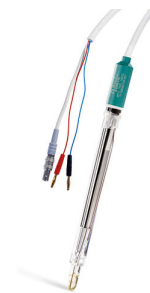

Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m). This electrode is suitable for routine pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.



Maintenance-free combined pH electrode (gel electrolyte) for piercing measurements of all types (e.g., with cheese, meat, dough) with integrated Pt1000 temperature sensor. The electrode is stored in saturated potassium chloride solution $c(\text{KCl}) = \text{sat.}$ (6.2308.000) and is not suitable for low-ion solutions. The ageing indicator gives early indication of when the electrode needs to be replaced.



6.0224.100

Biotrode

Combined pH electrode for measurements in very small sample volumes (>50 µL) and biological samples.

Idrolyte (6.2308.040) is used as reference electrolyte and storage solution.



6.2104.600

Electrode cable for plug in head U/plug F, 2x2 mm B, 1m

For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2313.000

Electrolyte 3 mol/L KCl (1000 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)

