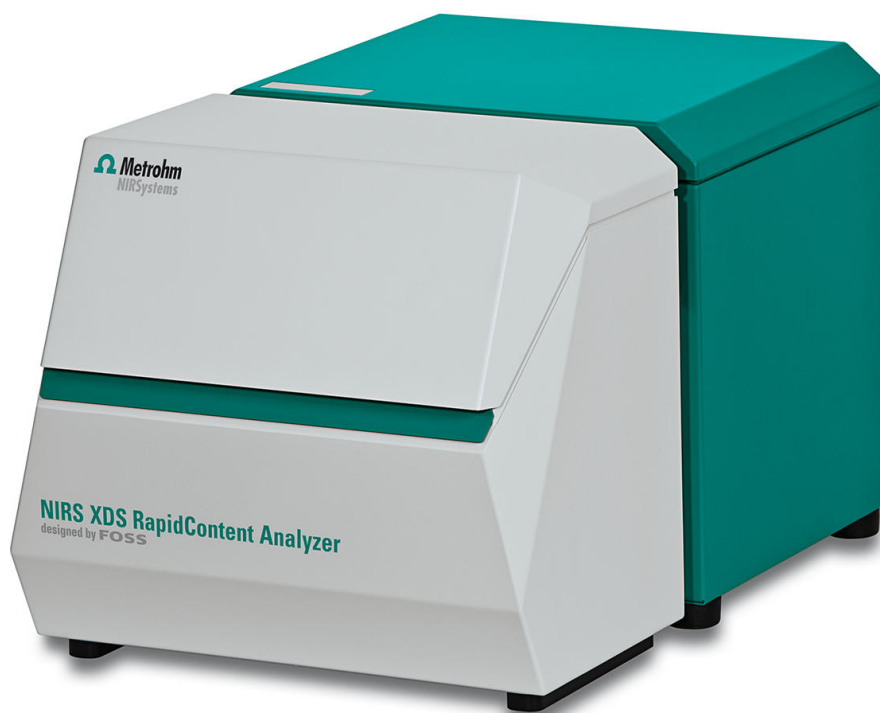




NIRS XDS RapidContent Analyzer – Solids

2.921.1120

Rapid, nondestructive analyses of solid and liquid substances

The NIRS XDS RapidContent Analyzer – Solids enables rapid, nondestructive analyses of solid powders, coarser granulates, pellets or flakes.

Composition or material identification testing is performed either in the laboratory or at-line, on samples contained in their original packaging, bags or vials. An optional variable spot-size feature enables the sample illumination to be adjusted based on the physical properties of your samples.

The measurement of liquid samples succeeds through the addition of the NIRS Transflection Kit for liquid samples (order number: 6.7400.010).

The sample is centered with the aid of the NIRS XDS iris (order number: 6.7425.000), thus enhancing the reproducibility of the analysis.

Thanks to the integrated mobile sample slide in X direction, the NIRS XDS RapidContent Analyzer enables the investigation of larger-sized nonhomogenous sample quantities.

Phạm vi giao hàng 2.921.1120

Qt.	Order no.	Description
-----	-----------	-------------

The result of 50 years experience in the area of spectroscopy was incorporated in the development of this NIRSystems spectrometer. The unique technology improves the signal-noise ratio at high resolution and enhances accuracy and precision. The patented monochromator and the modular design of the measuring module system make NIRSystems instruments indispensable for the spectral measurement of every imaginable sample type in the visible and near-infrared wavelength range (400 - 2,500 nm). The measuring modules can be replaced without difficulty, even during operation. This means that the user profits from the great flexibility in all areas of application.



Rapid, nondestructive analyses of solid and liquid substances

The NIRS XDS RapidContent Analyzer enables rapid, nondestructive analyses of solid and liquid formulations and substances. The sample is centered with the aid of an iris, thus enhancing the reproducibility of the analysis. A gold diffuse reflector enables the measurement of liquid samples. The Solids Module of the NIRS XDS RapidContent Analyzer extends the area of application to include practically every form of solid samples, whether it be fine powders, coarse granulates, pellets or flakes.



Accessory kit for NIRS XDS laboratory monochromators



1 PCS

6.7425.000

NIRS XDS iris

Iris for focusing mini sample vessels, reflection standard, various vials or a transfection vessel in the NIRS XDS RCA Analyzer, the NIRS XDS MasterLab Analyzer and the NIRS XDS MultiVial Analyzer



1 PCS

6.7402.040

NIRS XDS coarse granular sample cell

NIRS XDS coarse granular sample cell for the reflection measurement of samples with different surface properties.



1 PCS

6.7410.000

NIRS XDS tray for reflection standard

A reflection standard is required for the calibration of the monochromator with the sample module.

A holder in which the standard is positioned is used to fasten the calibration standard in place.

A reflection standard and its associated holder are required for:

- NIRS XDS MultiVial Analyzer (order number: 29211210)
- NIRS XDS MasterLab Analyzer (order number: 2.921.1310)



Phụ kiện tùy chọn

Order no.	Description	
6.7402.020	NIRS disposable backs for 6.7402.030, 100 pcs. The disposable backs are used for the secure sealing of the mini sample cups. The sample cups are used in the following: • NIRS DS2500 Analyzer (order number: 2.922.0010) • NIRS XDS MasterLab Analyzer (order number: 2.921.1310) • NIRS XDS RCA Analyzer (order number: 2.921.1110) • NIRS XDS MultiVial Analyzer (order number: 2.921.1210) The package is comprised of 100 disposable backs.	
6.7402.030	NIRS mini sample cups, 10 pcs. including 100 disposable backs Mini sample cup for the spectral recording of powders and granulates in reflection. The sample cup can be sealed with disposable backs in order to avoid sample loss and for the uniform distribution of the powder or granulate in the sample cup. The sample cups are used with the following instruments: • NIRS DS2500 Analyzer (order number: 2.922.0010) • NIRS XDS MasterLab Analyzer (order number: 2.921.1310) • NIRS XDS MultiVial Analyzer (order number: 2.921.1120) • NIRS XDS RapidContent Analyzer (order number: 2.921.1110) • NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)	
6.7430.010	NIRS XDS lab spare parts kit NIRS XDS spare parts kit of the lab instrument is comprised of two spare lamps, two fan filters and four fuses.	

6.7430.000

NIRS XDS lab spare lamp

Spare lamp with 12 V and 50 W for the monochromator of the NIRS XDS laboratory instruments



6.7450.000

NIRS reflection standard, set of 2

Set of wavelength standard and 80% reflection standard for the calibration of an NIRS XDS MasterLab Analyzer, NIRS XDS RCA Analyzer, NIRS XDS MultiVial Analyzer, NIRS XDS SmartProbe Analyzer, NIRS XDS Interactance OptiProbe Analyzer and NIRS DS2500 Analyzer.



6.7450.010

NIRS reflection standard, set of 7

Set comprised of six reflection standards and one wavelength standard for the calibration of an NIRS XDS MasterLab Analyzer, NIRS XDS RCA Analyzer, NIRS XDS MultiVial Analyzer, NIRS XDS SmartProbe Analyzer, NIRS XDS Interactance OptiProbe Analyzer and NIRS DS2500 Analyzer in regulated environments.



6.7400.010

NIRS Liquid Sample Kit Transflection

The NIRS Liquid Sample Kit Transflection supplements the following instruments with the function of the transflection measurement of liquids:

- NIRS DS2500 Analyzer (order number: 2.922.0010)
- NIRS XDS MasterLab Analyzer (order number: 2.921.1310)
- NIRS XDS MultiVial Analyzer (order number: 2.921.1120)
- NIRS XDS RapidContent Analyzer (order number: 2.921.1110)
- NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)



Vision Air - Universal spectroscopy software.

Vision Air Complete is a modern and simple-to-operate software solution for use in a regulated environment.



Overview of the advantages of Vision Air:

- Individual software applications with adapted user interfaces ensure intuitive and simple operation
- Simple creation and maintenance of operating procedures
- SQL database for secure and simple data management

The Vision Air Complete version (66072208) includes all applications for quality assurance using Vis-NIR spectroscopy:

- Application for instrument and data management
- Application for method development
- Application for routine analysis

Additional Vision Air Complete solutions:

- 66072207 (Vision Air Network Complete)
- 66072209 (Vision Air Pharma Complete)
- 66072210 (Vision Air Pharma Network Complete)

Vision Air - Universal spectroscopy software.

Vision Air Network is a modern and simple-to-operate software solution. The network version enables the structuring of a global server-client network with a central database.



Overview of the advantages of Vision Air:

- Individual software applications with adapted user interfaces ensure intuitive and simple operation
- Centralized data and instrument management
- Simple configuration and updating of several devices in just a few steps

The Vision Air Network Complete version (66072207) contains all applications for quality assurance using Vis-NIR spectroscopy:

- Application for a global instrument and data management
- Application for method development
- Application for routine analysis

Vision Air - Universal spectroscopy software.

Vision Air Pharma Complete is a modern and simple-to-operate software solution for use in regulated environments.



Overview of the advantages of Vision Air:

- Individual software applications with adapted user interfaces ensure intuitive and simple operation
- Simple creation and maintenance of operating procedures
- SQL database for safe and simple data management
- In compliance with all requirements pursuant to Title 21 CFR Part

11

The Vision Air Pharma Complete version (66072209) contains all applications for quality assurance using Vis-NIR spectroscopy:

- Application for instrument and data management
- Application for method development
- Application for routine analysis

Vision Air - Universal spectroscopy software.

Vision Air Pharma Network Complete is a modern and simple-to-operate software solution. The network version enables the setup of a global server-client network with a central database.



Overview of the advantages of Vision Air:

- Individual software applications with adapted user interfaces ensure intuitive and simple operation
- Centralized data and instrument management
- Simple configuration and updating of several instruments in just a few steps
- In compliance with all requirements pursuant to Title 21 CFR Part 11

The Vision Air Pharma Network Complete version (66072210) contains all applications for quality assurance using Vis-NIR spectroscopy:

- Application for a global instrument and data management
- Application for method development
- Application for routine analysis

6.7401.000

NIRS transfection vessel, optically flat

Optically flat transfection vessel for the spectral measurement of liquids. Can be used in combination with the following instruments:

- NIRS DS2500 Analyzer (order number: 2.922.0010)
- NIRS XDS MasterLab Analyzer (order number: 2.921.1310)
- NIRS XDS MultiVial Analyzer (order number: 2.921.1120)
- NIRS XDS RapidContent Analyzer (order number: 2.921.1110)
- NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)



6.7420.000

NIRS gold diffuse reflector, 1 mm total pathlength

Gold diffuse reflector for the transfection measurement of liquids. Can be used in combination with the following instruments:

- NIRS DS2500 Analyzer (order number: 2.922.0010)
- NIRS XDS MasterLab Analyzer (order number: 2.921.1310)
- NIRS XDS MultiVial Analyzer (order number: 2.921.1120)
- NIRS XDS RapidContent Analyzer (order number: 2.921.1110)
- NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)



6.7420.010

NIRS gold diffuse reflector, 2 mm total pathlength

Gold diffuse reflector for the transfection measurement of liquids. Can be used in combination with the following instruments:

- NIRS DS2500 Analyzer (order number: 2.922.0010)
- NIRS XDS MasterLab Analyzer (order number: 2.921.1310)
- NIRS XDS MultiVial Analyzer (order number: 2.921.1120)
- NIRS XDS RapidContent Analyzer (order number: 2.921.1110)
- NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)



6.7420.020

NIRS gold diffuse reflector, 4 mm total pathlength

Gold diffuse reflector for the transfection measurement of liquids. Can be used in combination with the following instruments:

NIRS DS2500 Analyzer (order number: 2.922.0010)
NIRS XDS MasterLab Analyzer (order number: 2.921.1310)
NIRS XDS MultiVial Analyzer (order number: 2.921.1120)
NIRS XDS RapidContent Analyzer (order number: 2.921.1110)
NIRS XDS RapidContent Analyzer - Solids (order number: 2.921.1210)



6.7430.060

Micro Fiber Cleaning Cloth

Cleaning cloth for monitors and optics.



6.7402.110

Disposable vials, 15 mm, reflection

123 lockable disposable glass vials with a diameter of 15 mm for analyses of solids in reflection. Suitable for NIRS RapidContent (Solid) Analyzers of the XDS, DS2500 and OMNIS product families.



6.7402.120

Disposable vials, 19 mm, reflection

225 lockable disposable glass vials with a diameter of 19 mm for analyses of solids in reflection. Suitable for NIRS RapidContent (Solid) Analyzers of the XDS, DS2500 and OMNIS product families.



6.7402.130

Disposable vials, 22 mm, reflection

172 lockable disposable glass vials with a diameter of 22 mm for analyses of solids in reflection. Suitable for NIRS RapidContent (Solid) Analyzers of the XDS, DS2500 and OMNIS product families.



6.7402.140

Disposable vials, 28 mm, reflection

216 lockable disposable glass vials with a diameter of 28 mm for analyses of solids in reflection. Suitable for the following analyzers:

- NIRS DS2500 Analyzer
- NIRS XDS RapidContent (Solid) Analyzer
- NIRS XDS MultiVial Analyzer
- NIRS XDS MasterLab Analyzer



6.7402.220

Stopper for 6.7402.130

172 stoppers for 22 mm Reflection disposable vials.



6.7402.230

Stopper for 6.7402.140

216 stoppers for 28 mm Reflection disposable vials.



6.7402.200

Stopper for 6.7402.110

123 stoppers for 15 mm Reflection disposable vials.



6.7402.210

Stopper for 6.7402.120

225 stoppers for 19 mm Reflection disposable vials.

