



MIRA P Flex

The MIRA P Flex package allows users to customize MIRA P to meet their needs. The Flex package includes the basic components needed to operate the MIRA P without sampling attachments. At least one sampling attachment is required for operation. The MIRA P Flex package includes the USP library, calibration/verification accessories, and a USB cable. Class 3B operation.

Phạm vi giao hàng

Qt.	Order no.	Description
1 PCS	6.06071.040	Calibrate/Verify Accessory (CVA)

50:50 (v:v) Toluene/Acetonitrile ASTM standard for calibration of the wave number axis of MIRA M-3/P spectrometers and polystyrene for wave number verification according to EP 2.2.48.



The Metrohm Instant Raman Analyzer (MIRA) P is a high-performance, handheld Raman spectrometer used for rapid, nondestructive determination and verification of different material types, such as Pharmaceutical APIs and excipients. Despite the small size of the instrument, the Mira P has a ruggedized design and features a high-efficiency spectrograph design equipped with our unique Orbital-Raster-Scan (ORS) technology. The Mira P is fully compliant with FDA 21 CFR Part 11 regulations.



The Advanced Package includes an attachment lens for analyzing materials directly or through containers (laser class 3b), as well as a vial holder attachment for analyzing samples contained in glass vials (laser class 1).

Phụ kiện tùy chọn

Order no.	Description	
6.02707.020	MIRA DS Robust Carrying Case	
	Robust carrying case for MIRA DS and all accessories.	
6.02707.400	MIRA P Protective Case	
	Protective case for the Mira P and all accessory parts.	
60607310	MIRA P Comprehensive USP Raman Library	
	Metrohm Comprehensive USP Raman Library for MIRA P containing ~ 900 traceable certified items.	
6.07501.010	ASTM Calibration Standard	
	50:50 (v:v) Toluene/Acetonitrile ASTM standard for calibration of the wavenumber axis of MIRA/MISA systems. Class 3B operation.	

6.07502.000 Vial Holder

Vial holder attachment for MIRA systems. Accommodates 15 x 26 mm glass vials.



Tablet Holder

Tablet holder attachment for MIRA systems that is used for analysis of formulated tablets of various sizes and shapes.



6.07505.000 LWD Attachment Lens

Long distance point-and-shoot Adapter (LWD) for MIRA systems with a focal length of 7.6 mm. Class 3B operation.



6.07505.010 SWD Attachment Lens

Short distance point-and-shoot Adapter (SWD) for MIRA systems with a focal length of 1.00 mm. Class 3B operation.



Extra-long distance point-and-shoot Adapter (XLWD) for MIRA systems with a focal length of 18 mm. Class 3B operation.



6.07506.000 Right-Angle Attachment

The right-angle attachment is designed to allow the user to collect data by placing the substance on a surface and then laying the MIRA down next to it with the right angle tip covering the substance. Ideal for a baggy on hood of a squad car. Eliminates many of the variables during testing. Class 3B operation.



6.07506.010 Universal attachment

The Universal attachment is designed to be a universal tip for use in three different positions. Position 3 is used for direct contact. Position 2 is used for thin plastic bags, for which the focal point is approximately 1.0 mm from the end of the adapter. Position 1 is used for focusing through bottles, for which the focal point is approximately 8 mm from the end of the metal tip. Class 3B operation.



6.07506.030 Contact Ball Probe Attachment

The Contact Ball Probe Attachment is designed to allow the user to collect data from a sample with no concern of proper focus. Simply contact the sample or immerse the probe to acquire the data. The stainless steel construction allows for easy cleaning. Class 3B operation.



6.07506.050 Protective Sheath for Contact Ball Probe

The disposable protective sheath is made of extremely thin LDPE and is used to prevent cross-contamination when using the ball probe in several containers. Quantity: 250 pieces.




6.07506.060 Intelligent Universal Attachment (iUA)

iUA is an intelligent 3-position collection optics accessory for MIRA systems. The iUA combines the flexibility of a universal attachment with the intelligent mode of operation of MIRA SmartTip Attachments. Each position indicates the ideal purpose of that position, i.e. surface, bag, or bottle.



6.07507.000 Mira Sampling Blackout Cloth

Optical Blackout cloth for sampling in sunlight conditions.

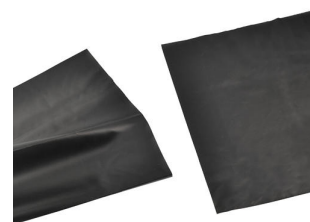

6.07507.100 MIRA DS Field Carry Pouch

The MIRA DS Field Carry Pouch is a utility pouch with padded sides for carrying the MIRA into the field. It has space for MIRA, four AA batteries, the right-angle attachment, iUA, calibration standard, and a few ID Kits.

The Field Carry Pouch can be attached to any of the current military pouch /backpack systems.



6.2133.020 Lithium battery AA 1.5V/ 2900 mAh

1.5V lithium AA batteries, 4 pieces.



6.7430.060 Micro Fiber Cleaning Cloth

Cleaning cloth for monitors and optics.



6.7502.010 Glass Vials for MIRA, 144 pieces


Glass vials for use with MIRA systems. 15 mm x 26 mm.

144 glass vials are included in the scope of delivery.




6.7503.100 MIRA Protective Boot

MIRA Protective Boot provides an extra layer of safety to MIRA devices. Made with a chemically resistant TPA material and specially designed to work with all SmartTip Attachments, the protective boot easily installs for a ruggedized feel and extra protection.

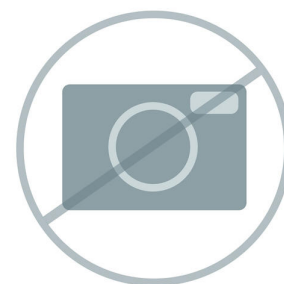

PowerPack Protective Boot provides an extra layer of safety to PowerPack. Made with a chemically resistant TPA material and specially designed to work with Mira, the protective boot easily installs for a ruggedized feel and extra protection.

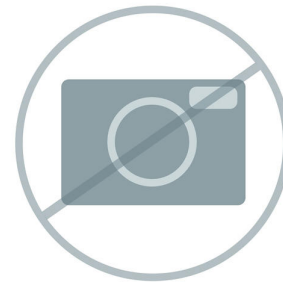

Proper laser safety precautions should be followed for the use of the MIRA products. Protection pursuant to EN 207 is offered with these protective glasses for the 785 nm laser excitation wavelength (also suitable for those who already wear eyeglasses). Work for long periods without difficulties, thanks to the lightweight construction and comfortable design of the protective laser glasses.



MIRA PowerPack Kit contains the MIRA PowerPack rechargeable external battery, a charger with USB-C socket, and a charger cable with USB-C plug. Designed for MIRA systems, the kit increases the period of usable charge to more than 8 hours.



MIRA (Metrohm Instant Raman Analyzer) library of illicit materials including drugs, TICs, TIMs, precursors, and general chemicals.



66071700

MIRA Comprehensive Raman Library

MIRA Metrohm Comprehensive RamanLibrary- More than 29,000 spectral entries of chemicals, minerals, organic chemicals, solvents, polymers, dyes and pigments, TICs, TIMs, inorganics and organometallics, personal care products, and illicit materials..



66071790

MIRA MCRL Organics - Sulfur Group


Metrohm Comprehensive Raman Library (MCRL) Organics - Sulfur Group for MIRA.



66071720

MIRA MCRL - Inorganics and Organometallics

Metrohm Comprehensive Raman Library (MCRL) Inorganics and Organometallics for MIRA.

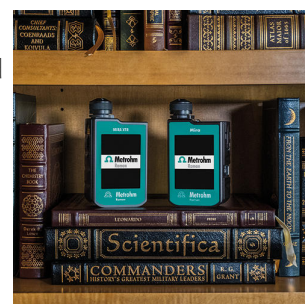

66071740

MIRA MCRL - Personal Care Products


Metrohm Comprehensive Raman Library (MCRL) Personal Care Products for MIRA.



66071710

MIRA MCRL - Polymers and Polymer Additives

Metrohm Comprehensive Raman Library (MCRL) Polymers and Polymer Additives for MIRA.

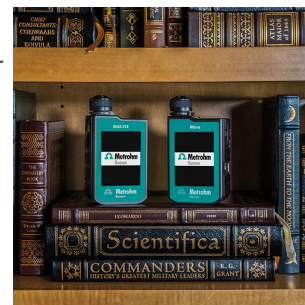

66071730

MIRA MCRL - Dyes, Pigments, and Stains

Metrohm Comprehensive Raman Library (MCRL) Dyes, Pigments, and Stains for MIRA.



66071810

MIRA MCRL - Pesticides

Metrohm Comprehensive Raman Library (MCRL) Pesticides for MIRA.



66071820

MIRA MCRL - Biochemicals


Metrohm Comprehensive Raman Library (MCRL) Biochemicals for MIRA.



66071830

MIRA MCRL - API

Metrohm Comprehensive Raman Library (MCRL) Active Pharmaceutical Ingredients (API) for MIRA.

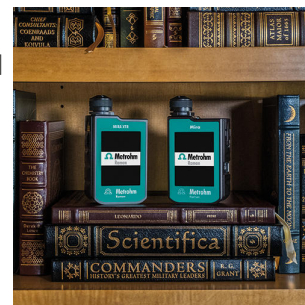

66071750

MIRA MCRL - Organics


Metrohm Comprehensive Raman Library (MCRL) Organics for MIRA.



66071760

MIRA MCRL Organics - Halogenated

Metrohm Comprehensive Raman Library (MCRL) Organics - Halogenated for MIRA.

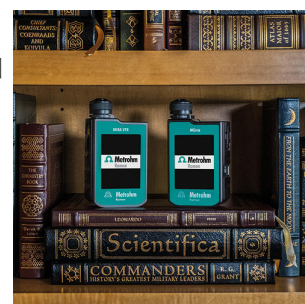

66071770

MIRA MCRL Organics - Nitrogen Group

Metrohm Comprehensive Raman Library (MCRL) Organics - Nitrogen Group for MIRA.



Metrohm Comprehensive Raman Library (MCRL) Organics - Oxygen Group for MIRA.

