



dPt ring electrode

6.00403.300

Digital, combined platinum ring electrode for OMNIS with fixed ground-joint diaphragm. This electrode is suitable for redox titrations when the pH value varies, e.g.:

- Oxygen content according to Winkler
- Determination of hydrogen peroxide with KMnO_4
- Diazotization titrations

The fixed ground-joint diaphragm is insensitive to contamination. $c(\text{KCl}) = 3 \text{ mol/L}$ is used as reference electrolyte and for storage.

dTrodes can be used on OMNIS Titrators.

Phạm vi giao hàng 6.00403.300

Qt.	Order no.	Description
-----	-----------	-------------

100 Stk

6.1236.050

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm, polyethylene.



100 Stk

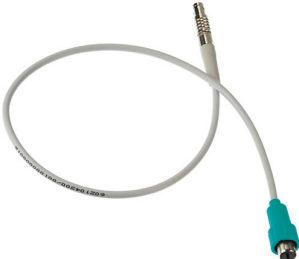
6.2008.040

Storage vessel

Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.



Phụ kiện tùy chọn

Order no.	Description	
6.2308.020	Electrolyte 3 mol/L KCl (250 mL) Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$, (for Ag/AgCl reference systems)	
6.02104.300	Electrode cable plug-in head Q / plug P, 0.55 m Cable for connecting a dTrobe (plug-in head Q) to the digital measuring module from OMNIS.	
6.02104.310	Electrode cable plug-in head Q / plug P, 1.5 m Cable for connecting a dTrobe (plug-in head Q) to the digital measuring module from OMNIS.	