



Ag Titrode with Ag₂S coating

6.00430.100S

Combined silver ring electrode with a pH glass membrane as reference electrode.

The silver ring is coated with sulfide (Ag₂S) for higher sensitivity and a better limit of detection.

This maintenance-free electrode is suitable for precipitation titrations when the pH value remains constant (titrant silver nitrate), for example of:

- chloride, bromide, iodide
- sulfides
- hydrogen sulfide
- mercaptans
- cyanides

This electrode is stored in distilled water.

Phạm vi giao hàng 6.00430.100S

Qt.	Order no.	Description
-----	-----------	-------------

100 Stk

6.1236.050

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm, polyethylene.



100 Stk

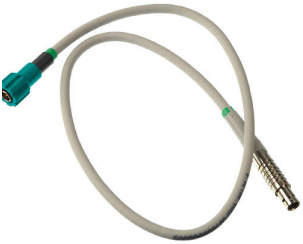
6.2008.040

Storage vessel

Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.



Phụ kiện tùy chọn

Order no.	Description	
6.2104.020	Electrode cable / 1 m / F For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).	
6.02104.000	Electrode cable plug-in head G / plug P, 0.55 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.02104.010	Electrode cable plug-in head G / plug P, 1.5 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	