


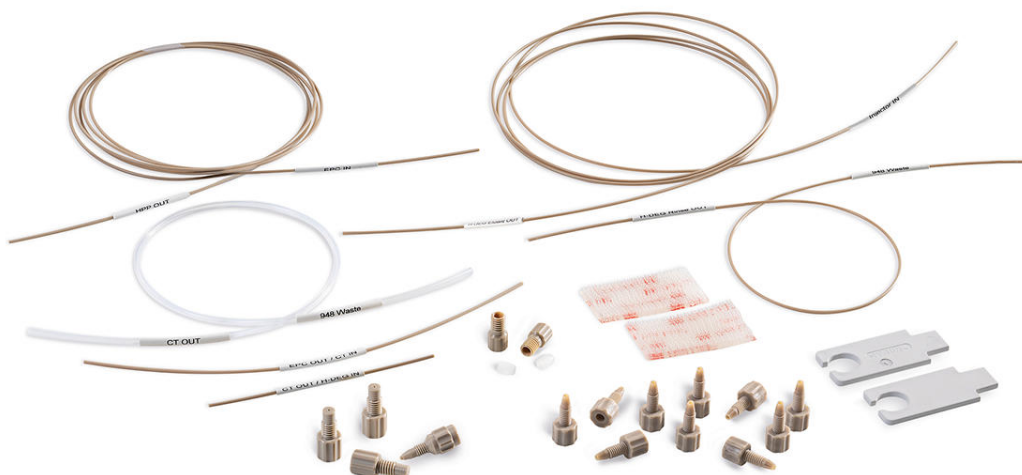

Accessory kit: CEP

6.05000.300

Accessory kit for installation of the 948 Continuous IC Module, CEP

Phạm vi giao hàng 6.05000.300

Qt.	Order no.	Description
1 PCS	6.01805.470	FEP tubing M6, 1/8" x 1.59 mm, 500 mm, CT OUT

An FEP tubing that connects the Continuous Anion Trap outlet with the waste collector of the ion chromatograph.




1 PCS

6.01831.010

PEEK capillary, 0.25 mm, HPP OUT

A PEEK capillary that connects the high-pressure pump of the ion chromatograph with the Eluent Producer Cartridge (EPC) of the 948 Continuous IC Module, CEP.

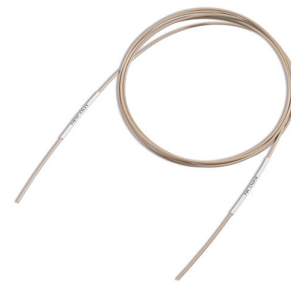

1 PCS

6.01831.020

PEEK capillary, 0.25 mm, EPC OUT

A PEEK capillary that connects the Eluent Producer Cartridge (EPC) with the Continuous Trap (CT) of the 948 Continuous IC Module, CEP.



1 PCS

6.01831.030

PEEK capillary, 0.25 mm, CT OUT

A PEEK capillary that connects the Continuous Trap (CT) with the High-Pressure Degasser (H-DEG) of the 948 Continuous IC Module, CEP.




1 PCS

6.01831.040

PEEK capillary, 0.25 mm, H-DEG OUT

A PEEK capillary that connects the High-Pressure Degasser (H-DEG) of the 948 Continuous IC Module, CEP with the injector of the ion chromatograph.

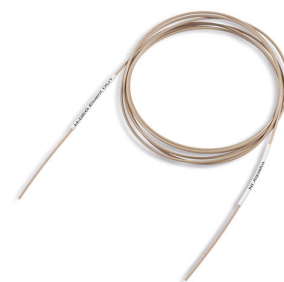

1 PCS

6.01831.050

PEEK capillary, 0.25 mm, H-DEG Rinse

A PEEK capillary that connects the High-Pressure Degasser (H-DEG) of the 948 Continuous IC Module, CEP with the waste collector of the ion chromatograph.



1 PCS

6.2053.030

Mounting ribbon for 6.2744.340, 25x50 mm

Mounting ribbon for Dosino gradient IC equipment



2 PCS

6.2744.070

Pressure screw short

Short version. With UNF 10/32 connection. 5 pieces. For the connection of PEEK capillaries



2 PCS

6.2744.080

M6 thread / UNF 10/32 adapter

With M6 outer thread. For IC instruments.



1 PCS

6.2744.200

Adapter UNF 10/32 outer / M6 inner

Connecting 1/16 in. PEEK capillaries and transfer tubing (M6)

