



797 VA Computrace for trace analysis

MVA-01

Analysis system for voltammetric trace analysis and education. For manual control. The system is suitable for quantitative analysis in routine and research laboratories. Comprised of 797 VA Computrace. Without PC.

Phụ kiện tùy chọn

Order no.	Description
2.800.0010	800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



2.843.0040

843 VA Membrane Pump Station for Compact VA systems

The 843 Membrane Pump Station has two built-in membrane pumps that can be controlled directly by the 797 VA Computrace via the remote signal. It is used for the automatic rinsing and emptying of the measuring vessel in partly automated and fully automated VA and CVS systems with the 797 VA Computrace.

The version for voltammetry is equipped with the complete accessories for the automatic rinsing and emptying of the measuring vessel.



2.843.0140

843 VA Peristaltic Pump Station for Compact VA systems

The 843 Peristaltic Pump Station has two built-in peristaltic pumps that can be controlled directly by the 797 VA Computrace via the remote signal. It is used for the automatic rinsing and emptying of the measuring vessel in partly automated and fully automated VA and CVS systems with the 797 VA Computrace.

The version for voltammetry is equipped with the complete accessories for the automatic rinsing and emptying of the measuring vessel.



2.846.0010

846 Dosing Interface


USB-capable control unit for connecting a maximum of four 800 Dosinos or 805 Dosimats for dosing and Liquid Handling tasks. A Touch Control or the connection to a PC with OMNIS Software, tiamo™, MagIC Net, viva, or 797 VA Computrace is required for operation.



6.2055.100

Bottle holder, double

Bottle holder for 2 x 1 L bottles. For 6.5617.000 aspiration equipment

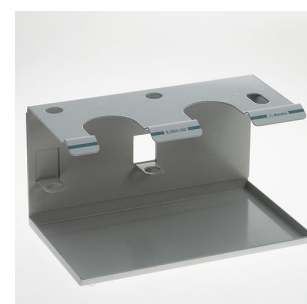

807 Dosing Unit with integrated data chip with 2 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.

